

## Single Drain Board DrainTops

DOMESTIC U.S. SIZES

Style No.	Exterior Dimensions (in)	
	Length	Width
454 (sink right)	60	26 - 31.9
455 (sink left)	60	26 - 31.9

## Dual Drain Board DrainTops

DOMESTIC U.S. SIZES

Style No.	Exterior Dimensions (in)	
	Length	Width
451	54	26 - 30
450	60	26 - 30
457	66	27.5 - 32
456	72	27.5 - 32

### Classic DrainTops


DrainTops have been specially designed for use with Classic and Classic MarineTops.

DrainTop worksurfaces gently slope from the containment rim to the sink area for improved drainage. The molded drain board sections offer curved raised ribs that facilitate glassware drying and improved cleanups. The sink area accommodates single or double Dropln Sinks.


Flexibly designed, DrainTops can be used at the end of islands and peninsulas or in wall runs as 54", 60", 66" or 72" [1372, 1524, 1676 or 1829 mm] sections.

### Drain Grooves


Drain grooves can be machined on one or both sides of undermount or Dropln sinks per specification. These recessed slots incline toward the sink basin, facilitating worksurface drainage while providing a temporary labware drying area.

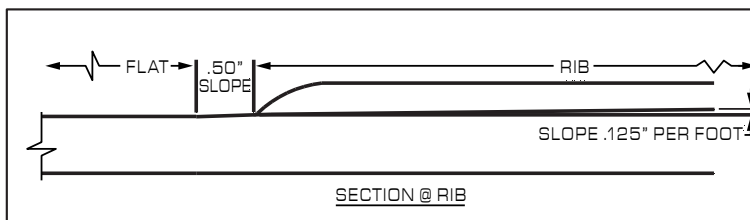


451 (18" x 30" Available Sink Area\*)

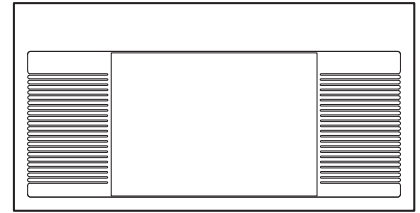

450 (18" x 30" Available Sink Area\*)

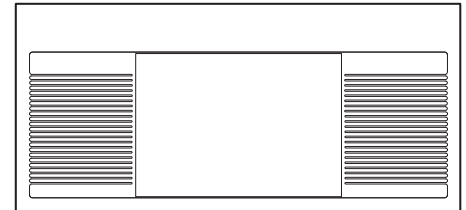


454 (20" x 38" Available Sink Area\*)

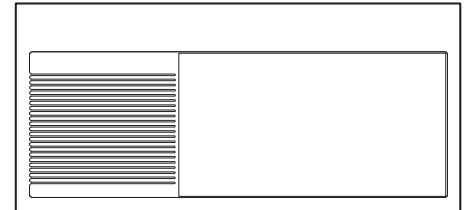


455 (20" x 38" Available Sink Area\*)



457 (20" x 38" Available Sink Area\*)

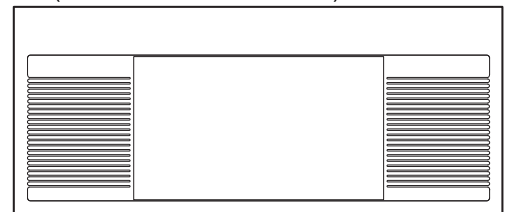

456 (20" x 38" Available Sink Area\*)

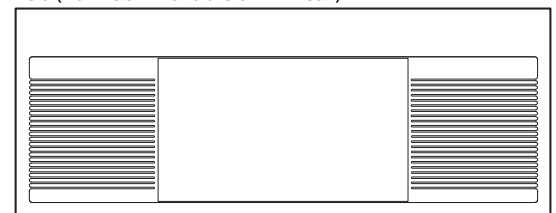

\*Sinks to be Measured from the Outer Sink Rim (OD) Plus Cutout Bevel