




## TRIFACIAL SERVICE CENTER

Save space, time and money without sacrificing versatility. The **Trifacial Service Center** provides an ample worksurface and central sink allowing the students to easily access the services. When the Trifacial is combined with Sheldon Student Tables, the classroom becomes a lecture and lab environment which reduces classroom disruption and provides more circulation area for the students.

### Features

- Three sides of triangular configuration
- Large, 34-gallon capacity sink for easy cleanup
- Center height adjustable at time of installation for elementary, middle school, and high school applications
- Equipped with three vandal-resistant unicast cold water and gas fixtures and three GFI outlets
- Support structure of fire-retardant fiberglass reinforced polyester

# TRIFACIAL SINK CENTER

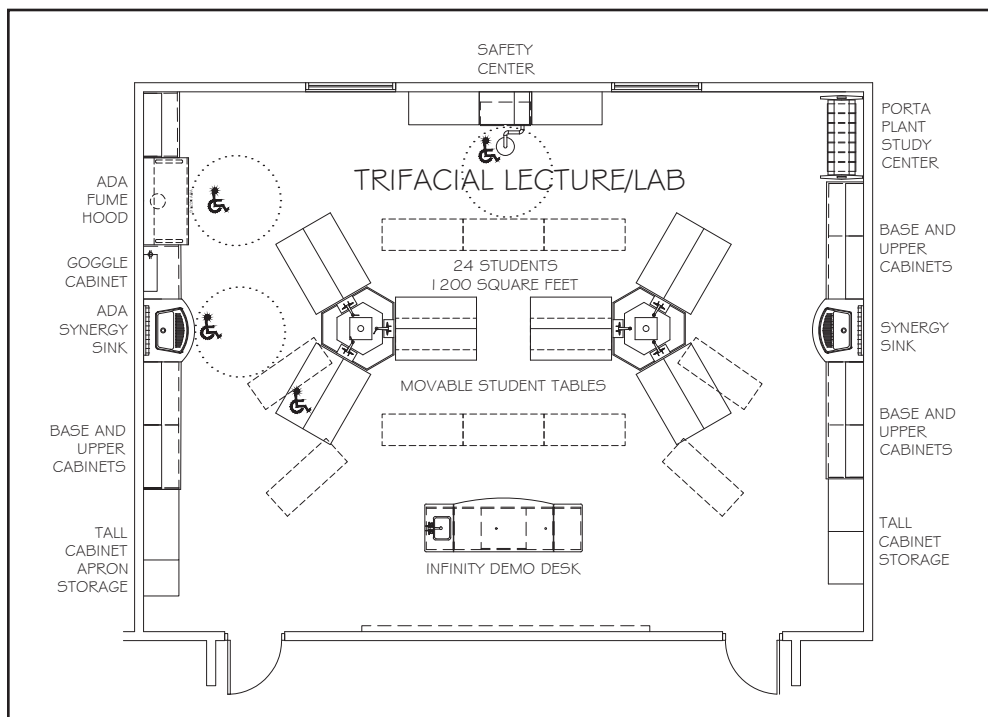

## Trifacial Lab

This lab is designed as lecture and lab layout. This allows the students to collaborate with each other while creating a sightline for the instructor to all locations. The student tables can also be arranged in various lecture layouts (see dotted lines).

## Worksurface

Molded fiberglass reinforced polyester top is standard. Epoxy resin top is optional. Both surfaces include integral raised perimeter edge to contain spillage.

## Sink Basin

One-piece molded fiberglass polyester sink is standard. Drop-in style (epoxy resin) is optional.

## Services

Unicast fixtures (cold water & double gas) are standard. Unimix fixtures (cold water, hot water and gas) are optional. Internal piping, wiring, conduit and drain lines are not included.

## Connectivity

Three (3) GFI Duplex receptacles

## Support Understructure

Molded one-piece fiberglass reinforced polyester understructure.

## Steel Pedestal Base

Height adjustable at time of installation. Includes removable steel panels for access to services. Standard color is Black.


## Accessibility

Accessibility worksurface available at ADA height (max. 34" high). Includes unicast fixtures with ADA wristblade handles.

## Color Options-Epoxy

Available in Black Onyx, Lunar White, Gray, Tan, and Dark Khaki.

## Color Options-Fiberglass\*

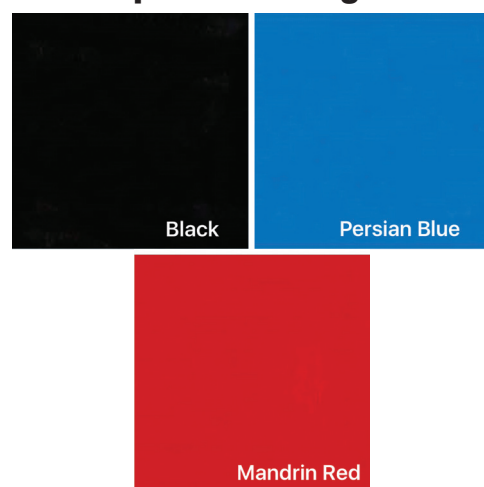

**\*If color match is critical, samples are available upon request.**