


# The Eliminator HP Series


The Eliminator HP hood comes standard with all the same features as other Eliminator models, along with some additional innovations to set this hood apart from its competitors.

## THE ELIMINATOR HP SERIES

**High Performance/Energy Efficient - A Green Alternative**

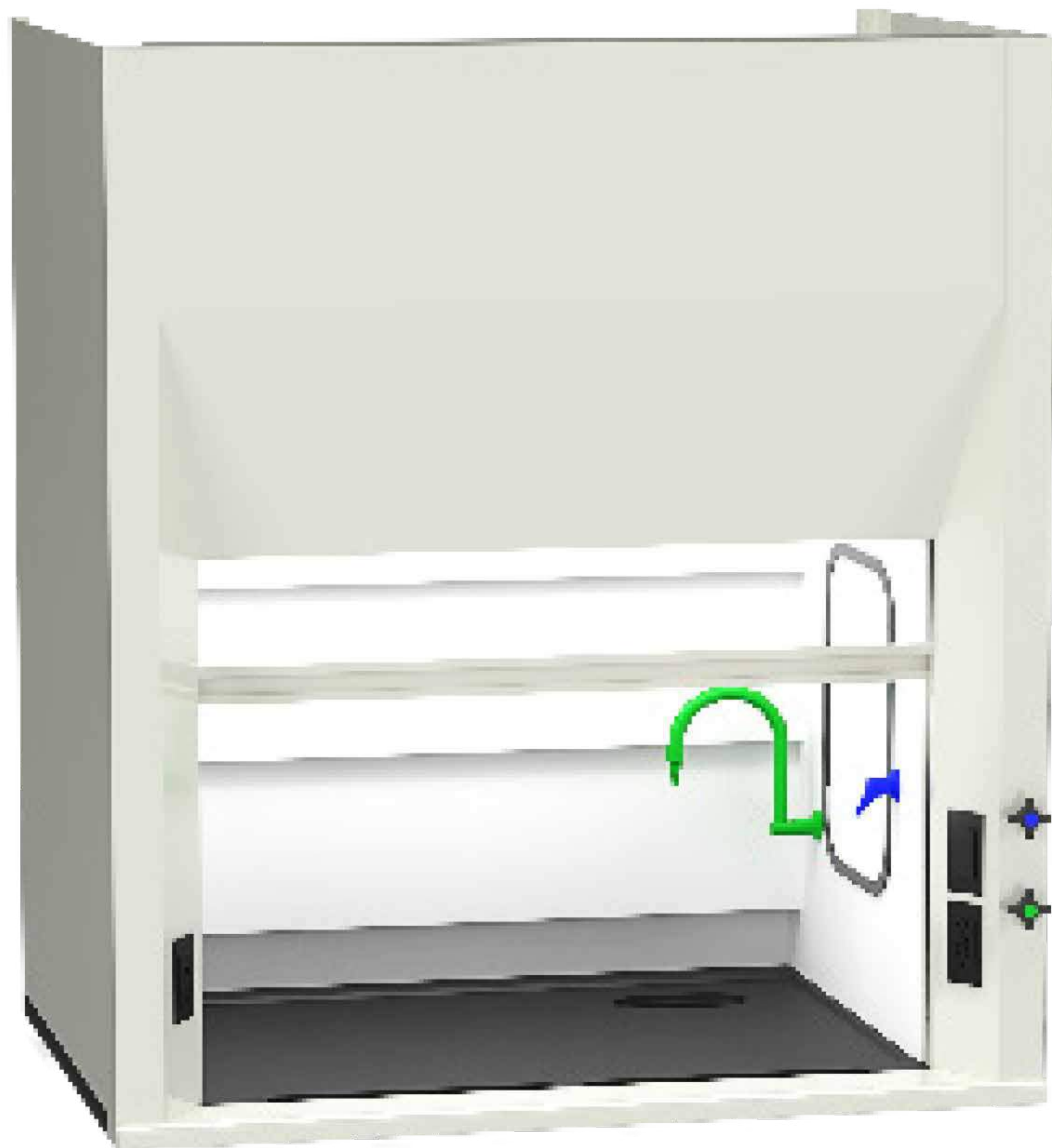

- Powder coat finish
- Full frame construction
- Shaft drive sash
- U.L. listed poly resin liner
- Flush-mount airfoil
- Knock down capability
- LED light with bulbs
- Full-length finger lift
- Sash sweep
- 25.5" interior working depth

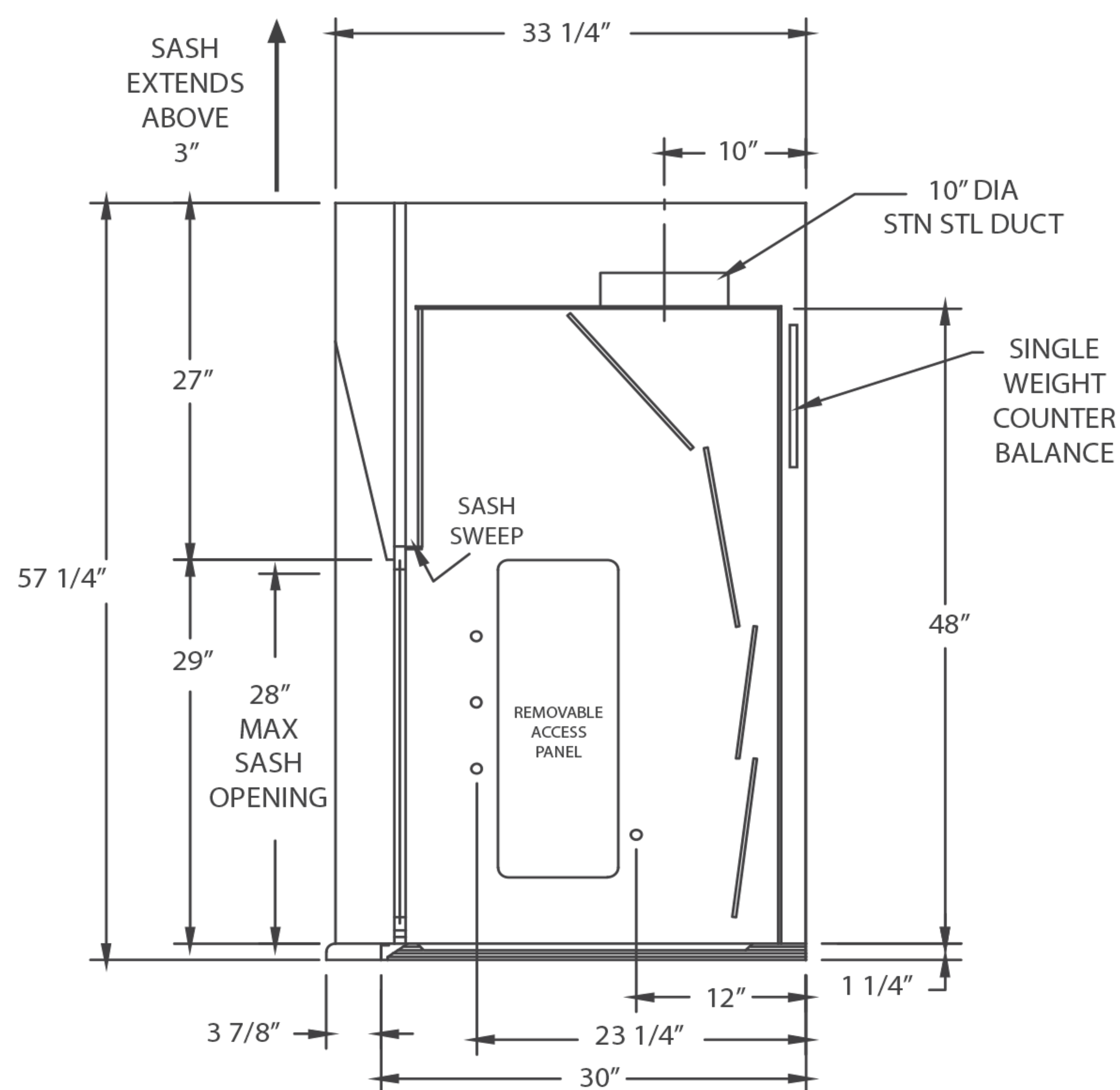

U.L. 1805  
Classified

- 60 FPM
- Fits through standard door
- 30" Countertop

The HP Eliminator is THE hood for remodeling, new construction, or retrofitting.

Part Number	Superstructure			Worktop			Hood Opening		CFM @ 28"*	CFM @ 18"*	SP	Blower Model #	Shipping Weight
	W	H	D	W	H	D	W	H					
EH-HP-48	48"	56"	33-1/4"	48"	1-1/4"	30"	38"	28"	469	313	.25	BL-HP-48	415
EH-HP-60	60"	56"	33-1/4"	60"	1-1/4"	30"	50"	28"	619	413	.34	BL-HP-60	490
EH-HP-72	72"	56"	33-1/4"	72"	1-1/4"	30"	62"	28"	769	513	.34	BL-HP-72	535
EH-HP-96	96"	56"	33-1/4"	96"	1-1/4"	30"	86"	28"	1069	713	.34	BL-HP-96	570

\* Height of open sash; CFM calculated at 60 FPM


For High Performance VAV Hoods, add -VAV to the part number

# The Eliminator HP Series

The HP model is designed to operate economically while creating a safe work environment without compromising containment. By operating at 60 FPM, air flow across the sash opening results in a 40% savings in fume hood operation costs. Features of the HP Series include:

- Passed ASHRAE 110 AM testing with sash in full open position.
- Standard 30" deep work surface for zero work space loss in your lab. Other depths available.

## THE ELIMINATOR HP SERIES

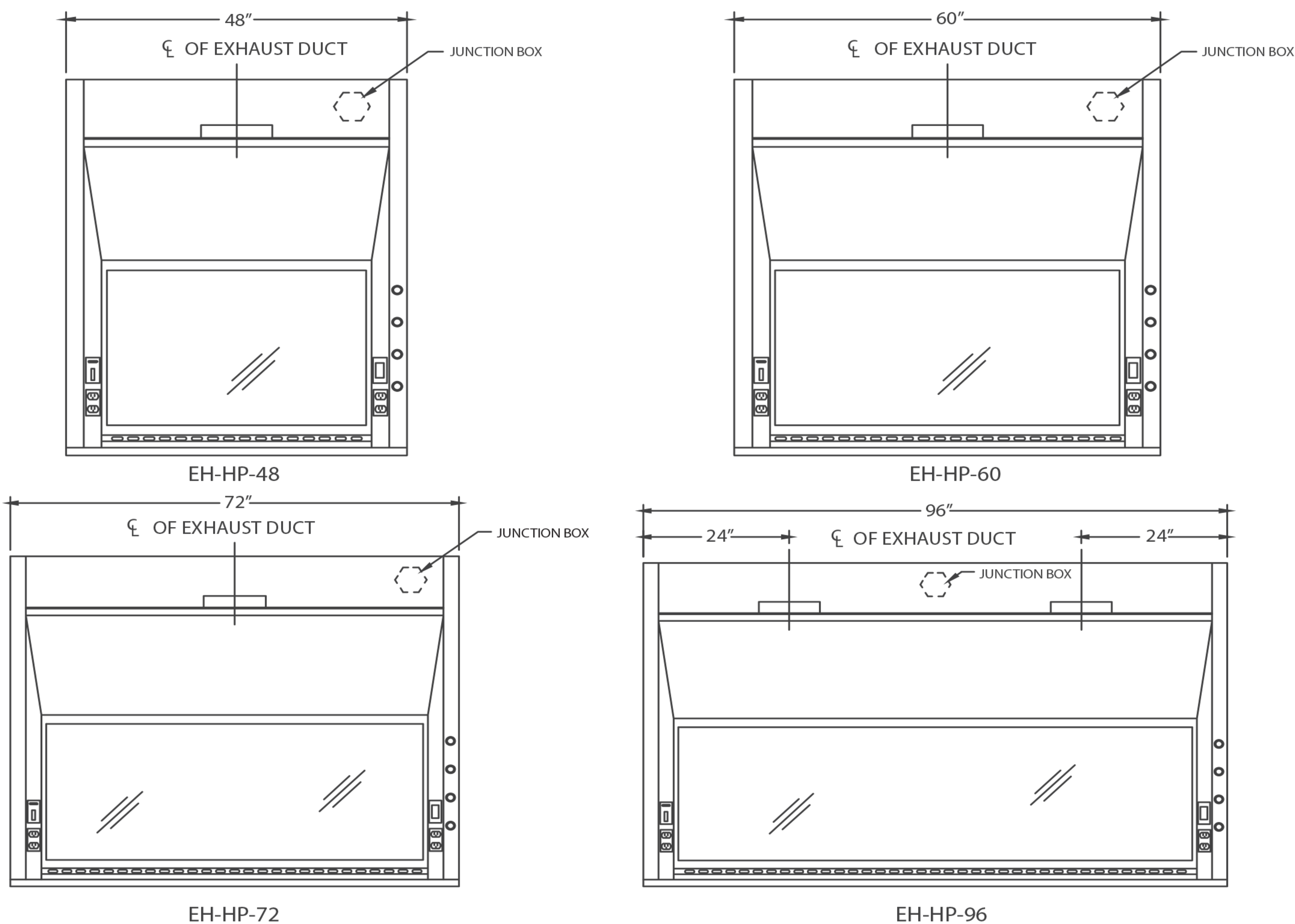

## ACCESSORIES

ITEM	AMS PT#	PAGE #	ITEM	AMS PT#	PAGE #
Ceiling Enclosures	CE	33	Sash Stops	SASHSTOP	38
Plumbing Accessories		34	Chain and Sprocket Drive	CHNSPK	CALL
Electrical Accessories		35	Stainless Steel Sash	SS SASH	38
Epoxy Tops	ET	39	Stainless Steel Airfoil	SS AF	
Alarms	ALARM	32	Push Button Sash	PUSH BUTTON	38
Base Cabinets	BC	26	Combination Sash	H&V	38
Blowers	BL	CALL	Lattice Assembly	LATTICE	CALL
Pre-Plumbing	PREP		Auto Sash Return	AUTO RETURN	38
Pre-Wiring	PREW				