

RESEARCH COLLECTION

Steel Laboratory Furniture



KEWAUNEE®

Research Collection

Steel Furniture

- Fully-welded, heavy-gauge steel construction
- Constructed to meet SEFA standards
- Powder coat, lab grade finish
- Stainless steel construction option
- Customize with multiple styles, pulls, and hardware





Contents

Steel Furniture Styles & Construction	4
Catalog Number Explanation	6
Steel Cabinet Options.....	8
Sitting Height Cabinets	13
ADA Height Cabinets.....	16
Standing Height Cabinets.....	24
Mobile Cabinets.....	36
Apron Assemblies.....	45
Wall Mounted Storage Cabinets	53
Full Height Storage Cabinets	60
Safety Cabinets.....	65
Fillers, Scribes & Panels.....	71
Sloping Tops & Enclosures	79
Fume Hoods & Ventilated Enclosures.....	81
Worktops.....	82
Center Bench Shelving	84
Wall Shelving	86
Reagent Racks & Shelving	89
Omega Shelving	90
Pegboards.....	92
Miscellaneous Accessories	94
General Plumbing Information.....	97
Service Fittings.....	98
Drain Accessories	106
Sinks & Troughs.....	107
Electrical Fittings	110
Technical Information	112
Typical Installations	118

Beautifully finished, smooth, clean, easy to maintain surfaces. Rugged, strong, heavy gauge steel, full-frame, four sided, flush front construction. That's Kewaunee's Steel Laboratory Furniture.

Kewaunee Steel Cabinet Construction

Kewaunee's fully-welded, heavy-gauge, steel cabinets are designed to SEFA 8 performance standards and built to rigorous benchmarks. The Kewaunee classic "unit" design allows doors and drawers to be reconfigured in place. The configurable design provides a rich offering of hardware and aesthetic options.




STANDARD DRAWER



Standard Drawers feature one-piece drawer body construction with full-extension, side-mounted, soft-close slides.

PULL STYLE OPTIONS

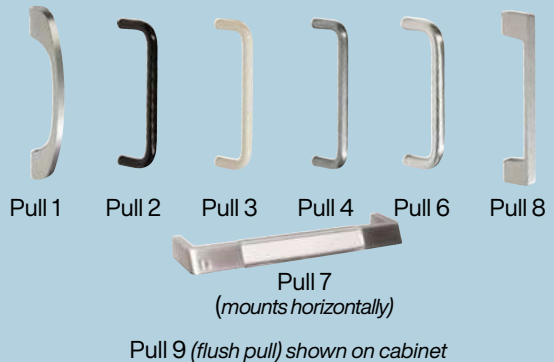


DOOR CATCH



Positive Catch

Catch furnished with 5-knuckle hinges only

DOOR HINGE OPTIONS

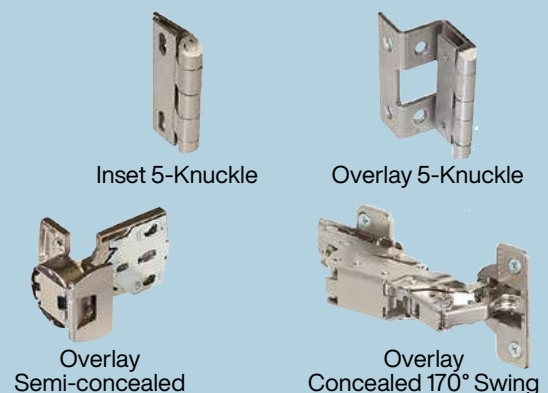

Additional Cabinet Features

- Four sided, flush front, frame construction
- Solid finished back with optional removable back
- Full-depth adjustable shelf
- Leveling bolts accessible through cabinet top

Kewaunee Steel Cabinet Styles



shown with five knuckle hinge and rounded aluminum pull style 1

Inset - Square Edge

Doors and drawer heads are flush with face of cabinet. Available with Positive Catch and a wide selection of pulls.



shown with five knuckle hinge and chrome plated wire pull style 6

Overlay - Square Edge

Doors and drawer heads overlay face of cabinet. Available with 5-knuckle, semi-concealed or concealed hinges and a wide selection of pulls.



shown with maple veneer door and drawer head, five knuckle hinge and square aluminum pull style 8

Overlay - Wood on Steel

Oak or maple veneer doors and drawer heads overlay face of cabinet. Available with 5-knuckle, semi-concealed or concealed hinges and a wide selection of pulls. *(Flush Pull 9 not available on wood)*



shown with concealed hinge and beveled label holder pull style 7

Overlay - Bevel Edge

Doors and drawer heads, with beveled top edge, overlay face of cabinet. Available with 5-knuckle, semi-concealed or concealed hinges and a wide selection of pulls.

Steel Cabinet Information

Catalog Number Explanation

E40C352236D-1611B11C...

Cabinet Type

A = Apron Assemblies and Legs
B = Open Cupboard Base Cabinet
C = All Cupboard Base Cabinet
D = All Drawer Base Cabinet
E = Drawer & Cupboard Base Cabinet
G = Special Purpose Base Cabinet
S = Full Height Storage Cabinet
W = Wall Storage Cabinet

Door & Drawer Configuration

Cabinet Construction & Design

Cabinet Height

Base Cabinets: 29" – 32" – 35"
Aprons: 04" Full Heights: 84"
Wall Cabinets: 24" – 30" – 36" – 48"
Mobile Cabinets: 26" – 29" – 31"
Suspended Cabinet: 19" – 25"

Cabinet Options

(See Cabinet Option list next page)

Style – Hardware – Material

(See Selection list below)

Special Configurations

Cabinet Length

Standard Lengths: 12" – 15" – 18" – 21"
24" – 30" – 36" – 42"
48" – 52" – 58" – 60"

Cabinet Depth

13" – 16" – 20" – 22" – 28"

actual dimensions can be found in the Specifications at the top of each section

Door & Drawer Head Style – Hardware – Material Selections

E40C352236D-1611B11C...

Cabinet Style - Hinge Type - Catch

0 = Inset - 5-knuckle Hinge - Positive Catch
1 = Full Overlay - 5-knuckle Hinge - Positive Catch
4 = Full Overlay - Semi-concealed Hinge
6 = Full Overlay - Concealed Hinge (170° Swing)

Door & Drawer Head Style

1 = Square Edge Steel
6 = Wood – available on Full Overlay Cabinet Styles only
7 = Beveled Edge Steel – available on Full Overlay Base Cabinet Styles only

Cabinet Material / D&D Material

0 = Cold Rolled Steel / Cold Rolled Steel
1 = Cold Rolled Steel / Red Oak
2 = Cold Rolled Steel / White Maple
7 = Cold Rolled Steel / Steamed Beech
*3 = Stainless Steel / Stainless Steel
5 = Cold Rolled Steel / FSC® Red Oak
6 = Cold Rolled Steel / FSC® White Maple
8 = Cold Rolled Steel / FSC® Steamed Beech

Pull Type


0 = No Option
1 = Rounded Aluminum
2 = Black Wire
3 = White Wire
4 = Stainless Steel Wire
6 = Chrome Plated Wire
‡7 = Beveled Label Holder
8 = Square Aluminum
†9 = Flush Aluminum

‡ mounts horizontally
† not available on wood doors & drawer heads

* Stainless Steel Cabinets available in D&D Style 1 only

Steel Cabinet Information

Drawer and Cupboard Size Coding

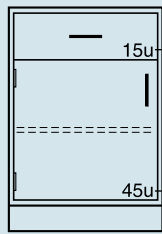

Drawer Code

Cupboard Code

E40C352215L- _____ 15"L
E40C352218L- _____ 18"L
E40C352221L- _____ 21"L
E40C352224L- _____ 24"L

Base cabinets are designed on a "unit" system to allow maximum flexibility in the selection of drawer and/or cupboard door sizes.

Standing height cabinets (35" high) have an opening 60 "units" high. Sitting height cabinets (29" high) have an opening 48 "units" high. ADA height cabinets (32" high) have an opening 54 "units" high. This means that any combination of drawer and cupboard door heights equaling 60, 54 or 48 units can be used in a base cabinet.

Mobile standing height cabinets (31 3/4" high) have an opening 53 "units" high. Mobile ADA height cabinets (29 1/4" high) have an opening 48 "units" high. Mobile sitting height cabinets (25 3/8" high) have an opening 40 "units" high.

The adjacent tables indicate the "unit" height and inside clearance of the Inset style base cabinets shown in this catalog. Simply match the unit number shown on the cabinet drawers and cupboard to the appropriate number in the tables.

Drawer Size Codes

"Units" High	Drawer Head Height	Inside Dimensions	
		Clear/Height	Length (Front-to-Back)
8	3 3/4"	2 3/8"	18"
9	4 1/4"	2 3/8"	18"
10	4 3/4"	3 7/16"	18"
11	5 3/16"	3 15/16"	18"
12	5 11/16"	4 7/16"	18"
14	6 5/8"	5 3/8"	18"
15	7 1/8"	5 7/8"	18"
16	7 5/8"	6 3/8"	18"
20	9 9/16"	8 5/16"	18"
24	11 1/2"	10 1/4"	18"
3	1 3/16"	† Pull Out Shelf	

† Pull Out Shelf mounted on full extension slides only

Door Opening Clearance Width I.D. (side-to-side)

Base Cabinet Length	Number of Drawers Across		
	1	2	3
12"	9 1/2"		
15"	12 1/2"		
18"	15 1/2"		
21"	18 1/2"		
24"	21 1/2"	9 1/2"	
30"	27 1/2"	12 1/2"	
36"	33 1/2"	15 1/2"	9 1/2"
42"	39 1/2"	18 1/2"	
48"	45 1/2"	21 1/2"	13 1/2"

Drawer Width I.D. (side-to-side)

Base Cabinet Length	Number of Drawers Across		
	1	2	3
12"	8 3/8"		
15"	11 3/8"		
18"	14 3/8"		
21"	17 3/8"		
24"	20 3/8"	8 3/8"	
30"	26 3/8"	11 3/8"	
36"	32 3/8"	14 3/8"	8 3/8"
42"	38 3/8"	17 3/8"	
48"	44 3/8"	20 3/8"	12 3/8"

Cupboard Size Codes

"Units" High	Door Height	Inside Dimensions	
		Height (Top-to-Bottom)	Depth (Front-to-Back)
30	14 9/16"	13 9/16"	20 1/2"
36	17 7/16"	16 7/16"	20 1/2"
45	21 13/16"	20 13/16"	20 1/2"
48	23 3/16"	22 3/16"	20 1/2"
54	26 1/16"	25 1/16"	20 1/2"
60	29"	28"	20 1/2"

Cupboard and Drawer sizes listed in these charts are for Inset style cabinets only. Sizes for Overlay styles are generally slightly larger. Contact your Kewaunee Dealer for more information.

Cabinet Options – (Add to End of Catalog Number)

B – Base Cabinet Caster Option

- B11 – 4" Swivel Caster w/ brakes & with counter balance
- B31 – 4" Swivel Caster w/ brakes at front & counter balance
- B41 – 3" Swivel Caster w/ brakes & with counter balance
- B61 – 3" Swivel Caster w/ brakes at front & counter balance

C – Toe space depth increased by 7/8" (for B00 & G07)

F – Removable Back (not available on Mobile Cabinets)

* H – Label Holders

I – 1/4" Laminated Safety Glass

* J – Disk Locks w/Lock Panels * J3 – Pin Locks w/Lock Panels

* K – Disk Locks w/o Lock Panels * K3 – Pin Locks w/o Lock Panels

O – Shelf Lips

P – Pullout Trays & Pullout Pans

- P1 – One Pullout Tray - at bottom of cupboard
- P2 – Two Pullout Trays - center & bottom
- P3 – One Pullout Pan - at bottom of cupboard
- P4 – Two Pullout Pans - center & bottom

R – 1/8" Tempered Safety Glass

U – Toe Space on Base Cabinet Ends

- U1 – Toe Space on Left End
- U2 – Toe Space on Right End
- U3 – Toe Space on Both Ends

V – Drawer Slides

- V1 – 200 lb Interlocking Full Extension
- V2 – 200 lb Full Extension
- V4 – 200 lb Soft Close Full Extension

* Furnished on all cupboards and drawers on cabinet. Lock keying and/or numbering must be specified.

Steel Cabinet Options

Base Cabinet Caster Option

Option B

Base Cabinet Caster Options add four (4) swivel casters to eligible base cabinets. Each option is furnished with counterbalance weight for improved stability and is assembled using the included self-tapping screw hardware. *(Casters ship loose.)*

Swivel Casters (4 per base cabinet)

B11 – *4" Swivel Casters (brakes on all)

B31 – *4" Swivel Casters (brakes at front only)


B41 – 3" Swivel Casters (brakes on all)

B61 – 3" Swivel Casters (brakes at front only)

* Adds 1" to height of cabinet

Load Rating **165 lbs** per caster

Load Rating **150 lbs** per caster

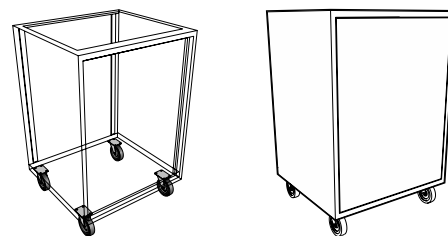

B11 – B31 – B41 – B61 Caster Option

To specify Casters, add **B__** to the end of the cabinet number.

Example: E40C352236D-0101**B31**
(Eligible base cabinets only.)

Note: Caster option *not* available on 12" length cabinets.

Load Rating Note: Load ratings shown are per caster; total capacity depends on load distribution and application.

Toe Space Depth Increase by 7/8"

Option C

Option **C** adds 7/8" to the depth of a base cabinet toe space. This option is only used on base cabinets without doors or drawers (B00 & G07) to allow the face of

the cabinet to be pulled forward to align with the face of the doors and drawers on adjacent cabinets.

To specify increased toe space depth add **C** to the end of the cabinet

number.
Example: B00C352236N-0000**C**

Removable Back


Option F

All cabinets are furnished with a solid back for improved cleanliness. A removable back is available when required for access to the pipe space behind the cabinet.

Note: Removable Backs are not available on Mobile Cabinets.

To specify Removable Back, add **F** to the end of the cabinet number.

Example: D30C352236N-0101**F**

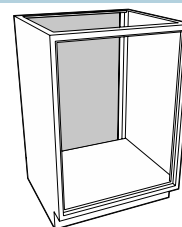

Label Holders

Option H

Option **H** adds the **F-0260-00** aluminum Label Holder to each drawer and door on the cabinet. Overall size is 3" x 1 1/8" and holds a 2 1/2" x 1 1/8" card.

Double doors receive only one Label Holder per pair.

Shipped loose for field installation.

To specify Label Holders, add **H** to the end of the cabinet number.

Example: E40C352236D-0101**H**
(Added to all cupboards and drawers.)



One-half Actual Size

Shelf Lips

Option O

Option **O** adds a 3/4" high steel lip to the front edge of all steel shelves in the cabinet. (**Note:** Shelf lips will not

be provided on the fixed shelf on full-overlay cabinets with a lock option selected.)

To specify Shelf Lips, add **O** to the end of the cabinet number.

Example: W00C361336N-0101**O**
(Added to all shelves in cabinet.)

Steel Cabinet Options

Pullout Trays & Pans

Option P

P1 – One Pullout Tray - at bottom of cupboard

P2 – Two Pullout Trays - center & bottom

P3 – One Pullout Pan - at bottom of cupboard

P4 – Two Pullout Pans - center & bottom

Pullout Trays and Pans are fabricated from 18 gauge cold rolled steel and are mounted on full-extension slides. Both have 2 $\frac{3}{8}$ " high sides and back. Trays are open at the front while pans are enclosed on all four sides. Each is available in two configurations: one tray or pan mounted at bottom of cupboard; or two trays or pans mounted at bottom

and in center of cupboard. When pullout trays or pans are added, the adjustable shelf is omitted.


Pullout Trays and Pans are available for base cabinets only.

To specify a Pullout Tray, add **P1** or **P2** to the end of the cabinet number.


Example: E40C352236D-0101P2
(Available for base cabinets only.)

To specify a Pullout Pan, add **P3** or **P4** to the end of the cabinet number.

Example: E40C352236D-0101P3
(Available for base cabinets only.)

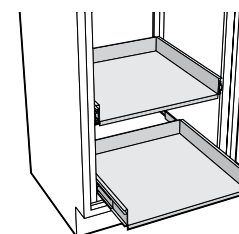

Pullout Tray

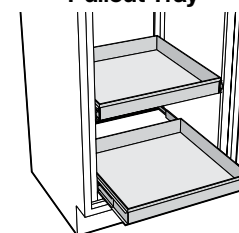

Pullout Pan

Laminated Safety Glass

Option I

Tempered Safety Glass

Option R

Option **I** replaces the standard float glass in plate glass or glazed door wall, counter mounted and full height cabinets with 1/4" Laminated Safety Glass.

Option **R** replaces the standard 1/8" float glass in glazed door wall, counter mounted and full height cabinets with 1/8" Tempered Safety Glass.

To specify Laminated Safety Glass, add **I** to the end of the cabinet number.

Example: W50C361336-0101I

To specify Tempered Glass, add **R** to the end of the cabinet number.

Example: W50C361336-0101R

Toe Space on Cabinet Ends

Option U

When cabinets are used at the end of center benches and in some other special applications, they may require a toe space on one or both ends of the

cabinets. These toe spaces match the toe space at the front of the cabinets and may be added to the left, right, or both ends.

U1 – Toe Space on Left End

U2 – Toe Space on Right End

U3 – Toe Space on Both Ends

To specify End Toe Spaces, add **U_** to the end of the cabinet number.

Example: E40C352236D-0101U2
(Available for base cabinets only.)

Optional Drawer Slides

Option V

Replaces standard 100 lb full-extension slides on all drawers in cabinet when **V** option is specified.


V1 – 200 lb Interlocking Full Extension

V2 – 200 lb Full Extension

V4 – 200 lb Soft Close Full Extension

To specify Drawer Slides, add **V_** to the end of the cabinet number.

Example: E40C352236D-1011V1

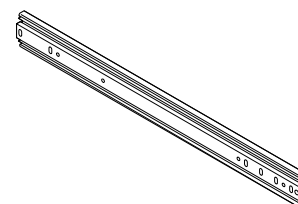

Steel Cabinet Options

Locks with Lock Panels

Locks without Lock Panels

Option J

Option K

All cabinets may be ordered with locks by specifying the K option. In addition, base cabinets may be specified with the J option to add lock panels. Lock panels are installed between drawers

and/or drawers and cupboards to provide additional security when the locks on a cabinet are keyed different.

Locks Types

Cabinet locks are available in two types. Disk Tumbler which offers medium security, and Pin Tumbler which offer a higher level of security. Pin Tumbler also allows more keying arrangements. Both are fabricated with die cast sleeves and cylinders with a satin nickel

face and zinc plated steel cams. Both locks use the same sleeve and cam assembly allowing the lock cylinder to be removed and switched after the lock has been installed. All locks are furnished with two keys unless otherwise specified.

Disk Tumbler Locks

Security: Medium

Keying: Different, Alike, Alike in Sets, Masterkeyed

Key Changes: Up to 2000

Master Key Changes: Up to 5

J – Disk Locks with Lock Panels

K – Disk Locks w/o Lock Panels

Pin Tumbler Locks

Security: Medium-High

Keying: Different, Alike, Alike in Sets, Masterkeyed, Grand-masterkeyed, Sub-masterkeyed, Mason Keyed

Key Changes: Up to 2000

Master Key Changes: Up to 52

J3 – Pin Locks with Lock Panels

K3 – Pin Locks w/o Lock Panels


Lock schedule with keying arrangement must be furnished when locks are specified. See following pages for keying arrangement diagrams; [Click here](#) for lock schedule examples and blank form.

To specify Locks without Lock Panels, add **K** or **K3** to the end of the cabinet number.


Example: E40C352236D-0101**K3**
(Added to all cupboards and drawers.)

To specify Locks with Lock Panels, add **J** or **J3** to the end of the cabinet number.


Example: E40C352236D-0101**J**
(Added to all cupboards and drawers.)

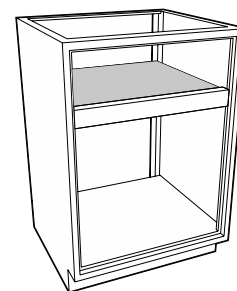

Lock Panel

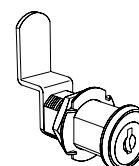

Disk Tumbler Lock
(Medium Security)

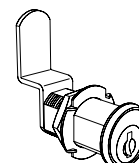

Pin Tumbler Lock
(Higher Security)


Lock Keying Options

Keyed Different


Disk Tumbler Locks: Available
Pin Tumbler Locks: Available

All locks are keyed different from each other with individual keys for each lock. No key will open more than one lock.

Individual Locks

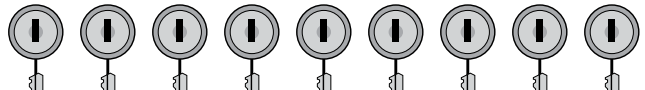

Individual Keys

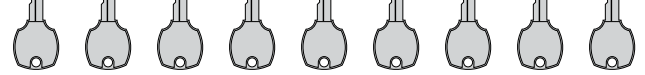

Keyed Alike

Disk Tumbler Locks: Available
Pin Tumbler Locks: Available

All locks are keyed alike. One key will open all locks.

Individual Locks



Single Key





Keyed Alike in Sets

Disk Tumbler Locks: Available
Pin Tumbler Locks: Available

Locks are grouped into sets, with as many locks as needed in each set. All locks in a set can be opened by a single key, but that key can open no other locks. No key will open all locks.

Locks Sets

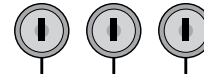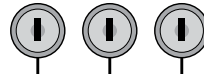

Single Key/Set



Master Keyed




Disk Tumbler Locks: Available
Pin Tumbler Locks: Available

All locks are keyed different from each other with individual keys for each lock. Master key provided that can open all locks in the set. Multiple sets are available, but no key can open all locks in all sets.




Master Keys



Individual Locks

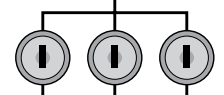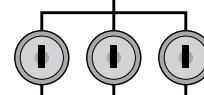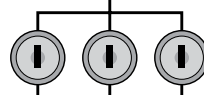

Individual Keys

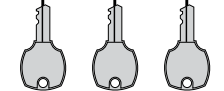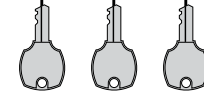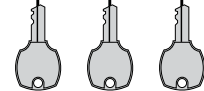


Lock Keying Options

Grand Master Keyed

Disk Tumbler Locks: N/A

Pin Tumbler Locks: Available


All locks are keyed different from each other with individual keys for each lock. Master key provided that can open all locks in the set. Multiple sets are available, with Grand master key that can open all locks in all sets.

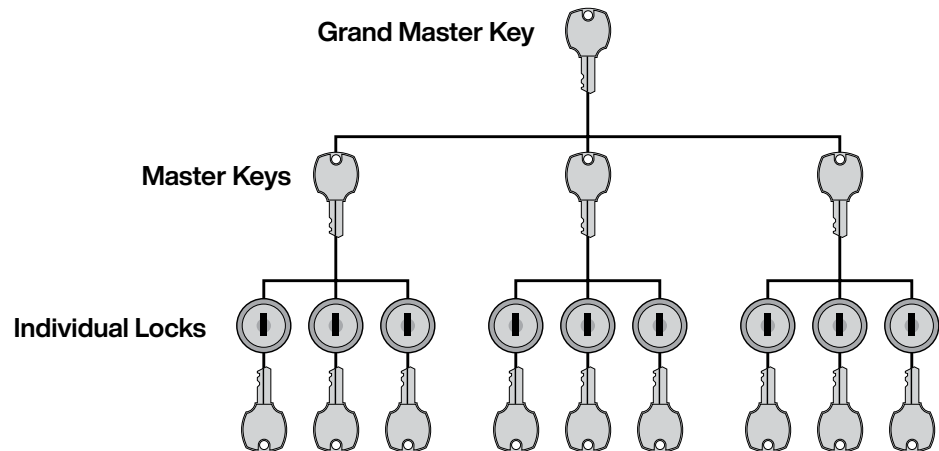

Sub-master Keyed & Mason Keyed

Disk Tumbler Locks: N/A

Pin Tumbler Locks: Available

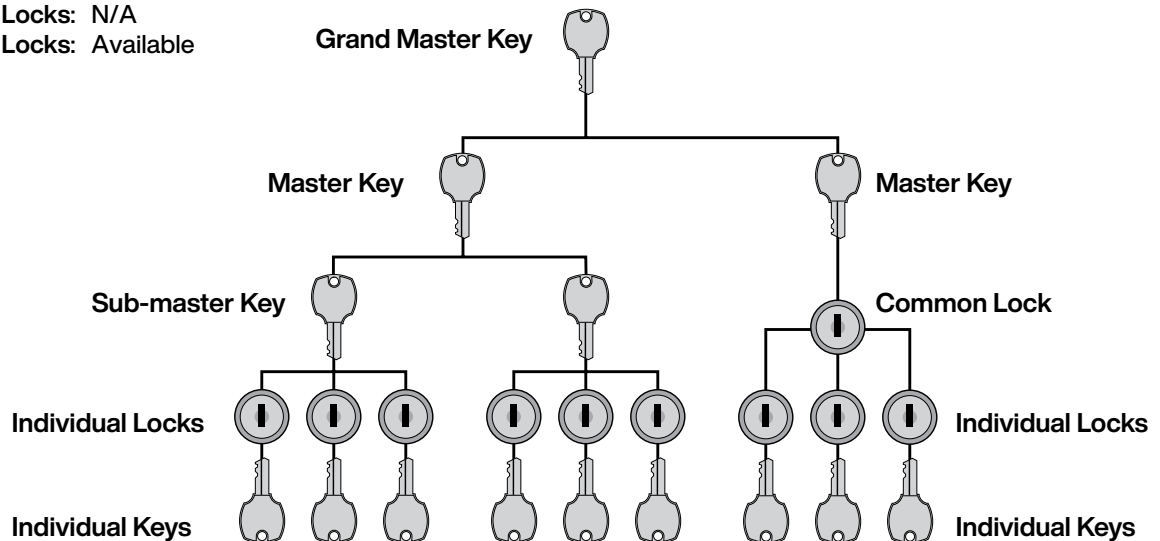


Sub-master Keyed

All locks are keyed different from each other with individual keys for each lock. Sub-master key provided that can open all locks in sub-set. Multiple sub-sets are available, with master key that can open all locks in all sub-sets. Multiple sets are available, with grand master key that can open all locks in all sets and sub-sets.

Mason Keyed

Locks are keyed different from each other with individual keys for each lock. Each individual key can also open all common locks in set. Master key provided that can open all locks in set. Multiple sets are available with grand master key that can open all locks in all sets.

Steel Sitting Height Base Cabinets

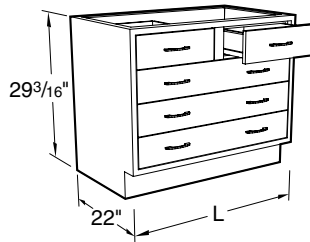



Specifications:

Sitting Height Cabinets are 29" high and 22" deep. Lengths as shown.

Open and Cupboard Cabinets are furnished with one adjustable shelf.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.

Open Cabinets

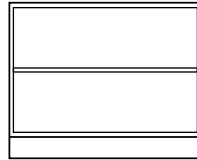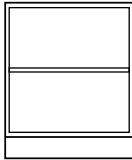




B00C292212N-00_0	12"L	B00C292230N-00_0	30"L
B00C292215N-00_0	15"L	B00C292236N-00_0	36"L
B00C292218N-00_0	18"L	B00C292242N-00_0	42"L
B00C292221N-00_0	21"L	B00C292248N-00_0	48"L
B00C292224N-00_0	24"L		

Open Cabinets may be ordered with the toespace $\frac{7}{8}$ " deeper so that cabinet fronts line up with the doors and drawer fronts on adjacent Full Overlay cabinets. To specify the deeper toe space add a **C** to the end of the part number.

Example:




B00C292236N-0000C

Cupboard Cabinets

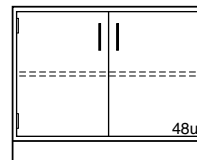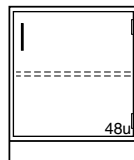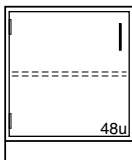

C03C292212L-_____	12"L	C03C292212R-_____	12"L	C03C292230D-_____	30"L
C03C292215L-_____	15"L	C03C292215R-_____	15"L	C03C292236D-_____	36"L
C03C292218L-_____	18"L	C03C292218R-_____	18"L	C03C292242D-_____	42"L
C03C292221L-_____	21"L	C03C292221R-_____	21"L	C03C292248D-_____	48"L
C03C292224L-_____	24"L	C03C292224R-_____	24"L		





Drawer Cabinets

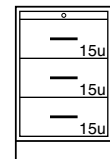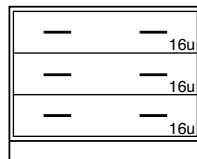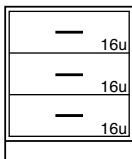

D00C292212N-_____	12"L	D00C292230N-_____	30"L
D00C292215N-_____	15"L	D00C292236N-_____	36"L
D00C292218N-_____	18"L	D00C292242N-_____	42"L
D00C292221N-_____	21"L	D00C292248N-_____	48"L
D00C292224N-_____	24"L		


D26C292218N-_____ 18"L

1 Pull Out Shelf

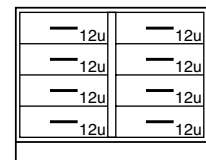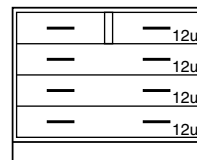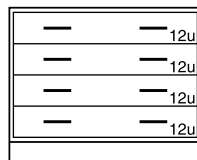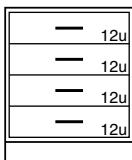

D30C292212N-_____	12"L	D30C292230N-_____	30"L	D31C292230N-_____	30"L	D32C292230N-_____	30"L
D30C292215N-_____	15"L	D30C292236N-_____	36"L	D31C292236N-_____	36"L	D32C292236N-_____	36"L
D30C292218N-_____	18"L	D30C292242N-_____	42"L	D31C292242N-_____	42"L	D32C292242N-_____	42"L
D30C292221N-_____	21"L	D30C292248N-_____	48"L	D31C292248N-_____	48"L	D32C292248N-_____	48"L
D30C292224N-_____	24"L						

Steel Sitting Height Base Cabinets

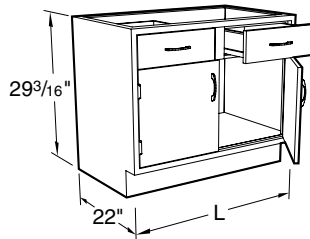


Specifications:

Sitting Height Cabinets are 29" high and 22" deep. Lengths as shown.


Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.

Drawer Cabinets

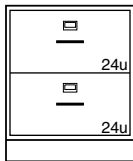


D81C292218N- 18"L
D81C292221N- 21"L
D81C292224N- 24"L

2 File Drawer with hangers

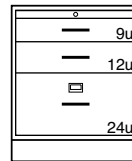

D84C292218N- 18"L
D84C292221N- 21"L
D84C292224N- 24"L


Pull Out Shelf,
1 File Drawer with hangers

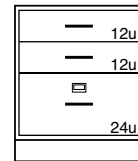

D85C292218N- 18"L
D85C292221N- 21"L
D85C292224N- 24"L

1 File Drawer with hangers


Cupboard and Drawer Cabinets

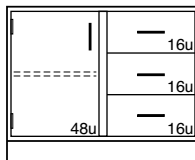


E05C292230L- 30"L
E05C292236L- 36"L
E05C292242L- 42"L
E05C292248L- 48"L




E05C292230R- 30"L
E05C292236R- 36"L
E05C292242R- 42"L
E05C292248R- 48"L

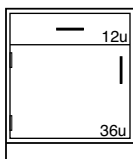

E40C292212L- 12"L
E40C292215L- 15"L
E40C292218L- 18"L
E40C292221L- 21"L
E40C292224L- 24"L

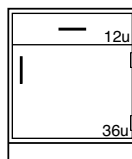


E40C292212R- 12"L
E40C292215R- 15"L
E40C292218R- 18"L
E40C292221R- 21"L
E40C292224R- 24"L

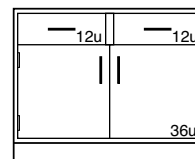

E41C292230D- 30"L
E41C292236D- 36"L
E41C292242D- 42"L
E41C292248D- 48"L

Steel Sitting Height Base Cabinets

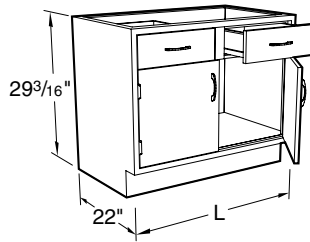



Specifications:

Sitting Height Cabinets are 29" high and 22" deep. Lengths as shown.

Sink cabinets are provided with lowered back for plumbing access. Sink support kits are available when needed to support under-mount sinks.




Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.

Cupboard and Drawer Cabinets

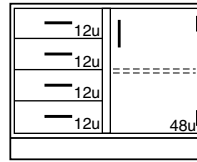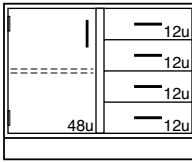


E45C292230L-_____	30"L	E45C292230R-_____	30"L
E45C292236L-_____	36"L	E45C292236R-_____	36"L
E45C292242L-_____	42"L	E45C292242R-_____	42"L
E45C292248L-_____	48"L	E45C292248R-_____	48"L

Sink Cabinets

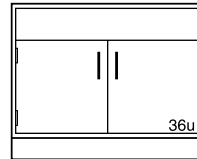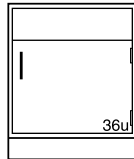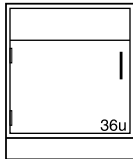


G00C292212L-_____	12"L	G00C292212R-_____	12"L	G00C292230D-_____	30"L
G00C292215L-_____	15"L	G00C292215R-_____	15"L	G00C292236D-_____	36"L
G00C292218L-_____	18"L	G00C292218R-_____	18"L	G00C292242D-_____	42"L
G00C292221L-_____	21"L	G00C292221R-_____	21"L	G00C292248D-_____	48"L
G00C292224L-_____	24"L	G00C292224R-_____	24"L		


Sink Support Kits

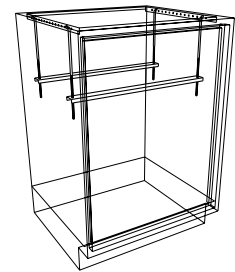

Diagonal Corner Cabinets

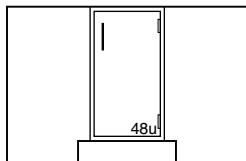

G14C292214R-_____

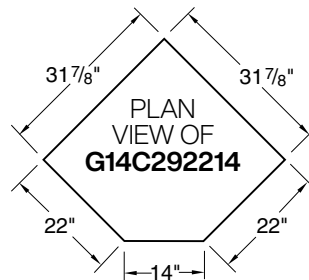


Part Number	Fits Cabinet
K01M000012N	12"
K01M000015N	15"
K01M000018N	18"
K01M000024N	24"
K01M000030N	30"
K01M000036N	36"
K01M000042N	42"
K01M000048N	48"

Sink support kit includes: four (4) threaded hanging rods and two (2) support rails and are designed to support under-mount tub sinks up to 14" deep. For sinks from 14" to 20" deep add a -D to the end of the part number.

Example: K01M000036N-D

Steel ADA Height Base Cabinets

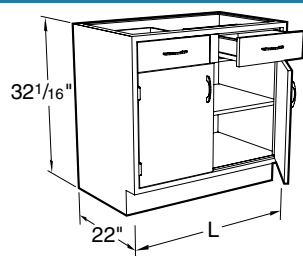



Specifications:

ADA Height Cabinets are 32" high and 22" deep. Lengths as shown.

Open and Cupboard Cabinets are furnished with one adjustable shelf.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.

Open Cabinets

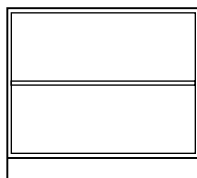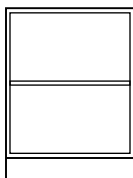

B00C322212N-00_0 12"L
B00C322215N-00_0 15"L
B00C322218N-00_0 18"L
B00C322221N-00_0 21"L
B00C322224N-00_0 24"L




B00C322230N-00_0 30"L
B00C322236N-00_0 36"L
B00C322242N-00_0 42"L
B00C322248N-00_0 48"L

Open Cabinets may be ordered with the toespace $\frac{7}{8}$ " deeper so that cabinet fronts line up with the doors and drawer fronts on adjacent Full Overlay cabinets. To specify the deeper toe space add a **C** to the end of the part number.

Example:

G00C322236D-0000C

Cupboard Cabinets

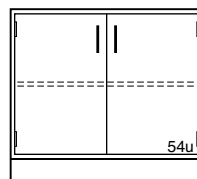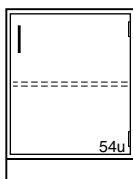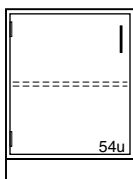


C03C322212L- 12"L
C03C322215L- 15"L
C03C322218L- 18"L
C03C322221L- 21"L
C03C322224L- 24"L

C03C322212R- 12"L
C03C322215R- 15"L
C03C322218R- 18"L
C03C322221R- 21"L
C03C322224R- 24"L

C03C322230D- 30"L
C03C322236D- 36"L
C03C322242D- 42"L
C03C322248D- 48"L

Steel ADA Height Base Cabinets

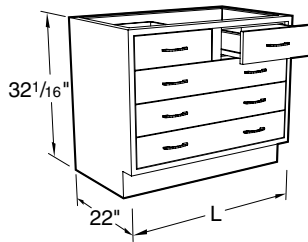


Specifications:


ADA Height Cabinets are 32" high and 22" deep. Lengths as shown.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.


Drawer Cabinets

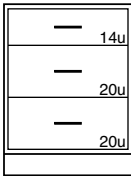


D00C322212N- 12"L
D00C322215N- 15"L
D00C322218N- 18"L
D00C322221N- 21"L
D00C322224N- 24"L

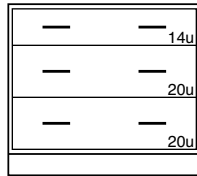


D00C322230N- 30"L
D00C322236N- 36"L
D00C322242N- 42"L
D00C322248N- 48"L

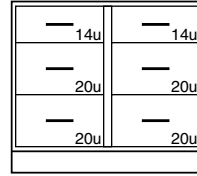


D02C322230N- 30"L
D02C322236N- 36"L
D02C322242N- 42"L
D02C322248N- 48"L

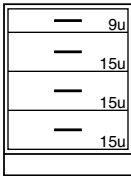


D30C322212N- 12"L
D30C322215N- 15"L
D30C322218N- 18"L
D30C322221N- 21"L
D30C322224N- 24"L

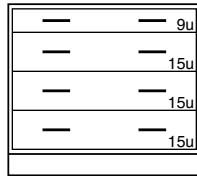


D30C322230N- 30"L
D30C322236N- 36"L
D30C322242N- 42"L
D30C322248N- 48"L

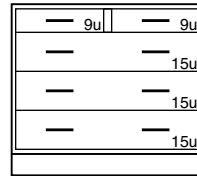


D31C322230N- 30"L
D31C322236N- 36"L
D31C322242N- 42"L
D31C322248N- 48"L

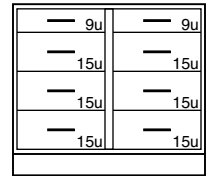


D32C322230N- 30"L
D32C322236N- 36"L
D32C322242N- 42"L
D32C322248N- 48"L

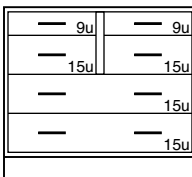


D36C322230N- 30"L
D36C322236N- 36"L
D36C322242N- 42"L
D36C322248N- 48"L

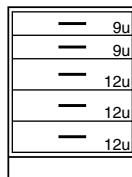


D50C322212N- 12"L
D50C322215N- 15"L
D50C322218N- 18"L
D50C322221N- 21"L
D50C322224N- 24"L

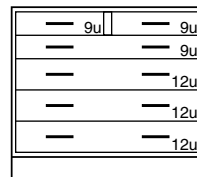


D51C322230N- 30"L
D51C322236N- 36"L
D51C322242N- 42"L
D51C322248N- 48"L

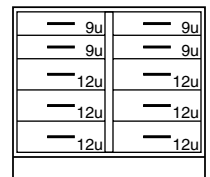

D52C322230N- 30"L
D52C322236N- 36"L
D52C322242N- 42"L
D52C322248N- 48"L

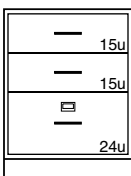


D85C322215N- 15"L
D85C322218N- 18"L
D85C322221N- 21"L
D85C322224N- 24"L

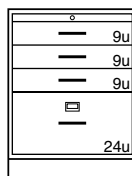

D84C322218N- 18"L
D84C322221N- 21"L
D84C322224N- 24"L

Steel ADA Height Base Cabinets

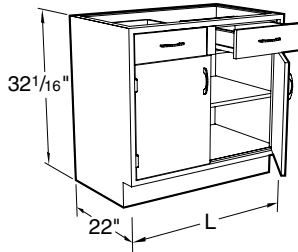


Specifications:


ADA Height Cabinets are 32" high and 22" deep. Lengths as shown.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.


Cupboard and Drawer Cabinets

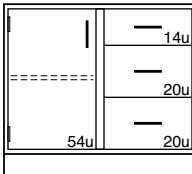


E05C322230L- 30"
E05C322236L- 36"
E05C322242L- 42"
E05C322248L- 48"

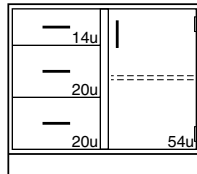


E05C322230R- 30"
E05C322236R- 36"
E05C322242R- 42"
E05C322248R- 48"

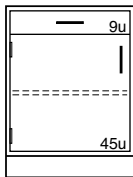


E40C322212L- 12"
E40C322215L- 15"
E40C322218L- 18"
E40C322221L- 21"
E40C322224L- 24"

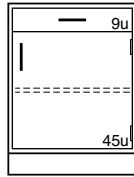


E40C322212R- 12"
E40C322215R- 15"
E40C322218R- 18"
E40C322221R- 21"
E40C322224R- 24"

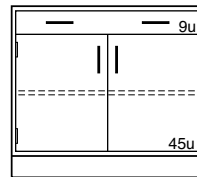


E40C322230D- 30"
E40C322236D- 36"
E40C322242D- 42"
E40C322248D- 48"

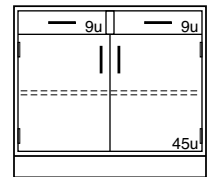


E41C322230D- 30"
E41C322236D- 36"
E41C322242D- 42"
E41C322248D- 48"

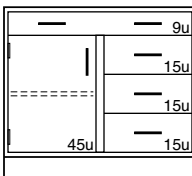


E42C322230L- 30"
E42C322236L- 36"
E42C322242L- 42"
E42C322248L- 48"

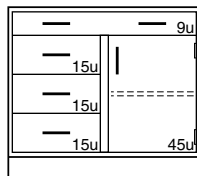


E42C322230R- 30"
E42C322236R- 36"
E42C322242R- 42"
E42C322248R- 48"

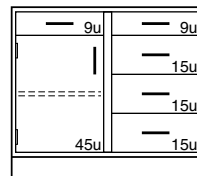


E43C322230L- 30"
E43C322236L- 36"
E43C322242L- 42"
E43C322248L- 48"

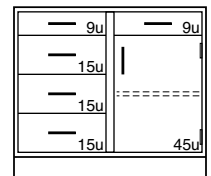


E43C322230R- 30"
E43C322236R- 36"
E43C322242R- 42"
E43C322248R- 48"

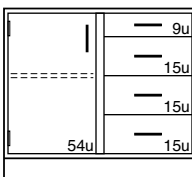

E45C322230L- 30"
E45C322236L- 36"
E45C322242L- 42"
E45C322248L- 48"

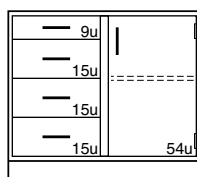


E45C322230R- 30"
E45C322236R- 36"
E45C322242R- 42"
E45C322248R- 48"

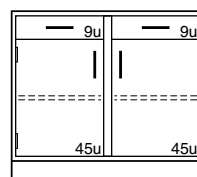

E56C322230D- 30"
E56C322236D- 36"
E56C322242D- 42"
E56C322248D- 48"

Steel ADA Height Base Cabinets

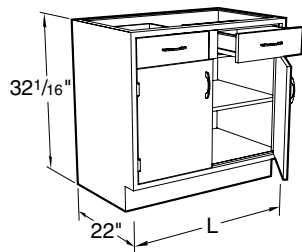


Specifications:


ADA Height Cabinets are 32" high and 22" deep. Lengths as shown.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.


Cupboard and Drawer Cabinets

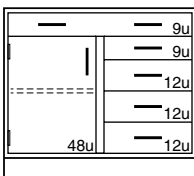


E62C322230L- 30"L
E62C322236L- 36"L
E62C322242L- 42"L
E62C322248L- 48"L

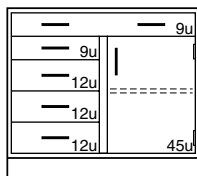


E62C322230R- 30"L
E62C322236R- 36"L
E62C322242R- 42"L
E62C322248R- 48"L

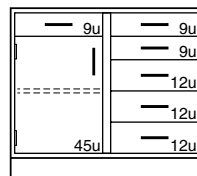


E63C322230L- 30"L
E63C322236L- 36"L
E63C322242L- 42"L
E63C322248L- 48"L

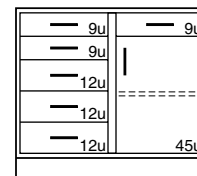

E63C322230R- 30"L
E63C322236R- 36"L
E63C322242R- 42"L
E63C322248R- 48"L

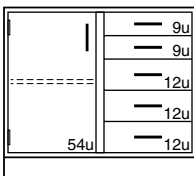


E65C322230L- 30"L
E65C322236L- 36"L
E65C322242L- 42"L
E65C322248L- 48"L

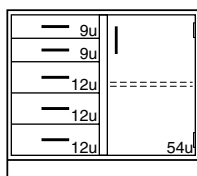

E65C322230R- 30"L
E65C322236R- 36"L
E65C322242R- 42"L
E65C322248R- 48"L

Steel ADA Height Base Cabinets

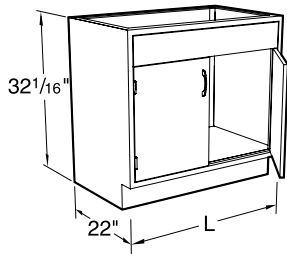


Specifications:


ADA Height Cabinets are 32" high and 22" deep. Lengths as shown.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.


Blind Corner Cabinets

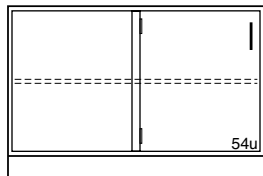

G12C322248L-____ 48"L

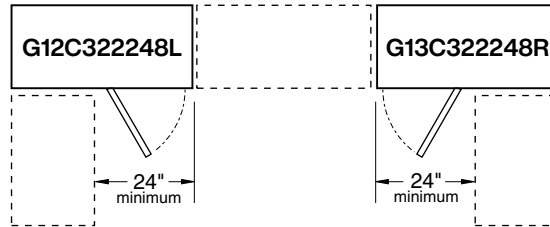


PLAN VIEW OF
G12C322248L & G13C322248R

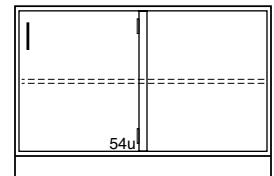


G13C322248R-____ 48"L

Acid Storage Cabinets

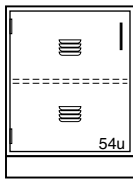


20" Deep

G21C322018L-____ 18"L
G21C322024L-____ 24"L

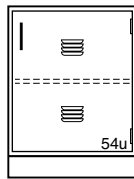

20" Deep

G21C322018R-____ 18"L
G21C322024R-____ 24"L

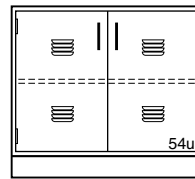

20" Deep

G21C322030D-____ 30"L
G21C322036D-____ 36"L
G21C322042D-____ 42"L
G21C322048D-____ 48"L

22" Deep

G21C322218L-____ 18"L
G21C322224L-____ 24"L

22" Deep


G21C322218R-____ 18"L
G21C322224R-____ 24"L

22" Deep


G21C322230D-____ 30"L
G21C322236D-____ 36"L
G21C322242D-____ 42"L
G21C322248D-____ 48"L

Note: G21 Acid Storage Cabinets are available with five-knuckle hinges only


Diagonal Corner Cabinet

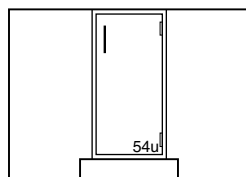

G14C322214R-____

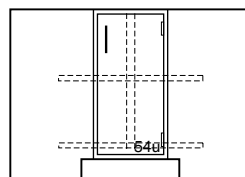


G15C322214R-____

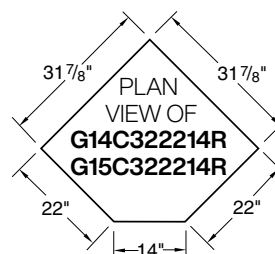

Note: G15 Cabinet contains a Lazy Susan

Steel ADA Height Base Cabinets

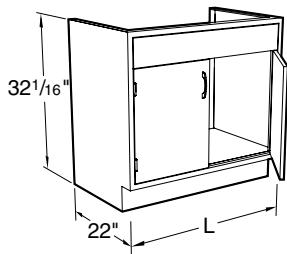


Specifications:

ADA Height Cabinets are 32" high and 22" deep. Lengths as shown.


Sink cabinets are provided with lowered back for plumbing access. Sink support kits are available when needed to support under-mount sinks.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.


Sink Cabinets

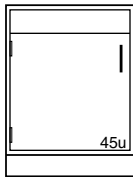

G00C322212L- 12"L
G00C322215L- 15"L
G00C322218L- 18"L
G00C322221L- 21"L
G00C322224L- 24"L

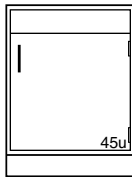


G00C322212R- 12"L
G00C322215R- 15"L
G00C322218R- 18"L
G00C322221R- 21"L
G00C322224R- 24"L

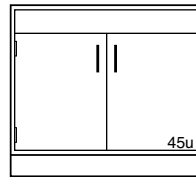

G00C322230D- 30"L
G00C322236D- 36"L
G00C322242D- 42"L
G00C322248D- 48"L


Bin Storage Cabinet

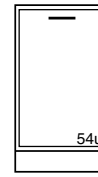

G23C322218N- 18"L

Furnished with 11" high drawer body with full height drawer front mounted on full extension slides.

Sink Support Kits

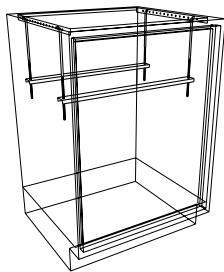


Part Number	Fits Cabinet
K01M000012N	12"
K01M000015N	15"
K01M000018N	18"
K01M000024N	24"
K01M000030N	30"
K01M000036N	36"
K01M000042N	42"
K01M000048N	48"

Sink support kit includes: four (4) threaded hanging rods and two (2) support rails and are designed to support under-mount tub sinks up to 14" deep. For sinks from 14" to 20" deep add a -D to the end of the part number.


Example: K01M000036N-D

Solvent Storage Cabinets

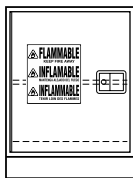


20" Deep

G63C322018L- 100 18"L
G63C322024L- 100 24"L

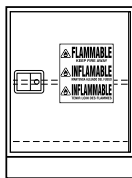

20" Deep

G63C322018R- 100 18"L
G63C322024R- 100 24"L

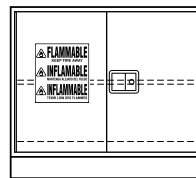

20" Deep

G63C322030D- 100 30"L
G63C322036D- 100 36"L
G63C322048D- 100 48"L

22" Deep

G63C322218L- 100 18"L
G63C322224L- 100 24"L

22" Deep

G63C322218R- 100 18"L
G63C322224R- 100 24"L

22" Deep


G63C322230D- 100 30"L
G63C322236D- 100 36"L
G63C322248D- 100 48"L

Note: Solvent Storage Cabinets are UL listed for the storage of flammable and combustible liquids, meeting UFC, OSHA and NFPA No. 30-1993 standards.

Note: G63 and G68 cabinets are provided with unique locks that will not match other keying options. The vendor provides one key-alike group.

Note: Cabinets are compatible with standard B10, B30, B40, B60 caster options.

Steel ADA Height Base Cabinets

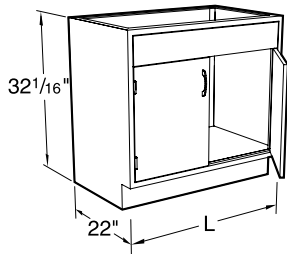


Specifications:

ADA Height Fume Hood Cabinets are 32" high and either 20" or 22" deep. Lengths as shown.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.

Fume Hood Base Cabinets

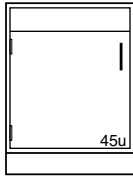


20" Deep

G08C322012L-	12"L
G08C322015L-	15"L
G08C322018L-	18"L
G08C322021L-	21"L
G08C322024L-	24"L

22" Deep

G08C322212L-	12"L
G08C322215L-	15"L
G08C322218L-	18"L
G08C322221L-	21"L
G08C322224L-	24"L

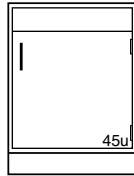


20" Deep

G08C322012R-	12"L
G08C322015R-	15"L
G08C322018R-	18"L
G08C322021R-	21"L
G08C322024R-	24"L

22" Deep

G08C322212R-	12"L
G08C322215R-	15"L
G08C322218R-	18"L
G08C322221R-	21"L
G08C322224R-	24"L

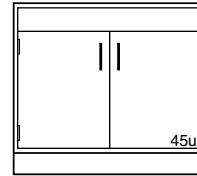


20" Deep

G08C322030D-	30"L
G08C322036D-	36"L
G08C322042D-	42"L
G08C322048D-	48"L

22" Deep

G08C322230D-	30"L
G08C322236D-	36"L
G08C322242D-	42"L
G08C322248D-	48"L

Fume Hood Solvent Storage Base Cabinets

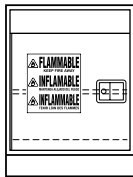


20" Deep

G68C322024L-100	24"L
-----------------	------

22" Deep

G68C322224L-100	24"L
-----------------	------

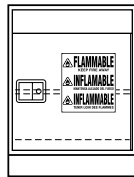


20" Deep

G68C322024R-100	24"L
-----------------	------

22" Deep

G68C322224R-100	24"L
-----------------	------

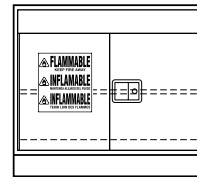

20" Deep

G68C322030D-100	30"L
G68C322036D-100	36"L
G68C322048D-100	48"L

22" Deep


G68C322230D-100	30"L
G68C322236D-100	36"L
G68C322248D-100	48"L

Note: Solvent Storage Cabinets are UL listed for the storage of flammable and combustible liquids, meeting UFC, OSHA and NFPA No. 30-1993 standards.

Note: G63 and G68 cabinets are provided with unique locks that will not match other keying options. The vendor provides one key-alike group.

Note: Cabinets are compatible with standard B10, B30, B40, B60 caster options.

Fume Hood Vacuum Storage Base Cabinets

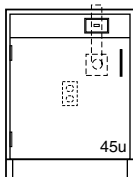


20" Deep

G35C322024L-	24"L
--------------	------

22" Deep

G35C322224L-	24"L
--------------	------

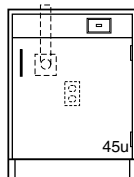


20" Deep

G35C322024R-	24"L
--------------	------

22" Deep

G35C322224R-	24"L
--------------	------

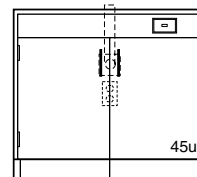

20" Deep


G35C322030D-	30"L
G35C322036D-	36"L
G35C322048D-	48"L

22" Deep

G35C322230D-	30"L
G35C322236D-	36"L
G35C322248D-	48"L


Steel ADA Height Base Cabinets

Fume Hood Acid Storage Base Cabinets

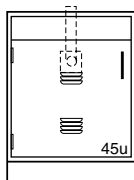


20" Deep

G80C322018L-____ 24"L
G80C322024L-____ 24"L

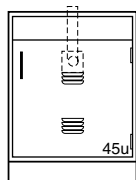

20" Deep

G80C322018R-____ 18"L
G80C322024R-____ 24"L

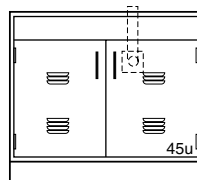

20" Deep

G80C322030D-____ 30"L
G80C322036D-____ 36"L
G80C322048D-____ 48"L

22" Deep

G80C322218L-____ 18"L
G80C322224L-____ 24"L

22" Deep


G80C322218R-____ 18"L
G80C322224R-____ 24"L

22" Deep

G80C322230D-____ 30"L
G80C322236D-____ 36"L
G80C322248D-____ 48"L


Note: G80 Acid Storage Cabinets are available with five-knuckle hinges only

Wheelchair Accessible Sink Cabinet with Doors

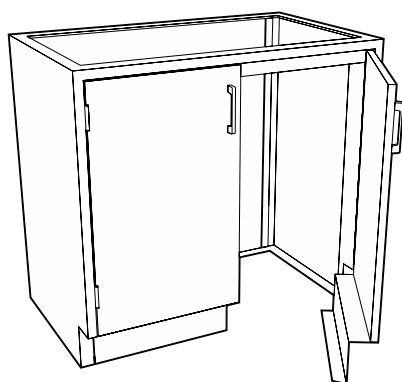


The Wheelchair Accessible Base Cabinet with Doors is designed to be used as a sink unit, but may be used anywhere wheelchair accessibility is needed. The cabinet is constructed without a bottom, and with the toe-

space rails attached to the doors so it will match adjacent cabinets when the doors are closed. The cabinet incorporates a 2¾" top rail to give a minimum of 29¼" of clearance from the floor to the bottom of the rail.

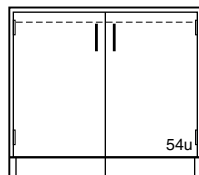

G03C322236D-____ 36"L
G03C322242D-____ 42"L
G03C322248D-____ 48"L


Wheelchair Accessible Base Cabinet without Doors

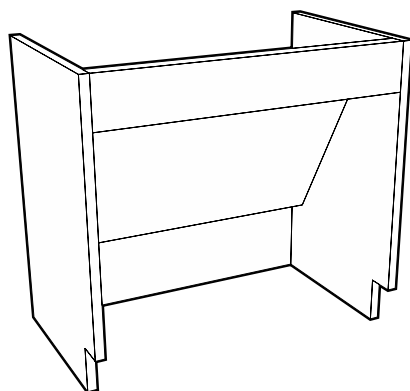

The Wheelchair Accessible Base Cabinet without Doors is designed to be used as a sink unit, but may be used anywhere wheelchair accessibility is needed. The cabinet is designed without a bottom, and without doors for easy accessibility. The cabinet incorporates a

4½" top rail to give a minimum of 27½" of clearance from the floor to the bottom of the rail and a removable sloping knee space panel to enclose piping and drain lines.

Delivered unassembled

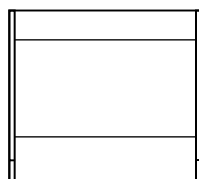


G07C322236N-___0 36"L
G07C322242N-___0 42"L
G07C322248N-___0 48"L

Wheelchair Accessible Base Cabinets without Doors may be ordered with the toespace ⅞" deeper so that cabinet fronts line up with the doors and drawer fronts on adjacent Full Overlay cabinets. To specify the deeper toe space, add a **C** to the end of the part number.

Example: G07C322236N-0000C

Steel Standing Height Base Cabinets

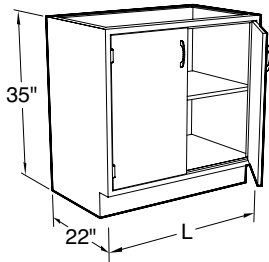



Specifications:

Standing Height Cabinets are 35" high and 22" deep. Lengths as shown.

Open and Cupboard Cabinets are furnished with one adjustable shelf.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.

Open Cabinets


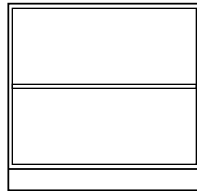
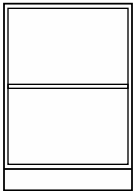


B00C352212N-00_0	12"L	B00C352230N-00_0	30"L
B00C352215N-00_0	15"L	B00C352236N-00_0	36"L
B00C352218N-00_0	18"L	B00C352242N-00_0	42"L
B00C352221N-00_0	21"L	B00C352248N-00_0	48"L
B00C352224N-00_0	24"L		

Open Cabinets may be ordered with the toespace $\frac{7}{8}$ " deeper so that cabinet fronts line up with the doors and drawer fronts on adjacent Full Overlay cabinets. To specify the deeper toe space add a **C** to the end of the part number.





Example: G00C352236N-0000C

Cupboard Cabinets

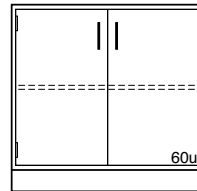
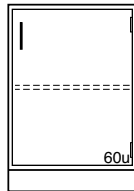
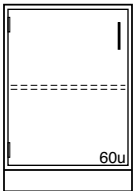


C03C352212L-_____	12"L	C03C352212R-_____	12"L	C03C352230D-_____	30"L
C03C352215L-_____	15"L	C03C352215R-_____	15"L	C03C352236D-_____	36"L
C03C352218L-_____	18"L	C03C352218R-_____	18"L	C03C352242D-_____	42"L
C03C352221L-_____	21"L	C03C352221R-_____	21"L	C03C352248D-_____	48"L
C03C352224L-_____	24"L	C03C352224R-_____	24"L		

Drawer Cabinets

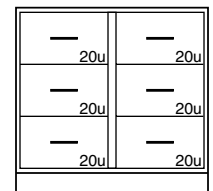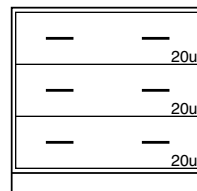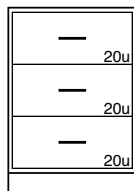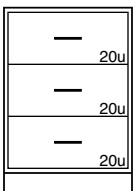


D00C352212N-_____	12"L	D00C352221N-_____	21"L	D00C352230N-_____	30"L	D02C352230N-_____	30"L
D00C352215N-_____	15"L	D00C352224N-_____	24"L	D00C352236N-_____	36"L	D02C352236N-_____	36"L
D00C352218N-_____	18"L			D00C352242N-_____	42"L	D02C352242N-_____	42"L
				D00C352248N-_____	48"L	D02C352248N-_____	48"L

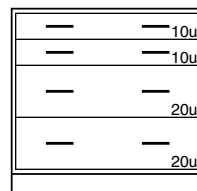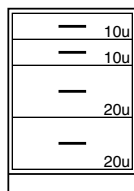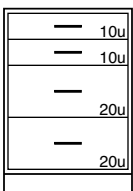

D08C352212N-_____	12"L	D08C352221N-_____	21"L	D08C352230N-_____	30"L
D08C352215N-_____	15"L	D08C352224N-_____	24"L	D08C352236N-_____	36"L
D08C352218N-_____	18"L			D08C352242N-_____	42"L
				D08C352248N-_____	48"L

Steel Standing Height Base Cabinets

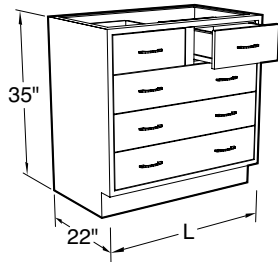


Specifications:


Standing Height Cabinets are 35" high and 22" deep. Lengths as shown.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.


Drawer Cabinets

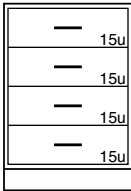


D30C352212N- 12"L
D30C352215N- 15"L
D30C352218N- 18"L
D30C352221N- 21"L
D30C352224N- 24"L

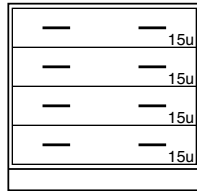


D30C352230N- 30"L
D30C352236N- 36"L
D30C352242N- 42"L
D30C352248N- 48"L

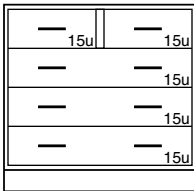


D31C352230N- 30"L
D31C352236N- 36"L
D31C352242N- 42"L
D31C352248N- 48"L

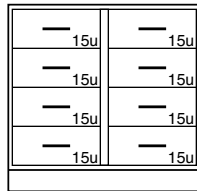


D32C352230N- 30"L
D32C352236N- 36"L
D32C352242N- 42"L
D32C352248N- 48"L

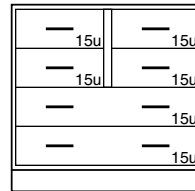


D36C352230N- 30"L
D36C352236N- 36"L
D36C352242N- 42"L
D36C352248N- 48"L

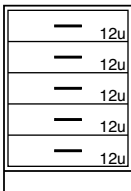


D50C352212N- 12"L
D50C352215N- 15"L
D50C352218N- 18"L
D50C352221N- 21"L
D50C352224N- 24"L

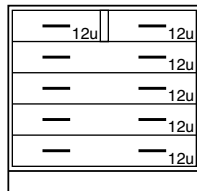

D51C352230N- 30"L
D51C352236N- 36"L
D51C352242N- 42"L
D51C352248N- 48"L

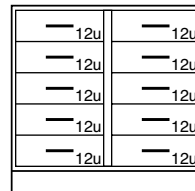


D52C352230N- 30"L
D52C352236N- 36"L
D52C352242N- 42"L
D52C352248N- 48"L

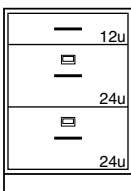

D81C352218N- 18"L
D81C352221N- 21"L
D81C352224N- 24"L


2 File Drawers with hangers

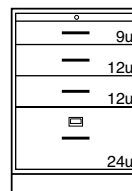

D84C352218N- 18"L
D84C352221N- 21"L
D84C352224N- 24"L

Pull Out Shelf
1 File Drawer with hangers

Steel Standing Height Base Cabinets

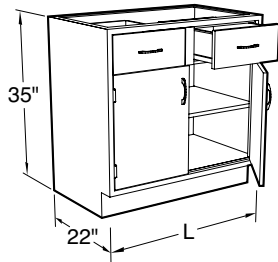


Specifications:

Standing Height Cabinets are 35" high and 22" deep. Lengths as shown.


Cupboards have one adjustable shelf.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.


Cupboard and Drawer Cabinets

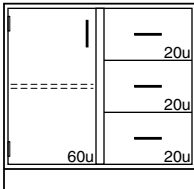


E05C352230L- 30"L
E05C352236L- 36"L
E05C352242L- 42"L
E05C352248L- 48"L

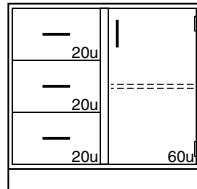


E05C352230R- 30"L
E05C352236R- 36"L
E05C352242R- 42"L
E05C352248R- 48"L

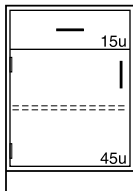


E40C352212L- 12"L
E40C352215L- 15"L
E40C352218L- 18"L
E40C352221L- 21"L
E40C352224L- 24"L

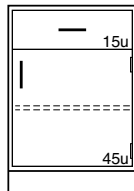


E40C352212R- 12"L
E40C352215R- 15"L
E40C352218R- 18"L
E40C352221R- 21"L
E40C352224R- 24"L

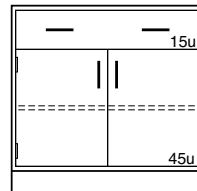


E40C352230D- 30"L
E40C352236D- 36"L
E40C352242D- 42"L
E40C352248D- 48"L

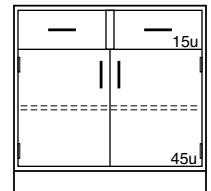


E41C352230D- 30"L
E41C352236D- 36"L
E41C352242D- 42"L
E41C352248D- 48"L

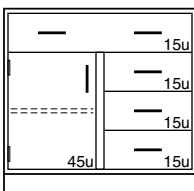


E42C352230L- 30"L
E42C352236L- 36"L
E42C352242L- 42"L
E42C352248L- 48"L

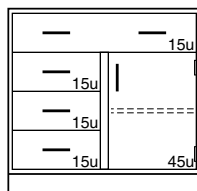


E42C352230R- 30"L
E42C352236R- 36"L
E42C352242R- 42"L
E42C352248R- 48"L

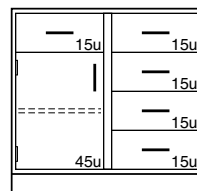


E43C352230L- 30"L
E43C352236L- 36"L
E43C352242L- 42"L
E43C352248L- 48"L

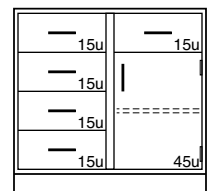


E43C352230R- 30"L
E43C352236R- 36"L
E43C352242R- 42"L
E43C352248R- 48"L

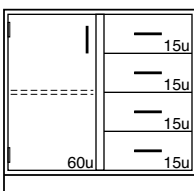

E45C352230L- 30"L
E45C352236L- 36"L
E45C352242L- 42"L
E45C352248L- 48"L

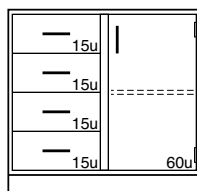


E45C352230R- 30"L
E45C352236R- 36"L
E45C352242R- 42"L
E45C352248R- 48"L

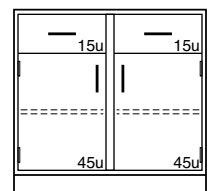

E56C352230D- 30"L
E56C352236D- 36"L
E56C352242D- 42"L
E56C352248D- 48"L

Steel Standing Height Base Cabinets

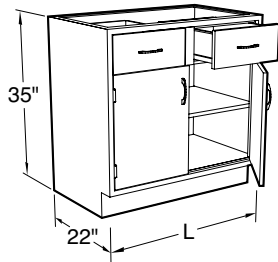


Specifications:

Standing Height Cabinets are 35" high and 22" deep. Lengths as shown.


Cupboards have one adjustable shelf.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.


Cupboard and Drawer Cabinets

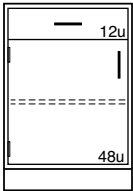


E60C352212L- 12"L
E60C352215L- 15"L
E60C352218L- 18"L
E60C352221L- 21"L
E60C352224L- 24"L

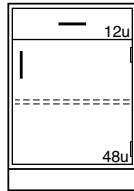


E60C352212R- 12"L
E60C352215R- 15"L
E60C352218R- 18"L
E60C352221R- 21"L
E60C352224R- 24"L

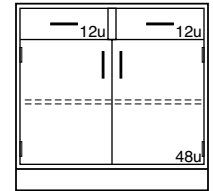


E61C352230D- 30"L
E61C352236D- 36"L
E61C352242D- 42"L
E61C352248D- 48"L

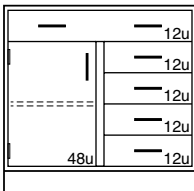


E62C352230L- 30"L
E62C352236L- 36"L
E62C352242L- 42"L
E62C352248L- 48"L

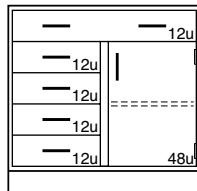


E62C352230R- 30"L
E62C352236R- 36"L
E62C352242R- 42"L
E62C352248R- 48"L

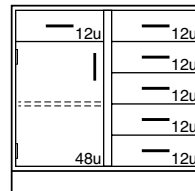


E63C352230L- 30"L
E63C352236L- 36"L
E63C352242L- 42"L
E63C352248L- 48"L

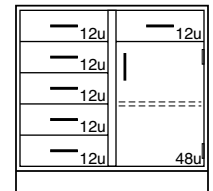

E63C352230R- 30"L
E63C352236R- 36"L
E63C352242R- 42"L
E63C352248R- 48"L

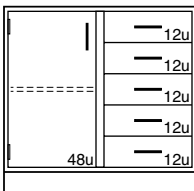


E65C352230L- 30"L
E65C352236L- 36"L
E65C352242L- 42"L
E65C352248L- 48"L

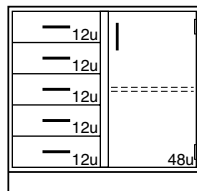

E65C352230R- 30"L
E65C352236R- 36"L
E65C352242R- 42"L
E65C352248R- 48"L

Steel Standing Height Base Cabinets

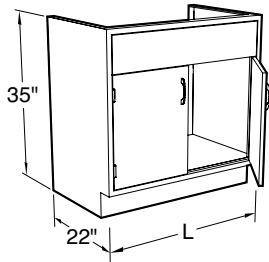

Specifications:


Standing Height Cabinets are 35" high and 22" deep. Lengths as shown.

Sink cabinets are provided with lowered back for plumbing access except where noted. Sink support kits are available when needed to support under-mount sinks.


Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.

Sink Cabinets


Sink Support Kits

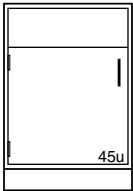



G00C352212L- 12"L
G00C352215L- 15"L
G00C352218L- 18"L
G00C352221L- 21"L
G00C352224L- 24"L

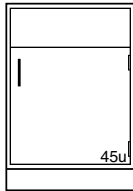


G00C352212R- 12"L
G00C352215R- 15"L
G00C352218R- 18"L
G00C352221R- 21"L
G00C352224R- 24"L

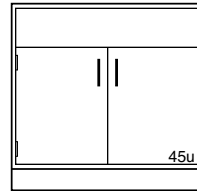

G00C352230D- 30"L
G00C352236D- 36"L
G00C352242D- 42"L
G00C352248D- 48"L

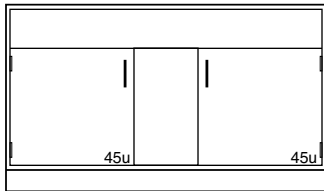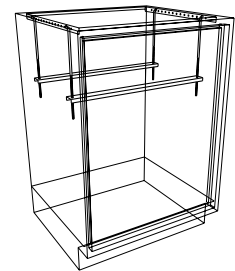


G00C352260D- 60"L

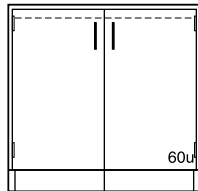

G03C352230D- 30"L
G03C352236D- 36"L
G03C352242D- 42"L
G03C352248D- 48"L

G03 Cabinets are constructed without a bottom and with the toe-space rail attached to the doors so it will match adjacent cabinets when the doors are closed.

Furnished with full height back

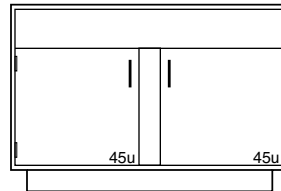

G04C352252D- 52"L
G04C352258D- 58"L
G04C352260D- 60"L

Center Bench End Sink Cabinet
Furnished with half height back

Part Number	Fits Cabinet
K01M000012N	12"
K01M000015N	15"
K01M000018N	18"
K01M000021N	21"
K01M000024N	24"
K01M000030N	30"
K01M000036N	36"
K01M000042N	42"
K01M000048N	48"
K01M000052N	52"
K01M000058N	58"
K01M000060N	60"

Sink support kit includes: four (4) threaded hanging rods and two (2) support rails and are designed to support under-mount tub sinks up to 14" deep. For sinks from 14" to 20" deep add a **-D** to the end of the part number.

Example: K01M000036N-D

Steel Standing Height Base Cabinets


Specifications

Acid Storage Base Cabinets are specifically designed for the storage of corrosive chemicals. These cabinets are lined with 1/8" thick corrosion resistant polyethylene sheets.

The cabinet doors are provided with louvers and are latched using a nylon roller catch.

An optional 1" deep polyethylene pan is available to contain spills.

Each cabinet is furnished with a PVC vent pipe for venting. (Requires a 2" hole in the counter top above.)

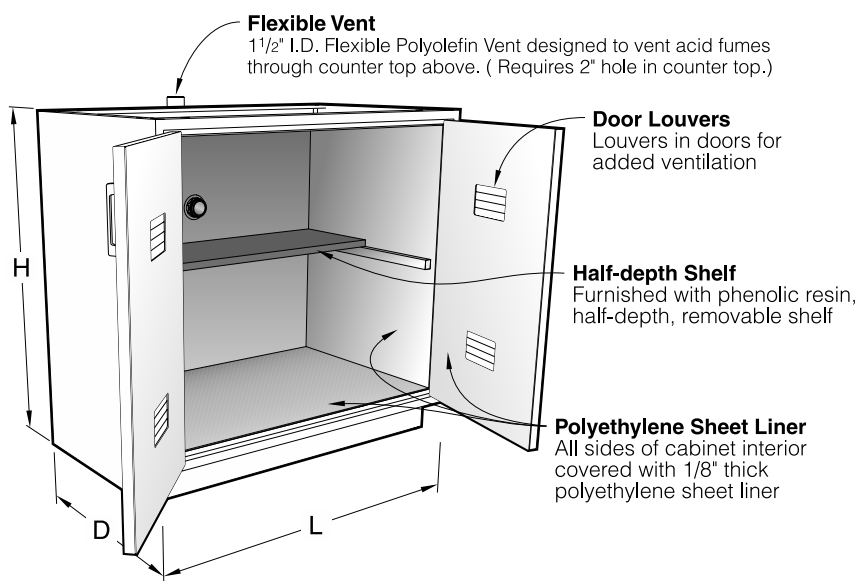


Optional Polyethylene Pan

Seamless 1" deep polyethylene pan designed to contain spills.

Part Number	Fits 20" Deep Cabinet
F-8440-15	15"L
F-8440-18	18"L
F-8440-21	21"L
F-8440-24	24"L
F-8440-30	30"L
F-8440-36	36"L
F-8440-42	42"L
F-8440-48	48"L

Part Number	Fits 22" Deep Cabinet
F-8441-15	15"L
F-8441-18	18"L
F-8441-21	21"L
F-8441-24	24"L
F-8441-30	30"L
F-8441-36	36"L
F-8441-42	42"L
F-8441-48	48"L

Acid Storage Cabinet

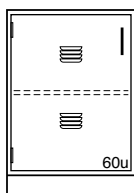


20" Deep

G21C352015L- ----	15"L
G21C352018L- ----	18"L
G21C352021L- ----	21"L
G21C352024L- ----	24"L

22" Deep

G21C352218L- ----	18"L
G21C352224L- ----	24"L

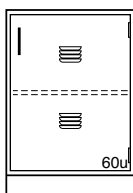


20" Deep

G21C352015R- ----	15"L
G21C352018R- ----	18"L
G21C352021R- ----	21"L
G21C352024R- ----	24"L

22" Deep

G21C352218R- ----	18"L
G21C352224R- ----	24"L

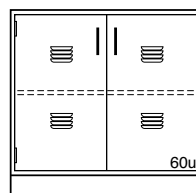

20" Deep

G21C352030D- ----	30"L
G21C352036D- ----	36"L
G21C352042D- ----	42"L
G21C352048D- ----	48"L

22" Deep

G21C352230D- ----	30"L
G21C352236D- ----	36"L
G21C352242D- ----	42"L
G21C352248D- ----	48"L

Note: G21 Acid Storage Cabinets are available with five-knuckle hinges only

Steel Standing Height Base Cabinets

Specifications


A 2" deep, steel, liquid-tight pan covers the entire bottom to contain liquid leaks and spills. The cabinet is provided with a full depth adjustable shelf and two diametrically opposed vents with spark screens in the back of the cabinet for cases when ventilation is required.

Steel Solvent Storage cabinets are 18 gauge steel, double panel construction with self-closing doors, synchronized so both doors always fully close. The right hand door is equipped with a three-point latching system that automatically engages the cabinet frame. Each door is equipped with a fusible-link hold-open


feature that ensures the doors close if the temperature outside the cabinet exceed 165°F. Cabinets are provided with a grounding screw at the rear.

All Solvent Storage Cabinets are labeled: CAUTION FLAMMABLE – KEEP FIRE AWAY in English, Spanish, and French.




Note: Solvent Storage Cabinets are UL listed for the storage of flammable and combustible liquids, meeting UFC, OSHA and NFPA No. 30-1993 standards.

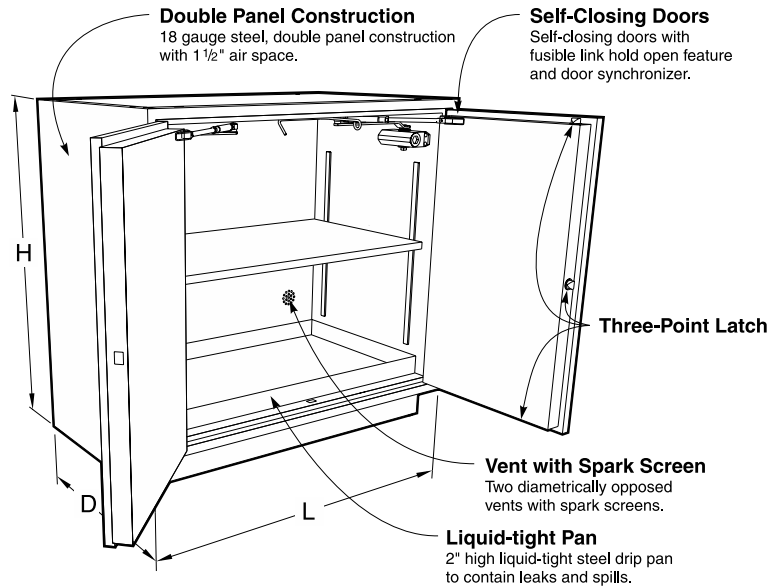

Vent Locations

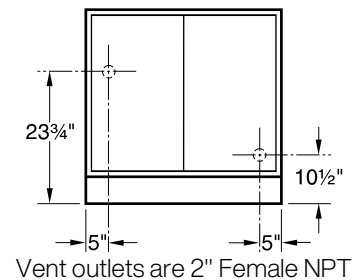

Solvent Storage Cabinets

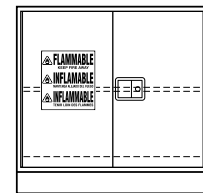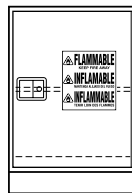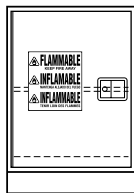

20" Deep	Capacity
G63C352024L-_100	24"L 12 gal

20" Deep	Capacity
G63C352024R-_100	24"L 12 gal

20" Deep	Capacity
G63C352030D-_100	30"L 12 gal
G63C352036D-_100	36"L 15 gal
G63C352048D-_100	48"L 20 gal

22" Deep	Capacity
G63C352224L-_100	24"L 12 gal


22" Deep	Capacity
G63C352224R-_100	24"L 12 gal

22" Deep	Capacity
G63C352230D-_100	30"L 12 gal
G63C352236D-_100	36"L 15 gal
G63C352248D-_100	48"L 20 gal

Note: G63 & G68 cabinets are provided with unique locks that will not match other keying options. The vendor provides one key-alike group.

Note: Cabinets are compatible with standard B10, B30, B40, B60 caster options.

Steel Standing Height Base Cabinets

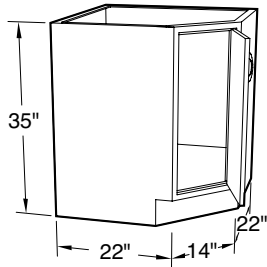


Specifications:

Standing Height Cabinets are 35" high and 22" deep. Lengths as shown.


Blind Corner Cabinets have one adjustable shelf.

Note: Last four digits of part number designate door and drawer head style, material, and pull selection. [Click here](#) for full explanation. Actual product may vary slightly from catalog illustration due to style and hardware selection.


Blind Corner Cabinets

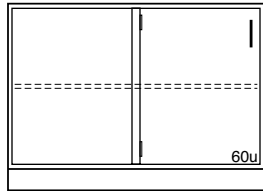

G12C352248L-____ 48"L

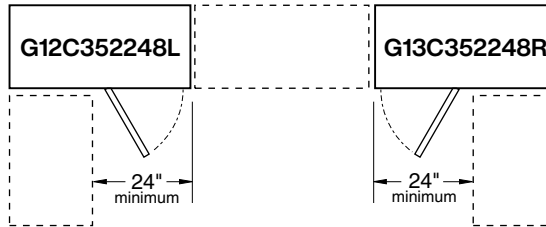


PLAN VIEW OF
G12C352248L & G13C352248R

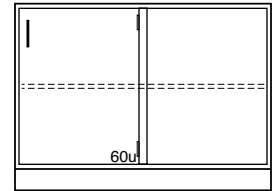


G13C352248R-____ 48"L

Diagonal Corner Cabinets

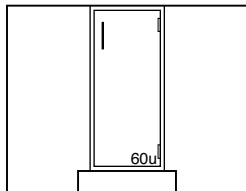


G14C352214R-____

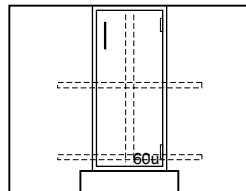


G15C352214R-____

Provided with Lazy Susan

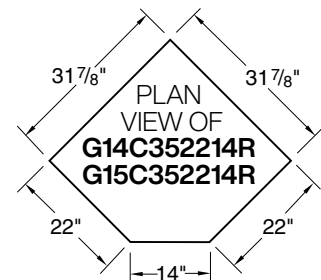

Waste Storage Cabinet

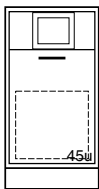

G20C352215N-____ 15"L
G20C352218N-____ 18"L
G20C352221N-____ 21"L
G20C352224N-____ 21"L

Cabinet furnished with 13" deep, cold rolled steel bin, mounted on full extension drawer slides for ease of access. Includes 6" x 5" (opening) stainless steel waste chute


Bin Storage Cabinet



G23C352218N-____ 18"L

Furnished with 11" high drawer body with full height drawer front mounted on full extension slides.

Fume Hood Base Cabinets

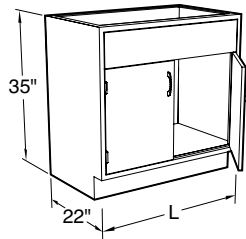


Specifications

Standing Height Steel Fume Hood Base Cabinets are 35" high and 22" deep; ADA Height Steel Fume Hood Base Cabinets are 32" high and 22" deep; lengths as shown.

Cabinet Style Option:


Specify cabinet door and hardware style when ordering by replacing blanks in last four digits of catalog number with style numbers.

Standing Height Fume Hood Base Cabinets

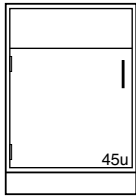


20" Deep

G08C352012L-_____	12"L
G08C352015L-_____	15"L
G08C352018L-_____	18"L
G08C352021L-_____	21"L
G08C352024L-_____	24"L

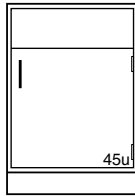


20" Deep

G08C352012R-_____	12"L
G08C352015R-_____	15"L
G08C352018R-_____	18"L
G08C352021R-_____	21"L
G08C352024R-_____	24"L

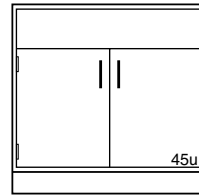

20" Deep

G08C352030D-_____	30"L
G08C352036D-_____	36"L
G08C352042D-_____	42"L
G08C352048D-_____	48"L

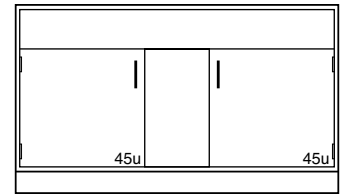

20" Deep

G08C352060D-_____	60"L
-------------------	------

22" Deep

G08C352212L-_____	12"L
G08C352215L-_____	15"L
G08C352218L-_____	18"L
G08C352221L-_____	21"L
G08C352224L-_____	24"L

22" Deep

G08C352212R-_____	12"L
G08C352215R-_____	15"L
G08C352218R-_____	18"L
G08C352221R-_____	21"L
G08C352224R-_____	24"L


22" Deep

G08C352230D-_____	30"L
G08C352236D-_____	36"L
G08C352242D-_____	42"L
G08C352248D-_____	48"L

22" Deep

G08C352260D-_____	60"L
-------------------	------

Distillation Hood Base Cabinets



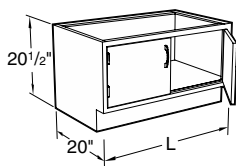


Specifications

Distillation Fume Hood Steel and Wood Base Cabinets are designed to be used under distillation fume hoods. Both are 20 1/2" high and are a nominal 20" deep. Lengths as shown.

Cabinet Style Option:

Specify cabinet door and hardware style when ordering by replacing blanks in last four digits of catalog number with style numbers.

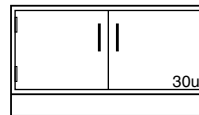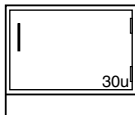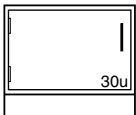


Distillation Fume Hood Base Cabinets (20" Deep)

G05C202012L-_____	12"L	G05C202012R-_____	12"L	G05C202030D-_____	30"L
G05C202015L-_____	15"L	G05C202015R-_____	15"L	G05C202036D-_____	36"L
G05C202018L-_____	18"L	G05C202018R-_____	18"L	G05C202042D-_____	42"L
G05C202024L-_____	24"L	G05C202024R-_____	24"L	G05C202048D-_____	48"L

Distillation Fume Hood Base Cabinets (22" Deep)

G05C202212L-_____	12"L	G05C202212R-_____	12"L	G05C202230D-_____	30"L
G05C202215L-_____	15"L	G05C202215R-_____	15"L	G05C202236D-_____	36"L
G05C202218L-_____	18"L	G05C202218R-_____	18"L	G05C202242D-_____	42"L
G05C202224L-_____	24"L	G05C202224R-_____	24"L	G05C202248D-_____	48"L

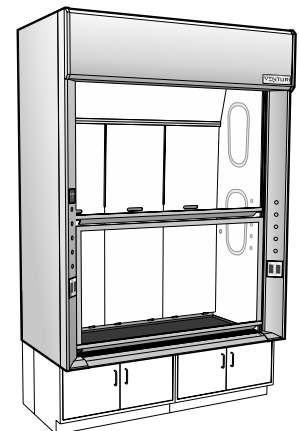


Acid Storage Fume Hood Base Cabinets

Specifications

Acid Storage Fume Hood Base Cabinets are specifically designed for the storage of corrosive chemicals. These cabinets are lined with a molded one piece linear low density polyethylene tub with coved corners and a 1" lip at the bottom front.

The cabinet doors are lined with 1/8" sheet polyethylene and the doors are latched using a nylon roller catch.


Each cabinet is furnished with a PVC vent pipe for venting to the fume hood above. (Requires a 2" hole in the hood work surface.)

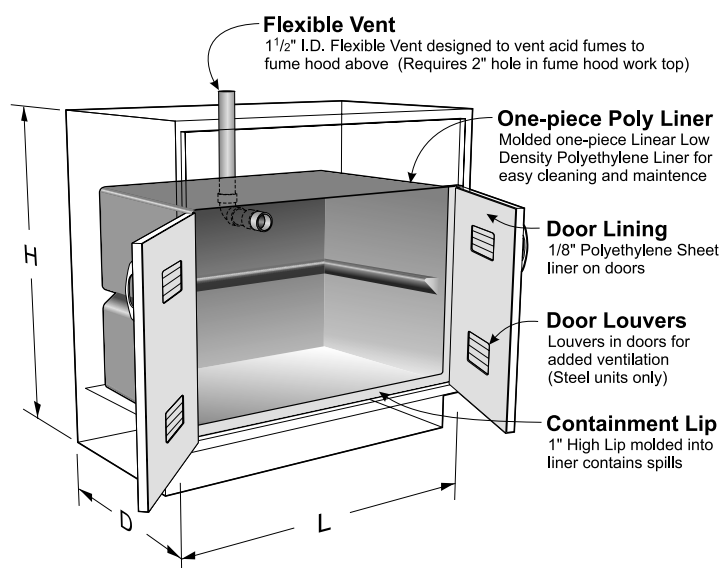

Optional Phenolic Resin Removable Half-Depth Shelf

Part Number	Shelf for Cabinet Sized
K08M010914N	18"L
K08M010919N	24"L
K08M010925N	30"L
K08M010931N	36"L
K08M010943N	48"L


Standing Height Acid Storage Fume Hood Cabinets

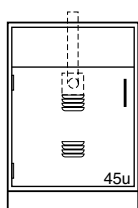

Vent Location
On Center
11 5/8" from Top of Cabinet

20" Deep


G80C352018L- ____ 18"L
G80C352024L- ____ 24"L

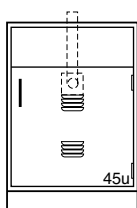

Vent Location
On Center
11 5/8" from Top of Cabinet

20" Deep


G80C352018R- ____ 18"L
G80C352024R- ____ 24"L

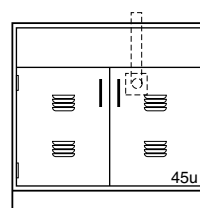

Vent Location
30" - 7 5/8" Right of Center
36" & 48" - 5" Right of Center
11 5/8" from Top of Cabinet

20" Deep

G80C352030D- ____ 30"L
G80C352036D- ____ 36"L
G80C352048D- ____ 48"L

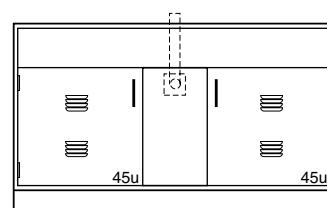

Vent Location
On Center
11 5/8" from Top of Cabinet

20" Deep

G80C352060D- ____ 60"L

22" Deep

G80C352218L- ____ 18"L
G80C352224L- ____ 24"L

22" Deep

G80C352218R- ____ 18"L
G80C352224R- ____ 24"L

22" Deep

G80C352230D- ____ 30"L
G80C352236D- ____ 36"L
G80C352248D- ____ 48"L

22" Deep

G80C352260D- ____ 60"L

Note: G80 Acid Storage Cabinets are available with five-knuckle hinges only

Solvent Storage Fume Hood Base Cabinets

Specifications


A 2" deep, steel, liquid-tight pan covers the entire bottom to contain liquid leaks and spills. The cabinet is provided with a full depth adjustable shelf and two diametrically opposed vents with spark screens in the back of the cabinet for cases when ventilation is required.

Steel Solvent Storage cabinets are 18 gauge steel, double panel construction with self-closing doors, synchronized so both doors will always fully close. The right hand door is equipped with a three-point latching system that automatically engages the cabinet frame. Each door is equipped with a fusible-link hold-open feature that ensures the doors close if the temperature outside the cabinet exceed


165°F. Cabinets are provided with a grounding screw at the rear.

All Solvent Storage Cabinets are labeled: CAUTION FLAMMABLE – KEEP FIRE AWAY in English, Spanish, and French.


Note: Solvent Storage Cabinets are UL listed for the storage of flammable and combustible liquids, meeting UFC, OSHA and NFPA No. 30-1993 standards.

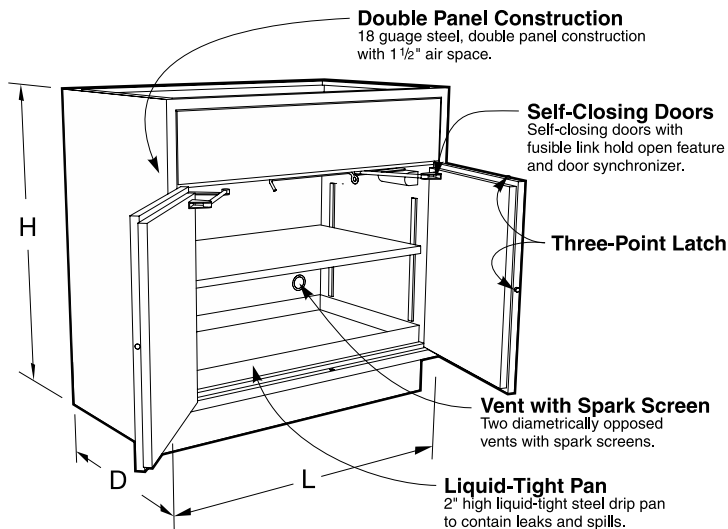


Vent Locations

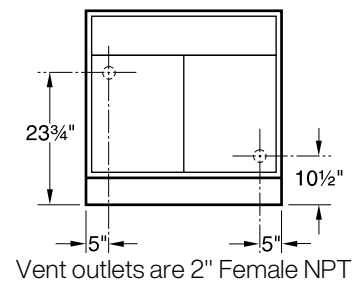


Fume Hood Solvent Storage Base Cabinets

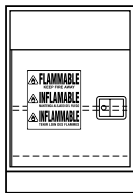

20" Deep	Capacity
G68C352024L-100	24"L 12 gal

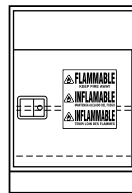

20" Deep	Capacity
G68C352024R-100	24"L 12 gal

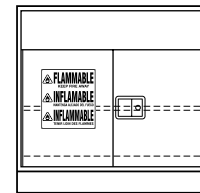

20" Deep	Capacity
G68C352030D-100	30"L 12 gal
G68C352036D-100	36"L 15 gal
G68C352048D-100	48"L 20 gal

22" Deep	Capacity
G68C352224L-100	24"L 12 gal

22" Deep	Capacity
G68C352224R-100	24"L 12 gal

22" Deep	Capacity
G68C352230D-100	30"L 12 gal
G68C352236D-100	36"L 15 gal
G68C352248D-100	48"L 20 gal

Note: G63 & G68 cabinets are provided with unique locks that will not match other keying options. The vendor provides one key-alike group.

Note: Cabinets are compatible with standard B10, B30, B40, B60 caster options.


Vacuum Pump Storage Base Cabinets

Specifications

Vacuum Pump Storage Fume Hood Cabinets are designed without a bottom to allow vacuum pumps and other equipment to be rolled in or out of the cabinets. The interior is lined with an easy to clean sound absorbent material. Each cabinet is furnished with a 120 VAC, 20

amp, duplex receptacle mounted on the inside cabinet back and a pilot lighted toggle switch mounted in the top front rail. (Wiring is not included.) Cabinets are also furnished with a 1½" diameter PVC vent pipe for venting or access to the hood above. (Requires a 2" hole in the

hood work surface.) The toespace rail is attached to the door to allow total access to the cabinet. Cabinet inside clearance at the floor is 14½" or 16½" front-to-back, 25" high, and 3" less than the overall cabinet length.

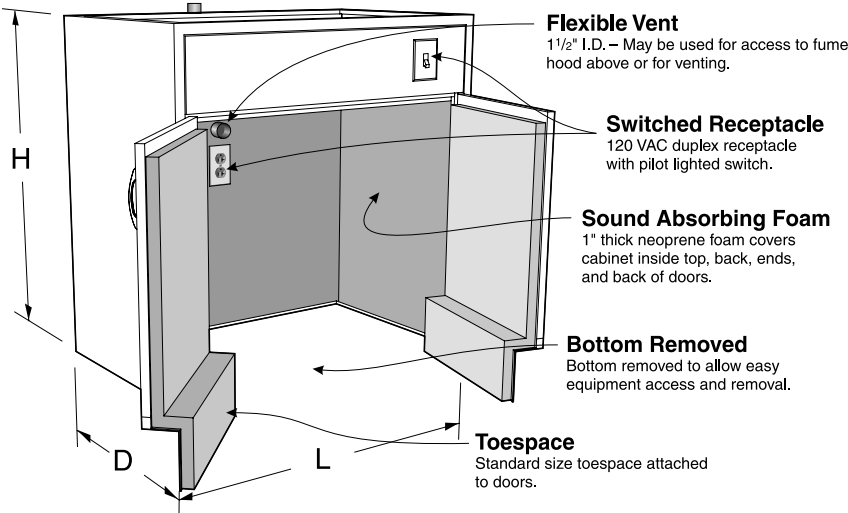

Fume Hood Vacuum Pump Storage Cabinets

20" Deep
G35C352024L-____ 24"L

20" Deep
G35C352024R-____ 24"L


20" Deep
G35C352030D-____ 30"L
G35C352036D-____ 36"L
G35C352048D-____ 48"L

22" Deep
G35C352224L-____ 24"L

22" Deep
G35C352224R-____ 24"L

22" Deep
G35C352230-____ 30"L
G35C352236-____ 36"L
G35C352248-____ 48"L

Vacuum Pump Carts

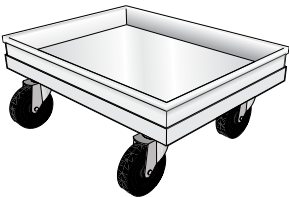

Vacuum Pump Carts are designed for use with Kewaunee's Vacuum Pump cabinets shown on this page. The cart consists of a 1" deep pan, fabricated from powder coated or stainless steel and mounted on 4" swivel casters.


Carts for 20" Deep Cabinets

Part Number	Length
VPCC061414N-00_0	14"
VPCC061420N-00_0	20"
VPCC061426N-00_0	26"
VPCC061432N-00_0	32"
VPCC061438N-00_0	38"
VPCC061444N-00_0	44"

Carts for 22" Deep Cabinets

Part Number	Length
VPCC061614N-00_0	14"
VPCC061620N-00_0	20"
VPCC061626N-00_0	26"
VPCC061632N-00_0	32"
VPCC061638N-00_0	38"
VPCC061644N-00_0	44"

Steel Sitting Height Mobile Cabinets

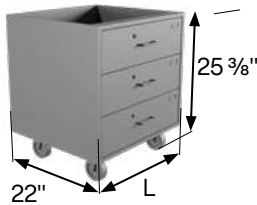

Specifications



Steel Sitting Height Mobile Cabinets are 25 3/8" high (with 3" casters) and 22" deep. Lengths as shown. Designed to fit under 30" high Workstation frames.

Drawer cabinets are furnished with counter weights. All mobile cabinets have finished backs.

Note: Last four digits of part number designates door and drawer head style, hinge, hardwood species and pull selection.

Actual product may vary slightly from catalog illustration due to style and hardware selection.

Cupboard Cabinets

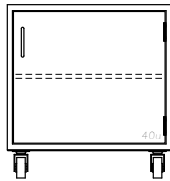
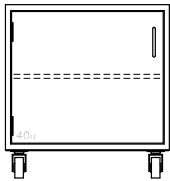


C03C262215L-__B__ 15"L
C03C262218L-__B__ 18"L
C03C262221L-__B__ 21"L
C03C262224L-__B__ 24"L

C03C262215R-__B__ 15"L
C03C262218R-__B__ 18"L
C03C262221R-__B__ 21"L
C03C262224R-__B__ 24"L

Note: Steel Mobile Cabinets are constructed with open top frames. Tops are required to close the top of the cabinets.



Drawer/Cupboard Cabinets

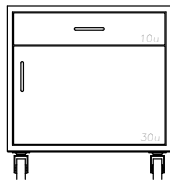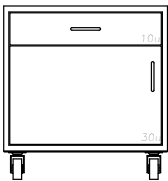

E40C262215L-__B__ 15"L
E40C262218L-__B__ 18"L
E40C262221L-__B__ 21"L
E40C262224L-__B__ 24"L


E40C262215R-__B__ 15"L
E40C262218R-__B__ 18"L
E40C262221R-__B__ 21"L
E40C262224R-__B__ 24"L

Drawer Cabinets

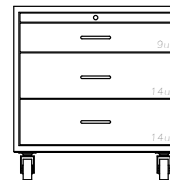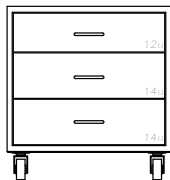

D00C262215N-__B__ 15"L
D00C262218N-__B__ 18"L
D00C262221N-__B__ 21"L
D00C262224N-__B__ 24"L

D26C262215N-__B__ 15"L
D26C262218N-__B__ 18"L
D26C262221N-__B__ 21"L
D26C262224N-__B__ 24"L

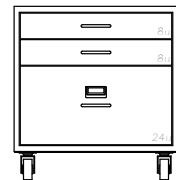

D85C262215N-__B__ 15"L
D85C262218N-__B__ 18"L
D85C262221N-__B__ 21"L
D85C262224N-__B__ 24"L

Steel ADA Height Mobile Cabinets



Specifications



Steel ADA Height Mobile Cabinets are 29¼" high (with 3" casters) and 22" deep. Lengths as shown. Designed to fit under 33" high Workstation frames.

Drawer cabinets are furnished with counter weights. All mobile cabinets have finished backs.

Note: Last four digits of part number designates door and drawer head style, hinge, hardwood species and pull selection.

Actual product may vary slightly from catalog illustration due to style and hardware selection.

Cupboard Cabinets

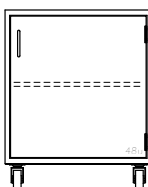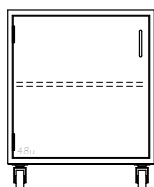



C03C292215L-__B_ 15"L
C03C292218L-__B_ 18"L
C03C292221L-__B_ 21"L
C03C292224L-__B_ 24"L

C03C292215R-__B_ 15"L
C03C292218R-__B_ 18"L
C03C292221R-__B_ 21"L
C03C292224R-__B_ 24"L

Note: Steel Mobile Cabinets are constructed with open top frames. Tops are required to close the top of the cabinets.

Additional ADA Height Mobile Cabinet configurations available by selecting the B41 or B61 option on Sitting Height Base Cabinets




Drawer/Cupboard Cabinets

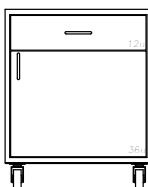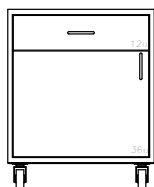

E40C292215L-__B_ 15"L
E40C292218L-__B_ 18"L
E40C292221L-__B_ 21"L
E40C292224L-__B_ 24"L

E40C292215R-__B_ 15"L
E40C292218R-__B_ 18"L
E40C292221R-__B_ 21"L
E40C292224R-__B_ 24"L

Drawer Cabinets

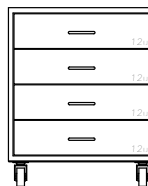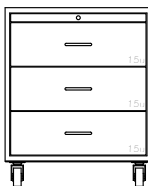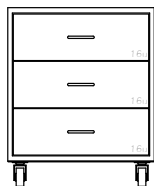


D00C292215N-__B_ 15"L
D00C292218N-__B_ 18"L
D00C292221N-__B_ 21"L
D00C292224N-__B_ 24"L

D26C292218N-__B_ 18"L

D30C292215N-__B_ 15"L
D30C292218N-__B_ 18"L
D30C292221N-__B_ 21"L
D30C292224N-__B_ 24"L

D85C292218N-__B_ 18"L
D85C292221N-__B_ 21"L
D85C292224N-__B_ 24"L

Steel Standing Height Mobile Cabinets

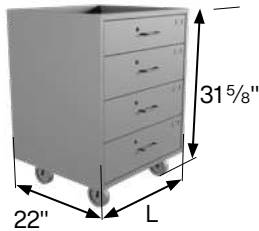

Specifications



Steel Standing Height Mobile Cabinets are 31 5/8" high (with 3" casters) and 22" deep. Lengths as shown. Designed to fit under 36" high Workstation frames.

Drawer cabinets are furnished with counter weights. All mobile cabinets have finished backs.

Note: Last four digits of part number designates door and drawer head style, hinge, hardwood species and pull selection.

Actual product may vary slightly from catalog illustration due to style and hardware selection.



Cupboard Cabinets

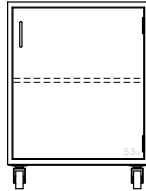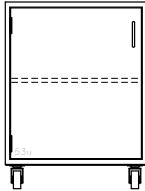

Note: Steel Mobile Cabinets are constructed with open top frames. Tops are required to close the top of the cabinets.





C03C312215L-__B_	15"L	C03C312215R-__B_	15"L
C03C312218L-__B_	18"L	C03C312218R-__B_	18"L
C03C312221L-__B_	21"L	C03C312221R-__B_	21"L
C03C312224L-__B_	24"L	C03C312224R-__B_	24"L

Drawer/Cupboard Cabinets

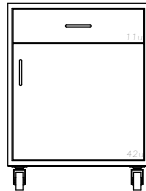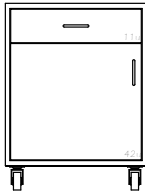

E40C312215L-__B_	15"L	E40C312215R-__B_	15"L
E40C312218L-__B_	18"L	E40C312218R-__B_	18"L
E40C312221L-__B_	21"L	E40C312221R-__B_	21"L
E40C312224L-__B_	24"L	E40C312224R-__B_	24"L

Drawer Cabinets

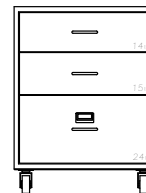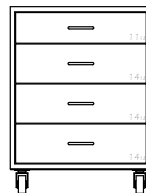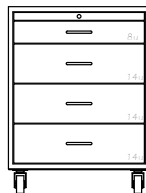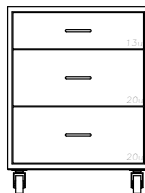

D00C312215N-__B_	15"L	D28C312215N-__B_	15"L	D30C312215N-__B_	15"L	D85C312215N-__B_	15"L
D00C312218N-__B_	18"L	D28C312218N-__B_	18"L	D30C312218N-__B_	18"L	D85C312218N-__B_	18"L
D00C312221N-__B_	21"L	D28C312221N-__B_	21"L	D30C312221N-__B_	21"L	D85C312221N-__B_	21"L
D00C312224N-__B_	24"L	D28C312224N-__B_	24"L	D30C312224N-__B_	24"L	D85C312224N-__B_	24"L

Mobile Cabinet Tops & Accessories

Specifications:

1" Thick Mobile Cabinet Tops are available in three depths, 22", 23", or 24" and three materials: Kemresin or Phenolic Resin. They may be used on any Mobile Cabinet and add 1" to the overall height of the cabinets. All edges of top are finished.

1" Thick — Mobile Cabinet Top (adds 1" to height of Mobile Cabinet)



22" Deep flush to front of case

Length	Kemresin®	Kemresin Colors	Black Phenolic Resin
15"	F94R012215N-__	BK = Black	F94T012215
18"	F94R012218N-__	GR = Grey	F94T012218
21"	F94R012221N-__	PT = Putty	F94T012221
24"	F94R012224N-__	SL = Slate	F94T012224
30"	F94R012230N-__		F94T012230
36"	F94R012236N-__		F94T012236
42"	F94R012242N-__		F94T012242
48"	F94R012248N-__		F94T012248



Inset Cabinet
End View



23" Deep flush to overlay door/drawer fronts - 1" overhang of case

Length	Kemresin®	Kemresin Colors	Black Phenolic Resin
15"	F94R012315N-__	BK = Black	F94T012315
18"	F94R012318N-__	GR = Grey	F94T012318
21"	F94R012321N-__	PT = Putty	F94T012321
24"	F94R012324N-__	SL = Slate	F94T012324
30"	F94R012330N-__		F94T012330
36"	F94R012336N-__		F94T012336
42"	F94R012342N-__		F94T012342
48"	F94R012348N-__		F94T012348



Inset Cabinet
End View



24" Deep 1" overhang of overlay door/drawer fronts

Length	Kemresin®	Kemresin Colors	Black Phenolic Resin
15"	F94R012415N-__	BK = Black	F94T012415
18"	F94R012418N-__	GR = Grey	F94T012418
21"	F94R012421N-__	PT = Putty	F94T012421
24"	F94R012424N-__	SL = Slate	F94T012424
30"	F94R012430N-__		F94T012430
36"	F94R012436N-__		F94T012436
42"	F94R012442N-__		F94T012442
48"	F94R012448N-__		F94T012448



Overlay Cabinet
End View

Mobile Cabinet Tops & Accessories

Specifications:

Steel Mobile Cabinet Tops are fabricated from 18 gauge steel and finished to match the cabinet. They are designed for Steel Mobile Cabinets only.

Mobile Cabinet Grab Bar allows easier movement and relocation of mobile cabinets. Fabricated from powder coated steel or stainless steel, they may

be used on any Mobile Cabinet, with or without Tops.

Bars must be attached in field with supplied screws and/or bolts.

22" Deep — Steel — Mobile Cabinet Top *(adds less than ¼" to height of Mobile Cabinet)*

22" Deep Top

for Steel Mobile Cabinets only



*Inset Cabinet
End View*

Lg	Painted Steel	Stainless Steel
15"	F94M012215N	F94S012215N
18"	F94M012218N	F94S012218N
21"	F94M012221N	F94S012221N
24"	F94M012224N	F94S012224N
30"	F94M012230N	F94S012230N
36"	F94M012236N	F94S012236N
42"	F94M012242N	F94S012242N
48"	F94M012248N	F94S012248N



*Overlay Cabinet
End View*

23" Deep — Steel — Mobile Cabinet Top *(adds ½" to height of Mobile Cabinet)*

23" Deep Top

for Steel Mobile Cabinets only



*Inset Cabinet
End View*

Lg	Painted Steel	Stainless Steel
15"	F94M012315N	F94S012315N
18"	F94M012318N	F94S012318N
21"	F94M012321N	F94S012321N
24"	F94M012324N	F94S012324N
30"	F94M012330N	F94S012330N
36"	F94M012336N	F94S012336N
42"	F94M012342N	F94S012342N
48"	F94M012348N	F94S012348N



*Overlay Cabinet
End View*


Mobile Cabinet Grab Bar



Length	Painted Steel	Stainless Steel
15"	F73M030815N	F73S030815N
18"	F73M030818N	F73S030818N
21"	F73M030821N	F73S030821N
24"	F73M030824N	F73S030824N



Suspended Cabinets - 19" High

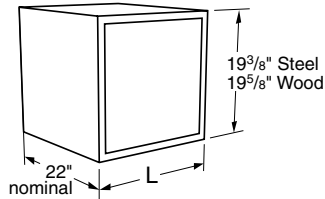


Specifications:

Suspended Base Cabinets are designed to hang beneath 30" or 36" deep work surface frames. They are available in either steel or wood construction.


Cabinets on this page are a nominal 19" high, and are available in lengths shown.

See start of catalog for designating cabinet style, door and drawer head material, and hardware selection.

Open Cabinets

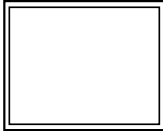

B00C192215N-____ 15"L
B00C192218N-____ 18"L
B00C192221N-____ 21"L
B00C192224N-____ 24"L

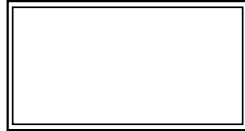


B00C192230N-____ 30"L
B00C192236N-____ 36"L
B00C192242N-____ 42"L
B00C192248N-____ 48"L


Cupboard Cabinets



C03C192215L-____ 15"L
C03C192218L-____ 18"L
C03C192221L-____ 21"L
C03C192224L-____ 24"L

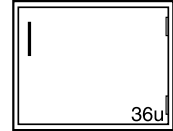


C03C192215R-____ 15"L
C03C192218R-____ 18"L
C03C192221R-____ 21"L
C03C192224R-____ 24"L

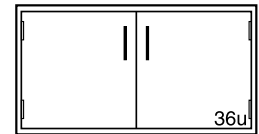


C03C192230D-____ 30"L
C03C192236D-____ 36"L
C03C192242D-____ 42"L
C03C192248D-____ 48"L


Drawer Cabinets

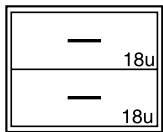


D00C192215N-_0_ 15"L
D00C192218N-_0_ 18"L
D00C192221N-_0_ 21"L
D00C192224N-_0_ 24"L

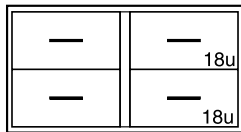


D02C192230N-_0_ 30"L
D02C192236N-_0_ 36"L
D02C192242N-_0_ 42"L
D02C192248N-_0_ 48"L

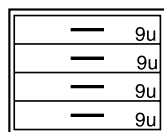


D30C192215N-_0_ 15"L
D30C192218N-_0_ 18"L
D30C192221N-_0_ 21"L
D30C192224N-_0_ 24"L

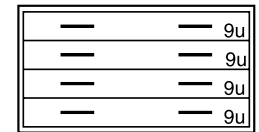


D30C192230N-_0_ 30"L
D30C192236N-_0_ 36"L
D30C192242N-_0_ 42"L
D30C192248N-_0_ 48"L

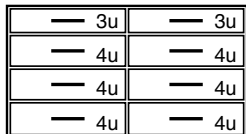


D32C192230N-_0_ 30"L
D32C192236N-_0_ 36"L
D32C192242N-_0_ 42"L
D32C192248N-_0_ 48"L

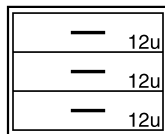

D50C192215N-_0_ 15"L
D50C192218N-_0_ 18"L
D50C192221N-_0_ 21"L
D50C192224N-_0_ 24"L

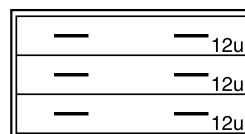


D50C192230N-_0_ 30"L
D50C192236N-_0_ 36"L
D50C192242N-_0_ 42"L
D50C192248N-_0_ 48"L

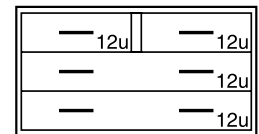

D51C192230N-_0_ 30"L
D51C192236N-_0_ 36"L
D51C192242N-_0_ 42"L
D51C192248N-_0_ 48"L

Suspended Cabinets - 19" High

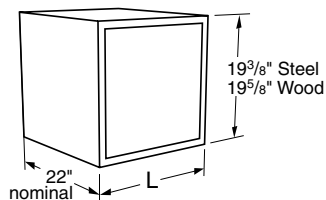


Specifications:

Suspended Base Cabinets are designed to hang beneath 30" or 36" deep work surface frames. They are available in either steel or wood construction.


Cabinets on this page are a nominal 19" high, and are available in lengths shown.

See start of catalog for designating cabinet style, door and drawer head material, and hardware selection.

Drawer Cabinets

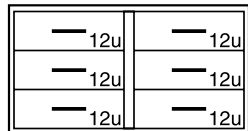


D52C192230N-0 30"L
D52C192236N-0 36"L
D52C192242N-0 42"L
D52C192248N-0 48"L

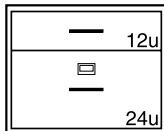


1 File Drawer with hangers

D81C192215N-0 15"L
D81C192218N-0 18"L
D81C192221N-0 21"L
D81C192224N-0 24"L

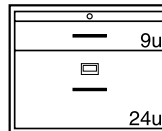


Pullboard
1 File Drawer with hangers


D85C192215N-0 15"L
D85C192218N-0 18"L
D85C192221N-0 21"L
D85C192224N-0 24"L

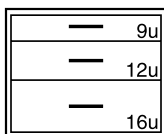

D86C192215N-0 15"L
D86C192218N-0 18"L
D86C192221N-0 21"L
D86C192224N-0 24"L

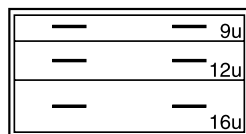


D86C192230N-0 30"L
D86C192236N-0 36"L
D86C192242N-0 42"L
D86C192248N-0 48"L

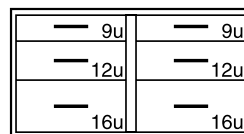


D87C192230N-0 30"L
D87C192236N-0 36"L
D87C192242N-0 42"L
D87C192248N-0 48"L

Door and Drawer Cabinets

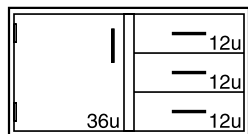


E43C192230L-0 30"L
E43C192236L-0 36"L
E43C192242L-0 42"L
E43C192248L-0 48"L

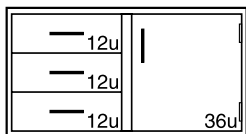

E43C192230R-0 30"L
E43C192236R-0 36"L
E43C192242R-0 42"L
E43C192248R-0 48"L

Suspended Base Cabinets – 25" High

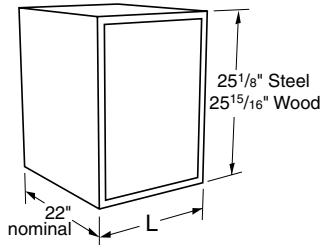


Specifications:

Suspended Base Cabinets are designed to hang beneath 30" or 36" deep work surface frames. They are available in either steel or wood construction.


Cabinets on this page are a nominal 25" high, and are available in lengths shown.

See start of catalog for designating cabinet style, door and drawer head material, and hardware selection.


Open Cabinets

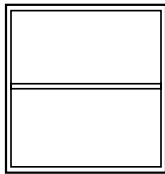


B00C252215N-00_0 15"L
B00C252218N-00_0 18"L
B00C252221N-00_0 21"L
B00C252224N-00_0 24"L




B00C252230N-00_0 30"L
B00C252236N-00_0 36"L
B00C252242N-00_0 42"L
B00C252248N-00_0 48"L

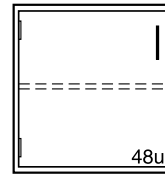

C03C252215L-____ 15"L
C03C252218L-____ 18"L
C03C252221L-____ 21"L
C03C252224L-____ 24"L

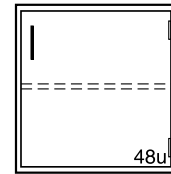


C03C252215R-____ 15"L
C03C252218R-____ 18"L
C03C252221R-____ 21"L
C03C252224R-____ 24"L

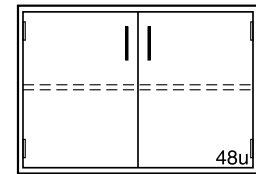


C03C252230D-____ 30"L
C03C252236D-____ 36"L
C03C252242D-____ 42"L
C03C252248D-____ 48"L


Drawer Cabinets

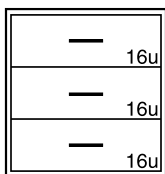


D00C252215N-__0__ 15"L
D00C252218N-__0__ 18"L
D00C252221N-__0__ 21"L
D00C252224N-__0__ 24"L

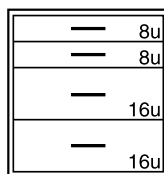


D08C252215N-__0__ 15"L
D08C252218N-__0__ 18"L
D08C252221N-__0__ 21"L
D08C252224N-__0__ 24"L

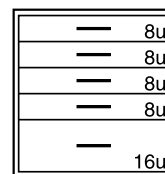


D10C252215N-__0__ 15"L
D10C252218N-__0__ 18"L
D10C252221N-__0__ 21"L
D10C252224N-__0__ 24"L

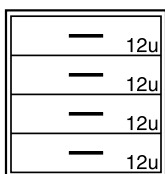


D30C252215N-__0__ 15"L
D30C252218N-__0__ 18"L
D30C252221N-__0__ 21"L
D30C252224N-__0__ 24"L

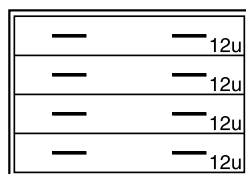

D30C252230N-__0__ 30"L
D30C252236N-__0__ 36"L
D30C252242N-__0__ 42"L
D30C252248N-__0__ 48"L

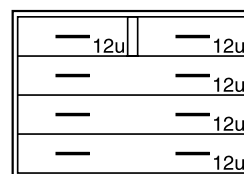


D31C252230N-__0__ 30"L
D31C252236N-__0__ 36"L
D31C252242N-__0__ 42"L
D31C252248N-__0__ 48"L

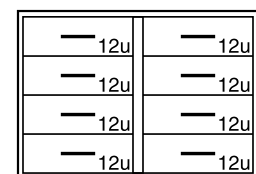

D32C252230N-__0__ 30"L
D32C252236N-__0__ 36"L
D32C252242N-__0__ 42"L
D32C252248N-__0__ 48"L

Suspended Base Cabinets – 25" High

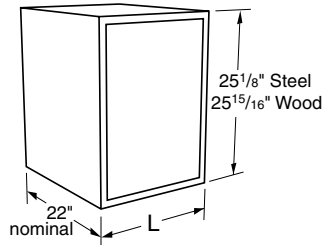


Specifications:

Suspended Base Cabinets are designed to hang beneath 30" or 36" deep work surface frames. They are available in either steel or wood construction.

Cabinets on this page are a nominal 25" high, and are available in lengths shown.


See start of catalog for designating cabinet style, door and drawer head material, and hardware selection.

Drawer Cabinets

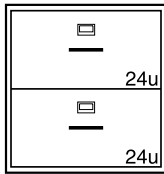

2 File Drawers with hangers


D81C252215N- _0	15"L
D81C252218N- _0	18"L
D81C252221N- _0	21"L
D81C252224N- _0	24"L

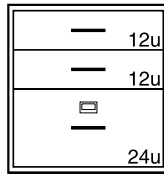


1 File Drawer with hangers

D85C252215N- _0	15"L
D85C252218N- _0	18"L
D85C252221N- _0	21"L
D85C252224N- _0	24"L


Door and Drawer Cabinets

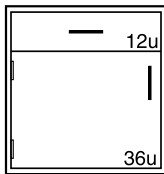


E40C252215L- _ _ _ _	15"L
E40C252218L- _ _ _ _	18"L
E40C252221L- _ _ _ _	21"L
E40C252224L- _ _ _ _	24"L

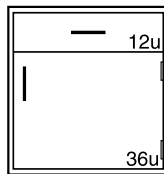


E40C252215R- _ _ _ _	15"L
E40C252218R- _ _ _ _	18"L
E40C252221R- _ _ _ _	21"L
E40C252224R- _ _ _ _	24"L

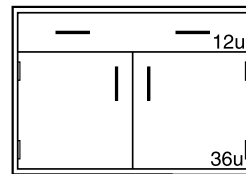


E40C252230D- _ _ _ _	30"L
E40C252236D- _ _ _ _	36"L
E40C252242D- _ _ _ _	42"L
E40C252248D- _ _ _ _	48"L

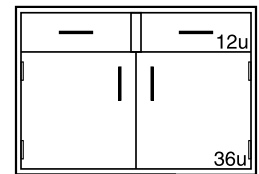

E41C252230D- _ _ _ _	30"L
E41C252236D- _ _ _ _	36"L
E41C252242D- _ _ _ _	42"L
E41C252248D- _ _ _ _	48"L

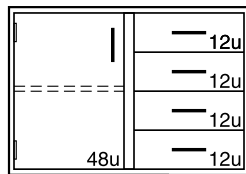


E44C252230L- _ _ _ _	30"L
E44C252236L- _ _ _ _	36"L
E44C252242L- _ _ _ _	42"L
E44C252248L- _ _ _ _	48"L

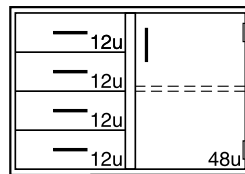


E44C252230R- _ _ _ _	30"L
E44C252236R- _ _ _ _	36"L
E44C252242R- _ _ _ _	42"L
E44C252248R- _ _ _ _	48"L


Sink Cabinets

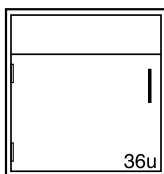

G00C252215L- _ _ _ _	15"L
G00C252218L- _ _ _ _	18"L
G00C252221L- _ _ _ _	21"L
G00C252224L- _ _ _ _	24"L

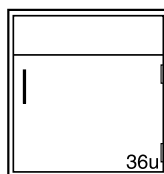


G00C252215R- _ _ _ _	15"L
G00C252218R- _ _ _ _	18"L
G00C252221R- _ _ _ _	21"L
G00C252224R- _ _ _ _	24"L

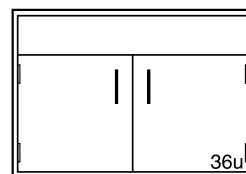

G00C252230D- _ _ _ _	30"L
G00C252236D- _ _ _ _	36"L
G00C252242D- _ _ _ _	42"L
G00C252248D- _ _ _ _	48"L

Steel Apron Assemblies

Components of an Apron Assembly

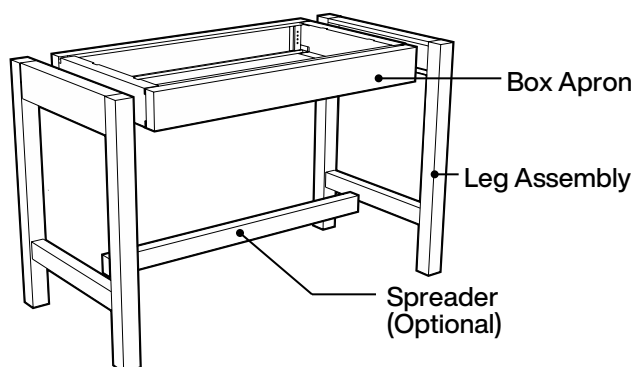

Specifications:



Box Aprons are 4" high in lengths and depths as shown. Dimensions shown are sizes of the apron without legs. Each leg or pedestal adds 2" to the overall length of the apron assembly. A free-standing table with legs on both ends will therefore be 4" longer than the apron length.

All aprons can be used between cabinets, with leg assemblies or with pedestal legs.

For aprons 42" long and longer, spreaders are available to connect the two leg assemblies. When a spreader is used, order the A08 or A28 leg assembly which are designed to accept the spreader.

Note: Caster option not available.

22" Deep Plain Rail Box Apron

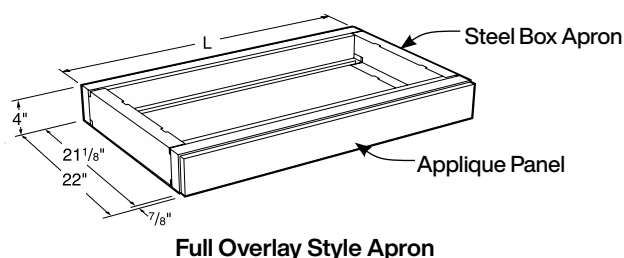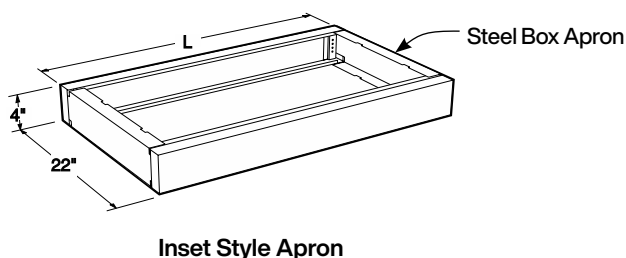

All Plain Rail Box Aprons available in either Inset or Full Overlay Style

22" Deep Apron		22" Deep Apron		22" Deep Apron	
	Apron Length		Apron Length		Apron Length
A00C042215N - ___0	15"	A00C042236N - ___0	36"	A00C042257N - ___0	57"
A00C042216N - ___0	16"	A00C042237N - ___0	37"	A00C042258N - ___0	58"
A00C042217N - ___0	17"	A00C042238N - ___0	38"	A00C042259N - ___0	59"
A00C042218N - ___0	18"	A00C042239N - ___0	39"	A00C042260N - ___0	60"
A00C042219N - ___0	19"	A00C042240N - ___0	40"	A00C042261N - ___0	61"
A00C042220N - ___0	20"	A00C042241N - ___0	41"	A00C042262N - ___0	62"
A00C042221N - ___0	21"	A00C042242N - ___0	42"	A00C042263N - ___0	63"
A00C042222N - ___0	22"	A00C042243N - ___0	43"	A00C042264N - ___0	64"
A00C042223N - ___0	23"	A00C042244N - ___0	44"	A00C042265N - ___0	65"
A00C042224N - ___0	24"	A00C042245N - ___0	45"	A00C042266N - ___0	66"
A00C042225N - ___0	25"	A00C042246N - ___0	46"	A00C042267N - ___0	67"
A00C042226N - ___0	26"	A00C042247N - ___0	47"	A00C042268N - ___0	68"
A00C042227N - ___0	27"	A00C042248N - ___0	48"	A00C042269N - ___0	69"
A00C042228N - ___0	28"	A00C042249N - ___0	49"	A00C042270N - ___0	70"
A00C042229N - ___0	29"	A00C042250N - ___0	50"	A00C042271N - ___0	71"
A00C042230N - ___0	30"	A00C042251N - ___0	51"	A00C042272N - ___0	72"
A00C042231N - ___0	31"	A00C042252N - ___0	52"		
A00C042232N - ___0	32"	A00C042253N - ___0	53"		
A00C042233N - ___0	33"	A00C042254N - ___0	54"		
A00C042234N - ___0	34"	A00C042255N - ___0	55"		
A00C042235N - ___0	35"	A00C042256N - ___0	56"		


Note: Box apron lengths available in 1" increments from 15" to 72".

Length

A00C0422   N-___0

Steel Apron Assemblies

Components of an Apron Assembly

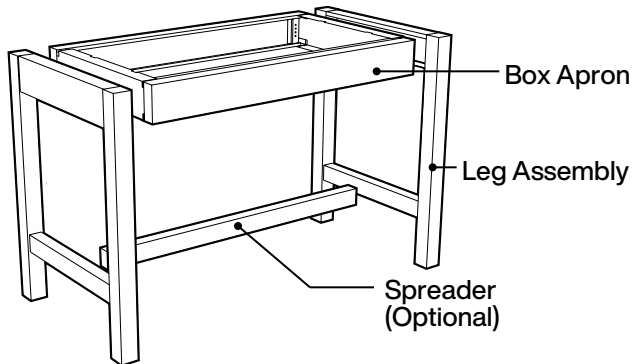

Specifications:



Box Aprons are 4" high in lengths and depths as shown. Dimensions shown are sizes of the apron without legs. Each leg or pedestal adds 2" to the overall length of the apron assembly. A free-standing table with legs on both ends will therefore be 4" longer than the apron length.

All aprons can be used between cabinets, with leg assemblies or with pedestal legs.

For aprons 42" long and longer, spreaders are available to connect the two leg assemblies. When a spreader is used, order the A08 or A28 leg assembly which are designed to accept the spreader.

Note: Caster option not available.

28" Deep Plain Rail Box Apron

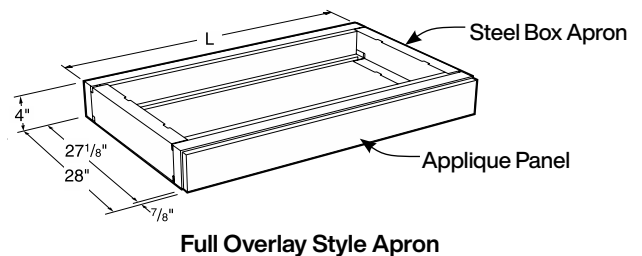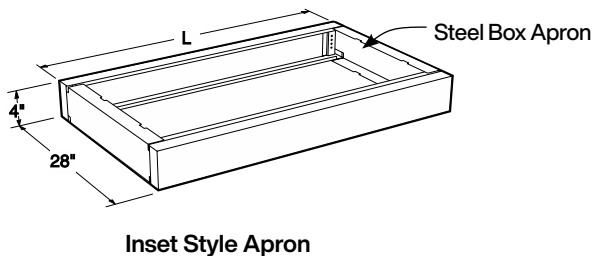

All Plain Rail Box Aprons available in either Inset or Full Overlay Style


28" Deep Apron		28" Deep Apron		28" Deep Apron	
	Apron Length		Apron Length		Apron Length
A00C042815N - ___0	15"	A00C042836N - ___0	36"	A00C042857N - ___0	57"
A00C042816N - ___0	16"	A00C042837N - ___0	37"	A00C042858N - ___0	58"
A00C042817N - ___0	17"	A00C042838N - ___0	38"	A00C042859N - ___0	59"
A00C042818N - ___0	18"	A00C042839N - ___0	39"	A00C042860N - ___0	60"
A00C042819N - ___0	19"	A00C042840N - ___0	40"	A00C042861N - ___0	61"
A00C042820N - ___0	20"	A00C042841N - ___0	41"	A00C042862N - ___0	62"
A00C042821N - ___0	21"	A00C042842N - ___0	42"	A00C042863N - ___0	63"
A00C042822N - ___0	22"	A00C042843N - ___0	43"	A00C042864N - ___0	64"
A00C042823N - ___0	23"	A00C042844N - ___0	44"	A00C042865N - ___0	65"
A00C042824N - ___0	24"	A00C042845N - ___0	45"	A00C042866N - ___0	66"
A00C042825N - ___0	25"	A00C042846N - ___0	46"	A00C042867N - ___0	67"
A00C042826N - ___0	26"	A00C042847N - ___0	47"	A00C042868N - ___0	68"
A00C042827N - ___0	27"	A00C042848N - ___0	48"	A00C042869N - ___0	69"
A00C042828N - ___0	28"	A00C042849N - ___0	49"	A00C042870N - ___0	70"
A00C042829N - ___0	29"	A00C042850N - ___0	50"	A00C042871N - ___0	71"
A00C042830N - ___0	30"	A00C042851N - ___0	51"	A00C042872N - ___0	72"
A00C042831N - ___0	31"	A00C042852N - ___0	52"		
A00C042832N - ___0	32"	A00C042853N - ___0	53"		
A00C042833N - ___0	33"	A00C042854N - ___0	54"		
A00C042834N - ___0	34"	A00C042855N - ___0	55"		
A00C042835N - ___0	35"	A00C042856N - ___0	56"		

Note: Box apron lengths available in 1" increments from 15" to 72".

Length
A00C0428   **N-___0**

Steel Apron Assemblies

Components of an Apron Assembly

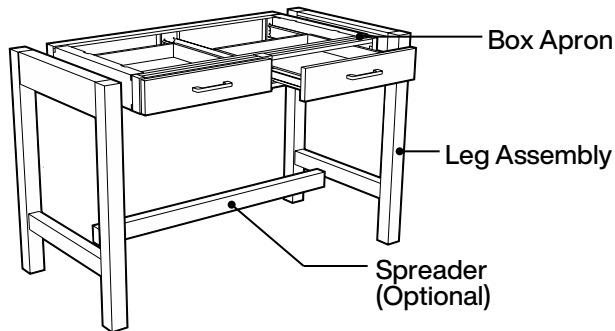

Specifications:


Box Aprons are 4" high in lengths and depths as shown. Dimensions shown are sizes of the apron without legs. Each leg or pedestal adds 2" to the overall length of the apron assembly. A free-standing table with legs on both ends will therefore be 4" longer than the apron length.

All aprons can be used between cabinets, with leg assemblies or with pedestal legs.

For aprons 42" long and longer, spreaders are available to connect the two leg assemblies. When a spreader is used, order the A08 or A28 leg assembly which are designed to accept the spreader.


Note: Caster option not available.

22" Deep Box Apron with Drawers

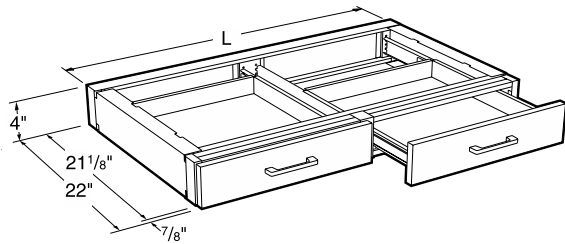

Note: Apron Drawers greater than 30"L receive two pulls.


1 Drawer

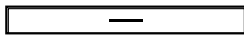

22" Deep Apron

	Apron Length	Dwr Size
A51C042220N -1	20"	18"
A51C042224N -1	24"	21"
A51C042226N -1	26"	24"
A51C042230N -1	30"	27"
A51C042232N -1	32"	30"
A51C042236N -1	36"	33"


2 Drawers

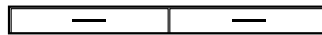

22" Deep Apron

	Apron Length	Dwr Size
A52C042226N -1	26"	12"
A52C042230N -1	30"	12"
A52C042232N -1	32"	15"
A52C042236N -1	36"	15"
A52C042238N -1	38"	18"
A52C042242N -1	42"	18"
A52C042244N -1	44"	21"
A52C042248N -1	48"	21"
A52C042250N -1	50"	24"
A52C042254N -1	54"	24"
A52C042256N -1	56"	27"
A52C042260N -1	60"	27"
A52C042262N -1	62"	30"
A52C042266N -1	66"	30"
A52C042268N -1	68"	33"
A52C042272N -1	72"	33"

3 Drawers

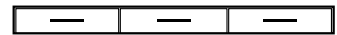


22" Deep Apron

	Apron Length	Dwr Size
A53C042250N -1	50"	15"
A53C042254N -1	54"	15"
A53C042256N -1	56"	17"
A53C042260N -1	60"	17"
A53C042262N -1	62"	19"
A53C042266N -1	66"	19"
A53C042268N -1	68"	21"
A53C042272N -1	72"	21"

Steel Apron Assemblies

Components of an Apron Assembly

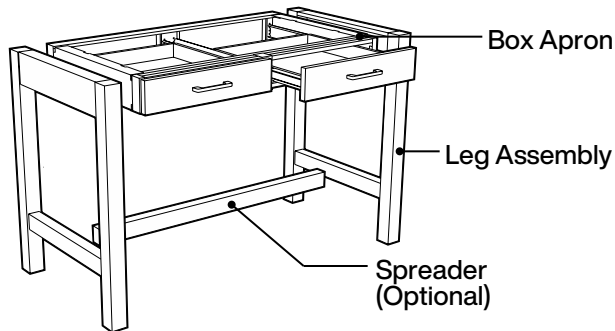

Specifications:


Box Aprons are 4" high in lengths and depths as shown. Dimensions shown are sizes of the apron without legs. Each leg or pedestal adds 2" to the overall length of the apron assembly. A free-standing table with legs on both ends will therefore be 4" longer than the apron length.

All aprons can be used between cabinets, with leg assemblies or with pedestal legs.

For aprons 42" long and longer, spreaders are available to connect the two leg assemblies. When a spreader is used, order the A08 or A28 leg assembly which are designed to accept the spreader.


Note: Caster option not available.

28" Deep Box Apron with Drawers

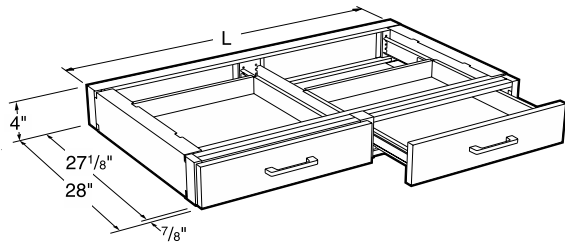

Note: Apron Drawers greater than 30"L receive two pulls.


1 Drawer

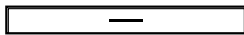

28" Deep Apron

	Apron Length	Dwr Size
A51C042820N -1	20"	18"
A51C042824N -1	24"	21"
A51C042826N -1	26"	24"
A51C042830N -1	30"	27"
A51C042832N -1	32"	30"
A51C042836N -1	36"	33"


2 Drawers

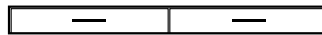

28" Deep Apron

	Apron Length	Dwr Size
A52C042826N -1	26"	12"
A52C042830N -1	30"	12"
A52C042832N -1	32"	15"
A52C042836N -1	36"	15"
A52C042838N -1	38"	18"
A52C042842N -1	42"	18"
A52C042844N -1	44"	21"
A52C042848N -1	48"	21"
A52C042850N -1	50"	24"
A52C042854N -1	54"	24"
A52C042856N -1	56"	27"
A52C042860N -1	60"	27"
A52C042862N -1	62"	30"
A52C042866N -1	66"	30"
A52C042868N -1	68"	33"
A52C042872N -1	72"	33"

3 Drawers

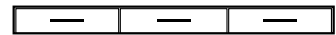

28" Deep Apron

	Apron Length	Dwr Size
A53C042850N -1	50"	15"
A53C042854N -1	54"	15"
A53C042856N -1	56"	17"
A53C042860N -1	60"	17"
A53C042862N -1	62"	19"
A53C042866N -1	66"	19"
A53C042868N -1	68"	21"
A53C042872N -1	72"	21"

Steel Apron Assemblies


Specifications:

Corner Box Aprons provide support for corner worktops as well as the extra storage space of an apron drawer. The blind end may either be attached to the wall or supported by a leg or pedestal assembly.


Corner Box Aprons are furnished with steel panels when specified with wood drawer heads.

2" High Box Aprons are designed for use when knee space height is critical. They are available 22" deep and in lengths from 15" to 72" on 1" increments. 2" High Box Aprons are not designed for use with legs.


Corner Box Apron

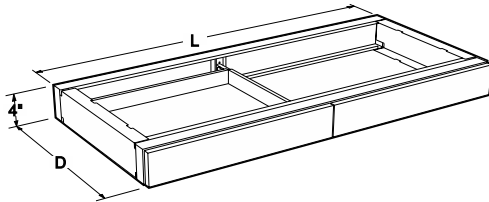

Furnished with steel panels when specified with wood drawer heads.

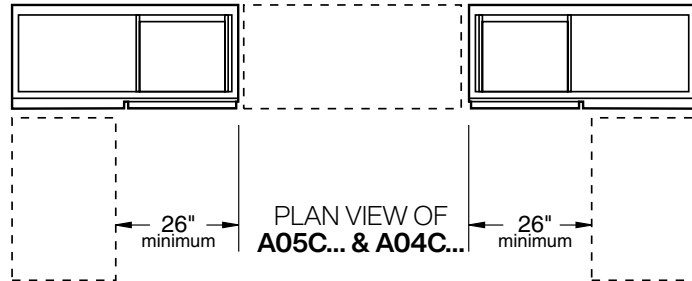

1 Drawer on Right

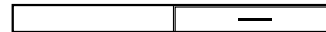

22" Deep

	Apron Length	Dwr Size
A05C042248R -1	48"	21"
A05C042254R -1	54"	21"
A05C042260R -1	60"	21"
A05C042266R -1	66"	21"
A05C042272R -1	72"	21"

28" Deep

	Apron Length	Dwr Size
A05C042848R -1	48"	21"
A05C042854R -1	54"	21"
A05C042860R -1	60"	21"
A05C042866R -1	66"	21"
A05C042872R -1	72"	21"

1 Drawer on Left



22" Deep

	Apron Length	Dwr Size
A04C042248L -1	48"	21"
A04C042254L -1	54"	21"
A04C042260L -1	60"	21"
A04C042266L -1	66"	21"
A04C042272L -1	72"	21"

28" Deep

	Apron Length	Dwr Size
A04C042848L -1	48"	21"
A04C042854L -1	54"	21"
A04C042860L -1	60"	21"
A04C042866L -1	66"	21"
A04C042872L -1	72"	21"

Steel Apron Assemblies

Specifications:


Apron Rails may be used in place of Box Aprons to provide work surface support between two units. They are available 2" or 4" high and from 15" to 72" long in 1" increments. They may be installed vertically to provide a 2" or 4" high rail, or

turned 90 degrees to provide a 1½" high rail. Apron rails should be ordered in pairs to provide both front and rear top support.

Worktop Support Rails and Assemblies may be used anywhere maximum knee space is wanted but top support is

required. Fabricated of heavy gauge steel channels, it is only ⅞" high. The Keyboard Tray may be added to any of the frames or may be used independently.

Inset 4" High Apron Rail

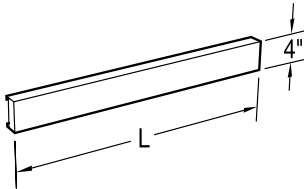


A06C0402□□N-00_0 4" High

Inset Apron Rails are available from 15" to 72" long in 1" increments. To specify apron rail, insert the length required into the part number.

Example: A06C040248N-0000 is a 48" long cold rolled steel apron rail.

Apron rails should be ordered in pairs to provide both front and rear work surface support.

Overlay 4" High Apron Rail

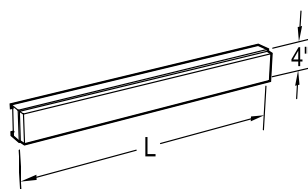


A06C0403□□N-1__0 4" High

Overlay Apron Rails are available from 15" to 72" long in 1" increments. To specify apron rail, insert the length required into the part number.

Example: A06C040348N-1610 is a 48" long apron rail with Oak overlay panel.

Apron rails should be ordered in pairs to provide both front and rear work surface support.

4" High Adjustable Rail

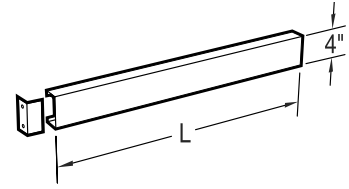


A56C040236N-00_0 36"L A56C040248N-00_0 48"L A56C040260N-00_0 60"L A56C040272N-00_0 72"L

Same as A06 Apron Rails except furnished with a loose angle at one end to allow the rail to be cut in the field to fit available space.

Apron rails should be ordered in pairs to provide both front and rear work surface support.

2" High Plain Rail Box Apron

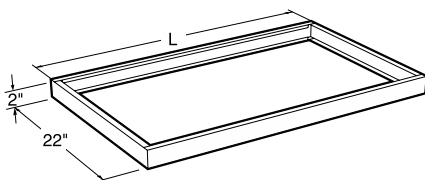


A00C0222□□N-00_0

2" high Box Aprons are available from 15" to 72" long in 1" increments. To specify apron, insert the length required into the part number.

Example: A00C022248N-0000 is a 48" long 2" high cold rolled steel box apron.

Not for use with legs

2" High Apron Rail

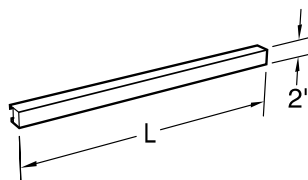


A06C0202□□N-00_0 2" High

Apron rails are available from 15" to 72" long in 1" increments. To specify apron rail, insert the length required into the part number.

Example: A06C020248N-0000 is a 48" long cold rolled steel apron rail.

Apron rails should be ordered in pairs to provide both front and rear work surface support.

Support Rails

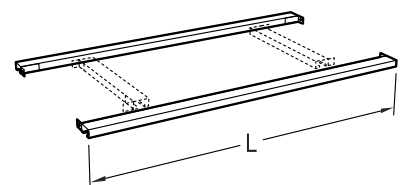

A81C0102□□N-0000

Support rails may be used individually or assembled to form a support framework. Rails are available from 12" to 72" long in 1" increments and are furnished with angles attached to each end. To specify rail, insert the length required in the part number.

Example: A81C010248N-0000 is a 48" long support rail.

Support rails should be ordered in pairs to provide both front and rear work surface support.

Steel Leg Assemblies & Pedestals

Specifications:



Legs are 2" square and all except the Top Mounted Leg are joined in pairs with a 2" high stretcher and a 4" high apron rail. Dimensions are as shown.

For aprons 42" long and longer, spreaders are available to connect the two leg assemblies. When a spreader is used, order the A08 or A28 leg assembly which are designed for use with a spreader.



Adjustable Height Legs are adjustable in 1" increments from 29" to 35" high.

End Leg Assemblies

Intermediate Leg Assemblies

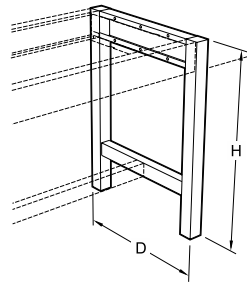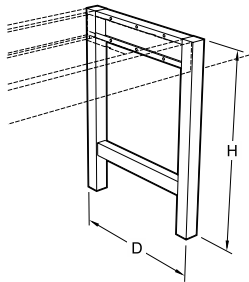

For use with Spreader

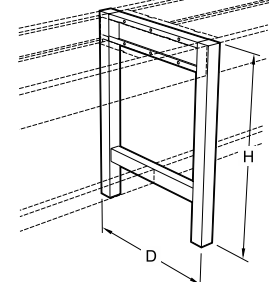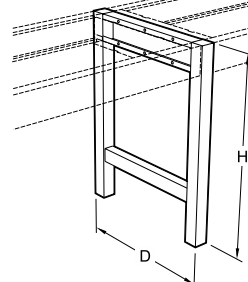

For use with Spreader

22³/₁₆" Deep

A07C292202N-00_0 29"H
A07C322202N-00_0 32"H
A07C352202N-00_0 35"H

22³/₁₆" Deep

A08C292202N-00_0 29"H
A08C322202N-00_0 32"H
A08C352202N-00_0 35"H

22³/₁₆" Deep

A27C292202N-00_0 29"H
A27C322202N-00_0 32"H
A27C352202N-00_0 35"H

22³/₁₆" Deep

A28C292202N-00_0 29"H
A28C322202N-00_0 32"H
A28C352202N-00_0 35"H

28³/₁₆" Deep

A07C292802N-00_0 29"H
A07C322802N-00_0 32"H
A07C352802N-00_0 35"H

28³/₁₆" Deep

A08C292802N-00_0 29"H
A08C322802N-00_0 32"H
A08C352802N-00_0 35"H

28³/₁₆" Deep

A27C292802N-00_0 29"H
A27C322802N-00_0 32"H
A27C352802N-00_0 35"H



28³/₁₆" Deep

A28C292802N-00_0 29"H
A28C322802N-00_0 32"H
A28C352802N-00_0 35"H



Each Leg Assembly adds 2" to the overall length of the apron assembly to which it is attached.

Adjustable Height End Leg Assemblies

Adjustable Height Intermediate Leg Assemblies

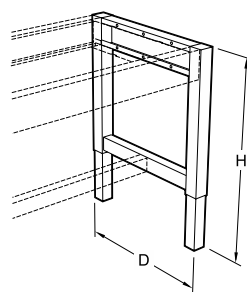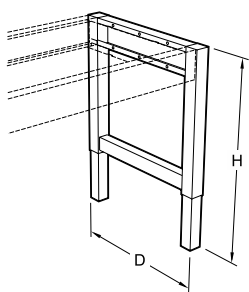

For use with Spreader

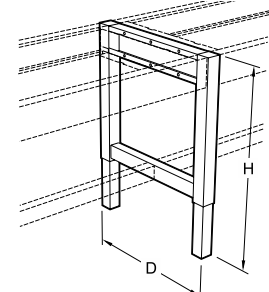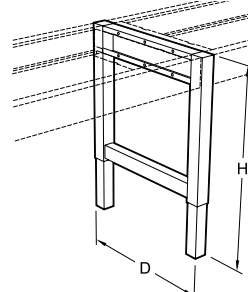

For use with Spreader

22³/₁₆" Deep

A57C352202N-00_0 Adj.
Adjustable Height 29"-35"

22³/₁₆" Deep

A58C352202N-00_0 Adj.
Adjustable Height 29"-35"

22³/₁₆" Deep

A77C352202N-00_0 Adj.
Adjustable Height 29"-35"

22³/₁₆" Deep

A78C352202N-00_0 Adj.
Adjustable Height 29"-35"

28³/₁₆" Deep

A57C352802N-00_0 Adj.
Adjustable Height 29"-35"

28³/₁₆" Deep

A58C352802N-00_0 Adj.
Adjustable Height 29"-35"

28³/₁₆" Deep

A77C352802N-00_0 Adj.
Adjustable Height 29"-35"

28³/₁₆" Deep

A78C352802N-00_0 Adj.
Adjustable Height 29"-35"

Steel Apron Rails & Panels

Specifications:


Knee Space panels are designed to close off the space between the floor and the bottom of the apron or countertop. They are sized for 29" high, 32" high, and 35" high cabinets. Each Knee Space Panel is made with a separate base rail, to allow the

removal of the panel without removing the base molding.

Knee Space Panels are furnished with two angles to attach to panel to the cabinets.

Pedestal Legs are 2" wide, depths and heights as shown. They are designed to be used with or without Box Aprons to form knee spaces at the end of assemblies.

Knee Space Panels

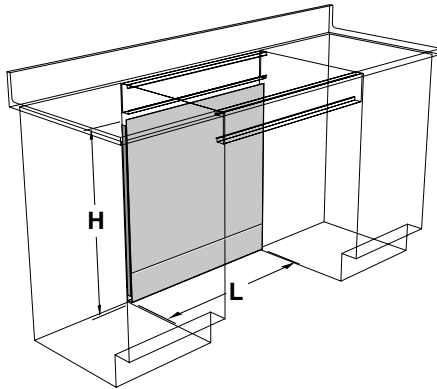

Knee space panels are available from 15" to 72" long in 1" increments. To specify knee space panel, insert the length required into the part number.

Example: A10C310048N-0000 is a 48" long knee space panel.

Sitting Height Knee Space Panels

Apron Part Number	Height	Length	for use with:
A10C2500 ☐ N-00_0	25"	15"-72"L	4" high apron
A10C2800 ☐ N-00_0	28"	15"-72"L	A80 support frame
A10C2900 ☐ N-00_0	29"	15"-72"L	no apron


ADA Height Knee Space Panels

Apron Part Number	Height	Length	for use with:
A10C2800 ☐ N-00_0	28"	15"-72"L	4" high apron
A10C3100 ☐ N-00_0	31"	15"-72"L	A80 support frame
A10C3200 ☐ N-00_0	32"	15"-72"L	no apron or frame

Standing Height Knee Space Panels

Apron Part Number	Height	Length	for use with:
A10C3100 ☐ N-00_0	31"	15"-72"L	4" high apron
A10C3400 ☐ N-00_0	34"	15"-72"L	A80 support frame
A10C3500 ☐ N-00_0	35"	15"-72"L	no apron or frame

Spreader

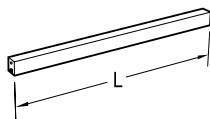

Spreader Part Number

A09C0202 ☐ N-00_0


Spreader lengths available in 1" increments from 15" to 72".

Example:
A09C020248N-00_0 48"L

Note:

Spreader to be used with leg assemblies only, not between cabinets.

Top Mounted Leg

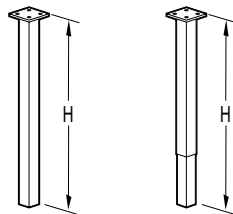


A15C290202N-00_0	29"H
A15C320202N-00_0	32"H
A15C350202N-00_0	35"H

A65C350202N-00_0 Adj.
Adjustable Height 29"-35"

2" square leg welded to 4" square steel plate

Wall Pedestal Leg

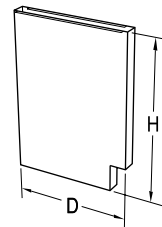

Left Side 22" Deep

A11C292202N-00_0	29"H
A11C322202N-00_0	32"H
A11C352202N-00_0	35"H

Left Side 28" Deep

A11C292802N-00_0	29"H
A11C322802N-00_0	32"H
A11C352802N-00_0	35"H


Right Side 22" Deep

A12C292202N-00_0	29"H
A12C322202N-00_0	32"H
A12C352202N-00_0	35"H

Right Side 28" Deep

A12C292802N-00_0	29"H
A12C322802N-00_0	32"H
A12C352802N-00_0	35"H

Center Pedestal Leg

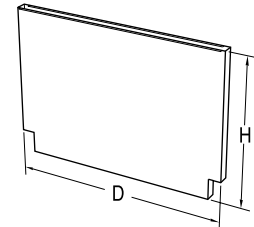

52" Deep


A13C295202N-00_0	29"H
A13C325202N-00_0	32"H
A13C355202N-00_0	35"H

Each Pedestal Leg adds 2" to the overall length of the assembly to which it is attached.

Pedestal Legs may be ordered with the toespace 7/8" deeper so that the front edge lines up with the doors and drawer fronts on adjacent Full Overlay cabinets. To specify the deeper toe space, add a **C** to the end of the part number.

Example: A11C352202N-0000C

Steel Wall Mounted Storage Cabinets

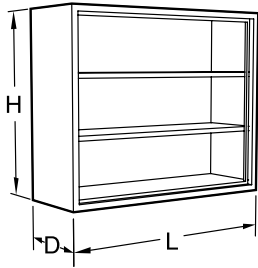


Specifications:

Open Wall Mounted Storage Cabinets are available in four heights, 24", 30", 36" and 48" and two depths; 12 $\frac{3}{4}$ " and 16". Lengths are as shown. All are furnished with one, two, or three adjustable shelves as shown.

Wall cabinets are designed to be mounted on hanging rails. [Click here for hanging rails.](#)

Open Storage

24" High

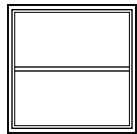


12 $\frac{3}{4}$ " Deep

W00C241315N-00_0	15"L
W00C241318N-00_0	18"L
W00C241321N-00_0	21"L
W00C241324N-00_0	24"L

16" Deep

W00C241615N-00_0	15"L
W00C241618N-00_0	18"L
W00C241621N-00_0	21"L
W00C241624N-00_0	24"L

24" High

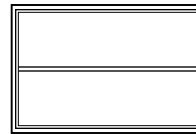


12 $\frac{3}{4}$ " Deep

W00C241330N-00_0	30"L
W00C241336N-00_0	36"L
W00C241342N-00_0	42"L
W00C241348N-00_0	48"L

16" Deep

W00C241630N-00_0	30"L
W00C241636N-00_0	36"L
W00C241642N-00_0	42"L
W00C241648N-00_0	48"L

30" High

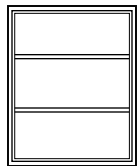

12 $\frac{3}{4}$ " Deep

W00C301315N-00_0	15"L
W00C301318N-00_0	18"L
W00C301321N-00_0	21"L
W00C301324N-00_0	24"L

16" Deep

W00C301615N-00_0	15"L
W00C301618N-00_0	18"L
W00C301621N-00_0	21"L
W00C301624N-00_0	24"L

30" High




12 $\frac{3}{4}$ " Deep

W00C301330N-00_0	30"L
W00C301336N-00_0	36"L
W00C301342N-00_0	42"L
W00C301348N-00_0	48"L

16" Deep

W00C301630N-00_0	30"L
W00C301636N-00_0	36"L
W00C301642N-00_0	42"L
W00C301648N-00_0	48"L

36" High

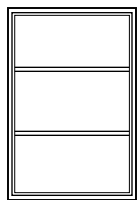


12 $\frac{3}{4}$ " Deep

W00C361315N-00_0	15"L
W00C361318N-00_0	18"L
W00C361321N-00_0	21"L
W00C361324N-00_0	24"L

16" Deep

W00C361615N-00_0	15"L
W00C361618N-00_0	18"L
W00C361621N-00_0	21"L
W00C361624N-00_0	24"L

36" High

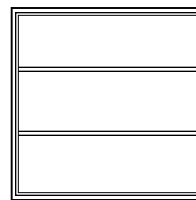


12 $\frac{3}{4}$ " Deep

W00C361330N-00_0	30"L
W00C361336N-00_0	36"L
W00C361342N-00_0	42"L
W00C361348N-00_0	48"L

16" Deep

W00C361630N-00_0	30"L
W00C361636N-00_0	36"L
W00C361642N-00_0	42"L
W00C361648N-00_0	48"L

48" High

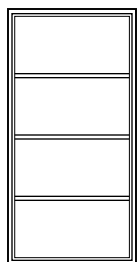


12 $\frac{3}{4}$ " Deep

W00C481315N-00_0	15"L
W00C481318N-00_0	18"L
W00C481321N-00_0	21"L
W00C481324N-00_0	24"L

16" Deep

W00C481615N-00_0	15"L
W00C481618N-00_0	18"L
W00C481621N-00_0	21"L
W00C481624N-00_0	24"L

48" High




12 $\frac{3}{4}$ " Deep

W00C481330N-00_0	30"L
W00C481336N-00_0	36"L
W00C481342N-00_0	42"L
W00C481348N-00_0	48"L

16" Deep

W00C481630N-00_0	30"L
W00C481636N-00_0	36"L
W00C481642N-00_0	42"L
W00C481648N-00_0	48"L

Steel Wall Mounted Storage Cabinets

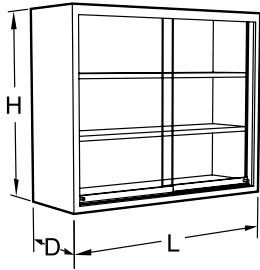

Specifications:

Plate Glass Door Wall Mounted Storage Cabinets are available in three heights, 24", 30", 36" and 48" and two depths; 12³/₄" and 16". Lengths are as shown. Cabinets are furnished with one, two, or three adjustable shelves as shown.


Doors are 1/4" float glass and are mounted on friction-free roller glides.

Wall cabinets are designed to be mounted on hanging rails. [See below.](#)

Sliding Plate Glass Doors — Clear Glass

Wall Mounted Cabinet Hanging Rails

24" High

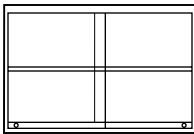


12³/₄" Deep

W10C241330N-00_0	30"L
W10C241336N-00_0	36"L
W10C241342N-00_0	42"L
W10C241348N-00_0	48"L

16" Deep

W10C241630N-00_0	30"L
W10C241636N-00_0	36"L
W10C241642N-00_0	42"L
W10C241648N-00_0	48"L

30" High

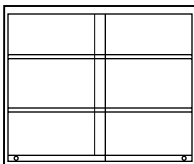


12³/₄" Deep

W10C301330N-00_0	30"L
W10C301336N-00_0	36"L
W10C301342N-00_0	42"L
W10C301348N-00_0	48"L

16" Deep

W10C301630N-00_0	30"L
W10C301636N-00_0	36"L
W10C301642N-00_0	42"L
W10C301648N-00_0	48"L

36" High

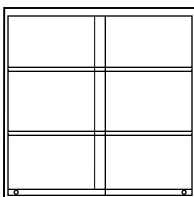


12³/₄" Deep

W10C361330N-00_0	30"L
W10C361336N-00_0	36"L
W10C361342N-00_0	42"L
W10C361348N-00_0	48"L

16" Deep

W10C361630N-00_0	30"L
W10C361636N-00_0	36"L
W10C361642N-00_0	42"L
W10C361648N-00_0	48"L

48" High

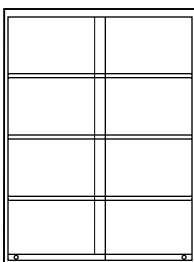


12³/₄" Deep

W10C481330N-00_0	30"L
W10C481336N-00_0	36"L
W10C481342N-00_0	42"L
W10C481348N-00_0	48"L

16" Deep


W10C481630N-00_0	30"L
W10C481636N-00_0	36"L
W10C481642N-00_0	42"L
W10C481648N-00_0	48"L

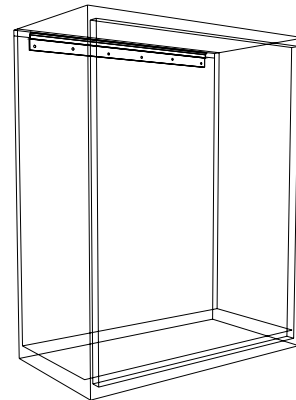

L01M020109N	for 15" wall cabinet
L01M020112N	for 18" wall cabinet
L01M020115N	for 21" wall cabinet
L01M020118N	for 24" wall cabinet
L01M020124N	for 30" wall cabinet
L01M020130N	for 36" wall cabinet
L01M020136N	for 42" wall cabinet
L01M020142N	for 48" wall cabinet

Formed of heavy gauge cold-rolled steel.
Designed to capture lip on top rear of wall cabinets for easy installation.

Steel Wall Mounted Storage Cabinets

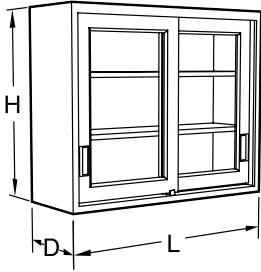

Specifications:

Sliding Door Wall Mounted Storage Cabinets are available in four heights, 24", 30", 36" and 48" and two depths; 12 $\frac{3}{4}$ " and 16". Lengths are as shown. Cabinets are glazed with $\frac{1}{8}$ " float glass unless otherwise specified and are furnished with one, two or three adjustable shelves as shown. Doors are mounted on friction-free roller glides.


Sliding Door Wall Mounted Storage Cabinets are available in Door Style 1 and 6 only.

Wall cabinets are designed to be mounted on hanging rails. [Click here for hanging rails.](#)

Sliding Glazed Doors

Sliding Panel Doors

24" High

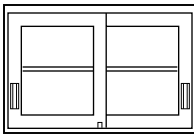


12 $\frac{3}{4}$ " Deep

W41C241330N-00_0	30"L
W41C241336N-00_0	36"L
W41C241342N-00_0	42"L
W41C241348N-00_0	48"L

16" Deep

W41C241630N-00_0	30"L
W41C241636N-00_0	36"L
W41C241642N-00_0	42"L
W41C241648N-00_0	48"L

24" High

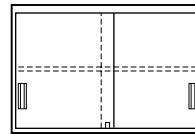


12 $\frac{3}{4}$ " Deep

W12C241330N-00_0	30"L
W12C241336N-00_0	36"L
W12C241342N-00_0	42"L
W12C241348N-00_0	48"L

16" Deep

W12C241630N-00_0	30"L
W12C241636N-00_0	36"L
W12C241642N-00_0	42"L
W12C241648N-00_0	48"L

30" High

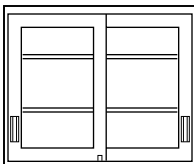


12 $\frac{3}{4}$ " Deep

W41C301330N-00_0	30"L
W41C301336N-00_0	36"L
W41C301342N-00_0	42"L
W41C301348N-00_0	48"L

16" Deep

W41C301630N-00_0	30"L
W41C301636N-00_0	36"L
W41C301642N-00_0	42"L
W41C301648N-00_0	48"L

30" High

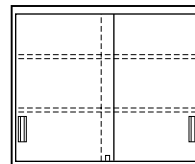


12 $\frac{3}{4}$ " Deep

W12C301330N-00_0	30"L
W12C301336N-00_0	36"L
W12C301342N-00_0	42"L
W12C301348N-00_0	48"L

16" Deep

W12C301630N-00_0	30"L
W12C301636N-00_0	36"L
W12C301642N-00_0	42"L
W12C301648N-00_0	48"L

36" High

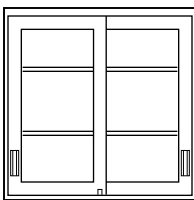


12 $\frac{3}{4}$ " Deep

W41C361330N-00_0	30"L
W41C361336N-00_0	36"L
W41C361342N-00_0	42"L
W41C361348N-00_0	48"L

16" Deep

W41C361630N-00_0	30"L
W41C361636N-00_0	36"L
W41C361642N-00_0	42"L
W41C361648N-00_0	48"L

36" High

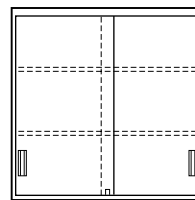


12 $\frac{3}{4}$ " Deep

W12C361330N-00_0	30"L
W12C361336N-00_0	36"L
W12C361342N-00_0	42"L
W12C361348N-00_0	48"L

16" Deep

W12C361630N-00_0	30"L
W12C361636N-00_0	36"L
W12C361642N-00_0	42"L
W12C361648N-00_0	48"L

48" High

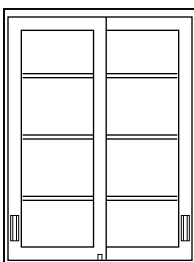


12 $\frac{3}{4}$ " Deep

W41C481330N-00_0	30"L
W41C481336N-00_0	36"L
W41C481342N-00_0	42"L
W41C481348N-00_0	48"L

16" Deep

W41C481630N-00_0	30"L
W41C481636N-00_0	36"L
W41C481642N-00_0	42"L
W41C481648N-00_0	48"L

48" High

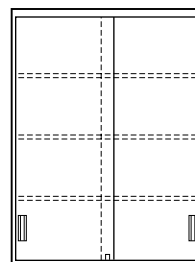


12 $\frac{3}{4}$ " Deep

W12C481330N-00_0	30"L
W12C481336N-00_0	36"L
W12C481342N-00_0	42"L
W12C481348N-00_0	48"L

16" Deep

W12C481630N-00_0	30"L
W12C481636N-00_0	36"L
W12C481642N-00_0	42"L
W12C481648N-00_0	48"L

Steel Wall Mounted Storage Cabinets

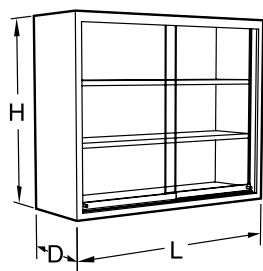

Specifications:


Double Sided Plate Glass Pass-Thru Cabinets are available 24", 30", and 36" high, 12³/₄" and 16" deep and in lengths shown. Double Sided Pass-Thru Cabinets have doors on both sides and are designed to be mounted on center bench reagent racks.

Doors are 1/4" float glass and are mounted on friction-free roller glides.

Wall cabinets are designed to be mounted on hanging rails. [Click here for hanging rails.](#)

Double Sided Pass-Thru Cabinet – Clear Glass

24" High

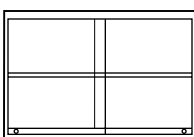


12³/₄" Deep

W13C241330N-00_0	30"L
W13C241336N-00_0	36"L
W13C241342N-00_0	42"L
W13C241348N-00_0	48"L

16" Deep

W13C241630N-00_0	30"L
W13C241636N-00_0	36"L
W13C241642N-00_0	42"L
W13C241648N-00_0	48"L

30" High

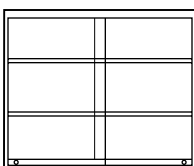


12³/₄" Deep

W13C301330N-00_0	30"L
W13C301336N-00_0	36"L
W13C301342N-00_0	42"L
W13C301348N-00_0	48"L

16" Deep

W13C301630N-00_0	30"L
W13C301636N-00_0	36"L
W13C301642N-00_0	42"L
W13C301648N-00_0	48"L

36" High

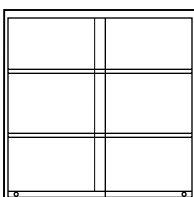


12³/₄" Deep

W13C361330N-00_0	30"L
W13C361336N-00_0	36"L
W13C361342N-00_0	42"L
W13C361348N-00_0	48"L

16" Deep

W13C361630N-00_0	30"L
W13C361636N-00_0	36"L
W13C361642N-00_0	42"L
W13C361648N-00_0	48"L

Steel Wall Mounted Storage Cabinets

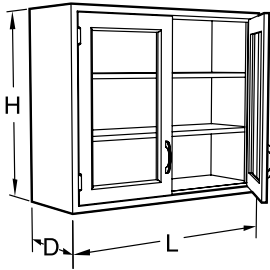

Specifications:


Swinging Glazed Door Wall Mounted Storage Cabinets are available in four heights, 24", 30", 36" and 48" and two depths; 12 $\frac{3}{4}$ " and 16". Lengths are as shown. Cabinets are glazed with $\frac{1}{8}$ " float glass unless otherwise specified and are furnished with one, two or three adjustable shelves as shown.

Note: Last four digits of part number designate door style, material, and pull selection. [Click here](#) for full explanation. Beveled Edge, Door Style 7 is not available for wall cabinets.

Wall cabinets are designed to be mounted on hanging rails. [Click here for hanging rails.](#)


Swinging Glazed Doors

24" High

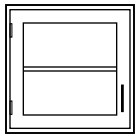


Left Hand	Right Hand	12 $\frac{3}{4}$ " D
W50C241315L- ----	W50C241315R- ----	15"L
W50C241318L- ----	W50C241318R- ----	18"L
W50C241321L- ----	W50C241321R- ----	21"L
W50C241324L- ----	W50C241324R- ----	24"L
16" D		
W50C241615L- ----	W50C241615R- ----	15"L
W50C241618L- ----	W50C241618R- ----	18"L
W50C241621L- ----	W50C241621R- ----	21"L
W50C241624L- ----	W50C241624R- ----	24"L

24" High

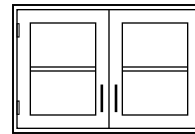


12 $\frac{3}{4}$ " Deep	
W50C241330D- ----	30"L
W50C241336D- ----	36"L
W50C241342D- ----	42"L
W50C241348D- ----	48"L
16" Deep	
W50C241630D- ----	30"L
W50C241636D- ----	36"L
W50C241642D- ----	42"L
W50C241648D- ----	48"L

30" High

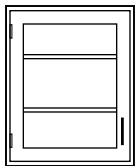


Left Hand	Right Hand	12 $\frac{3}{4}$ " D
W50C301315L- ----	W50C301315R- ----	15"L
W50C301318L- ----	W50C301318R- ----	18"L
W50C301321L- ----	W50C301321R- ----	21"L
W50C301324L- ----	W50C301324R- ----	24"L
16" D		
W50C301615L- ----	W50C301615R- ----	15"L
W50C301618L- ----	W50C301618R- ----	18"L
W50C301621L- ----	W50C301621R- ----	21"L
W50C301624L- ----	W50C301624R- ----	24"L

30" High

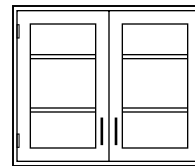


12 $\frac{3}{4}$ " Deep	
W50C301330D- ----	30"L
W50C301336D- ----	36"L
W50C301342D- ----	42"L
W50C301348D- ----	48"L
16" Deep	
W50C301630D- ----	30"L
W50C301636D- ----	36"L
W50C301642D- ----	42"L
W50C301648D- ----	48"L

36" High

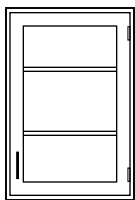


Left Hand	Right Hand	12 $\frac{3}{4}$ " D
W50C361315L- ----	W50C361315R- ----	15"L
W50C361318L- ----	W50C361318R- ----	18"L
W50C361321L- ----	W50C361321R- ----	21"L
W50C361324L- ----	W50C361324R- ----	24"L
16" D		
W50C361615L- ----	W50C361615R- ----	15"L
W50C361618L- ----	W50C361618R- ----	18"L
W50C361621L- ----	W50C361621R- ----	21"L
W50C361624L- ----	W50C361624R- ----	24"L

36" High

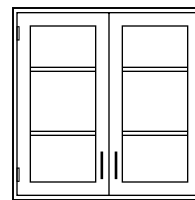


12 $\frac{3}{4}$ " Deep	
W50C361330D- ----	30"L
W50C361336D- ----	36"L
W50C361342D- ----	42"L
W50C361348D- ----	48"L
16" Deep	
W50C361630D- ----	30"L
W50C361636D- ----	36"L
W50C361642D- ----	42"L
W50C361648D- ----	48"L

48" High

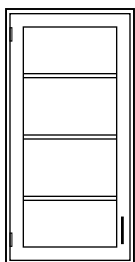


Left Hand	Right Hand	12 $\frac{3}{4}$ " D
W50C481315L- ----	W50C481315R- ----	15"L
W50C481318L- ----	W50C481318R- ----	18"L
W50C481321L- ----	W50C481321R- ----	21"L
W50C481324L- ----	W50C481324R- ----	24"L
16" D		
W50C481615L- ----	W50C481615R- ----	15"L
W50C481618L- ----	W50C481618R- ----	18"L
W50C481621L- ----	W50C481621R- ----	21"L
W50C481624L- ----	W50C481624R- ----	24"L

48" High

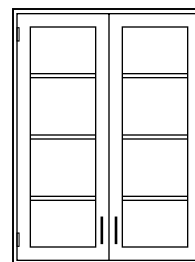

12 $\frac{3}{4}$ " Deep	
W50C481330D- ----	30"L
W50C481336D- ----	36"L
W50C481342D- ----	42"L
W50C481348D- ----	48"L
16" Deep	
W50C481630D- ----	30"L
W50C481636D- ----	36"L
W50C481642D- ----	42"L
W50C481648D- ----	48"L

Steel Wall Mounted Storage Cabinets

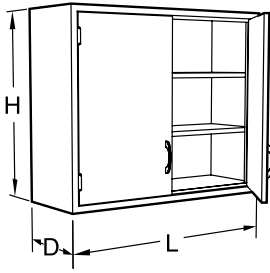



Specifications:



Swinging Panel Door Wall Mounted Storage Cabinets are available in four heights, 24", 30", 36" and 48" and two depths; 12 $\frac{3}{4}$ " and 16". Lengths are as shown. Cabinets are furnished with one, two or three adjustable shelves as shown.



Note: Last four digits of part number designate door style, material, and pull selection. [Click here](#) for full explanation. Beveled Edge, Door Style 7 is not available for wall cabinets.



Wall cabinets are designed to be mounted on hanging rails. [Click here for hanging rails.](#)

Swinging Panel Doors

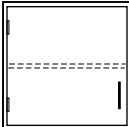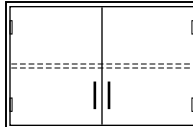
	Left Hand		Right Hand		12¾" D		12¾" Deep	
	W25C241315L- _____		W25C241315R- _____		15"L		W25C241330D- _____ 30"L	
	W25C241318L- _____		W25C241318R- _____		18"L		W25C241336D- _____ 36"L	
	W25C241321L- _____		W25C241321R- _____		21"L		W25C241342D- _____ 42"L	
	W25C241324L- _____		W25C241324R- _____		24"L		W25C241348D- _____ 48"L	
						16" D		
	W25C241615L- _____		W25C241615R- _____		15"L		16" Deep	
	W25C241618L- _____		W25C241618R- _____		18"L		W25C241630D- _____ 30"L	
	W25C241621L- _____		W25C241621R- _____		21"L		W25C241636D- _____ 36"L	
	W25C241624L- _____		W25C241624R- _____		24"L		W25C241642D- _____ 42"L	
							W25C241648D- _____ 48"L	

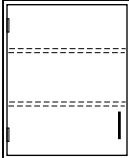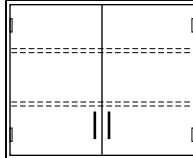	Left Hand		Right Hand		12¾" D		12¾" Deep	
	W25C301315L- _____		W25C301315R- _____		15"L		W25C301330D- _____ 30"L	
	W25C301318L- _____		W25C301318R- _____		18"L		W25C301336D- _____ 36"L	
	W25C301321L- _____		W25C301321R- _____		21"L		W25C301342D- _____ 42"L	
	W25C301324L- _____		W25C301324R- _____		24"L		W25C301348D- _____ 48"L	
						16" D		
	W25C301615L- _____		W25C301615R- _____		15"L		16" Deep	
	W25C301618L- _____		W25C301618R- _____		18"L		W25C301630D- _____ 30"L	
	W25C301621L- _____		W25C301621R- _____		21"L		W25C301636D- _____ 36"L	
	W25C301624L- _____		W25C301624R- _____		24"L		W25C301642D- _____ 42"L	
							W25C301648D- _____ 48"L	

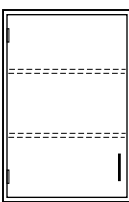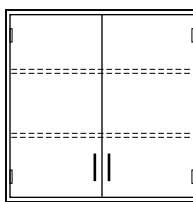	Left Hand		Right Hand		12¾" D		12¾" Deep	
	W25C361315L- _____		W25C361315R- _____		15"L		W25C361330D- _____ 30"L	
	W25C361318L- _____		W25C361318R- _____		18"L		W25C361336D- _____ 36"L	
	W25C361321L- _____		W25C361321R- _____		21"L		W25C361342D- _____ 42"L	
	W25C361324L- _____		W25C361324R- _____		24"L		W25C361348D- _____ 48"L	
						16" D		
	W25C361615L- _____		W25C361615R- _____		15"L		16" Deep	
	W25C361618L- _____		W25C361618R- _____		18"L		W25C361630D- _____ 30"L	
	W25C361621L- _____		W25C361621R- _____		21"L		W25C361636D- _____ 36"L	
	W25C361624L- _____		W25C361624R- _____		24"L		W25C361642D- _____ 42"L	
							W25C361648D- _____ 48"L	

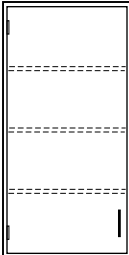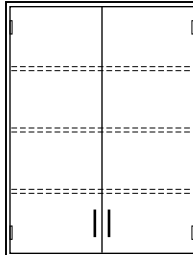	Left Hand		Right Hand		12¾" D		12¾" Deep	
	W25C481315L- _____		W25C481315R- _____		15"L		W25C481330D- _____ 30"L	
	W25C481318L- _____		W25C481318R- _____		18"L		W25C481336D- _____ 36"L	
	W25C481321L- _____		W25C481321R- _____		21"L		W25C481342D- _____ 42"L	
	W25C481324L- _____		W25C481324R- _____		24"L		W25C481348D- _____ 48"L	
						16" D		
	W25C481615L- _____		W25C481615R- _____		15"L		16" Deep	
	W25C481618L- _____		W25C481618R- _____		18"L		W25C481630D- _____ 30"L	
	W25C481621L- _____		W25C481621R- _____		21"L		W25C481636D- _____ 36"L	
	W25C481624L- _____		W25C481624R- _____		24"L		W25C481642D- _____ 42"L	
							W25C481648D- _____ 48"L	

Steel Wall Mounted Storage Cabinets

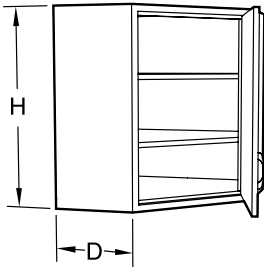

















Specifications:


Corner Wall Mounted Storage Cabinets are available in four heights, 24", 30", 36" and 48" and two depths, 12 $\frac{3}{4}$ " and 16" and with either glazed or panel swinging doors. Cabinets are furnished with one, two or three adjustable shelves as shown.

Note: Last four digits of part number designate door style, material, and pull selection. [Click here](#) for full explanation. Beveled Edge, Door Style 7 is not available for wall cabinets.


Wall cabinets are designed to be mounted on hanging rails. [Click here for hanging rails.](#)

Swinging Door Corner Cabinets

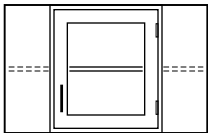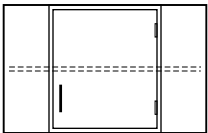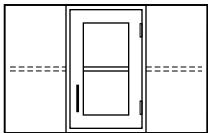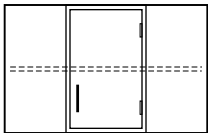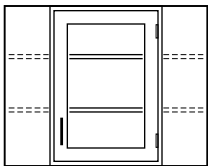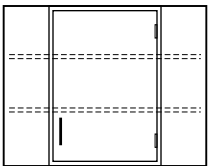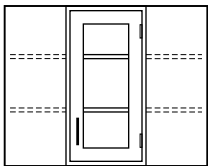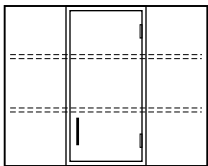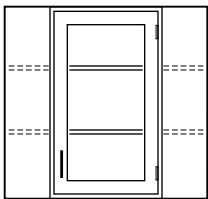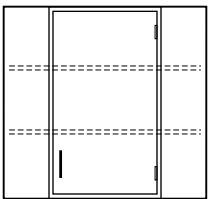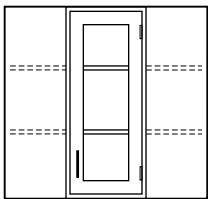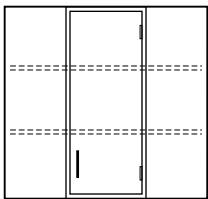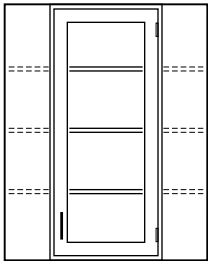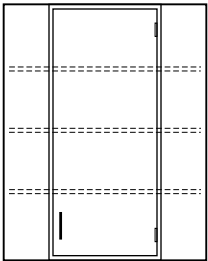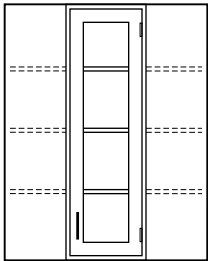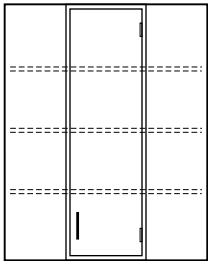Glazed Door	Panel Door	Glazed Door	Panel Door
24" High	12$\frac{3}{4}$" D	24" High	16" D
			
W51C241320R-____	W26C241320R-____	W51C241615R-____	W26C241615R-____
30" High	12$\frac{3}{4}$" D	30" High	16" D
			
W51C301320R-____	W26C301320R-____	W51C301615R-____	W26C301615R-____
36" High	12$\frac{3}{4}$" D	36" High	16" D
			
W51C361320R-____	W26C361320R-____	W51C361615R-____	W26C361615R-____
48" High	12$\frac{3}{4}$" D	48" High	16" D
			
W51C481320R-____	W26C481320R-____	W51C481615R-____	W26C481615R-____

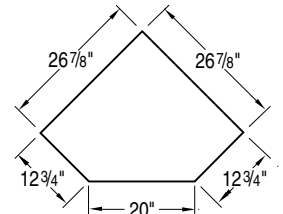


PLAN
VIEW OF
W51C241320R
W51C301320R
W51C361320R
W51C481320R
W26C241320R
W26C301320R
W26C361320R
W26C481320R

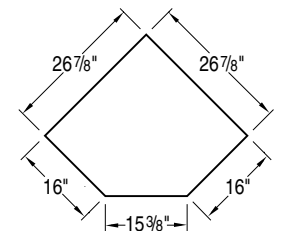

PLAN
VIEW OF
W51C241615R
W51C301615R
W51C361615R
W51C481615R
W26C241615R
W26C301615R
W26C361615R
W26C481615R

Steel Full Height Storage Cabinets

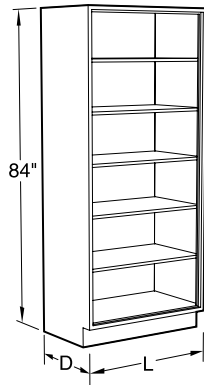


Specifications:

Full Height Storage Cabinets are 84" high and 12¾", 16", or 22" deep. Lengths as shown. They are furnished with five adjustable shelves.

CAUTION: Full height cabinets may present a potential tipping over problem. Proper installation requires these cabinets to be anchored to a wall or other cabinets to improve stability and ensure smooth door movement.

NOTE: Full height storage cabinets are not available in Overlay – Beveled Edge Door Style 7.

Open Storage



12¾" Deep


S00C841318N-00_0	18"L
S00C841324N-00_0	24"L

16" Deep

S00C841618N-00_0	18"L
S00C841624N-00_0	24"L

22" Deep

S00C842218N-00_0	18"L
S00C842224N-00_0	24"L



12¾" Deep

S00C841330N-00_0	30"L
S00C841336N-00_0	36"L
S00C841342N-00_0	42"L
S00C841348N-00_0	48"L


16" Deep

S00C841630N-00_0	30"L
S00C841636N-00_0	36"L
S00C841642N-00_0	42"L
S00C841648N-00_0	48"L

22" Deep

S00C842230N-00_0	30"L
S00C842236N-00_0	36"L
S00C842242N-00_0	42"L
S00C842248N-00_0	48"L

Sliding Glazed Doors

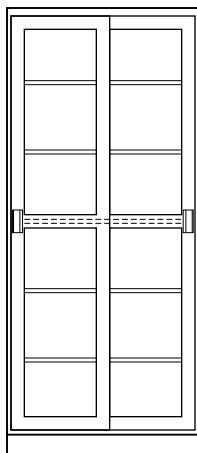

12¾" Deep

S41C841330N-00_0	30"L
S41C841336N-00_0	36"L
S41C841342N-00_0	42"L
S41C841348N-00_0	48"L


16" Deep

S41C841630N-00_0	30"L
S41C841636N-00_0	36"L
S41C841642N-00_0	42"L
S41C841648N-00_0	48"L

22" Deep

S41C842230N-00_0	30"L
S41C842236N-00_0	36"L
S41C842242N-00_0	42"L
S41C842248N-00_0	48"L

Sliding Panel Doors

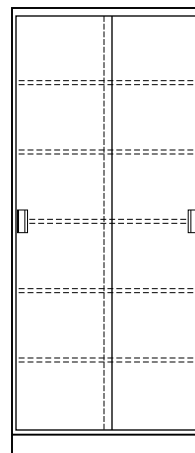

12¾" Deep

S12C841330N-00_0	30"L
S12C841336N-00_0	36"L
S12C841342N-00_0	42"L
S12C841348N-00_0	48"L


16" Deep

S12C841630N-00_0	30"L
S12C841636N-00_0	36"L
S12C841642N-00_0	42"L
S12C841648N-00_0	48"L

22" Deep

S12C842230N-00_0	30"L
S12C842236N-00_0	36"L
S12C842242N-00_0	42"L
S12C842248N-00_0	48"L

Steel Full Height Storage Cabinets

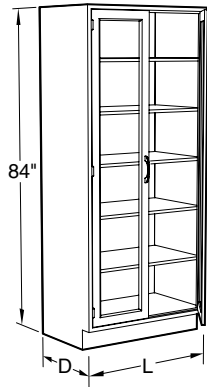


Specifications:

Full Height Storage Cabinets are 84" high and 12 $\frac{3}{4}$ ", 16", or 22" deep. Lengths as shown. They are furnished with five adjustable shelves.


CAUTION: Full height cabinets may present a potential tipping over problem. Proper installation requires these cabinets to be anchored to a wall or other cabinets to improve stability and ensure smooth door movement.

NOTE: Full height storage cabinets are not available in Overlay – Beveled Edge Door Style 7.

Swinging Glazed Doors

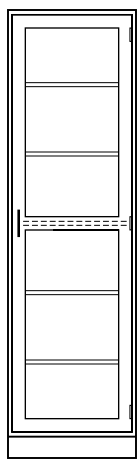


Left Hand	Right Hand	12 $\frac{3}{4}$ " D
S50C841318L-_____	S50C841318R-_____	18"L
S50C841324L-_____	S50C841324R-_____	24"L
		16" D
S50C841618L-_____	S50C841618R-_____	18"L
S50C841624L-_____	S50C841624R-_____	24"L
		22" D
S50C842218L-_____	S50C842218R-_____	18"L
S50C842224L-_____	S50C842224R-_____	24"L

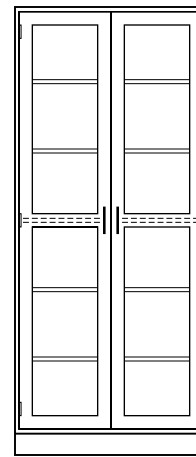


12 $\frac{3}{4}$ " Deep	
S50C841330D-_____	30"L
S50C841336D-_____	36"L
S50C841342D-_____	42"L
S50C841348D-_____	48"L
16" Deep	
S50C841630D-_____	30"L
S50C841636D-_____	36"L
S50C841642D-_____	42"L
S50C841648D-_____	48"L
22" Deep	
S50C842230D-_____	30"L
S50C842236D-_____	36"L
S50C842242D-_____	42"L
S50C842248D-_____	48"L

Swinging Panel Doors

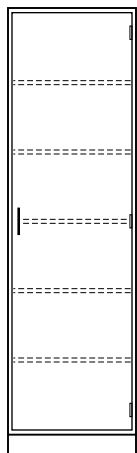


Left Hand	Right Hand	12 $\frac{3}{4}$ " D
S25C841318L-_____	S25C841318R-_____	18"L
S25C841324L-_____	S25C841324R-_____	24"L
		16" D
S25C841618L-_____	S25C841618R-_____	18"L
S25C841624L-_____	S25C841624R-_____	24"L
		22" D
S25C842218L-_____	S25C842218R-_____	18"L
S25C842224L-_____	S25C842224R-_____	24"L

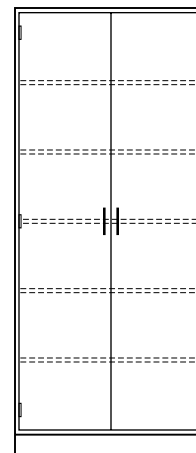

12 $\frac{3}{4}$ " Deep	
S25C841330D-_____	30"L
S25C841336D-_____	36"L
S25C841342D-_____	42"L
S25C841348D-_____	48"L
16" Deep	
S25C841630D-_____	30"L
S25C841636D-_____	36"L
S25C841642D-_____	42"L
S25C841648D-_____	48"L
22" Deep	
S25C842230D-_____	30"L
S25C842236D-_____	36"L
S25C842242D-_____	42"L
S25C842248D-_____	48"L

Steel Full Height Storage Cabinets

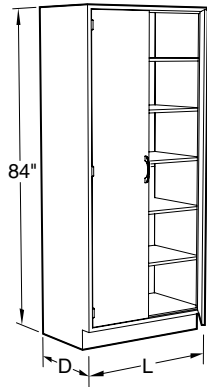

Specifications:


Microscope Storage Cabinets are 84" high and 16" or 22" deep. Lengths as shown. They are furnished with five adjustable shelves.

Full height swinging doors are hung on three hinges per door.

CAUTION: Full height cabinets may present a potential tipping over problem. Proper installation requires these cabinets to be anchored to a wall or other cabinets to improve stability and ensure smooth door movement.

NOTE: Full height storage cabinets are not available in Overlay – Beveled Edge Door Style 7.

Microscope Storage Cabinets

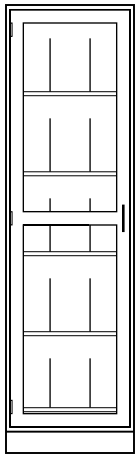


16" Deep

S51C841624L-____ 24"L

22" Deep

S51C842224L-____ 24"L

Provides 15 storage compartments

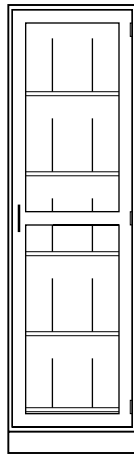


16" Deep

S51C841624R-____ 24"L

22" Deep

S51C842224R-____ 24"L

Provides 15 storage compartments

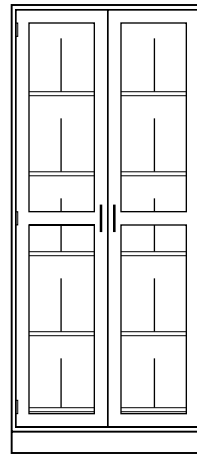


16" Deep

S51C841636D-____ 36"L

22" Deep

S51C842236D-____ 36"L

Provides 20 storage compartments

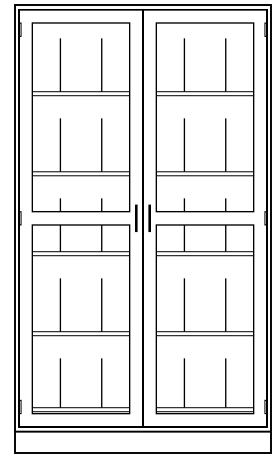

16" Deep


S51C841648D-____ 48"L

22" Deep

S51C842248D-____ 48"L

Provides 30 storage compartments

Steel Full Height Storage Cabinets

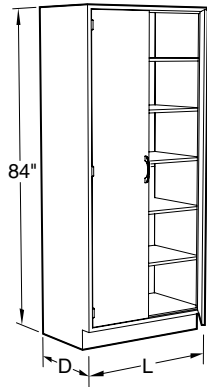

Specifications:


Full Height Storage Cabinets are 84" high and 12 $\frac{3}{4}$ ", 16", or 22" deep. Lengths as shown. Split door cabinets are furnished with a fixed center shelf and four adjustable shelves.

Wardrobe Storage Cabinets are 84" high and 22" deep. Lengths as shown. They are provided with a hat shelf and a horizontal bar for hanging garments.

CAUTION: Full height cabinets may present a potential tipping over problem. Proper installation requires these cabinets to be anchored to a wall or other cabinets to improve stability and ensure smooth door movement.

NOTE: Full height storage cabinets are not available in Overlay – Beveled Edge Door Style 7.

Split Swinging Glazed Doors

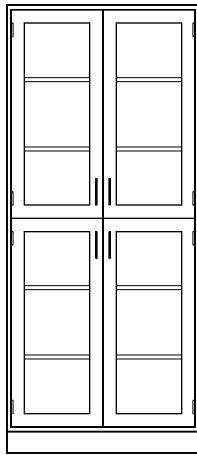

12 $\frac{3}{4}$ " Deep

S60C841330N-	30"L
S60C841336N-	36"L
S60C841342N-	42"L
S60C841348N-	48"L


16" Deep

S60C841630N-	30"L
S60C841636N-	36"L
S60C841642N-	42"L
S60C841648N-	48"L

22" Deep

S60C842230N-	30"L
S60C842236N-	36"L
S60C842242N-	42"L
S60C842248N-	48"L

Split Swinging Panel Doors

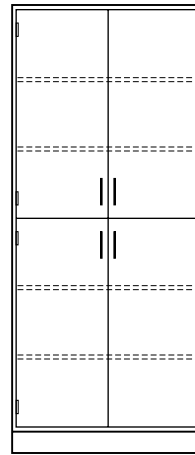

12 $\frac{3}{4}$ " Deep

S35C841330N-	30"L
S35C841336N-	36"L
S35C841342N-	42"L
S35C841348N-	48"L


16" Deep

S35C841630N-	30"L
S35C841636N-	36"L
S35C841642N-	42"L
S35C841648N-	48"L

22" Deep

S35C842230N-	30"L
S35C842236N-	36"L
S35C842242N-	42"L
S35C842248N-	48"L

Wardrobe Storage Cabinets

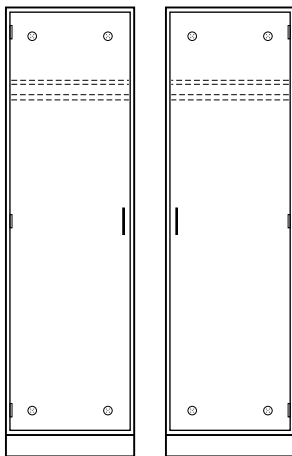


Left Hand — 22" Deep

S40C842218L-	18"L
S40C842224L-	24"L

Right Hand — 22" Deep


S40C842218R-	18"L
S40C842224R-	24"L

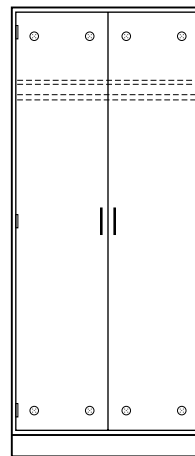

22" Deep

S40C842236D-	36"L
S40C842248D-	48"L

Steel Full Height Storage Cabinets

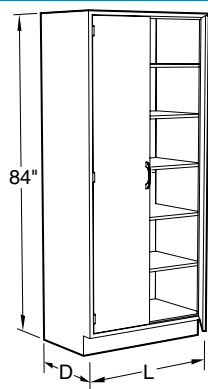

Specifications:

Full Height Storage Cabinets are 84" high. Lengths as shown.


Acid Storage Cabinets are furnished with a polyethylene corrosion resistant lining and five full-width corrosion resistant shelves. Ventilation is provided by louvers at the bottom of the door and a 3½" O.D. vent outlet at the top of the cabinet.

Full height swinging doors are hung on three hinges per door.

CAUTION: Full height cabinets may present a potential tipping over problem. Proper installation requires these cabinets to be anchored to a wall or other cabinets to improve stability and ensure smooth door movement.

NOTE: Full height storage cabinets are not available in Overlay – Beveled Edge Door Style 7.

Acid Storage Cabinets

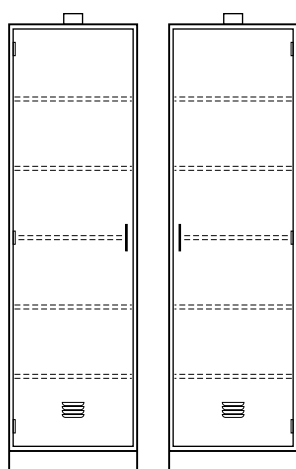

Left Hand — 16" Deep

S70C841624L-____ 24"L

Left Hand — 22" Deep


S70C842224L-____ 24"L

Right Hand — 16" Deep

S70C841624R-____ 24"L

Right Hand — 22" Deep

S70C842224R-____ 24"L

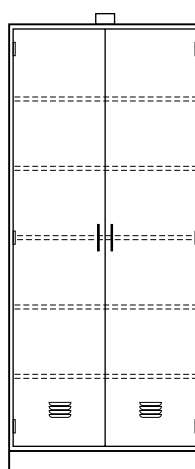

16" Deep

S70C841636D-____ 36"L

S70C841648D-____ 48"L

22" Deep

S70C842236D-____ 36"L

S70C842248D-____ 48"L

Safety Cabinets

We understand your need for stringent safety measures.

For over 100 years, we have built the world's finest laboratory furniture including steel, wood, and stainless steel casework, the Alpha family of flexible and mobile casework, Supreme Air fume hoods, bio-safety cabinets, Kemresin work surfaces, Alpha Columns and Enterprise. Kewaunee provides total Laboratory Solutions, Innovative Design, and Technological Expertise for labs throughout the world.



Dual vents with 2" threaded fittings placed one on the top left and one on the bottom right with built-in flame arrestor

Fully welded construction, increasing fire protection by reducing gapping

Continuous full height piano hinges

Environmentally friendly, electrostatically applied powder coat finish that meets SEFA 8 chemical requirement

Shelves are easily adjustable in 2 1/4" increments for versatile storage. The shelves are made of heavy gauge galvanized steel for supporting the storage capacity up to 90 gallons


18 gauge double walled steel construction with 1 1/2" air space

Adjustable leveling glide for uneven flooring

For convenience of loading the cabinet, each door can be latched open with a fusible link. In case of fire, the link melts at 165° closing the doors automatically

Self-closing, self latching, synchronized doors, ensure that the left door closes prior to the right door (For self-closing models)

Trilingual warning labels

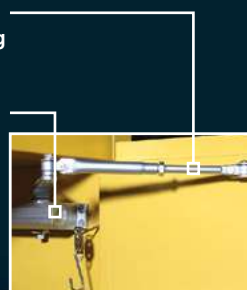

Removable panel to service three point latching mechanism



2" high rim on cabinet bottom to contain leaks and spills



Additional security provided by keyed paddle latch



Non-sparking, three-point latching mechanism

Safety Cabinet Features & Specifications

Safety Cabinet Features

- Kewaunee Safety Cabinets are produced in accordance with NFPA 30 Chapter 9.5.3 and OSHA 29 CFR 1910.106 regulations
- Kewaunee Safety Cabinets 12 gallons and greater are FM Global tested and approved
- Kewaunee Safety Cabinets are the first FM approved safety cabinets in India
- 2" deep liquid tight containment sump covering the entire floor of the cabinet
- 18-gauge double walled steel construction, along with 1½" air spaces between the walls, store flammable liquids and waste securely and safely
- Fully welded construction, increasing fire protection by reducing gapping
- Integrally welded, interlocking shelf hangers increase stability while reducing hardware
- Durable powder coat paint will resist solvents, acids and alkalis
- Two grounding attachment points at the rear of all cabinets for easy grounding
- Continuous full door height piano hinges
- FM & UL approved fusible-link hold open feature ensures the door closes should the temperature outside the cabinet exceed 165 degrees Fahrenheit
- Manual and self-closing door options are available
- Self-closing, self-latching synchronized doors, ensure that the left door closes prior to the right door
- Equipped with a lockable, non-sparking, three-point latching mechanism for additional security
- Adjustable leveling screws lower or raise the cabinets by 1.6 cm (5/8") at each corner
- Polyethylene trays for galvanized steel shelves offer additional protection and have a slightly raised curved edge to contain the spills (Corrosive units only)
- Polyethylene trays in the bottom sill area catch incidental spills and leaks (Corrosive units only)
- Safety labels are provided in three languages: English, Spanish and French
- Upper and lower diametrically opposed side vent with flame arrestor for proper venting when specified
- Laboratory grade formulation meets SEFA 8 chemical requirement

Flammable & Corrosive Cabinet Specifications

NFPA 30 Chapter 9.5.3 (1)

States that "storage cabinets are designed and constructed to limit the internal temperature at the center of the cabinet and 1 in. (25 mm) from the top of the cabinet to not more than 325° F (163°C), when subjected to a 10-minute fire test that simulates the fire exposure of the standard time-temperature curve specified in NFPA 251, Standard Methods of Fire Resistance of Building Construction and Materials, shall be acceptable. All joints and seams shall remain tight and the door shall remain securely closed during the test".

NFPA 30 Chapter 9.5.3 (2) and OSHA 29 CFR 1910.106 (d)(3)(ii)(a)

States that "metal storage cabinets constructed in the following manner shall be acceptable:

- a. The bottom, top, door and sides of the cabinet shall be at least No. 18 gauge sheet steel and shall be double-walled, with 1½ in. (38 mm) air space.
- b. Joints shall be riveted, welded, or made tight by some equally effective means.
- c. The door shall be provided with a three-point latch arrangement, and the door sill shall be raised at least 2 in. (50 mm) above the bottom of the cabinet to retain spilled liquid within the cabinet."

OSHA 29 CFR 1910.106 (e)(2)(ii)(b)

States that "the quantity of liquid that may be located outside of an inside storage room or storage cabinet in a building or in any one fire area of a building shall not exceed:

- a. 25 gallons of Class 1A liquids in containers
- b. 120 gallons of Class 1B, 1C, or II or III liquids in containers
- c. 660 gallons of Class 1B, 1C, II or III liquids in a single portable tank"



Flammable Storage Cabinets



Bi-Fold
Self-Closing
Doors



Each door
latches open
with a fusible link

Flammable Storage Cabinets

Cabinets ensure maximum protection for both personnel and products

They are specifically designed for the storage of flammable and combustible liquids. The cabinets are designed in accordance with NFPA (National Fire Protection Association), OSHA (Occupational Safety and Health Administration) and FM (Factory Mutual International) regulations. Available in a full range of sizes and designs including manual close, self-closing and bi-fold self-closing doors.

- 15 different capacities
- Three door options
- Factory Mutual Approved
- Upper and lower diametrically opposed side vents. [View vent location diagrams.](#)

Model	Type Doors	Shelves	Gallons Capacity	H x D x L in inches
G65M151731N G65M151731V	Manual Self-Close	1	4	15" x 17" x 31"
G65M181843N G65M181843V	Manual Self-Close	1	12	18" x 18" x 43"
G65M221843N G65M221843V	Manual Self-Close	1	15	22" x 18" x 43"
G65M181847N G65M181847N	Manual Self-Close	1	15	18" x 18" x 47"
G65M241843N G65M241843V	Manual Self-Close	1	17	24" x 18" x 43"
G65M441243N	Manual	2	20	44" x 12" x 43"
G65M352235N G65M352235V	Manual Self-Close	1	23	35" x 22" x 35"
G65M352436N G65M352436V	Manual Self-Close	1	23	35" x 24" x 36"
G69M362235V	Bi-Fold Self-Close	1	24	36" x 22" x 35"
G65M441843N G65M441843V G69M441843V	Manual Self-Close Bi-Fold Self-Close	1	30	44" x 18" x 43"
G65M461843N G65M461843V	Manual Self-Close	2	30	46" x 18" x 43"
G69M362247V	Bi-Fold Self-Close	1	36	36" x 22" x 47"
G65M651843N G65M651843V G69M651843V	Manual Self-Close Bi-Fold Self-Close	2	45	65" x 18" x 43"
G65M653131N G65M653131V	Manual Self-Close	2	60	65" x 31" x 31"
G65M653434N G65M653434V G69M653434V	Manual Self-Close Bi-Fold Self-Close	2	60	65" x 34" x 34"
G65M653443N G65M653443V G69M653443V	Manual Self-Close Bi-Fold Self-Close	2	90	65" x 34" x 43"

Acid/Corrosive Storage Cabinets



Removable panel to service three point latching mechanism



Polyethylene trays for spill containment



Acid/Corrosive Storage Cabinets

Cabinets are designed and built for the safe storage of corrosive chemicals


These cabinets adhere to Kewaunee's rigid quality construction standards.

- 15 different capacities
- Three door options
- Factory Mutual Approved
- Upper and lower diametrically opposed side vents.
[View vent location diagrams.](#)

Model	Type Doors	Shelves	Gallons Capacity	H x D x L in inches
G85M151731N G85M151731V	Manual Self-Close	1	4	15" x 17" x 31"
G85M181843N G85M181843V	Manual Self-Close	1	12	18" x 18" x 43"
G85M221843N G85M221843V	Manual Self-Close	1	15	22" x 18" x 43"
G85M181847N G85M181847V	Manual Self-Close	1	15	18" x 18" x 47"
G85M241843N G85M241843V	Manual Self-Close	1	17	24" x 18" x 43"
G85M441243N	Manual	2	20	44" x 12" x 43"
G85M352235N G85M352235V	Manual Self-Close	1	23	35" x 22" x 35"
G85M352436N G85M352436V	Manual Self-Close	1	23	35" x 24" x 36"
G89M362235V	Bi-Fold Self-Close	1	24	36" x 22" x 35"
G85M441843N G85M441843V	Manual Self-Close	1	30	44" x 18" x 43"
G85M461843N G85M461843V	Manual Self-Close	2	30	46" x 18" x 43"
G89M441843V	Bi-Fold Self-Close	1	30	44" x 18" x 43"
G89M362247V	Bi-Fold Self-Close	1	36	36" x 22" x 47"
G85M651843N G85M651843V G89M651843V	Manual Self-Close Bi-Fold Self-Close	2	45	65" x 18" x 43"
G85M653131N G85M653131V	Manual Self-Close	2	60	65" x 31" x 31"
G85M653434N G85M653434V G89M653434V	Manual Self-Close Bi-Fold Self-Close	2	60	65" x 34" x 34"
G85M653443N G85M653443V G89M653443V	Manual Self-Close Bi-Fold Self-Close	2	90	65" x 34" x 43"

Base Cabinet Scribes & Fillers

End Scribes

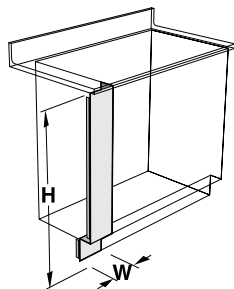

To close opening between wall and front of adjacent base cabinets. Includes detached toe space scribe and wall angle.

Available in widths from 1" - 24" in one inch increments

Sitting Height — 29" High

BESC2900 ☐ ☐ N-00_0


ADA Height — 32" High

BESC3200 ☐ ☐ N-00_0

Standing Height — 35" High

BESC3500 ☐ ☐ N-00_0

End Scribe Applique Panels

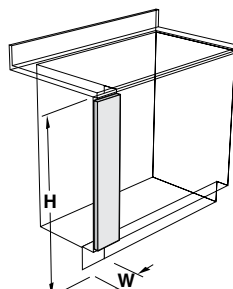

End Applique Panels may be used with End Scribes to match the face of doors and drawer heads on full overlay steel cabinets.

End Scribe must be ordered separately

Steel Square Edge

Widths from 3" - 24" in one inch increments

Sitting Height — 29" High

BEAC2900 ☐ ☐ N-11_0

ADA Height — 32" High — Steel

BEAC3200 ☐ ☐ N-11_0

Standing Height — 35" High — Steel

BEAC3500 ☐ ☐ N-11_0

Steel Beveled Edge

Widths from 3" - 24" in one inch increments

BEAC2900 ☐ ☐ N-1700

BEAC3200 ☐ ☐ N-1700

BEAC3500 ☐ ☐ N-1700

Wood Square Edge

Widths from 3" - 24" in one inch increments

BEAC2900 ☐ ☐ N-16_0


BEAC3200 ☐ ☐ N-16_0

BEAC3500 ☐ ☐ N-16_0

Insert width required before "N-" in part number to specify width of End Scribe or Applique Panel.

Example: BESC3500**06**N-0000 for 6" wide End Scribe or BEAC3500**06**N-1100 for 6" wide Steel Square Edge Applique Panel

Intermediate Fillers

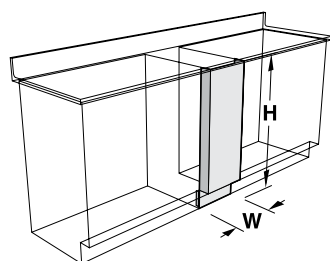

To close opening between two adjacent base cabinets. Includes detached toe space filler.

Available in widths from 6" - 48" in one inch increments

Sitting Height — 29" High

BIFC2900 ☐ ☐ N-00_0


ADA Height — 32" High

BIFC3200 ☐ ☐ N-00_0

Standing Height — 35" High

BIFC3500 ☐ ☐ N-00_0

Intermediate Applique Panels

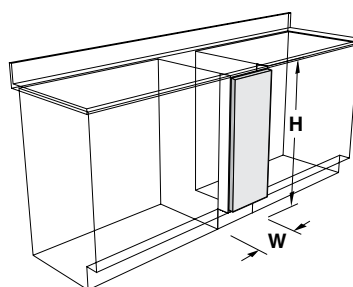

Intermediate Applique Panels may be used with Intermediate Fillers to match the face of doors and drawer heads on full overlay steel cabinets.

Intermediate Filler must be ordered separately

Steel Square Edge

Widths from 6" - 48" in one inch increments

Sitting Height — 29" High

BIAC2900 ☐ ☐ N-11_0

ADA Height — 32" High

BIAC3200 ☐ ☐ N-11_0

Standing Height — 35" High

BIAC3500 ☐ ☐ N-11_0

Steel Beveled Edge

Widths from 6" - 48" in one inch increments

BIAC2900 ☐ ☐ N-1700

BIAC3200 ☐ ☐ N-1700

BIAC3500 ☐ ☐ N-1700

Wood Square Edge

Widths from 6" - 48" in one inch increments

BIAC2900 ☐ ☐ N-16_0

BIAC3200 ☐ ☐ N-16_0


BIAC3500 ☐ ☐ N-16_0

Insert width required before "N-" in part number to specify width of Intermediate Filler or Applique Panel.

Example: BIFC3500**12**N-0000 for 12" wide Intermediate Filler or BIAC3500**12**N-1100 for 12" wide Steel Square Edge Applique Panel

Base Cabinet Scribes & Fillers

Corner Fillers

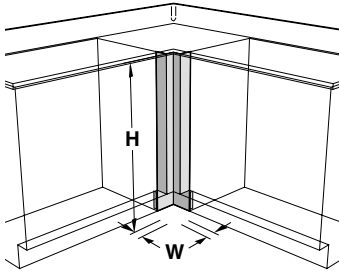

To close opening between two adjacent base cabinets joined at right angles. Includes detached toe space filler.

Available in widths from 1" - 12" in one inch increments

Sitting Height — 29" High

BCFC2900 N-00_0


ADA Height — 32" High

BCFC3200 N-00_0

Standing Height — 35" High

BCFC3500 N-00_0

Corner Applique Panels

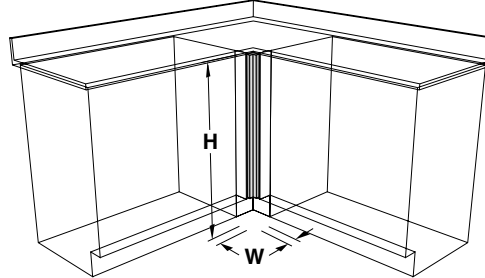

Corner Applique Panels may be used with Corner Fillers to match the face of doors and drawer heads on full overlay steel cabinets.

Corner Filler must be ordered separately

Steel Square Edge

Widths from 3" - 12" in one inch increments

Sitting Height — 29" High

BCAC2900 N-11_0

ADA Height — 32" High

BCAC3200 N-11_0

Standing Height — 35" High

BCAC3500 N-11_0

Steel Beveled Edge

Widths from 3" - 12" in one inch increments

BCAC2900 N-1700

BCAC3200 N-1700

BCAC3500 N-1700

Wood Square Edge

Widths from 3" - 12" in one inch increments

BCAC2900 N-16_0


BCAC3200 N-16_0

BCAC3500 N-16_0

Insert width required before "N-" in part number to specify size of Corner Filler or Applique Panel.

Example: BCFC350004N-0000 for 4" Corner Filler or BCAC350004N-1100 for 4" Steel Square Edge Applique Panel

Rear Scribes

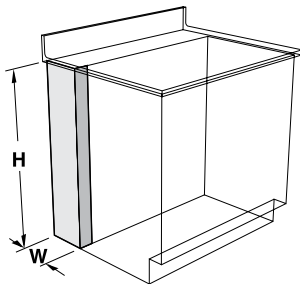

To close opening between wall and rear of adjacent base cabinet. Includes wall mounting angle.

Widths from 1" - 24" in one inch increments

Sitting Height — 29" High

BRSC2900 N-00_0


ADA Height — 32" High

BRSC3200 N-00_0

Standing Height — 35" High

BRSC3500 N-00_0

Pipe Chase Fillers

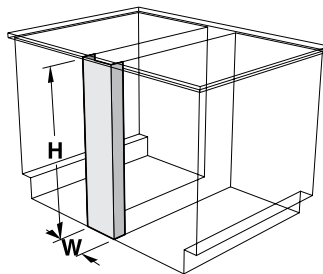

To close opening between two adjacent base cabinets installed back-to-back or a cabinet and a finished back.

Widths from 1" - 24" in one inch increments

Sitting Height — 29" High

BPCC2900 N-00_0


ADA Height — 32" High

BPCC3200 N-00_0

Standing Height — 35" High

BPCC3500 N-00_0

Cabinet End Applique Panel

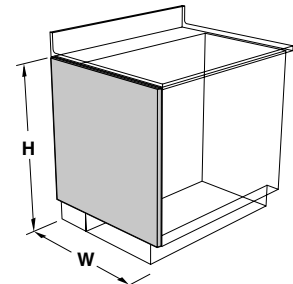

Wood Applique Panels to match exposed ends of base cabinets to cabinet faces.

Wood Square Edge

Sitting Height — 29" High

BPAC290023N-16_0 23" W

BPAC290029N-16_0 29" W

ADA Height — 32" High

BPAC320023N-16_0 23" W

BPAC320029N-16_0 29" W

Standing Height — 35" High


BPAC350023N-16_0 23" W

BPAC350029N-16_0 29" W

Insert width required before "N-" in part number to specify width of Rear Scribe or Pipe Chase Filler. **Example:** BRSC350008N-0000 for 8" wide Rear Scribe or BPCC350008N-0000 for 8" wide Pipe Chase Filler.

Wall Cabinet Scribes & Fillers

End Scribes

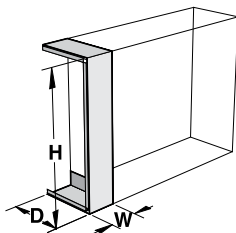

To close opening between wall and front of adjacent wall cabinet. Includes wall mounting angles.

12¾" Deep

Available in widths from 1" - 24" in one inch increments

16" Deep

24" High

WESC2413 ☐ ☐ N-00_0 WESC2416 ☐ ☐ N-00_0

30" High

WESC3013 ☐ ☐ N-00_0 WESC3016 ☐ ☐ N-00_0


36" High

WESC3613 ☐ ☐ N-00_0 WESC3616 ☐ ☐ N-00_0

48" High

WESC4813 ☐ ☐ N-00_0 WESC4816 ☐ ☐ N-00_0

End Scribe Applique Panels

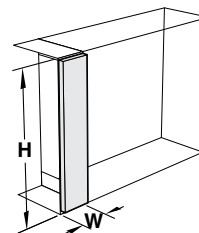

End Applique Panels may be used with End Scribes to match the face of doors and drawer heads on full overlay steel cabinets.

End Scribe must be ordered separately

Steel Square Edge

Widths from 3" - 24" in one inch increments

Steel Beveled Edge

Wood Square Edge

Widths from 3" - 24" in one inch increments

24" High

WEAC2400 ☐ ☐ N-11_0 not available WEAC2400 ☐ ☐ N-16_0

30" High

WEAC3000 ☐ ☐ N-11_0 not available WEAC3000 ☐ ☐ N-16_0

36" High

WEAC3600 ☐ ☐ N-11_0 not available WEAC3600 ☐ ☐ N-16_0


48" High

WEAC4800 ☐ ☐ N-11_0 not available WEAC4800 ☐ ☐ N-16_0

Insert width required before "N-" in part number to specify width of End Scribe or Applique Panel.

Example: WESC3613**06**N-0000 for 6" wide End Scribe or WEAC3600**06**N-1100 for 6" wide Steel Square Edge Applique Panel

Intermediate Fillers

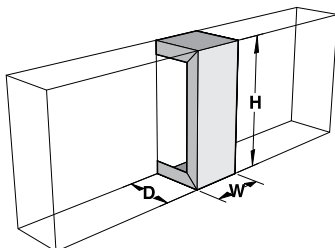

To close opening between two adjacent wall cabinets.

12¾" Deep

Available in widths from 6" - 48" in one inch increments

16" Deep

24" High

WIFC2413 ☐ ☐ N-00_0 WIFC2416 ☐ ☐ N-00_0

30" High

WIFC3013 ☐ ☐ N-00_0 WIFC3016 ☐ ☐ N-00_0


36" High

WIFC3613 ☐ ☐ N-00_0 WIFC3616 ☐ ☐ N-00_0

48" High

WIFC4813 ☐ ☐ N-00_0 WIFC4816 ☐ ☐ N-00_0

Intermediate Applique Panels

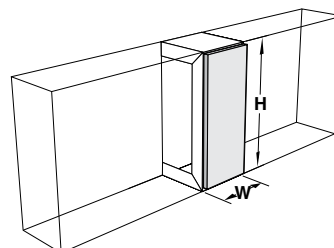

Intermediate Applique Panels may be used with Intermediate Fillers to match the face of doors and drawer heads on full overlay steel cabinets.

Intermediate Filler must be ordered separately

Steel Square Edge

Widths from 6" - 48" in one inch increments

Steel Beveled Edge

Wood Square Edge

Widths from 6" - 48" in one inch increments

24" High

WIAC2400 ☐ ☐ N-11_0 not available WIAC2400 ☐ ☐ N-16_0

30" High

WIAC3000 ☐ ☐ N-11_0 not available WIAC3000 ☐ ☐ N-16_0

36" High

WIAC3600 ☐ ☐ N-11_0 not available WIAC3600 ☐ ☐ N-16_0

48" High


WIAC4800 ☐ ☐ N-11_0 not available WIAC4800 ☐ ☐ N-16_0

Insert width required before "N-" in part number to specify width of Intermediate Filler or Applique Panel.

Example: WIFC3613**12**N-0000 for 12" wide Intermediate Filler or WIAC3600**12**N-1100 for 12" wide Steel Square Edge Applique Panel

Wall Cabinet Scribes & Fillers

Corner Fillers

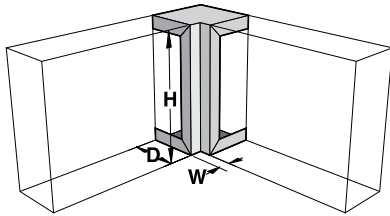

To close opening between two adjacent wall cabinets joined at right angles.

12¾" Deep

Available in depth and widths from 1" - 12" in one inch increments.

16" Deep

24" High

WCFC2413 N-00_0 WCFC2416 N-00_0

30" High

WCFC3013 N-00_0 WCFC3016 N-00_0


36" High

WCFC3613 N-00_0 WCFC3616 N-00_0

48" High

WCFC4813 N-00_0 WCFC4816 N-00_0

Corner Applique Panels

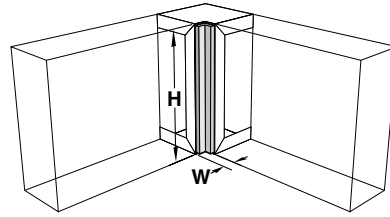

Corner Applique Panels may be used with Corner Fillers to match the face of doors and drawer heads on full overlay steel cabinets.

Corner Filler must be ordered separately

Steel Square Edge

Widths from 3" - 12" in one inch increments

Steel Beveled Edge

Wood Square Edge

Widths from 3" - 12" in one inch increments

24" High

WCAC2400 N-11_0 not available WCAC2400 N-16_0

30" High

WCAC3000 N-11_0 not available WCAC3000 N-16_0

36" High

WCAC3600 N-11_0 not available WCAC3600 N-16_0


48" High

WCAC4800 N-11_0 not available WCAC4800 N-16_0

Insert width required before "N-" in part number to specify width of Corner Filler or Applique Panel.

Example: WCFC3613**06**N-0000 for 6" wide Corner Filler or WCAC3600**06**N-1100 for 6" wide Steel Square Edge Applique Panel

Cabinet End Applique Panel

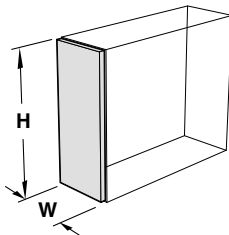

Wood Applique Panels to match exposed ends of wall cabinets to cabinet faces.

Wood Square Edge

24" High

WPAC240013N-16_0 12¾" W

WPAC240016N-16_0 16" W

30" High

WPAC300013N-16_0 12¾" W

WPAC300016N-16_0 16" W

36" High

WPAC360013N-16_0 12¾" W

WPAC360016N-16_0 16" W


48" High

WPAC480013N-16_0 12¾" W

WPAC480016N-16_0 16" W

Full Height Cabinet Scribes & Fillers

Full Height End Scribes

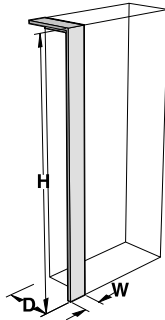

To close opening between wall and front of adjacent full height cabinets. Includes detached toe space scribe and wall mounting angles.

Available in widths from 1" - 24" in one inch increments

12¾" Deep

FESC8413 N-00_0


16" Deep

FESC8416 N-00_0

22" Deep

FESC8422 N-00_0

Full Height End Applique Panel

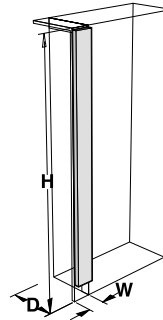

End Applique Panels may be used with End Scribes to match the face of doors and drawer heads on full overlay steel cabinets.

End Scribe must be ordered separately

Steel Square Edge

Widths from 3" - 24"
in one inch increments

Steel Beveled Edge

not available

Wood Square Edge

Widths from 3" - 24"
in one inch increments


FEAC8400 N-11_0

FEAC8400 N-16_0

Insert width required before "N-" in part number to specify width of End Scribe or Applique Panel.

Example: FESC841306N-0000 for 6" wide End Scribe or FEAC840006N-1100 for 6" wide Steel Square Edge Applique Panel

Full Height Intermediate Fillers

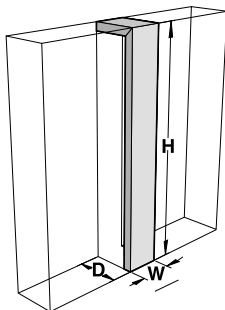

To close opening between two adjacent full height cabinets. Includes detached toe space filler.

Available in widths from 6" - 48" in one inch increments

12¾" Deep

FIFC8413 N-00_0


16" Deep

FIFC8416 N-00_0

22" Deep

FIFC8422 N-00_0

Full Height Intermediate Applique Panel

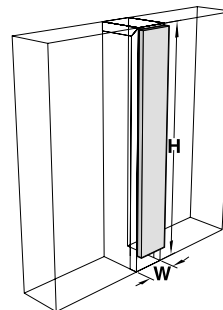

Intermediate Applique Panels may be used with Intermediate Fillers to match the face of doors and drawer heads on full overlay steel cabinets.

Intermediate Fillers must be ordered separately

Steel Square Edge

Widths from 6" - 48"
in one inch increments

Steel Beveled Edge

not available

Wood Square Edge

Widths from 6" - 48"
in one inch increments

FIAC8400 N-11_0


FIAC8400 N-16_0

Insert width required before "N-" in part number to specify width of Intermediate Filler or Applique Panel.

Example: FIFC841306N-0000 for 6" wide Intermediate Filler or FIAC840006N-1100 for 6" wide Steel Square Edge Applique Panel

Full Height Cabinet Scribes & Fillers

Corner Filler

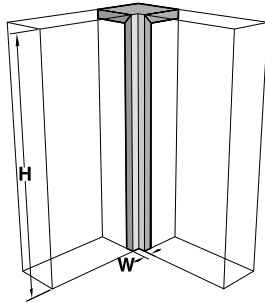

To close opening between two adjacent full height cabinets joined at right angles. Includes detached toe space filler.

Available in widths from 1" - 12" in one inch increments.

12¾" Deep

FCFC8413 ☐☐ **N-00_0**


16" Deep

FCFC8416 ☐☐ **N-00_0**

22" Deep

FCFC8422 ☐☐ **N-00_0**

Corner Fillers Applique Panel

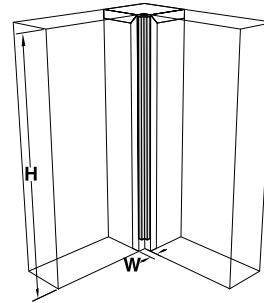

Corner Applique Panels may be used with Corner Fillers to match the face of doors and drawer heads on full overlay steel cabinets.

Corner Filler must be ordered separately

Steel Square Edge

Widths from 3" - 12"
in one inch increments

Steel Beveled Edge

not available

Wood Square Edge

Widths from 3" - 12"
in one inch increments


FCAC8400 ☐☐ **N-11_0**

FCAC8400 ☐☐ **N-16_0**

Insert width required before "N-" in part number to specify width of Corner Filler or Applique Panel.

Example: FCFC8413**06**N-0000 for 6" wide Corner Filler or FCAC8400**06**N-1100 for 6" wide Steel Square Edge Applique Panel

Full Height Rear Scribe

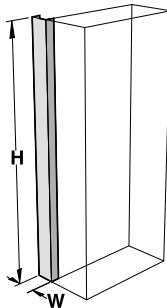

To close opening between wall and rear of adjacent full height cabinets.

Available in widths from 1" - 24" in one inch increments

56" High

FRSC5600 ☐☐ **N-00_0**


84" High

FRSC8400 ☐☐ **N-00_0**

Insert width required before "N-" in part number to specify width of Rear Scribe.

Example: FRSC8400**06**N-0000 for 6" wide Rear Scribe.

Cabinet End Applique Panel

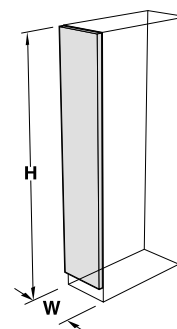

Wood Applique Panels to match exposed ends of full height cabinets to cabinet faces.

Wood Square Edge

84" High

FPAC840013N-16_0

12¾" W

FPAC840016N-16_0


16" W

FPAC840022N-16_0

22" W

Steel Cabinet Finished Backs

Base Cabinet Finished Back

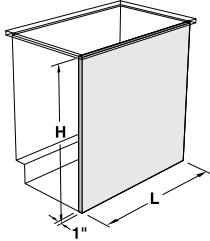

To cover backs of base cabinets exposed to view. Back adds 1" to the depth of the base cabinet.

Sitting Height — 29" High

BFBC290112N-00_0	12"W
BFBC290115N-00_0	18"W
BFBC290118N-00_0	18"W
BFBC290121N-00_0	21"W
BFBC290124N-00_0	24"W
BFBC290130N-00_0	30"W
BFBC290136N-00_0	36"W
BFBC290142N-00_0	42"W
BFBC290148N-00_0	48"W
BFBC290160N-00_0	60"W

ADA Height — 32" High


BFBC320112N-00_0	12"W
BFBC320115N-00_0	15"W
BFBC320118N-00_0	18"W
BFBC320121N-00_0	21"W
BFBC320124N-00_0	24"W
BFBC320130N-00_0	30"W
BFBC320136N-00_0	36"W
BFBC320142N-00_0	42"W
BFBC320148N-00_0	48"W
BFBC320160N-00_0	60"W

Standing Height — 35" High

BFBC350112N-00_0	12"W
BFBC350115N-00_0	15"W
BFBC350118N-00_0	18"W
BFBC350121N-00_0	21"W
BFBC350124N-00_0	24"W
BFBC350130N-00_0	30"W
BFBC350136N-00_0	36"W
BFBC350142N-00_0	42"W
BFBC350148N-00_0	48"W
BFBC350152N-00_0U3*	52"W
BFBC350158N-00_0U3*	58"W
BFBC350160N-00_0U3*	60"W
BFBC350160N-00_0	60"W

* Includes toespace notches for End Sink Cabinets

Wall Cabinet Finished Back

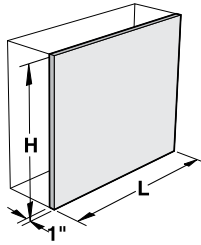

To cover backs of wall cabinets exposed to view. Adds 1" to depth of cabinet.

24" High

WFBC240115N-00_0	15"W
WFBC240118N-00_0	18"W
WFBC240121N-00_0	21"W
WFBC240124N-00_0	24"W
WFBC240130N-00_0	30"W
WFBC240136N-00_0	36"W
WFBC240142N-00_0	42"W
WFBC240148N-00_0	48"W

30" High

WFBC300115N-00_0	15"W
WFBC300118N-00_0	18"W
WFBC300121N-00_0	21"W
WFBC300124N-00_0	24"W
WFBC300130N-00_0	30"W
WFBC300136N-00_0	36"W
WFBC300142N-00_0	42"W
WFBC300148N-00_0	48"W


36" High

WFBC360115N-00_0	15"W
WFBC360118N-00_0	18"W
WFBC360121N-00_0	21"W
WFBC360124N-00_0	24"W
WFBC360130N-00_0	30"W
WFBC360136N-00_0	36"W
WFBC360142N-00_0	42"W
WFBC360148N-00_0	48"W

48" High

WFBC480115N-00_0	15"W
WFBC480118N-00_0	18"W
WFBC480121N-00_0	21"W
WFBC480124N-00_0	24"W
WFBC480130N-00_0	30"W
WFBC480136N-00_0	36"W
WFBC480142N-00_0	42"W
WFBC480148N-00_0	48"W

Full Height Finished Back

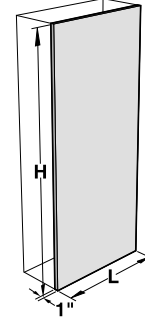


To cover backs of full height cabinets exposed to view. Adds 1" to depth of cabinet.

84" High

FFBC840118N-00_0	18"W
FFBC840124N-00_0	24"W
FFBC840130N-00_0	30"W
FFBC840136N-00_0	36"W
FFBC840142N-00_0	42"W
FFBC840148N-00_0	48"W

Steel Adjustable Slip Scribes

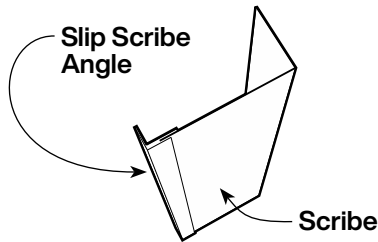


Specifications:

Steel Adjustable Slip Scribes are designed the same as standard scribes with the addition of a Slip Scribe Angle at the wall edge to dress-up the edge of the scribe.

The Steel Adjustable Slip Scribe is designed to be cut to width in the field to more easily fit the space between the wall and the end of the furniture assembly.

Adjustable Base Cabinet Scribe

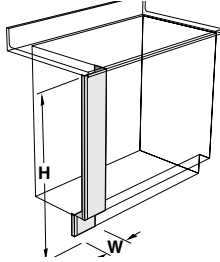

To close opening between wall and front of adjacent base cabinets. Includes detached toe space scribe, end scribe and slip scribe angle.

Sitting Height — 29" High x 12"W

ABSC290012N-00_0


ADA Height — 32" High

ABSC320012N-00_0

Standing Height — 35" High

ABSC350012N-00_0

Adjustable Wall Cabinet Scribe

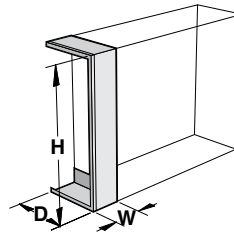

To close opening between wall and front of adjacent wall cabinet. Includes end scribe and slip scribe angle.

24" High

AWSC241312N-00_0

AWSC241612N-00_0

12"W

12³/₄"D

16"D

30" High

AWSC301312N-00_0

AWSC301612N-00_0

12³/₄"D

16"D

36" High

AWSC361312N-00_0

AWSC361612N-00_0

12³/₄"D

16"D

48" High


AWSC481312N-00_0

AWSC481612N-00_0

12³/₄"D

16"D

Adjustable Full Height Scribe

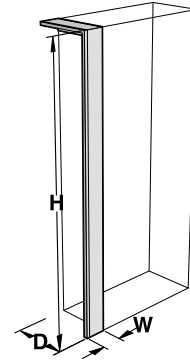

To close opening between wall and front of adjacent full height cabinets. Includes detached toe space scribe, end scribe, and slip scribe angle.

84" High

AFSC841312-00_0

AFSC841612-00_0

AFSC842212-00_0

12"W

12³/₄"D

16"D

22"D

Slip Scribe Angle




May be used with any standard scribe or filler to make it adjustable. Available in lengths from 3" to 96" in one inch increments.

SCAC0101□□N-00_0

Insert length required before "N-" in part number to specify length of angle.

Spring Clips

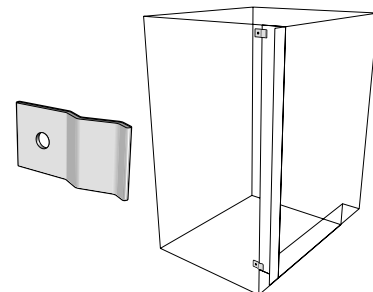


18 gauge spring steel clip, 3/4" x 1 1/2". May be used with any scribe or filler to make them easily removable.

F-4995-PK package of 10

Steel Sloping Tops


Sloping Tops

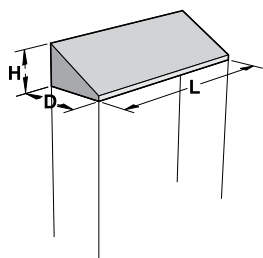

May be used on wall cabinets and full height cabinets to keep the tops dust free.

8 $\frac{3}{8}$ " High 12 $\frac{3}{4}$ " Deep	10 $\frac{3}{4}$ " High 16" Deep	13 $\frac{3}{4}$ " High 22" Deep	
SLTC081315N-00_0	SLTC111615N-00_0	SLTC142215N-00_0	15"L
SLTC081318N-00_0	SLTC111618N-00_0	SLTC142218N-00_0	18"L
SLTC081321N-00_0	SLTC111621N-00_0	SLTC142221N-00_0	21"L
SLTC081324N-00_0	SLTC111624N-00_0	SLTC142224N-00_0	24"L
SLTC081330N-00_0	SLTC111630N-00_0	SLTC142230N-00_0	30"L
SLTC081336N-00_0	SLTC111636N-00_0	SLTC142236N-00_0	36"L
SLTC081342N-00_0	SLTC111642N-00_0	SLTC142242N-00_0	42"L
SLTC081348N-00_0	SLTC111648N-00_0	SLTC142248N-00_0	48"L
SLTC081360N-00_0	SLTC111660N-00_0	SLTC142260N-00_0	60"L
SLTC081372N-00_0	SLTC111672N-00_0	SLTC142272N-00_0	72"L
SLTC081384N-00_0	SLTC111684N-00_0	SLTC142284N-00_0	84"L
SLTC081396N-00_0	SLTC111696N-00_0	SLTC142296N-00_0	96"L


Sloping Top Fillers

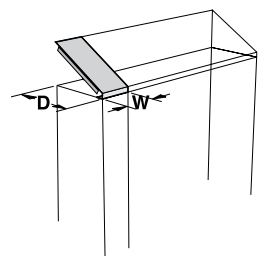

To close opening between wall and front of adjacent sloping top. Includes wall mounting angle.

8 $\frac{3}{8}$ " High 12 $\frac{3}{4}$ " Deep	10 $\frac{3}{4}$ " High 16" Deep	13 $\frac{3}{4}$ " High 22" Deep	
STFC081301N-00_0	STFC111601N-00_0	STFC142201N-00_0	1"W
STFC081306N-00_0	STFC111606N-00_0	STFC142206N-00_0	6"W
STFC081312N-00_0	STFC111612N-00_0	STFC142212N-00_0	12"W


Corner Cabinet Sloping Tops

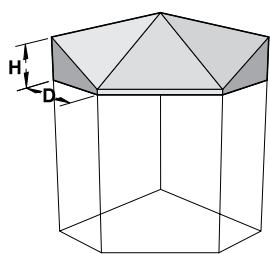

To be used on corner wall cabinets.

8 $\frac{3}{8}$ " High 12 $\frac{3}{4}$ " Deep	10 $\frac{3}{4}$ " High 16" Deep	
CSTC081327N-00_0	CSTC111627N-00_0	26 $\frac{7}{8}$ "L

Sloping Top Corner Fillers

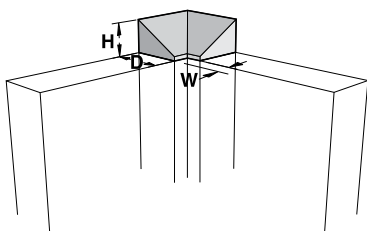


May be used on wall cabinets and full height cabinets when cabinets are joined at right angles.

8 $\frac{3}{8}$ " High 12 $\frac{3}{4}$ " Deep	10 $\frac{3}{4}$ " High 16" Deep	13 $\frac{3}{4}$ " High 22" Deep	
SCFC081301N-00_0	SCFC111601N-00_0	SCFC142201N-00_0	1"W
SCFC081306N-00_0	SCFC111606N-00_0	SCFC142206N-00_0	6"W

Steel Ceiling & Pipe Enclosures

Ceiling Enclosures

Corner Ceiling Enclosures

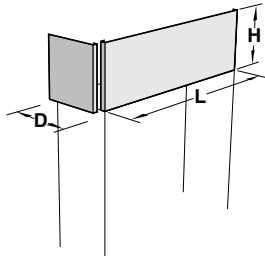


To close opening between ceiling and top of wall or full height cabinets.

Front Ceiling Enclosure

12" High	24" High	
FCEC120018N-00_0	FCEC240018N-00_0	18"L
FCEC120024N-00_0	FCEC240024N-00_0	24"L
FCEC120030N-00_0	FCEC240030N-00_0	30"L
FCEC120036N-00_0	FCEC240036N-00_0	36"L
FCEC120042N-00_0	FCEC240042N-00_0	42"L
FCEC120048N-00_0	FCEC240048N-00_0	48"L
FCEC120060N-00_0	FCEC240060N-00_0	60"L
FCEC120072N-00_0	FCEC240072N-00_0	72"L
FCEC120084N-00_0	FCEC240084N-00_0	84"L
FCEC120096N-00_0	FCEC240096N-00_0	96"L

End Ceiling Enclosure


12" High	24" High	
ECEC120013N-00_0	ECEC240013N-00_0	12¾"D
ECEC120016N-00_0	ECEC240016N-00_0	16"D
ECEC120022N-00_0	ECEC240022N-00_0	22"D

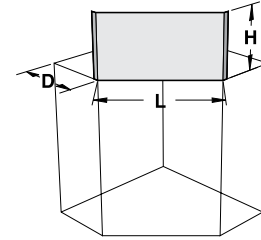

To close opening between ceiling and top of corner wall cabinets.

Pipe Enclosures

Wall Bench Enclosure

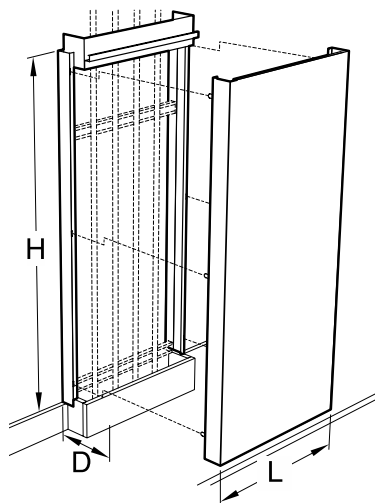


Pipe Enclosures conceal vertical service lines coming down from the ceiling to the work surface. One side is removable for easy access with no exposed fasteners. Designed to sit on top of curb on work surface and penetrate drop ceiling. Available in powder coated steel or stainless steel.

U01C ☐ ☐ ☐ ☐ N-00_0 Wall Bench Enclosure
 U02C ☐ ☐ ☐ ☐ N-00_0 Center Bench Enclosure

Available Heights	Available Depths	Available Lengths	Available Materials
60 = 60"	05 = 5"	12 = 12"	0 = Powder Coated Steel
66 = 66"	06 = 6"	18 = 18"	3 = Stainless Steel
72 = 72"	07 = 7"	24 = 24"	
78 = 78"	08 = 8"	30 = 30"	
84 = 84"	09 = 9"	36 = 36"	
90 = 90"	10 = 10"	42 = 42"	
96 = 96"		48 = 48"	

Center Bench Enclosure

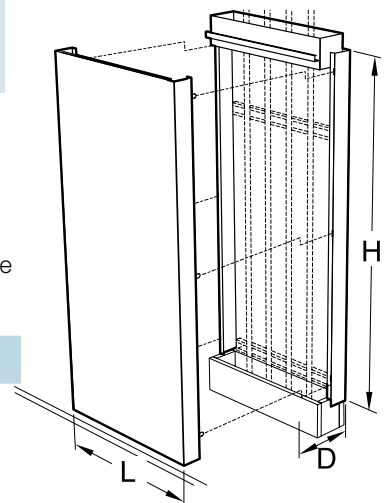

Fume Hoods & Ventilated Enclosures

Ducted Ventilated Products


Kewaunee offers a full line of Laboratory Fume Hoods to meet the needs of every lab application. In addition to our newest

line of **Supreme Air Venturi Fume Hoods**, we also offer uniquely designed hoods for education and specialty

hoods for custom applications.

Supreme Air – General Purpose Fume Hoods


The design is based on rigid frame construction that assures a solid installation and low vibration and sound levels. Access panels are easily removable to access service lines and fittings. Radius corner posts and airfoils, plus, easy operator control of interior baffle settings assure a high level of comfort, safety and efficiency.

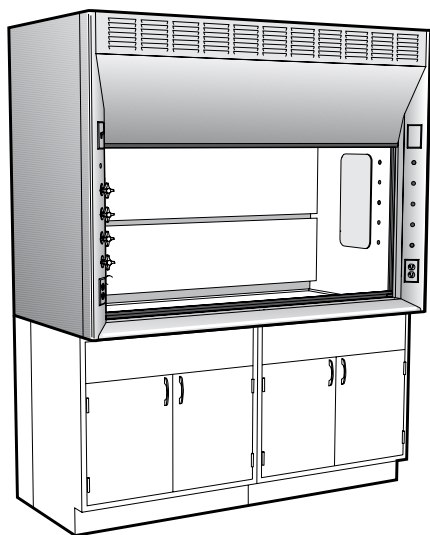

For more information on Kewaunee's **Supreme Air Fume Hoods** contact your local Kewaunee sales representative or visit the Kewaunee website.

Supreme Air Venturi – Next Generation Fume Hoods

The Supreme Air Venturi utilizes proven containment enhancing techniques coupled with innovative laboratory convenience features, including: adjustable LED lights, electromechanical sash stop, ultra-smooth sash system, and push-button sash management. The Supreme Air Venturi provides the safest user friendly fume hood manufactured today.

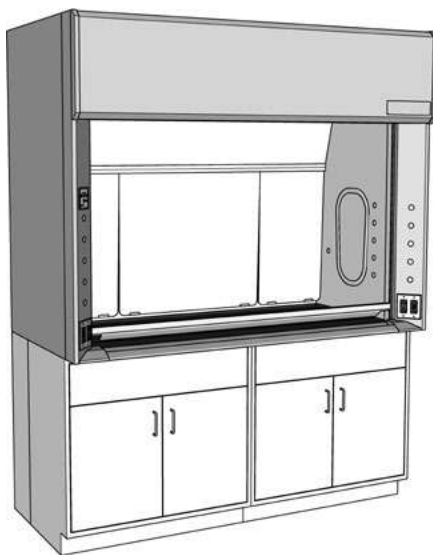

For more information on Kewaunee's **Supreme Air Venturi Fume Hoods** contact your local Kewaunee sales representative or visit the Kewaunee website.

TruView – Teaching Fume Hoods

Designed for today's classroom allowing clear line-of-sight through multiple hoods to allow easier supervision of students by the instructor. May be used as a stand-alone hood or in assemblies, both single-sided for use against the wall, or double-sided for use in the center of the laboratory.



For more information on Kewaunee's **TruView Teaching Fume Hoods** contact your local Kewaunee sales representative or visit the Kewaunee website.

Worktops


Worktop material selections are among the most critical of the laboratory furniture decisions. The intended applications for your worktops dictate material choice. Surfaces could be exposed

to reagents and caustic substances. Heavy equipment may need to be supported or only light duty tasks may be performed.

Kewaunee can recommend and furnish materials fabricated to your individual needs, complete with appropriate sinks, drains and accessories.

Kemresin® Worktops

manufactured exclusively by Kewaunee

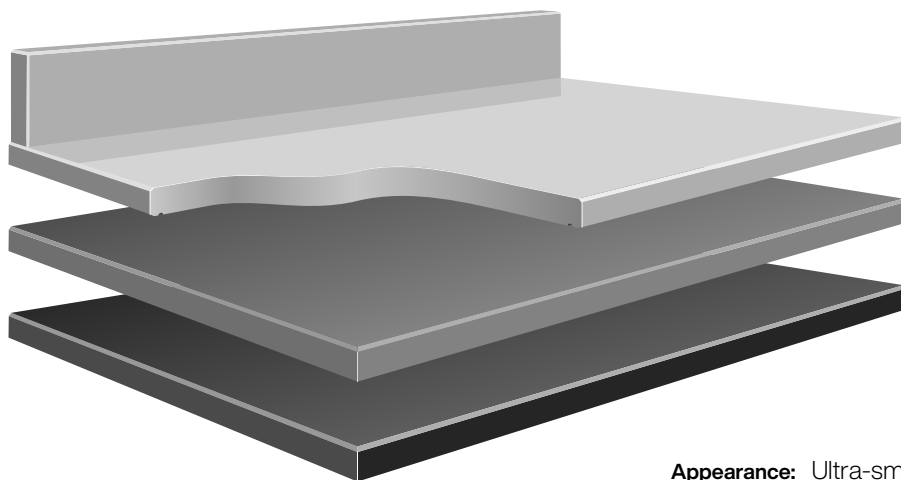

Appearance: Ultra-smooth non-porous surface.

Available Colors: Black, Grey, Putty, and Slate

Composition: 1" thick molded epoxy resin.

Physical Properties: Excellent resistance to practically all acids, alkalis and solvents. Highly resistant to heat and normal physical abuse.

Applications: Use wherever optimum physical and chemical resisting properties are required.

Other Worktop Materials

MATERIAL	APPEARANCE	COMPOSITION	CHEMICAL RESISTANCE	APPLICATIONS
STAINLESS STEEL	No.4 Finish	Type 304 furnished as standard, but other types are available when specified. Type 316 supplied for Perchloric Acid fume hoods.	Good chemical resistance to a wide range of reagents.	Recommended for sterile and radiological laboratories.
HARDWOOD, OAK OR MAPLE	Natural varnish, or penetrating oil finish	Solid edge grain Maple or Oak strips, glued together electronically, and finished with baked-on clear synthetic varnish or penetrating oil.	Fair chemical resistance, but not recommended for use in chemical laboratories.	Excellent top material for Physics tables, Electronic benches, and other tables where mechanical or physical work levels are prevalent.
COMPOSITE RESIN	Hard, smooth non-glare surface	Cellulose fiber reinforced phenolic resin core with a polyurethane copolymer surface. Available in multiple colors and patterns.	Good chemical resistance to wide range of acids, alkalis and solvents.	Use wherever a high degree of chemical resistance is required. (Not recommended for fume hood work surfaces.)

Worktops

Specifications:

The configurations shown are for Kewaunee base cabinets. Standard top thickness is 1", standard curb height is 4".

To order, select configuration and top material desired. Specify number of linear or square feet per product number.

Dimensions are nominal. Not all configurations are available in all materials. Indented tops available on special order.


Wall Bench Tops



Flat Top no curb		Style B
Standard	18"	
Sizes:	23" 24" 30" 31"	



Top with Applied Curb		Style C
Standard	18"	
Sizes:	23" 24" 30" 31"	

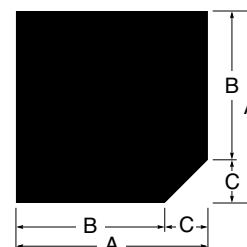

Corner Top no curbs		Style K		
	A	B	C	
Standard	33"	24"	9"	
Sizes:	39" 40"	30" 31"	9" 9"	

Flat Top & Ledge no curb

Style E

	A	B	C
Standard	30"	23"	7"
Sizes:	31"	23"	8"

The diagram shows a cross-section of a top and ledge assembly. It consists of a main rectangular block with a smaller rectangular section attached to its right side. Dimension A is the total width of the assembly. Dimension B is the width of the main block. Dimension C is the width of the smaller section. The assembly is shown in a perspective view with a top surface and a side surface.

Top & Ledge with Applied curb

Style F

	A	B	C
Standard	30"	23"	7"
Sizes:	31"	23"	8"

Corner Top with Applied Curbs		Style L		
	A	B	C	
Standard	33"	24"	9"	
Sizes:	39" 40"	30" 31"	9" 9"	

Center Bench Tops

Ledges and Shelves



Top without Center Joint		Style H
Standard	54"	
Sizes:	60" 62"	

The diagram illustrates the dimensions for a 'Top with Center Joint' in 'Style I'. It shows a rectangular top with a central joint. The total width is labeled 'A'. The width of the left side is labeled 'B', and the width of the right side is labeled 'C'. The joint is indicated by a vertical line in the center, with a double vertical line on each side of the joint, suggesting a seam or fold.

Top with Center Joint

Style I

	A	B	C
Standard	54"	27"	27"
Sizes:	54"	23"	31"
	60"	30"	30"
	62"	31"	31"

The diagram illustrates the dimensions of a rectangular top with a center ledge. The total width is labeled 'A'. The width is divided into three sections: 'B' (left section), 'C' (center ledge), and 'D' (right section). The sections are separated by vertical lines, and the dimensions are indicated by horizontal arrows below the top.

Top with Center Ledge

Style J

	A	B	C	D
Standard	54"	23"	8"	23"
Sizes:	60"	23"	14"	23"
	62"	23"	16"	23"



Ledge Unfinished Edges		Style A
Standard	8"	
Sizes:	10" 12" 14" 16"	



Shelf Finished Edges		Style N
Standard	8"	
Sizes:	10" 12" 14" 16"	

Center Bench Shelving Assemblies

Specifications:


Shelving Support Assemblies are designed to add shelving to center benches. Each assembly consists of four uprights with slots on 1" centers and a top frame. They are 12" deep in lengths as shown.

Each assembly is furnished with mounting brackets for attaching to the counter top and cabinets below but no fasteners. Fasteners must be furnished in field depending on top material and cabinet configuration.


Shelving Assembly Support Struts are especially designed to provide added support to the Shelving Support Assemblies while allowing space for connection of plumbing or electrical connections.

12" Deep Shelving Support Assembly

Shelving Assembly Support Strut

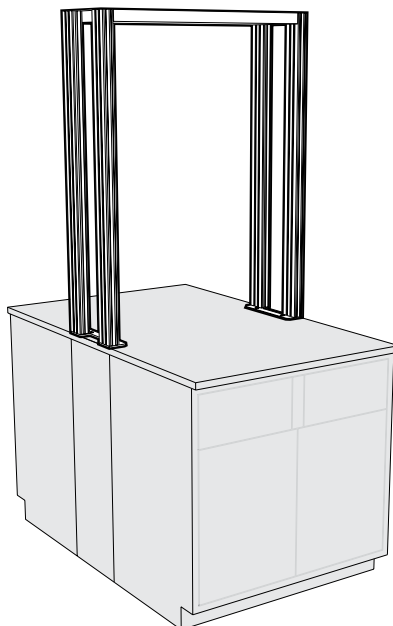

Length	Part Number
24"	F17N481224N
30"	F17N481230N
36"	F17N481236N
42"	F17N481242N
48"	F17N481248N
54"	F17N481254N
60"	F17N481260N
66"	F17N481266N
72"	F17N481272N

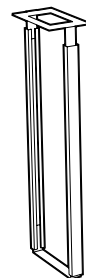

Width	Part Number
12"	K17M360112N
14"	K17M360114N

Shelving Uprights



Powder Coated Steel Uprights

Length	Single Slot Upright	Double Slot Upright	
24"	F13M240102N	F12M240102N	each
30"	F13M300102N	F12M300102N	each
36"	F13M360102N	F12M360102N	each
42"	F13M420102N	F12M420102N	each
48"	F13M480102N	F12M480102N	each
72"	F13M720102N	F12M720102N	each
79"	F13M790102N	F12M790102N	each
84"	F13M840102N	F12M840102N	each

Stainless Steel Uprights

Length	Single Slot Upright	Double Slot Upright	
24"	F13S240102N	F12S240102N	each
30"	F13S300102N	F12S300102N	each
36"	F13S360102N	F12S360102N	each
42"	F13S420102N	F12S420102N	each
48"	F13S480102N	F12S480102N	each
72"	F13S720102N	F12S720102N	each
79"	F13S790102N	F12S790102N	each
84"	F13S840102N	F12S840102N	each

Designed to be used with Adjustable Module Shelving Assemblies to create wall mounted shelving.

Center Bench Shelving Assemblies

Specifications:

Pass-Thru Shelf Assemblies for 12" Island Modules are six inches deeper than the Shelf Bracket to allow the Shelf to pass through the module core and form a continuous surface. These shelves may also be moved forward and mounted to form a 2" cord gap between

the shelves. Shelf Brackets are 11 gauge steel, adjustable on 1" centers.


Painted Steel Shelves – 1" thick, formed from 16 gauge cold rolled steel.

Veneer Plywood Shelves – 1" thick hardwood plywood with matching


veneer on all exposed surfaces. Available in Oak, Maple or Beech.

Phenolic Resin Shelves – 1" thick phenolic resin.

Pass-Thru Shelf Assemblies

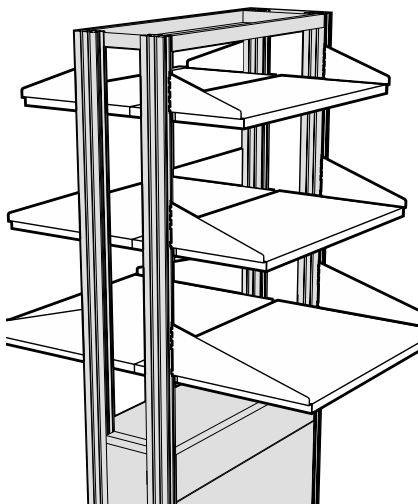

Pass-Thru Shelves

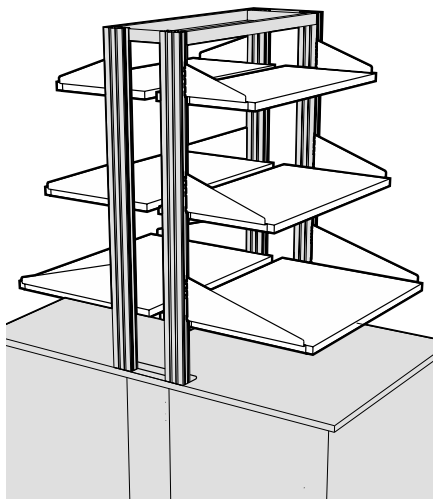

Pass-Thru Shelves with Cord Gap

	Length	Painted Steel	Veneer Plywood	Phenolic Resin
15" DEEP SHELF 9" Deep Bracket	24"	F21M040924X	F21W040924X	F21T040924X
	30"	F21M040930X	F21W040930X	F21T040930X
	36"	F21M040936X	F21W040936X	F21T040936X
	42"	F21M040942X	F21W040942X	F21T040942X
	48"	F21M040948X	F21W040948X	F21T040948X
	54"	F21M040954X		
	60"	F21M040960X		
	66"	F21M040966X		
18" DEEP SHELF 12" Deep Bracket	72"	F21M040972X		
	24"	F21M051224X	F21W051224X	F21T051224X
	30"	F21M051230X	F21W051230X	F21T051230X
	36"	F21M051236X	F21W051236X	F21T051236X
	42"	F21M051242X	F21W051242X	F21T051242X
	48"	F21M051248X	F21W051248X	F21T051248X
	54"	F21M051254X		
	60"	F21M051260X		
24" DEEP SHELF 18" Deep Bracket	66"	F21M051266X		
	72"	F21M051272X		
	24"	F21M061824X	F21W061824X	F21T061824X
	30"	F21M061830X	F21W061830X	F21T061830X
	36"	F21M061836X	F21W061836X	F21T061836X
	42"	F21M061842X	F21W061842X	F21T061842X
	48"	F21M061848X	F21W061848X	F21T061848X
	54"	F21M061854X		
	60"	F21M061860X		


To specify Veneer Plywood Species replace blank at end of part number with number that indicates material required.

Veneer Plywood Material Options

- | | |
|-------------------|------------------------|
| 1 = Red Oak | 5 = FSC® Red Oak |
| 2 = White Maple | 6 = FSC® White Maple |
| 7 = Steamed Beech | 8 = FSC® Steamed Beech |

Shelf Retaining Lips

Length	Painted Steel	Stainless Steel
24"	F63M020124N	F63S020124N
30"	F63M020130N	F63S020130N
36"	F63M020136N	F63S020136N
42"	F63M020148N	F63S020148N
54"	F63M020154N	F63S020154N
60"	F63M020160N	F63S020160N
66"	F63M020166N	F63S020166N
72"	F63M020172N	F63S020172N

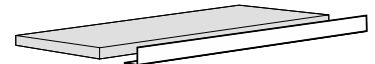

Designed to be added to Pass-thru and Wall Shelf Assemblies. Steel or Stainless Steel Shelf Lips may be mounted to steel, veneer plywood, or phenolic resin shelves with supplied screws.

Wall Shelving Assemblies

Specifications:

Vertical Uprights are 1¼" deep and either 79" or 42" high. They are manufactured with a double-slotted face with slots on 1" centers. They are designed to be mounted between two Support Rails to hold accurate spacing.


Support Rails are available in 8 foot

lengths and may be cut to the required length. Support Rails are notched on 6" centers to accept the Vertical Uprights.

Support Rail Covers are required to finish the assembly and must be ordered in pairs to match the Vertical Upright spacing.

Optional Center Rails may be ordered to provide a horizontal member at counter top height. They should be sized to match the Vertical Upright spacing.

Full Height Wall Suspension System

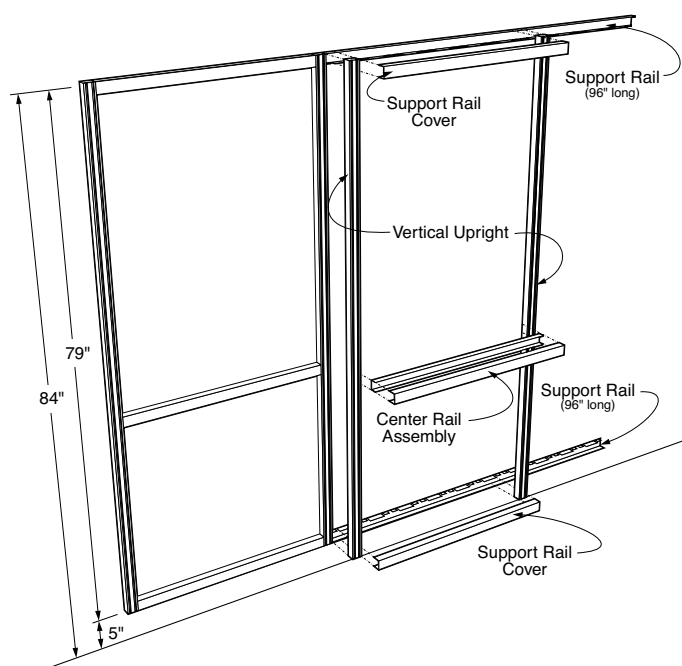

Vertical Upright

F10M790102N	79"H	pair
-------------	------	------

Support Rail

F11M020196S	96"L	each
-------------	------	------


Support Rail Covers

F11M030124P	24"L	pair
F11M030130P	30"L	pair
F11M030136P	36"L	pair
F11M030142P	42"L	pair
F11M030148P	48"L	pair
F11M030154P	54"L	pair
F11M030160P	60"L	pair
F11M030166P	66"L	pair
F11M030172P	72"L	pair

Optional Center Rail Assembly

F11M020124C	24"L	each
F11M020130C	30"L	each
F11M020136C	36"L	each
F11M020142C	42"L	each
F11M020148C	48"L	each
F11M020154C	54"L	each
F11M020160C	60"L	each
F11M020166C	66"L	each
F11M020172C	72"L	each

Above Countertop Wall Suspension System

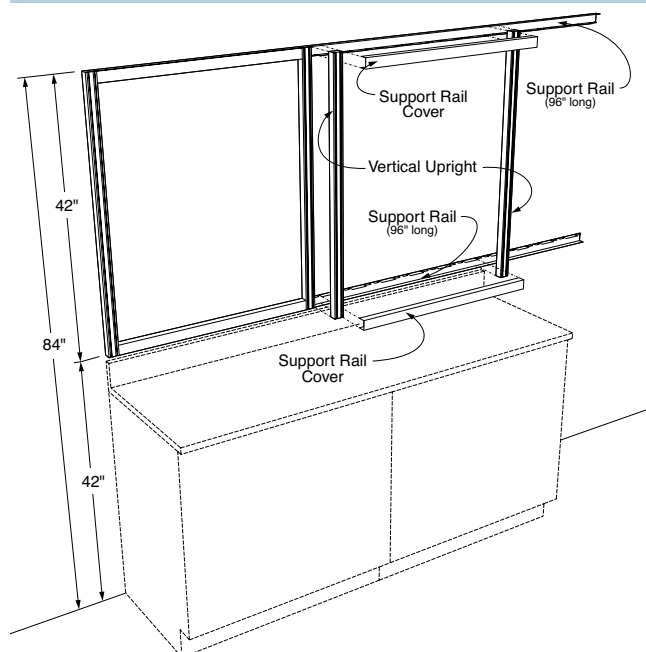

Vertical Upright

F10M420102N	42"H	pair
-------------	------	------

Support Rail

F11M020196S	96"L	each
-------------	------	------

Support Rail Covers

F11M030124P	24"L	pair
F11M030130P	30"L	pair
F11M030136P	36"L	pair
F11M030142P	42"L	pair
F11M030148P	48"L	pair
F11M030154P	54"L	pair
F11M030160P	60"L	pair
F11M030166P	66"L	pair
F11M030172P	72"L	pair

Wall Shelving Assemblies

Specifications:

Adjustable Wall Shelf Assemblies consist of a pair of 11 gauge steel Shelf Brackets and a Shelf. The Shelf Bracket engages the inner slot of the double slotted Wall Vertical Upright and is adjustable on 1" centers. Shelves are



available in Painted Steel, Wood veneer plywood, or Phenolic Resin.

Painted Steel Shelves – 1" thick, formed from 16 gauge cold rolled steel.

Veneer Plywood Shelves – 1" thick hardwood plywood with matching veneer on all exposed surfaces. Available in Oak, Maple or Beech.

Phenolic Resin Shelves – 1" thick phenolic resin.

Adjustable Wall Shelf Assemblies

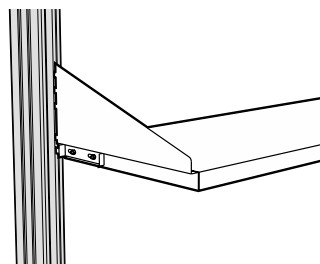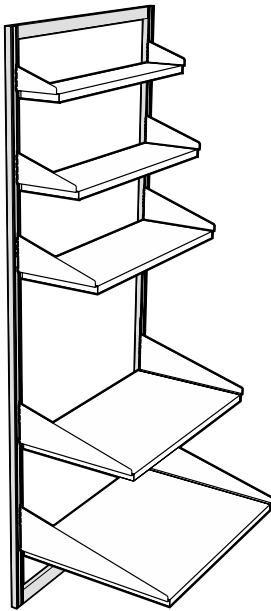

Optional Shelf Locking Tab

Sliding tab mounted on Shelf Bracket locks Shelf into upright to keep it from being accidentally removed.

To specify, add an **K** to the end of the Shelf Assembly part number.

Example: F21M051236NK

	Length	Painted Steel	Veneer Plywood	Phenolic Resin
6" DEEP SHELF 6" Deep Bracket	24"	F21M030624N	F21W030624N	F21T030624N
	30"	F21M030630N	F21W030630N	F21T030630N
	36"	F21M030636N	F21W030636N	F21T030636N
	42"	F21M030642N	F21W030642N	F21T030642N
	48"	F21M030648N	F21W030648N	F21T030648N
	54"	F21M030654N		
	60"	F21M030660N		
	66"	F21M030666N		
9" DEEP SHELF 9" Deep Bracket	24"	F21M040924N	F21W040924N	F21T040924N
	30"	F21M040930N	F21W040930N	F21T040930N
	36"	F21M040936N	F21W040936N	F21T040936N
	42"	F21M040942N	F21W040942N	F21T040942N
	48"	F21M040948N	F21W040948N	F21T040948N
	54"	F21M040954N		
	60"	F21M040960N		
	66"	F21M040966N		
12" DEEP SHELF 12" Deep Bracket	24"	F21M051224N	F21W051224N	F21T051224N
	30"	F21M051230N	F21W051230N	F21T051230N
	36"	F21M051236N	F21W051236N	F21T051236N
	42"	F21M051242N	F21W051242N	F21T051242N
	48"	F21M051248N	F21W051248N	F21T051248N
	54"	F21M051254N		
	60"	F21M051260N		
	66"	F21M051266N		
18" DEEP SHELF 18" Deep Bracket	24"	F21M061824N	F21W061824N	F21T061824N
	30"	F21M061830N	F21W061830N	F21T061830N
	36"	F21M061836N	F21W061836N	F21T061836N
	42"	F21M061842N	F21W061842N	F21T061842N
	48"	F21M061848N	F21W061848N	F21T061848N
	54"	F21M061854N		
	60"	F21M061860N		
24" DEEP SHELF 24" Deep Bracket	24"	F21M072424N	F21W072424N	F21T072424N
	30"	F21M072430N	F21W072430N	F21T072430N
	36"	F21M072436N	F21W072436N	F21T072436N
	42"	F21M072442N	F21W072442N	F21T072442N
	48"	F21M072448N	F21W072448N	F21T072448N
	54"	F21M072454N		
	60"	F21M072460N		

To specify Veneer Plywood Species replace blank at end of part number with number that indicates material required.

Veneer Plywood Material Options

- | | |
|-------------------|----------------------|
| 1 = Red Oak | 5 = FSC® Red Oak |
| 2 = White Maple | 6 = FSC® White Maple |
| 7 = Steamed Beech | 8 = FSC® St'md Beech |

Reagent Rack & Shelving Systems

Reagent Rack Uprights are available in two heights and two depths. All are 1½" thick and available in either powder coated steel to match or complement the casework color or in type 304 stainless steel. 12" deep uprights are designed for use on Center Benches


with a 14" deep shelf, 10" deep uprights are designed for use on Wall Benches with a 12" deep shelf, but both can be used on either bench.

Installation hardware is included with each upright.


Spacing between upright should not exceed 48".

Plastic Laminate Shelves are 1" thick particle board faced on top and all four edges with black plastic laminate. A backer sheet covers the bottom.


Typical Reagent Rack Assemblies

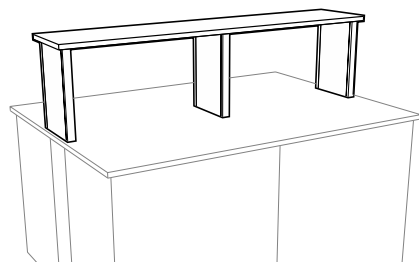

Single Tier Center Bench Reagent Rack
using 18" high by 12" deep uprights and 14" deep shelving.

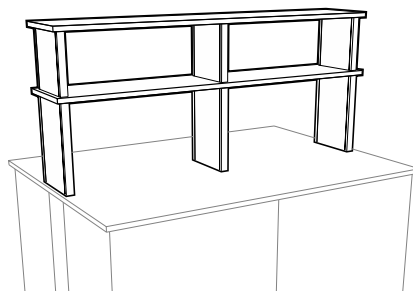

Double Tier Center Bench Reagent Rack
using 18" high by 12" deep and 12" by 12" deep uprights and 14" deep shelving.

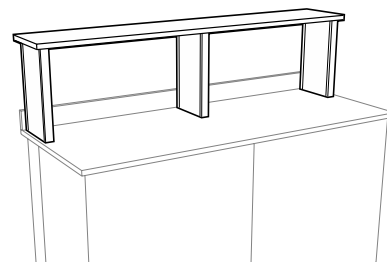

Single Tier Wall Bench Reagent Rack
using 18" high by 10" deep uprights and 12" deep shelving.

Reagent Rack Uprights


18" High x 12" Deep



18" H x 12" D x 1½" T

U03C181202N-0000
U03C181202N-0030

12" High x 12" Deep

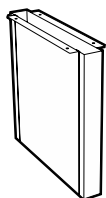

12" H x 12" D x 1½" T

U03C121202N-0000
U03C121202N-0030


18" High x 10" Deep



18" H x 10" D x 1½" T

U03C181002N-0000
U03C181002N-0030

12" High x 10" Deep

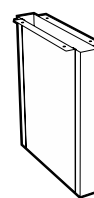


12" H x 10" D x 1½" T

U03C121002N-0000
U03C121002N-0030

P.C. Steel
Stainless Steel


Plastic Laminate Shelves

14" Deep Shelf

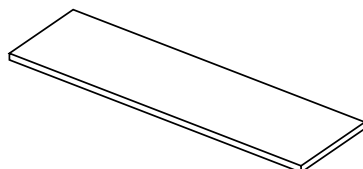

X54P011448N 14" Deep x 48" Long
X54P011472N 14" Deep x 72" Long

12" Deep Shelf

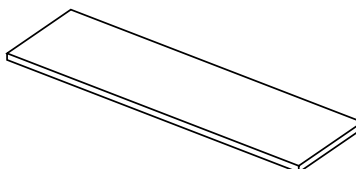

X54P011248N 12" Deep x 48" Long
X54P011272N 12" Deep x 72" Long

Reagent Rack & Shelving Systems

Fixed Wall Mounted Shelf Brackets are designed to mount directly to the wall and furnished without mounting hardware.

Adjustable Wall Mounted Shelf Brackets are designed to mount to


KF-0170-04 Shelving Uprights and include mounting bolts and nuts.

KF-0170-04 Shelving Uprights are designed to hold a maximum of three shelves.


Brackets and Uprights are finished

in Kewaunee's unique corrosion resistant powder coat finish to match or complement the casework color.


Typical Reagent Rack Assemblies

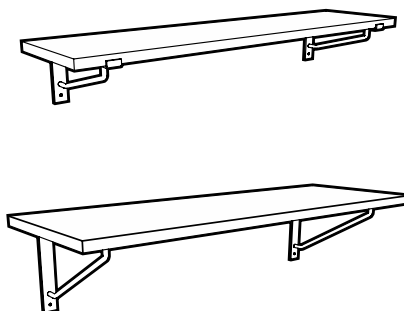

Fixed Wall Shelving

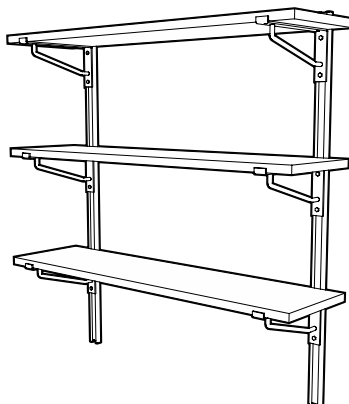


Adjustable Wall Shelving

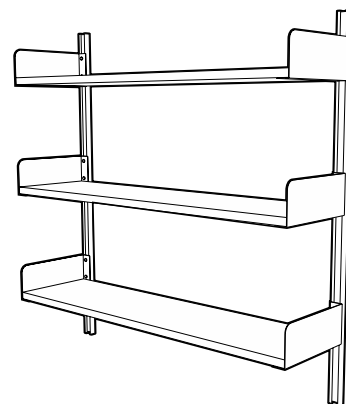


Adjustable Wall Shelving


Flat Top Shelf Bracket

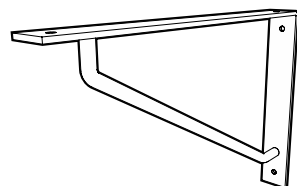

Lipped Shelf Bracket

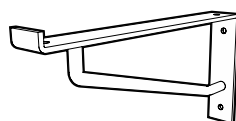

Book End Bracket

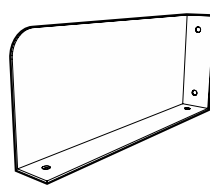

Shelving Upright



Fixed Shelf Brackets

	bracket size	shelf size
K22M081102N	11" x 8"	12"
K22M101502N	15" x 10"	16"
K22M121702N	17" x 12"	18"

part number is for single bracket

Fixed Shelf Brackets

	bracket size	shelf size
K23M061302N	12 1/8" x 6"	12"
K23M061702N	16 1/8" x 6"	16"
K23M061902N	18 1/8" x 6"	18"

part number is for single bracket


Fixed Shelf Brackets

	bracket size	shelf size
K24M061202N	12" x 6"	12"
K24M061602N	16" x 6"	16"
K24M061802N	18" x 6"	18"

part number is for pair of brackets

K21M480102N 48" L

High-strength, cold-rolled steel shelf support upright allows adjustable shelf bracket to be positioned along 48" length for maximum use of space. A pair of uprights holds a maximum of three shelf brackets.

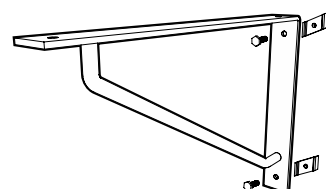

Adjustable Shelf Brackets

	bracket size	shelf size
K22M081102A	11" x 8"	12"
K22M101502A	15" x 10"	16"
K22M121702A	17" x 12"	18"

part number is for single bracket

Adjustable Shelf Brackets

	bracket size	shelf size
K23M061302A	12 1/8" x 6"	12"
K23M061702A	16 1/8" x 6"	16"
K23M061902A	18 1/8" x 6"	18"

part number is for single bracket

Adjustable Shelf Brackets

	bracket size	shelf size
K24M061202A	12" x 6"	12"
K24M061602A	16" x 6"	16"
K24M061802A	18" x 6"	18"

part number is for pair of brackets

Omega Shelving

Wall Rail Assemblies can be positioned end to end to extend assemblies. They can be used with shelves, hanger bars and other accessories.





Verticals are made of 11 gauge formed steel, with slots spaced on 1" increments.

Horizontal Rails are designed to facilitate verticals on 6" increments.


Rail Covers fit between verticals to close off exposed area.

Note: Mounting hardware not Included.

Above-Counter Wall Rail Assemblies

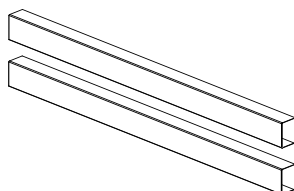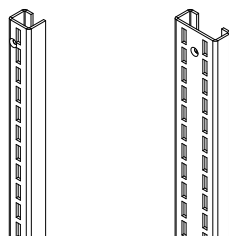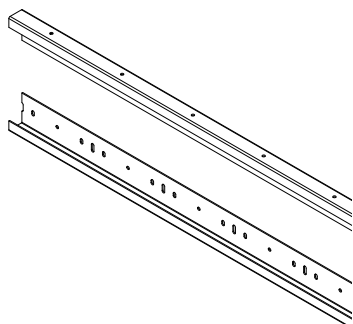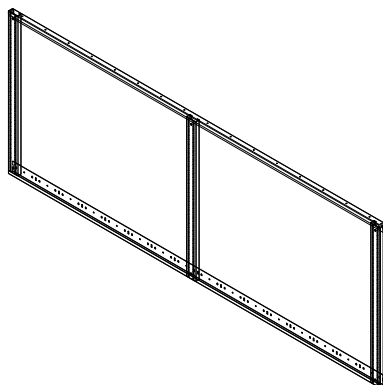

Full-Height Wall Rail Assemblies

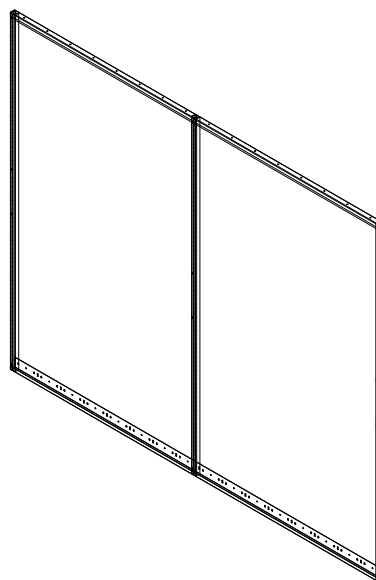

Horizontal Rails *(Set of 2)*

K11M020296N

Length

96"

Description

Upper & Lower

- Upper and lower rails are identical
- Rails can be positioned end to end to expand assemblies
- Rails can be field-cut for shorter assemblies

Full-Height Verticals

K13M800202N

Height

80"

Description

Single-slotted

K12M800202N

80"

Double-slotted

Above-Counter Verticals

K13M360202N

Height

36"

Description

Single-slotted

K12M360202N

36"

Double-slotted

Rail Covers *(Set of 2)*

K11M030224N

Length

24"

K11M030230N

30"

K11M030236N

36"

K11M030242N

42"

K11M030248N

48"

K11M030260N

60"

K11M030272N

72"

- For use with upper and lower horizontal rails


Omega Shelving

Outside Shelves are made of 16 gauge formed steel 1" thick.. Shelf lips are optional and field installed.


End brackets are made of 11 gauge formed steel, with 3 teeth spaced on 1" increments.

End brackets attach to shelf without fasteners using a unique tab and slot system.

Outside Shelves

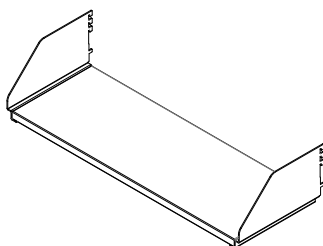


Product Number	Width	Depth	End Bracket Depth
K21M081224N	24"	12"	12"
K21M081230N	30"	12"	12"
K21M081236N	36"	12"	12"
K21M081242N	42"	12"	12"
K21M081248N	48"	12"	12"
K21M081260N	60"	12"	12"
K21M081272N	72"	12"	12"

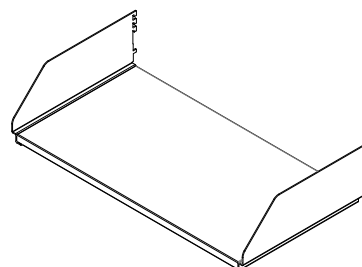


Product Number	Width	Depth	End Bracket Depth
K21M081824N	24"	18"	18"
K21M081830N	30"	18"	18"
K21M081836N	36"	18"	18"
K21M081842N	42"	18"	18"
K21M081848N	48"	18"	18"
K21M081860N	60"	18"	18"
K21M081872N	72"	18"	18"

Top Inverted Shelves

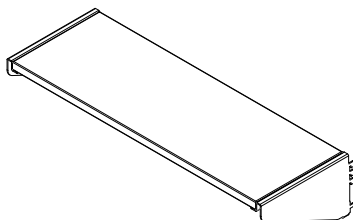

Product Number	Width	Depth	End Bracket Depth
K28M081224N	24"	12"	12"
K28M081230N	30"	12"	12"
K28M081236N	36"	12"	12"
K28M081242N	42"	12"	12"
K28M081248N	48"	12"	12"
K28M081260N	60"	12"	12"
K28M081272N	72"	12"	12"

- Designed to fit flush with the top of the core

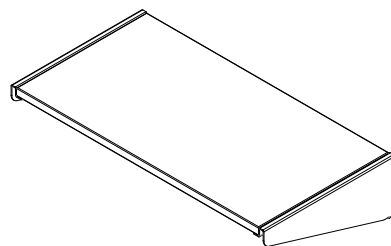


Product Number	Width	Depth	End Bracket Depth
K28M081824N	24"	18"	18"
K28M081830N	30"	18"	18"
K28M081836N	36"	18"	18"
K28M081842N	42"	18"	18"
K28M081848N	48"	18"	18"
K28M081860N	60"	18"	18"
K28M081872N	72"	18"	18"

- Designed to fit flush with the top of the core

Omega Shelving


Inverted Shelf Bracket Set

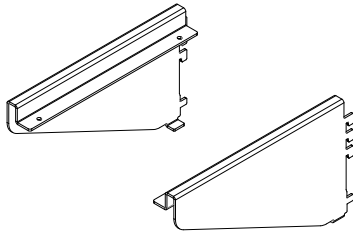

Product Number	Depth
K27M081202N	12"
K27M081802N	18"

- Designed to fit most 1" thick shelf material
- Mounting hardware not Included


Outer Shelf Bracket Set

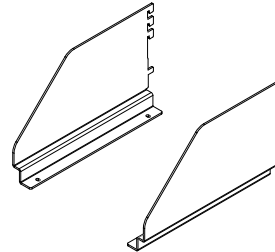

Product Number	Depth
K20M081202N	12"
K20M081802N	18"

- Designed to fit most 1" thick shelf material
- Mounting hardware not Included

Shelf Lips *(Quantity of 1)*

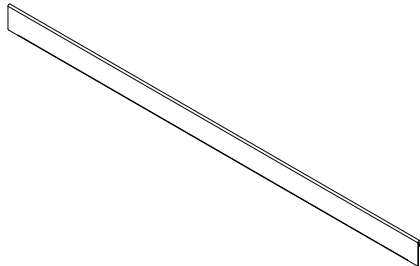

- $\frac{3}{4}$ " high steel lip
- Attaches to all shelves on the front or back of the shelf
- Mounting hardware not Included

Product Number	Width
K63M020124N	24"
K63M020130N	30"
K63M020136N	36"
K63M020142N	42"
K63M020148N	48"
K63M020160N	60"
K63M020172N	72"

Laboratory Pegboards


Kemresin Laboratory Pegboards

Kemresin Pegboards are 1" thick and available in the sizes listed below. The number of supplied pegs is listed per size. Kemresin pegboards are supplied exclusively with standard polypropylene pegs measuring ½" dia. x 6" L.

Peg Specification Note: Custom or non-standard peg lengths, styles, or attachment methods are not available and cannot be accommodated through special requests or pricing.


- At the end of the base part number, peg color will be selected by choosing **-PEGWH** = white or **-PEGBK** = black.
- After peg color, Kemresin color will be selected by choosing **BK** = Black, **GR** = Grey, **PT** = Putty, **SL** = Slate.
- Kemresin pegboards can also be ordered with Kemresin Drip Troughs and/or Finished Backs.

Example: PBDR180124-PEGBK BK is a 18" x 24" Black Kemresin Pegboard with Black Pegs

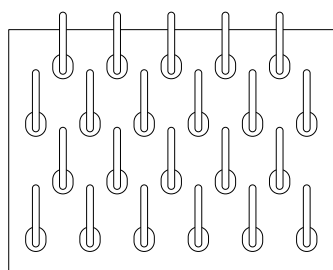

18" High	Length	Qty Pegs
PBDR180124-PEG_ _ _ _	24"L	22
24" High		
PBDR240124-PEG_ _ _ _	24"L	28
PBDR240130-PEG_ _ _ _	30"L	35
PBDR240136-PEG_ _ _ _	36"L	43
30" High		
PBDR300124-PEG_ _ _ _	24"L	39
PBDR300130-PEG_ _ _ _	30"L	49
PBDR300136-PEG_ _ _ _	36"L	60
PBDR300142-PEG_ _ _ _	42"L	70
PBDR300148-PEG_ _ _ _	48"L	81
36" High		
PBDR360124-PEG_ _ _ _	24"L	44
PBDR360130-PEG_ _ _ _	30"L	56
PBDR360136-PEG_ _ _ _	36"L	68
PBDR360142-PEG_ _ _ _	42"L	80
PBDR360148-PEG_ _ _ _	48"L	92

Type 304 Stainless Steel Laboratory Pegboards

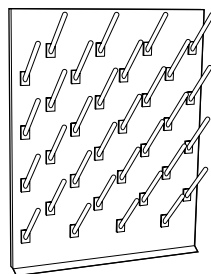


Pegboard without Drip Trough
30" High x ¾" Thick x 24" Long

X-020008 with Open Back

X-020010 with Finished Back

Peg Size: ½" diameter x 6" L
Peg Material: Polypropylene
Peg Color: White
Number of Pegs: 32

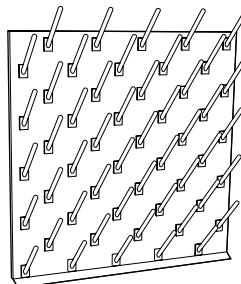


Pegboard without Drip Trough
30" High x ¾" Thick x 30" Long

X-020009 with Open Back

X-020011 with Finished Back

Peg Size: ½" diameter x 6" L
Peg Material: Polypropylene
Peg Color: White
Number of Pegs: 50

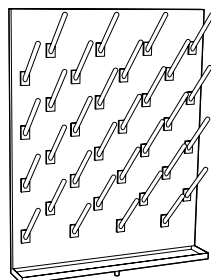


Pegboard with Drip Trough
30" High x ¾" Thick x 24" Long

X-020004 with Open Back

X-020006 with Finished Back

Peg Size: ½" diameter x 6" L
Peg Material: Polypropylene
Peg Color: White
Number of Pegs: 32

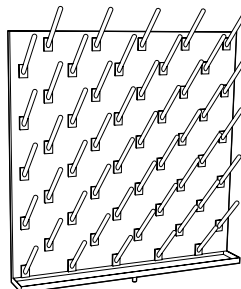

Pegboard with Drip Trough
30" High x ¾" Thick x 30" Long


X-020005 with Open Back

X-020007 with Finished Back

Peg Size: ½" diameter x 6" L
Peg Material: Polypropylene
Peg Color: White
Number of Pegs: 50

Miscellaneous Accessories


Writing Boards & Tackboards

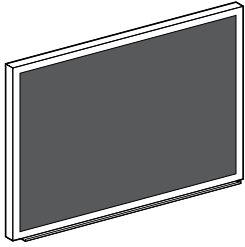

Black Chalkboard

F-5336-00	18¼"H x 24¼"L
F-5337-00	24¼"H x 36¼"L
F-5338-00	36¼"H x 48¼"L


Surface:	washable Duracite
Frame Material:	anodized aluminum
Trough Material:	anodized aluminum

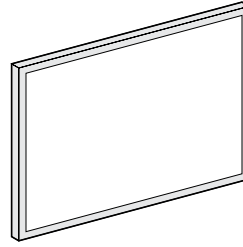

White Magnetic Markerboard

F-5342-00	18¼"H x 24¼"L
F-5343-00	24¼"H x 36¼"L
F-5344-00	36¼"H x 48¼"L


Surface:	¼" white magnetic
Frame Material:	grey mylar vinyl

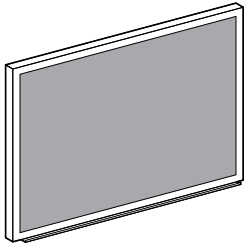

Green Chalkboard

F-5339-00	18¼"H x 24¼"L
F-5340-00	24¼"H x 36¼"L
F-5341-00	36¼"H x 48¼"L

Surface:	washable Duracite
Frame Material:	anodized aluminum
Trough Material:	anodized aluminum

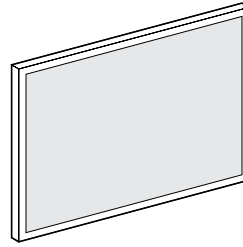


Cork Tackboard

F-5418-24	18"H x 24"L
F-5424-36	24"H x 36"L
F-5436-48	36"H x 48"L

Surface:	Resilient natural cork
Frame Material:	anodized aluminum

Kemresin Balance Table

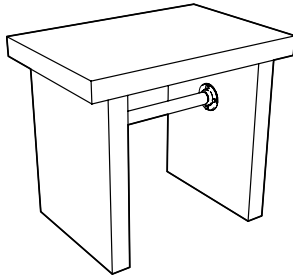


30" High x 24" Deep x 35" Long

K7R302435-BK	Black
K7R302435-GR	Grey
K7R302435-PT	Putty
K7R302435-SL	Slate

Worktop:	2" thick Kemresin
Worktop Size:	2"T x 24"D x 35"L
Pedestal Legs:	2" thick Kemresin
Total Weight:	300 pounds

Lund Key Cases

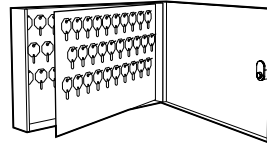

F-5332-00

Height:	11"
Depth:	2½"
Length:	17 ⁵ / ₁₆ "
Key Capacity:	90 Keys

F-5332-01

Height:	24¾"
Depth:	2½"
Length:	17 ⁵ / ₁₆ "
Key Capacity:	200 Keys

Gillis Stainless Steel Cart




F-5349-00

Overall Size:	33¾"H x 19"D x 31"L
Shelf Size:	(3) 18"D x 27"L (1" high lip on 3 sides of top & middle shelf)
Capacity:	400 pounds
Casters:	(4) Black rubber with 4" Diameter – Swivel

Miscellaneous Accessories

Pipe Space Kemstruts

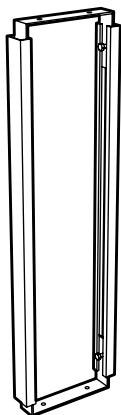

Floor mounted support to provide stability to rear overhang of fume hood tops or work surfaces, and for mounting piping and drain lines. Adjustable to match variances in height.

K12M **01** **A**

Available Heights

adjustable

20 = 19" - 23"

29 = 27" - 31"

32 = 31" - 35"

35 = 34" - 38"

Available Depths

depth

05 = 5"

06 = 6"

07 = 7"

08 = 8"

09 = 9"

10 = 10"

11 = 11"

12 = 12"

depth

13 = 13"

14 = 14"

15 = 15"

16 = 16"


17 = 17"

18 = 18"

19 = 19"

20 = 20"

Corner Kemstruts

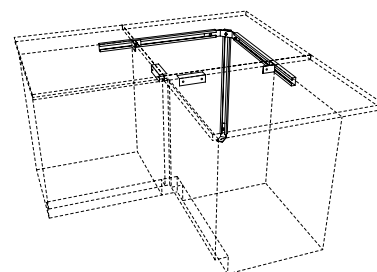

K20M293030A

K16M353030A


28" - 30" H

33" - 36" H

Floor and wall mounted support system for adding stability to corner tops in L-shaped furniture assemblies. Adjustable to match variances in height. Includes two support angles to attach to adjoining base cabinets.

NOTE: Kemstrut Corner Top Support Assemblies are not required when a corner base cabinet is used.

Kemstrut Pipe Hangers

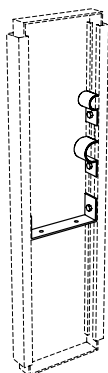


K05M030212N Drain line Support

K06M020102N Pipe Clamp
for 3/8" - 1/2" IPS

K07M020102N Pipe Clamp
for 3/4" - 1" IPS

Support Angles

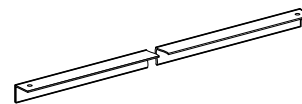

K25M020206N 6" L

K25M020212N 12" L

K25M020218N 18" L


K25M020224N 24" L

K25M020248N 48" L

K25M020296N 96" L

Formed from 16 gauge cold-rolled steel, with 1 1/2" sides.

LED Task Lighting

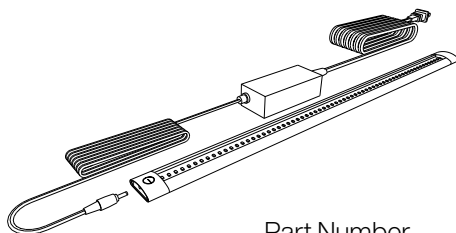


Part Number	Length
X-030001-00	300mm - 11 7/8"
X-030002-00	500mm - 19 5/8"
X-030003-00	1000mm - 39 3/8"

LED Task Lighting consists of a thin-profile, white, LED strip, 9.5mm high by 33mm deep with a touch dimmer switch. It is UL listed and furnished with 120 VAC to 24 VDC power supply with ten foot power cord. The light includes a magnetic mounting strip for mounting to steel shelving and integral tabs for mounting to wood and plastic laminate.

Miscellaneous Accessories


Rod Assembly

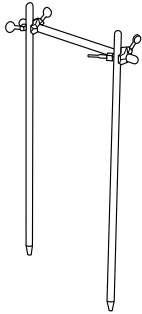

4740-00 Rod Assembly

Consists of two **0553-00** upright rods, one **0557-00** cross rod and two **0550-00** rod clamps.


Support Rods

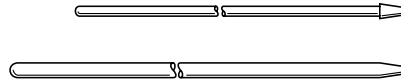

0549-18 1/2" Dia. x 18" H
0549-24 1/2" Dia. x 24" H
0549-36 1/2" Dia. x 36" H
0553-00 3/4" Dia. x 36" H

Duraluminum support rods are available in 1/2" and 3/4" diameter and have a tapered end to fit **0548-00** rod sockets.


Cross Rods

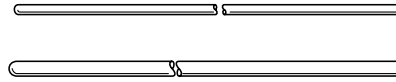

0554-00 1/2" Dia. x 39" L
0556-00 1/2" Dia. x 18" L
0557-00 3/4" Dia. x 39" L

Duraluminum cross rods are available in 1/2" and 3/4" diameter and have rounded ends


Rod Clamp

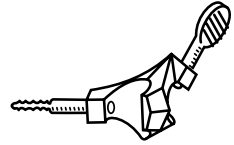

0550-00 Rod Clamp

Made of cast aluminum. Designed to clamp horizontal cross rods on vertical upright rods. Can be used for either 1/2" or 3/4" diameter rods


Rod Sockets

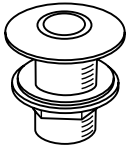

0548-00 Rod Socket

Made of chrome-plated brass. Use with 1/2" and 3/4" diameter upright rods. 1" diameter x 1 1/8" long threaded shank. Complete with locknut and washer.


Label Holder

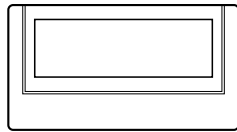

F-0260-00 Label Holder

Aluminum with adhesive back for easy installation. Holds 1 1/8" H x 2 1/2" L identification card.

Base Molding

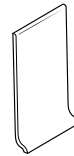


F-0862-04P 48" L

4" high black base molding. Made of pliable material to match variances in floor surface. Comes with high-strength adhesive applied to back.

Note: Specify number of 4' sections required.


Corner Guard

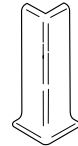

F-1368-00 Corner Guard

4" high black corner guard to finish off corner joints of base molding.

Wall Cabinet Hanger Rail

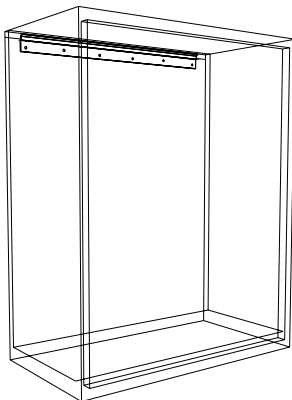

L01M020109N for 15" wall cabinet
L01M020112N for 18" wall cabinet
L01M020115N for 21" wall cabinet
L01M020118N for 24" wall cabinet
L01M020124N for 30" wall cabinet
L01M020130N for 36" wall cabinet
L01M020136N for 42" wall cabinet
L01M020142N for 48" wall cabinet

Formed of heavy gauge cold-rolled steel. Designed to capture lip on top rear of wall cabinets for easy installation.

General Plumbing Information

Recommended Supply Line Pressures

Gas.....	0.216 psi
Air (Ordinary).....	15 psi
Air (Air Powered Stirrer Operation).....	40 psi
Vacuum (Ordinary).....	18"-22" of Mercury
Vacuum (High).....	28.5" of Mercury
Hot and Cold Water.....	Local Pressure
Steam (Low Pressure).....	0-15 psi as required
Steam (High Pressure).....	15-30 psi as required

Service & Rough-in Abbreviations

A - Air
E - Electric - 120 volt AC
G - Gas
N - Nitrogen
O - Oxygen
S - Steam
V - Vacuum
W - Waste
CW - Cold Water
DW - Distilled Water
HV - Hood Vent
HW - Hot Water
WV - Waste Vent

Service Fitting Index Color Coding

When ordering service fittings, specify service required, so that the proper index button(s) will be supplied.

Air - Orange
Cold Water - Dark Green
Distilled Water - White
Gas - Dark Blue
Hot Water - Red
Hydrogen - Pink
Nitrogen - Brown
Oxygen - Light Green
Special Gases - Light Blue
Steam - Black
Vacuum - Yellow

Recommended Pipe & Drain Line Sizes

Service	Number of Outlets	Pipe Size	Pipe Material
Steam	1-8	3/8"	In accordance with governing codes
	9-16	1/2"	
	17-24	3/4"	
*Water	1-10	1/2"	In accordance with governing codes
	11-18	3/4"	
	19-24	1"	
Gas	1-12	3/8"	In accordance with governing codes
	13-18	1/2"	
	19-24	3/4"	
	25-30	1"	
Air or Vacuum	1-8	3/8"	In accordance with governing codes
	9-16	1/2"	
	17-24	3/4"	
	Fittings or Filter Pumps	Pipe Size	Pipe Material
Drain Lines	1-6 Fittings or 4 Pumps	1 1/2"	In accordance with governing codes
	7-10 Fittings or 8 Pumps	2"	


* If table is equipped with a series of filter pumps (aspiration), next larger pipe size should be used.

NOTE: For information not supplied, or if additional details are required, contact your local Kewaunee Representative.

Open Reagent Rack Horizontal Pipe Spacing


Water pipe lines should always be the lowest and a minimum of 9" above the working surface. Electric conduit should always be the uppermost service line. (Insulated pipe will require a greater spacing than indicated below.)

Pipe Diameters	Minimum Vertical Spacing Between Centerlines of 2 Service Pipe Lines
3/8" and 3/8"	2"
3/8" and 1/2"	2"
1/2" and 1/2"	2"
1/2" and 3/4"	2 1/2"
3/4" and 3/4"	2 1/2"
3/4" and 1"	3"
1" and 1"	3"

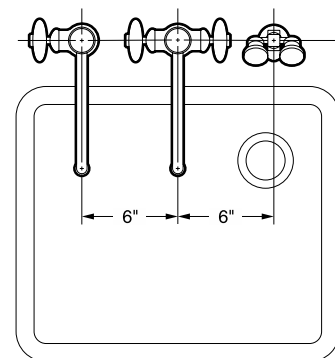

Minimum Spacing Between Fittings

These drawings show the allowable minimum spacing between service and electrical fittings on laboratory bench tops and ledges.

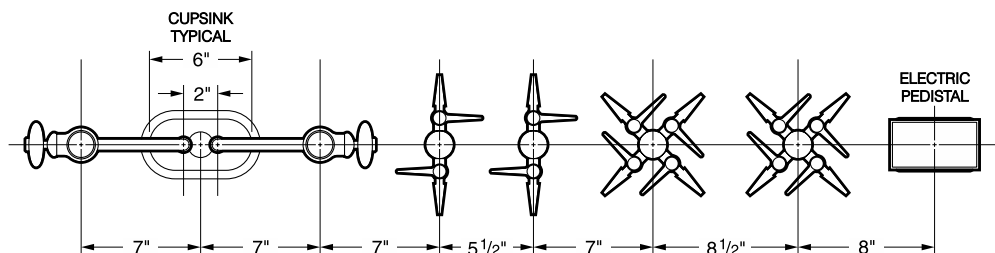

Laboratory Service Fittings


Most Laboratory Service Fittings require service selection upon ordering. This is represented by the blanks at the end of the part number. Example: **W-0242-__** in the blank spaces you would select from the available services for that fixture. The final part number would be **W-0242-0A** for Air or **W-0242-0V** for Vacuum.

For the Service Fittings on this page, selectable services are:


0A
0G
0N
0R
0V

Laboratory Ball Valves


(for use with air, vacuum or gases – not to be used for liquids, steam or oxygen)

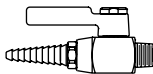


W-0242-__ Gas, Air or Vacuum Valve is chrome plated brass. Valve provides on/off control of gas, air or vacuum at pressures up to 75 PSI. $\frac{3}{8}$ " I.P.S. male inlet.

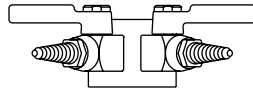


W-0286-02-__ Double Gas, Air or Vacuum Fitting (curb or panel mounted) is chrome plated brass. Consists of two-way 70° outlet panel flange and two laboratory ball valves. Complete with $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

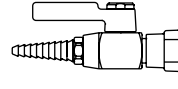


W-0287-__ Gas, Air or Vacuum Fitting (curb or panel mounted) is chrome plated brass. Consists of a single outlet panel flange and a laboratory ball valve. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

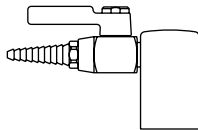


W-0260-__ Gas, Air or Vacuum Deck Mounted Fitting is chrome plated brass. Consists of a single outlet turret with a laboratory ball valve. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

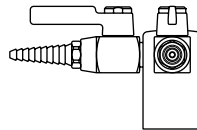

W-0263-__ Double Gas, Air or Vacuum Deck Mounted Fitting Same as No. W-0260 except consists of a 2-way 90° outlet turret and two laboratory ball valves. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

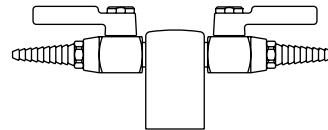


W-0266-__ Double Gas, Air or Vacuum Deck Mounted Fitting Same as No. W-0260, except consists of a 2-way 180° outlet turret and two laboratory ball valves. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

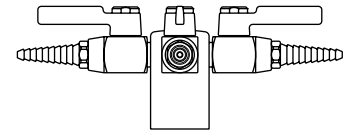


W-0272-__ Quadruple Gas, Air or Vacuum Deck Mounted Fitting Same as W-0260, except consists of a 4-way outlet turret and four valves. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

Needle Point Cock Valves


(for use with air, vacuum or gases – not to be used for liquids or steam)

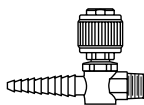

W-0243-__ Gas, Air or Vacuum Valve is chrome plated brass. Single needle point valve with integral straight serrated hose ends. Stainless steel floating cone and removable brass seat. Valve provides metering control of laboratory gases at pressures up to 75 PSI. $\frac{3}{8}$ " I.P.S. male inlet.

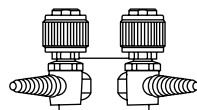

W-0400-__ Double Gas, Air or Vacuum Fitting (curb or panel mounted) with two needle point valves with integral serrated hose ends. Stainless steel floating cone and removable brass seat. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

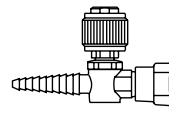


W-0401-__ Gas, Air or Vacuum Fitting (curb or panel mounted) Same as No. W-0400 except single needle point valve. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.


Laboratory Service Fittings

Needle Point Cock Valves (for use with air, vacuum or gases – not to be used for liquids or steam)


Selectable services are: **OA, OG, ON, OR, OV**

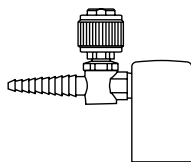


W-0259-__
Single Needle Cock Deck Mounted Fitting Turret Base is chrome plated brass. Single needle point valve with integral straight serrated hose end. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

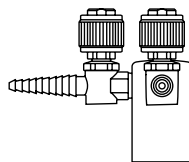

W-0402-__
Double Needle Cock Deck Mounted Fitting Turret Base is chrome plated brass. Two needle point valves with integral straight serrated hose ends at 90°. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

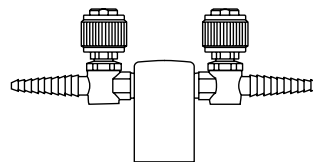


W-0403-__
Double Needle Cock Deck Mounted Fitting Turret Base is chrome plated brass. Same as No. W-0402 except integral straight serrated hose ends are at 180°. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

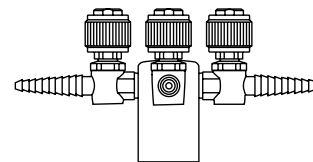

W-0404-__
Quadruple Needle Cock Deck Mounted Fitting Turret Base is chrome plated brass. Four needle point valves with integral straight serrated hose ends. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

Water Valve

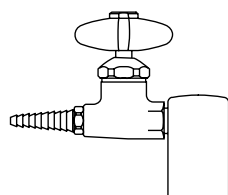


Selectable services are:


HW
CW

W-0306-__
Water Fitting (deck or panel mounted) Made of chrome plated brass. Consists of a single outlet turret and a straight water valve with removable straight serrated hose connector. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

Steam Valves

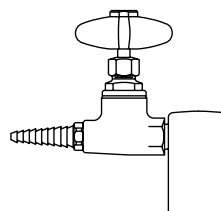


W-0300-STM
Steam Fitting (deck or panel mounted) Made of chrome plated brass. Consists of a single outlet turret and a straight steam valve. Designed for maximum steam pressure of 20 psi. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

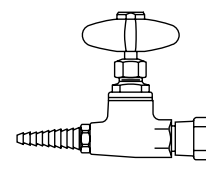

W-0301-STM
Steam Fitting (curb or panel mounted) Made of chrome plated brass. Consists of a single outlet panel flange and a straight steam valve. Designed for maximum steam pressure of 20 psi. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

Multiple Service Fitting

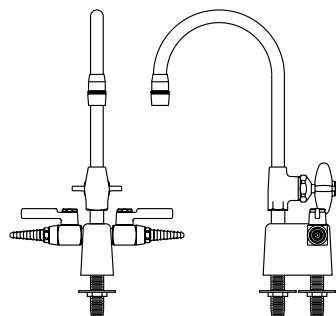


Selectable services are:

CWA
CWG
CWN
CWR
CWV

W-0273-__
Deck Mounted Combination Water/Gas Fitting Made of chrome plated brass. Water faucet has a 6" spread gooseneck with removable aerator outlet and disc type valve. Includes a two-way turret with two laboratory ball valves set at 180°. Complete with two $\frac{3}{8}$ " I.P.S. tank nipples, locknuts and washers.

W-0273-VB-__
Deck Mounted Combination Water/Gas Fitting Same as W-0273 except furnished with vacuum breaker type gooseneck and removable straight serrated hose connector.

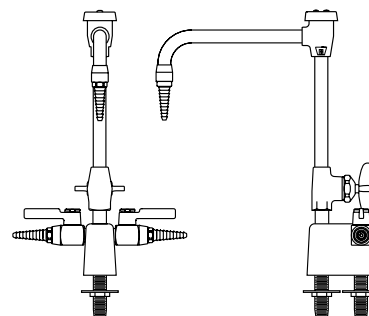


W-0274-VB-__
Deck Mounted Combination Water/Gas Fitting Made of chrome plated brass. Water faucet has a 8" spread vacuum breaker type gooseneck with removable straight serrated hose connector and disc type valve. Includes a two-way turret with two laboratory ball valves set at 180°. Complete with two $\frac{3}{8}$ " I.P.S. tank nipples, locknuts and washers.

W-0274-__
Deck Mounted Combination Water/Gas Fitting Same as W-0274 except furnished with standard gooseneck and removable aerator outlet.


Laboratory Service Fittings

Single Water Faucet


Selectable services are: **HW, CW**

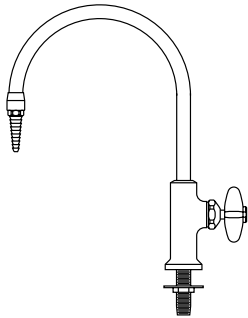


W-0337-8G-__
Deck Mounted Gooseneck Water Faucet Chrome plated brass with an 8" spread rigid gooseneck. Removable straight serrated hose connector. Valve is disc type with a renewable internal assembly. Faucet body is complete with $\frac{3}{8}$ " I.P.S. threaded mounting shank with locknut and washer.

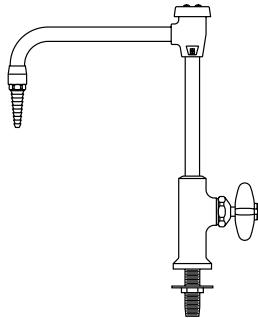


W-0337-8VB-__
Deck Mounted Vacuum Breaker Gooseneck Water Faucet Same as W-0337-8G except furnished with 8" spread rigid vacuum breaker type gooseneck.

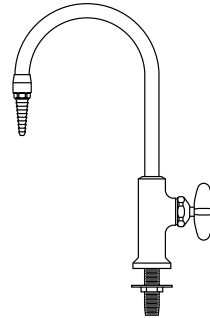


W-0337-6G-__
Deck Mounted Gooseneck Water Faucet Same as W-0337-8G except furnished with 6" spread rigid gooseneck

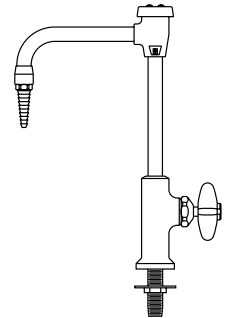


W-0337-6VB-__
Deck Mounted Vacuum Breaker Gooseneck Water Faucet Same as W-0337-8VB except furnished with 6" spread rigid vacuum breaker type gooseneck.

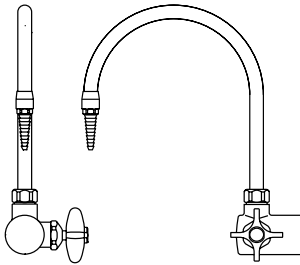


W-0350-8G-__
Right Hand Curb Mounted Gooseneck Water Faucet Made of chrome plated brass. Has an 8" spread rigid gooseneck with removable straight serrated hose connector. Valve is disc type with renewable internal assembly. Flanged inlet is complete with $\frac{1}{2}$ " I.P.S. tank nipple, locknut and washer.

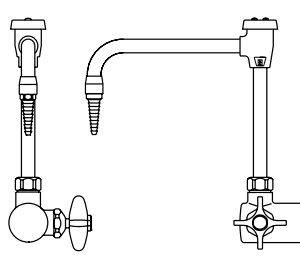


W-0350-8VB-__
Right Hand Curb Mounted Vacuum Breaker Gooseneck Water Faucet Same as W-0350-8G except supplied with vacuum breaker type gooseneck.

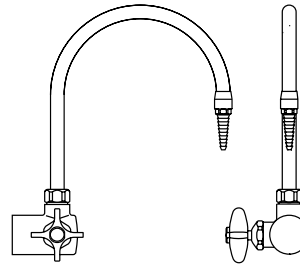

W-0351-8G-__
Left Hand Curb Mounted Gooseneck Water Faucet Same as W-0350-8G except supplied with valve handle on left.

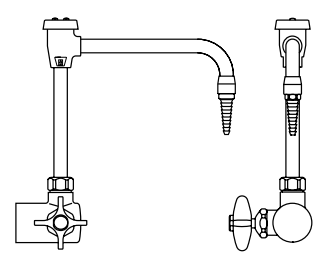


W-0351-8VB-__
Left Hand Curb Mounted Vacuum Breaker Gooseneck Water Faucet Same as W-0350-8VB except supplied with valve handle on left.

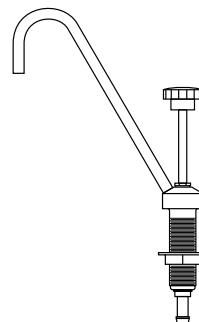


W-0405-00
Portable Pump Gooseneck Faucet for laboratory demonstration table. Made of chrome plated brass with 7 $\frac{1}{8}$ " spread, rigid, $\frac{1}{4}$ " I.P.S. gooseneck. Faucet/pump body is complete with an integral threaded mounting shank and flanged locknut. $\frac{5}{8}$ " O.D. hose connector inlet.

Laboratory Service Fittings


Mixing Faucet

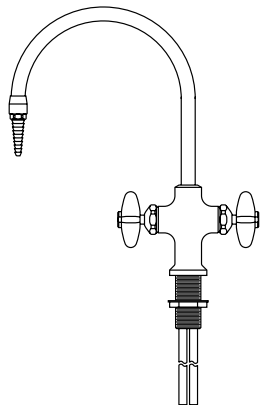


W-0340-8G-HCW
Deck Mounted Hot and Cold Water Mixing Gooseneck Faucet Made of chrome plated brass with an 8" spread rigid gooseneck. Removable straight serrated hose connector. Valves are disc type with renewable internal assemblies. Faucet body is complete with an integral threaded mounting shank and a flanged locknut. Furnished with $\frac{1}{2}$ " O.D. slip joint coupling nuts for $\frac{3}{8}$ " tubing.

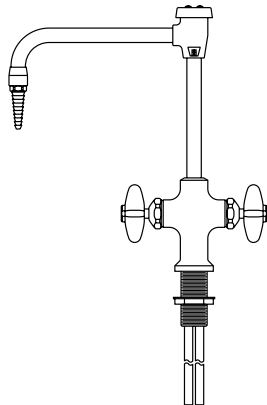

W-0340-8VB-HCW
Deck Mounted Hot and Cold Water Mixing Vacuum Breaker Gooseneck Faucet Same as W-0340-8G except supplied with vacuum breaker type gooseneck.

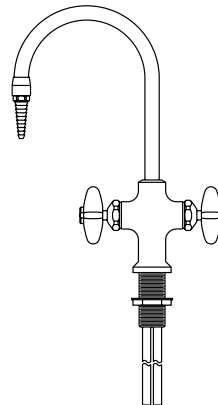


W-0340-6G-HCW
Deck Mounted Hot and Cold Water Mixing Vacuum Breaker Gooseneck Faucet Same as W-0340-8G except furnished with 6" spread rigid gooseneck

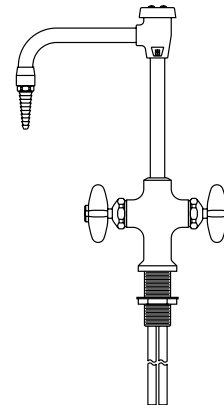

W-0340-6VB-HCW
Deck Mounted Hot and Cold Water Mixing Vacuum Breaker Gooseneck Faucet Same as W-0340-8VB except furnished with 6" spread rigid vacuum breaker type gooseneck.

Laboratory Service Fittings

Mixing Faucet

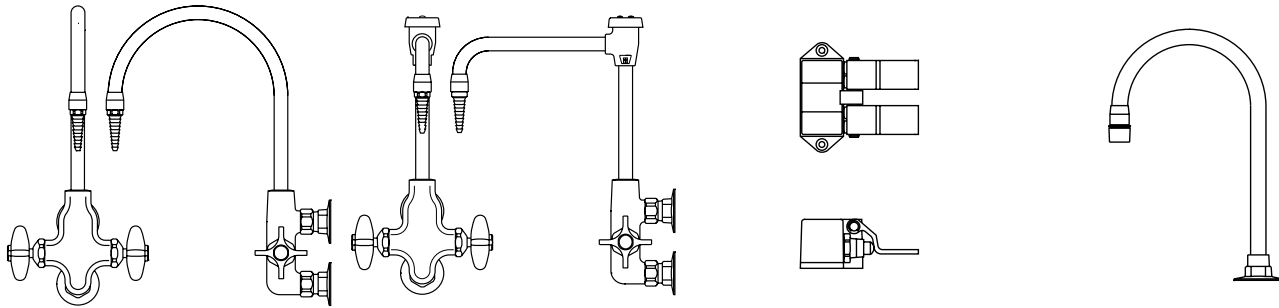

W-0361-8G-HCW

Panel Mounted Hot and Cold Water Mixing Gooseneck Faucet

Chrome plated brass. Designed for mounting on facia type reagent rack plumbing box uprights. It has an 8" spread rigid gooseneck equipped with a removable straight serrated hose connector. Valves are disc type with renewable internal assemblies. Flanged inlets on 3" centers are complete with $\frac{3}{8}$ " I.P.S. tank nipples, locknuts and washers.

W-0361-8VB-HCW

Panel Mounted Hot and Cold Water Mixing Vacuum Breaker Gooseneck Faucet

Same as W-0361-8G except supplied with vacuum breaker type gooseneck.

W-0825-HCW

Floor Mounted Foot Control Double Pedal Mixing Valve


is chrome plated brass. Has twin pedals with self-closing valves for mixing hot and cold water. $\frac{1}{2}$ " I.P.S. female thread inlets and outlet at the rear. Inlets are on 2" centers, $\frac{7}{8}$ " above floor. Outlet is on center $1\frac{1}{2}$ " above floor.

W-0833-HCW

Deck Mounted Water Gooseneck Fitting

is chrome plated brass. Consisting of a 8" gooseneck with mounting flange and W-0532-00 aerator. $\frac{3}{8}$ " I.P.S. female thread inlet complete with tank nipple, locknut and washer. **Intended for use as an outlet for W-0825-HCW.**

Faucet Accessories/Options

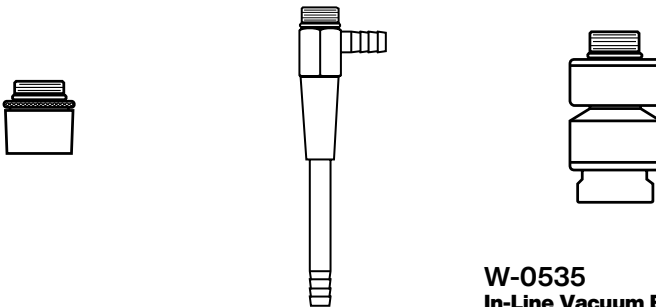

W-0532-00

Aerator Anti-splash aerator can be used in lieu of the removable straight serrated hose connectors supplied with the gooseneck water faucet. $\frac{3}{8}$ " I.P.S. male thread inlet.

W-0510-00

Filter Pump (Aspirator) 5 $\frac{1}{2}$ " high overall. Has serrated hose connector side and bottom outlets. $\frac{3}{8}$ " I.P.S. male thread inlet.


W-0535

In-Line Vacuum Breaker

Chrome plated. Mounts between the outlet end of a fitting (gooseneck, flange, etc.) and a removable serrated hose connector. $\frac{3}{8}$ " I.P.S. male thread inlet and female thread outlet. (Check plumbing codes for required location of vacuum breakers before ordering.)

Laboratory Service Fittings

Distilled Water Faucet

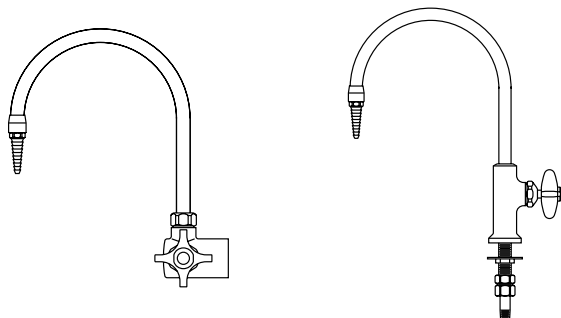

W-0367-00

Right Hand Distilled Water (Curb or Panel Mounted) Faucet


Chrome plated with a hot tin lined interior. It has a self-closing handle and a disc type valve with a renewable internal assembly. The 8" spread rigid gooseneck has a removable serrated hose end. 1/2" I.P.S. flanged female thread inlet complete with tank nipple, locknut and washer.

W-0368-00

Deck Mounted Distilled Water Gooseneck Faucet

Made of chrome plated brass. Has a hot tin lined interior. It has a self-closing handle and a disc type valve with renewable internal assembly. The 8" spread rigid gooseneck has a removable serrated hose end. Faucet body is complete with an integral threaded mounting shank and a flanged locknut. 1/2" I.P.S. male thread inlet with coupling nut and 1/4" I.P.S. tailpiece.

Ultra-Pure Deionized Water Faucet

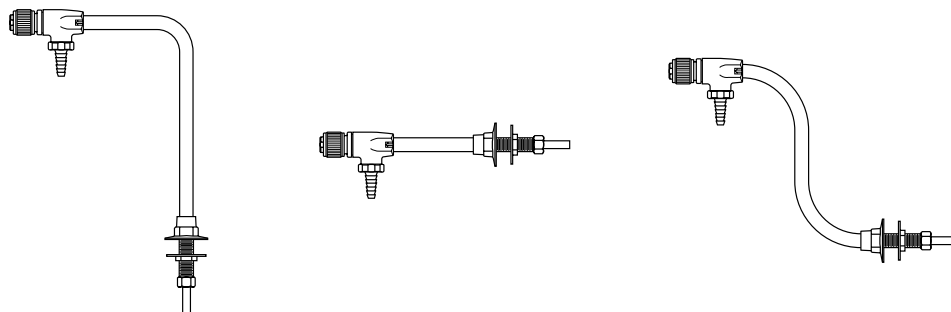

W-0407-00

Distilled, Deionized and Ultra-Pure Water Faucet (deck mounted)

made of chrome plated brass. All polypropylene internal construction with 6" spread gooseneck. Furnished with 1/2" I.P.S. shank with locknut and washer. Inlet is 3/8" O.D. compression.

W-0408-00

Distilled, Deionized and Ultra-Pure Water Faucet (panel mounted)

made of chrome plated brass. All polypropylene internal construction with 6" spread gooseneck. Furnished with 1/2" I.P.S. shank with locknut and washer. Inlet is 3/8" O.D. compression.


W-0409-00

Distilled Deionized and Ultra-Pure Water Faucet (panel mounted)


made of chrome plated brass. All polypropylene internal construction with 8" spread gooseneck. Furnished with 1/2" I.P.S. shank with locknut and washer. Inlet is 3/8" O.D. compression.

Laboratory Service Fittings


Vandal-Resistant Fittings

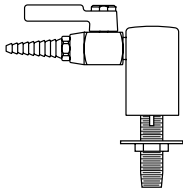

W-0410-VR-__
Vandal Resistant Deck Mounted Gas, Air or Vacuum Fitting with single vandal-resistant laboratory ball valve. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

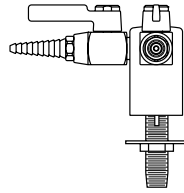


W-0411-VR-__
Vandal Resistant Deck Mounted Double Gas, Air or Vacuum Fitting with two vandal resistant laboratory ball valves with integral straight serrated hose ends at 90°. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

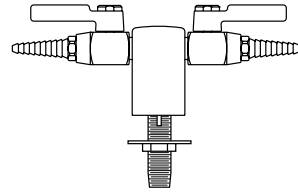


W-0412-VR-__
Vandal Resistant Deck Mounted Double Gas, Air or Vacuum Fitting Same as W-0411-VR except at 180°. Complete with a $\frac{3}{8}$ " I.P.S. tank nipple, locknut and washer.

Selectable services are:


OA
OG
ON
OR
OV

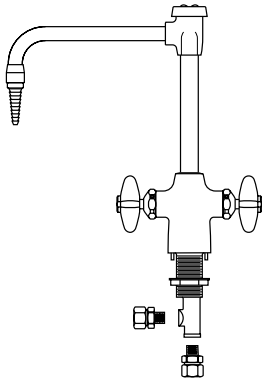

W-0413-VR-8VB
Vandal Resistant Deck Mounted Hot and Cold Mixing Faucet with 8" rigid vacuum breaker gooseneck. Valves are disc type with renewable internal assemblies. Faucet body is complete with an integral threaded mounting shank and a flanged locknut. $\frac{1}{2}$ " I.P.S. male thread inlets with coupling nuts for slip joint connections.

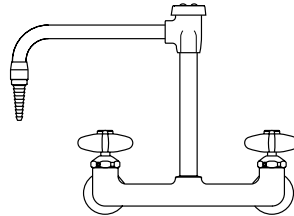


W-0414-VR-8VB
Vandal Resistant Panel Mounted Hot and Cold Mixing Faucet with 8" rigid vacuum breaker gooseneck and mounting shanks.

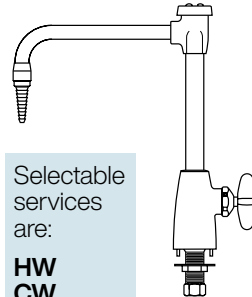


Selectable services are:
HW
CW


W-0415-VR-8VB-__
Vandal Resistant Deck Mounted Water Faucet with 8" rigid vacuum breaker gooseneck. Furnished with $\frac{3}{8}$ " IPS mounting shank with locknut and washer.

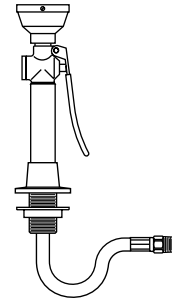

W-0417-VR
Vandal Resistant Deck Mounted Hand-Held Drench Hose with single spray type outlet heads, squeeze-lever operated valve and flexible stainless steel hose. $\frac{3}{8}$ " I.P.S. male inlet.

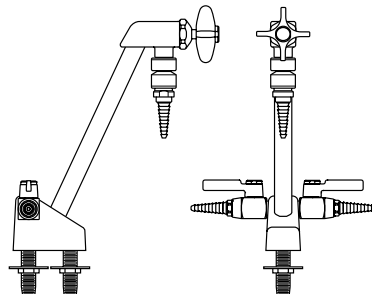

W-0416-VR-VB-__
Vandal Resistant Deck Mounted Combination Vacuum Breaker Water/Gas Fitting Made of chrome plated base. Rigid gooseneck with disc type valve and in-line vacuum breaker. Includes two-way turret with two laboratory ball valves. Complete with two $\frac{3}{8}$ " I.P.S. tank nipples, locknuts and washers.

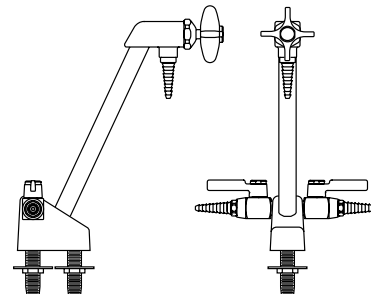


W-0416-VR-__
Vandal Resistant Deck Mounted Combination Water/Gas Fitting Same as W-0416-VR-VB except furnished without in-line vacuum breakers.

Selectable services are:

CWA
CWG
CWN
CWR
CWV

Laboratory Service Fittings


Emergency Eye Wash/Shower Equipment

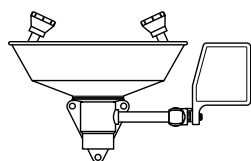

W-0927-00

Panel Mounted Eye Wash


Fixture Consists of stainless steel bowl, two Gentle Spray outlet heads, a "PUSH TO OPERATE" ball valve and mounting bracket. $\frac{1}{2}$ " I.P.S. female thread inlet and $\frac{1}{4}$ " I.P.S. female thread outlet.

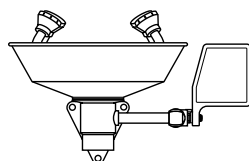

W-0933-00

Eye/Face Wash (wall mounted)


the same as No. W-0927-00, except with four Gentle Spray outlet heads for complete coverage of eyes and face.

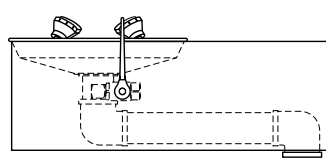

W-0930-00


BARRIER-FREE (ADA) Eye/Face Wash (wall mounted)

with four Gentle Spray outlet heads, stainless steel skirt and bowl and stay-open ball valve. $\frac{1}{2}$ " I.P.S. female thread inlet and $\frac{1}{4}$ " I.P.S. female thread outlet.

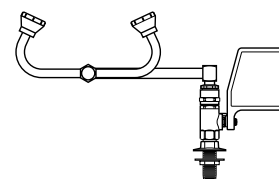


W-0931-00

Eye Wash (deck mounted) with two Gentle Spray outlet heads, stay-open ball valve and 90° swivel feature (right hand mounting).

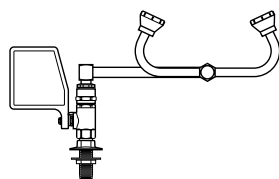

W-0932-00


Eye Wash (deck mounted) the same as No. W-0931-00 but (left hand mounting).

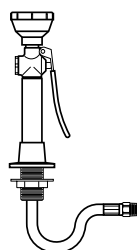

W-0923-00

Deck Mounted Hose Spray Fitting


Consists of 8 ft. reinforced PVC hose with a squeeze-lever operated valve and spray-type outlet head with a protective rubber collar. Flange, shank and locknut for attaching to table top are included. $\frac{3}{8}$ " I.P.S. male thread swivel inlet.

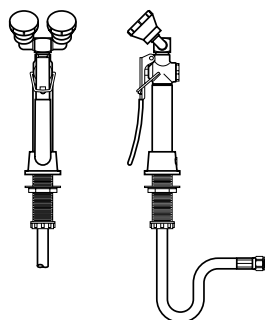

W-0934-00


Deck Mounted Dual Purpose Eye Wash/Drench Hose

with two Gentle Spray outlet heads angled at 45° and squeeze valve with locking clip. Flange, shank and locknut for attaching to table top are included. $\frac{3}{8}$ " I.P.S. male thread swivel inlet.

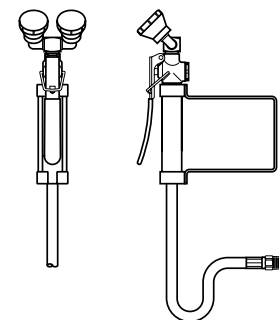

W-0935-00


Wall Mounted Dual Purpose Eye Wash/Drench Hose the same as No. W-0934-00 except wall mounted.

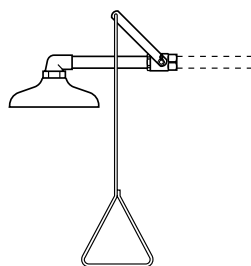

W-0374-01

Emergency Shower Head Assembly (horizontally mounted)


Includes 10" diameter orange ABS plastic shower head. Operated by a rigid pull-rod and easy to grip handle. Complete with an elbow, connecting pipe and ball valve. 1" I.P.S. female thread inlet.

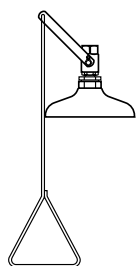

W-0928-00

Emergency Shower Head Assembly (vertically mounted)

The same as No. W-0374-01 except vertically mounted.

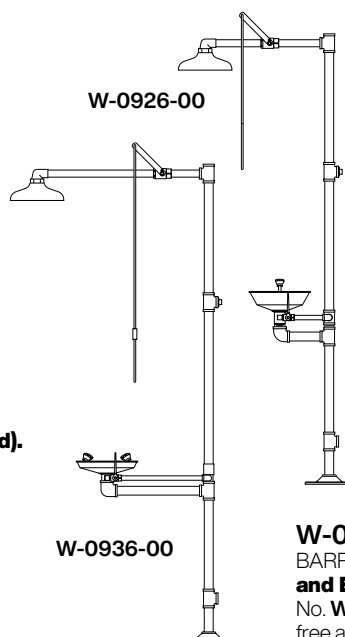

W-0926-00

Drench Shower and Eye Wash Combination

A necessity for chemical stockrooms and other areas handling a quantity of corrosive chemicals. Free-standing 10" dia. drench shower made of orange ABS plastic. Has a 1" ball valve operated by a rigid pull-rod and easy to grip handle. Eye wash fixture with stainless steel bowl, two Gentle Spray heads and a "PUSH TO OPERATE" ball valve. Eye wash and shower operate independently. Overall height is 7'6". $\frac{1}{4}$ " I.P.S. female thread inlet and outlet.

W-0936-00

BARRIER-FREE (ADA) Drench Shower and Eye Wash Combination

Similar to No. W-0926-00 except this unit is barrier-free and eye wash furnished with four Gentle Spray heads.

Cupsinks & Drain Accessories

Cupsinks



SCR0490N-BE Black Epoxy
SCR0490N-GE Grey Epoxy
SCR0490N-PE Putty Epoxy
SCR0490N-SE Slate Epoxy
 Complete with removable strainer.
 6"x 3" inside dimension. 1/2" I.P.S.
 male straight thread outlet.
 (Overall height is 7 1/2")



0491-BP Black Poly
0491-GP Grey Poly
 Complete with removable
 strainer. 6"x 3" inside dimen-
 sion. 1/2" I.P.S. male straight
 thread outlet.
 (Overall height is 8 1/2")



0492-BP Black Poly
0492-GP Grey Poly
 Complete with removable
 strainer. 9"x 3" inside dimen-
 sion. 1/2" I.P.S. male straight
 thread outlet.
 (Overall height is 7")




SCR0493N-BE Black Epoxy
SCR0493N-GE Grey Epoxy
SCR0493N-PE Putty Epoxy
SCR0493N-SE Slate Epoxy
 Complete with removable type
 316 stainless steel wire screen.
 13 3/8" x 4 1/8" x 5 1/4" I.D. 1/2" I.P.S. male
 straight thread outlet.
 (Overall height is 8")

Stainless Steel Cupsinks




0950-0F Stainless Steel
 Made of type 304 stainless
 steel and has integral cross
 bars. 5" I.D. 1/2" I.P.S. male
 straight thread outlet.

Drain Accessories

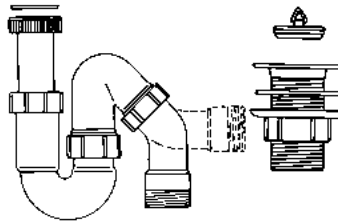

0464-00 Polyethylene
 Trap and Sink Outlet
 Adjustable riser (tailpiece)
 can be cut to the desired
 length. Sink outlet is complete
 with plug strainer, gasket and
 locknut. Has a 1/2" I.P.S. male
 straight thread outlet. Run-
 ning type trap can be used
 for either "P" or "S" type drain
 connections. 1/2" I.P.S. male
 straight thread outlet.

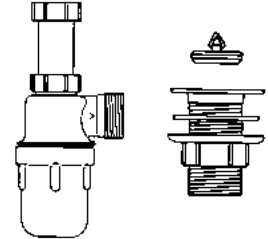

0469-00 Polyethylene
 Trap and Sink Outlet
 Adjustable riser (tailpiece)
 with 4 1/2" adjustment that
 does not require cutting. The
 sink outlet is complete with
 plug, strainer, gasket and
 locknut. Has a 1/2" I.P.S. male
 straight thread outlet. Trap
 is the bottle "P" type, with
 removable bottom for clea-
 nout. 1/2" I.P.S. male straight
 thread outlet.

Sink Outlet




0482-BP Black Poly
0482-GP Grey Poly
 Complete with locknut. 1/2"
 I.P.S. male straight thread out-
 let. For use with epoxy resin
 tub sinks and troughs.

Sink Overflow



0494BP-04 4" High
0494BP-06 6" High
0494BP-08 8" High
0494BP-10 10" High
 Sink Overflows are designed
 to be used with **0482-BP**
 sink outlets. 1/2" diameter.

St. Steel Sink Drain

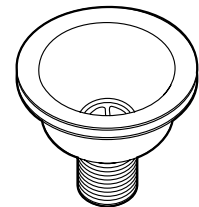

0995-00 Stainless Steel
 Crumb Cup Strainer
 Complete with gaskets, lock-
 nut and basket strainer. 1/2"
 male straight thread outlet.

Kemresin Tub Sinks & Drain Troughs

Kemresin® Tub Sinks and Drain Troughs are made from a combination of modified epoxy resins, upgraded for maximum properties with carefully selected curing resins and inert fillers to produce a thermosetting material. They are formed in metal molds to provide one-piece units with coved corners and bottoms pitched to the drain outlet. After being removed from the molds, they are cured at high temperatures for

several hours to complete the chemical reaction, and ensure maximum chemical, impact and thermal shock resistance. Drain outlets are sized for 1½" sink outlets unless otherwise specified.

Kemresin Drain Troughs are molded with no longitudinal joints. A slope of ⅛" per foot assures complete drainage, and a ⅞" cove radius in the bottom corners is provided for ease of cleaning. Five trough sections, or any part of

them, are available to be cemented together to provide Open End Troughs up to 21' long, Closed End Troughs with an end drain outlet up to 21'-8" long and Closed End Troughs with a center drain outlet up to 42'-8" long. Section TR-4 is provided with a center drain recess for a standard 1½" sink outlet.

Epoxy Resin Drain Troughs

Epoxy Resin Tub Sinks





Closed End Trough

Open End Trough

Double Tub Sink

Single Tub Sink

Center Bench
End Sink
with Trough Cutout

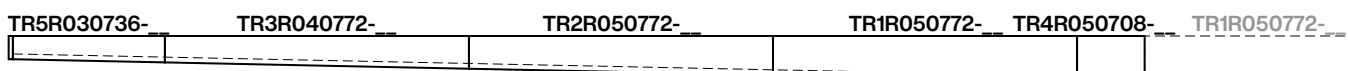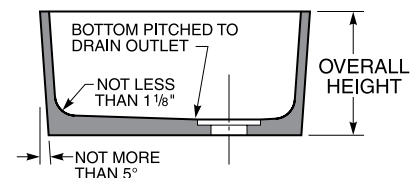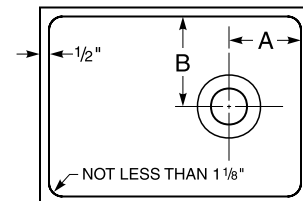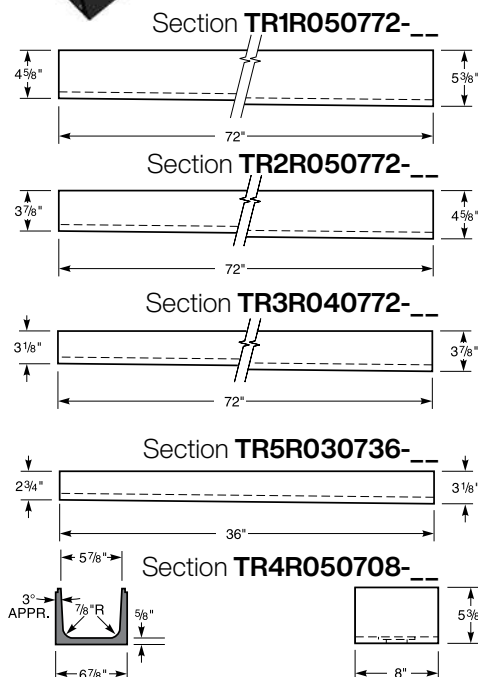

Kemresin Tub Sinks & Drain Troughs

Kemresin Colors

Kemresin® Tub Sinks and Drain Troughs are available in three standard colors: Black, Grey, and Slate.

To specify color, add color code to the end of the part number.

Example: SUR1005N-BK

Color Codes

BK = Black

GR = Grey

PT = Putty

SL = Slate

Kemresin Tub Sinks

Catalog Number	Inside Dimensions (at Top of Sink)			Drain Location		Overall Sink Height
	Length	Width	Height	A	B	
SUR1000N-__	14"	10"	6"	4"	5"	6 ⁷ / ₈ "
SUR1001N-__	16"	8"	7"	4"	4"	8"
SUR1003N-__	18"	15"	11 ³ / ₁₆ "	4"	4"	12 ⁷ / ₈ "
(1) SUR1003D-__	37"	15"	11"	4 ¹ / ₂ "	4 ¹ / ₂ "	12"
SUR1005N-__	25"	15"	10"	4"	4"	11 ¹ / ₄ "
SUR1006N-__	18"	15"	8"	4"	4"	9"
(2) SUR1006D-__	36 ³ / ₄ "	15"	8 ¹ / ₄ "	4 ¹ / ₂ "	4 ¹ / ₂ "	9"
SUR1007N-__	24"	16"	8"	4"	8"	9 ¹ / ₈ "
SUR1008N-__	16"	12"	8"	4"	6"	9 ¹ / ₈ "
SUR1011N-__	18"	6 ¹ / ₂ "	6"	9"	3 ¹ / ₄ "	7"
(3) SUR1016N-__	23 ³ / ₈ "	13 ³ / ₄ "	13 ³ / ₄ "	11 ¹ / ₈ "	9 ¹ / ₂ "	17 ³ / ₄ "
SUR1021N-__	16"	16"	7 ¹ / ₂ "	4"	4"	8 ³ / ₄ "
SUR1022N-__	21"	18"	10"	4 ¹ / ₂ "	9"	11"
SUR1023N-__	21 ¹ / ₂ "	15 ¹ / ₂ "	7"	4"	4"	8 ³ / ₈ "
SUR1024N-__	21 ¹ / ₂ "	15 ¹ / ₂ "	11"	4"	4"	12 ¹ / ₈ "
(4) SUR1024D-__	44"	15 ¹ / ₂ "	11"	4"	4"	12 ¹ / ₂ "
SUR1026N-__	35 ¹ / ₂ "	19 ¹ / ₂ "	10"	17 ³ / ₄ "	9 ³ / ₄ "	11 ¹ / ₈ "
SUR1027N-__	42"	20"	10"	21"	10"	11 ¹ / ₄ "
SUR1030N-__	16"	12"	8"	8"	6"	9"
SUR1031N-__	24"	16"	8"	12"	8"	9"
SUR1032N-__	29 ³ / ₄ "	14 ³ / ₄ "	6"	14 ⁷ / ₈ "	7 ³ / ₈ "	7"
SUR1034N-__	21 ¹ / ₂ "	19 ¹ / ₄ "	8"	10 ³ / ₄ "	9 ⁵ / ₈ "	9 ¹ / ₈ "
SUR1035N-__	45 ³ / ₄ "	17 ³ / ₈ "	12 ¹ / ₂ "	22 ⁷ / ₈ "	8 ³ / ₄ "	13 ³ / ₈ "
SUR1036N-__	48 ⁵ / ₈ "	24 ¹ / ₂ "	5 ³ / ₈ "	3 ³ / ₄ "	12 ¹ / ₄ "	6"
(5) SUR1038N-__	12" diameter		7 ⁵ / ₈ "	center		8 ³ / ₈ "
(6) SUR1039N-__	18" diameter		8 ¹ / ₈ "	center		8 ⁷ / ₈ "
SURSM01N-__	12"	6"	6"	6"	3"	7"
SURSM02N-__	19 ¹ / ₂ "	16 ¹ / ₂ "	8"	9 ³ / ₄ "	8 ¹ / ₄ "	9"
SURSM03N-__	24"	21"	8"	6"	4 ¹ / ₂ "	9 ¹ / ₈ "
SURSM04N-__	25 ¹ / ₂ "	16 ¹ / ₂ "	19"	4 ¹ / ₂ "	4 ¹ / ₂ "	20 ¹ / ₈ "
SURSM05N-__	28"	18"	10"	4"	4"	11 ¹ / ₈ "
SURSM06N-__	28"	18"	10"	24"	4"	11 ¹ / ₈ "
SURSM07N-__	30"	16"	10"	4"	4"	11 ¹ / ₈ "
SURSM08N-__	30"	16"	10"	26"	4"	11 ¹ / ₈ "
SURSM09N-__	30"	18"	18"	4 ¹ / ₂ "	4 ³ / ₄ "	18 ³ / ₄ "
SURSM12N-__	8"	8"	8"	4"	4"	8 ⁷ / ₈ "
SURSM13N-__	14"	10"	8"	4"	5"	9"
SURSM14N-__	23 ¹ / ₄ "	14 ¹ / ₄ "	12 ³ / ₄ "	11 ¹ / ₈ "	9 ¹ / ₄ "	13 ⁵ / ₈ "
SURSM15N-__	19 ³ / ₄ "	18"	24"	9 ⁷ / ₈ "	9"	25"
SURSM16N-__	40"	18"	18 ¹ / ₄ "	20"	9"	19 ¹ / ₄ "
SURSM17N-__	29 ³ / ₄ "	23 ¹ / ₂ "	6"	14 ⁷ / ₈ "	11 ³ / ₄ "	7"
SURSM18N-__	18"	15"	15 ¹ / ₄ "	4 ³ / ₄ "	4 ³ / ₄ "	16 ¹ / ₂ "
SURSM19N-__	35 ¹ / ₂ "	19 ¹ / ₂ "	15 ³ / ₈ "	3 ¹ / ₂ "	4 ¹ / ₂ "	16 ¹ / ₂ "
SURSM20N-__	21"	18"	15 ¹ / ₄ "	4 ¹ / ₂ "	4 ¹ / ₂ "	16 ¹ / ₂ "
SURSM21N-__	60 ¹ / ₂ "	22"	7 ¹ / ₄ "	6"	6"	8 ³ / ₈ "
SURSM22N-__	24"	16"	15 ³ / ₈ "	12"	8"	16 ¹ / ₈ "
SURSM23N-__	28"	18"	15 ¹ / ₂ "	5 ¹ / ₂ "	4"	16 ¹ / ₂ "
SURSM25N-__	21"	18"	15 ¹ / ₈ "	10 ¹ / ₂ "	9"	16 ¹ / ₈ "
SURSM26N-__	21"	11 ³ / ₈ "	15 ¹ / ₈ "	10 ¹ / ₂ "	5 ¹ / ₁₆ "	16 ¹ / ₈ "
SURSM29N-__	36"	24"	12 ¹ / ₄ "	18"	12"	13 ¹ / ₄ "
SURSM31N-__	24"	15 ³ / ₄ "	12 ¹ / ₈ "	4 ³ / ₄ "	8"	13"
SURSM33N-__	24"	18"	11"	4"	9"	11 ⁷ / ₈ "
(7) SURSM34N-__	15"	14"	10 ¹ / ₂ "	7 ¹ / ₂ "	7"	11 ¹ / ₄ "
SURSM36N-__	14"	10"	6"	7"	5"	6 ⁷ / ₈ "
(8) SURSM37N-__	24"	20"	10 ⁵ / ₈ "	6 ³ / ₄ "	10"	11 ¹ / ₄ "
SURSM39N-__	35"	6 ⁵ / ₈ "	4 ³ / ₄ "	17 ¹ / ₂ "	3 ⁵ / ₁₆ "	5 ¹ / ₂ "
SURSM41N-__	18"	15"	8 ¹ / ₄ "	9"	7 ¹ / ₂ "	9"
SURSM42N-__	25"	15"	10 ¹ / ₂ "	12 ¹ / ₂ "	7 ¹ / ₂ "	11 ¹ / ₄ "
SURSM43N-__	21"	18"	10 ³ / ₈ "	10 ¹ / ₂ "	9"	11 ¹ / ₈ "
SURSM44N-__	36"	16"	11"	12"	4"	12"
SURSM45N-__	27 ³ / ₄ "	14 ³ / ₄ "	8 ¹ / ₂ "	10 ¹ / ₂ "	7 ¹ / ₂ "	9 ³ / ₈ "
SURSM47N-__	43"	15 ¹ / ₂ "	11 ¹ / ₈ "	21 ¹ / ₂ "	7 ³ / ₄ "	12"
SURSM49N-__	36"	16"	8 ³ / ₈ "	6"	8"	9"
SURSM50N-__	24"	24"	12 ¹ / ₂ "	12"	12"	13 ¹ / ₄ "
SURSM51N-__	30"	15 ⁷ / ₈ "	16 ³ / ₈ "	6"	7 ¹⁵ / ₁₆ "	17 ¹ / ₈ "
SURSM52N-__	45 ¹ / ₈ "	11 ¹ / ₂ "	8 ¹ / ₈ "	22 ⁹ / ₁₆ "	5 ³ / ₄ "	9"
SURSM53N-__	31 ¹ / ₂ "	16"	16 ¹ / ₂ "	7 ¹ / ₈ "	8"	17 ¹ / ₄ "
SURSM55N-__	16"	16"	12 ³ / ₈ "	4 ¹ / ₄ "	4 ¹ / ₄ "	13"
SURSM56N-__	42 ¹ / ₄ "	18 ¹ / ₈ "	12 ⁵ / ₈ "	21 ¹ / ₈ "	9 ¹ / ₁₆ "	13 ³ / ₈ "

- (1) Double tub sink. Distance between centerlines of drain outlets is 19³/₈".
- (2) Double tub sink. Distance between centerlines of drain outlets is 19".
- (3) Center table end sink has 24¹/₂" long by 3³/₄" high back splash without trough cutout. Table side of sink is not tapered.
- (4) Double tub sink. Distance between centerlines of drain outlets is 22¹/₂".
- (5) Cylindrical sink -12" inside diameter at top of sink with center drain outlet.
- (6) Cylindrical sink -18" inside diameter at top of sink with center drain outlet.
- (7) Sides have 1¹/₂" taper per foot in height.
- (8) Exterior of one 21" long side is not tapered.

General Notes:

- Sinks with center drain outlets, except **SUR1011N**, can be furnished with a drain recess sized for a crumb cup strainer or garbage disposal when specified.
- Any sink may be furnished with trough cutout. Specify location and size. (Trough location must be at least 1¹/₂" from inside corner and 2" from bottom.)
- Most sinks can be joined together to form multiple tub sinks provided adjoining sides of sinks are the same size. Sinks longer than 30", wider than 24" or deeper than 12" cannot be factory joined, but they can be factory prepared for field assembly. Assemblies of three or more tub sinks should be no longer than 60" overall. For longer assemblies, one factory joint can be made with a second joint prepared for field assembly.

Kemresin Drop-in Sinks



Kemresin Drop-in Sinks are formed in metal molds to provide a lipped one-piece tub with coved corners and bottoms pitched to the drain outlet. They are made from a combination of modified epoxy resins, upgraded for maximum properties with carefully selected curing resins and inert fillers to produce a thermosetting material.

After being removed from the molds, they are cured at high temperatures for several hours to complete the chemical reaction, and ensure maximum chemical, impact and thermal shock resistance. Drain outlets are sized for 1½" sink outlets unless otherwise specified.

They are designed to be mounted in

Kemresin tops in a specially rabbeted cutout so that no support from below is required.

Drop-in Kemresin Sinks

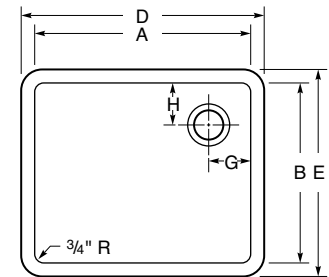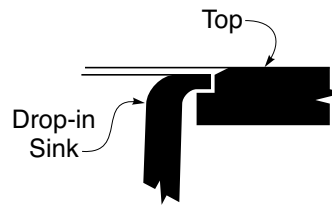


Part Number	I.D.		Depth	O.D.		Drain Location	
	A	B	C	D	E	G	H
SDR1000N-__	13 ⁷ / ₈ "	9 ⁷ / ₈ "	6"	16 ¹ / ₈ "	12 ¹ / ₈ "	7"	5"
SDR1001N-__	8"	16"	7"	10 ¹ / ₄ "	18 ¹ / ₄ "	4"	4"
SDR1003N-__	17 ¹⁵ / ₁₆ "	14 ¹⁵ / ₁₆ "	11"	20 ³ / ₁₆ "	17 ³ / ₁₆ "	3 ¹ / ₂ "	3 ¹ / ₂ "
SDR1005N-__	25"	15"	10"	27 ¹ / ₄ "	17 ¹ / ₄ "	4"	4"
SDR1006N-__	17 ⁷ / ₈ "	14 ⁷ / ₈ "	8"	20 ³ / ₈ "	17 ³ / ₈ "	3 ¹ / ₂ "	3 ¹ / ₂ "
SDR1008N-__	16"	12"	8"	18 ¹ / ₄ "	14 ¹ / ₄ "	4"	6"
SDR1011N-__	18"	6 ¹ / ₂ "	6"	20 ¹ / ₄ "	8 ³ / ₄ "	9"	3 ¹ / ₄ "
SDR1021N-__	16"	16"	7 ¹ / ₂ "	18 ¹ / ₄ "	18 ¹ / ₄ "	4"	4"
(1) SDR1038N-__	12"	diameter	7 ⁵ / ₈ "	14 ¹ / ₄ "	diameter	centered	
(2) SDR1039N-__	18"	diameter	8 ¹ / ₈ "	20 ¹ / ₄ "	diameter	centered	
SDR1040N-__	15 ³ / ₈ "	11 ³ / ₈ "	8"	17 ⁶ / ₈ "	13 ⁶ / ₈ "	4"	4"
SDR1041N-__	23 ³ / ₈ "	15 ³ / ₈ "	8.3"	25 ⁶ / ₈ "	17 ⁶ / ₈ "	4"	4"
SDR1042N-__	28"	15"	11.8"	30 ¹ / ₄ "	17 ¹ / ₄ "	5"	5"
SDR1043N-__	30"	16"	17.8"	32 ¹ / ₄ "	18 ¹ / ₄ "	5"	8"
SDR1044N-__	12"	8"	6"	14 ¹ / ₄ "	10 ¹ / ₄ "	6"	4"
SDR1045N-__	21"	17"	10"	23 ¹ / ₄ "	19 ¹ / ₄ "	4 ¹ / ₂ "	8 ¹ / ₂ "
SDR1046N-__	21 ¹ / ₂ "	15 ¹ / ₂ "	11"	23 ³ / ₄ "	17 ³ / ₄ "	4"	4"
SDR1047N-__	24"	18"	11"	26 ¹ / ₄ "	20 ¹ / ₄ "	4"	4"
(3) SDR1048N-__	17 ¹ / ₂ "	12 ³ / ₈ "	6 ¹ / ₂ "	19 ¹ / ₂ "	14 ⁵ / ₈ "	8 ³ / ₄ "	6 ³ / ₁₆ "
(4) SDR1049N-__	18"	15"	11"	20 ¹ / ₄ "	17 ¹ / ₄ "	4 ¹ / ₂ "	5"
SDR1070N-__	24"	16"	15 ¹ / ₂ "	26 ¹ / ₄ "	18 ¹ / ₄ "	12"	8"

ADA Compliant


SDR1000A-__	13 ⁷ / ₈ "	9 ⁷ / ₈ "	5"	16 ¹ / ₈ "	12 ¹ / ₈ "	4"	4"
SDR1005A-__	25"	15"	4.8"	27 ¹ / ₄ "	17 ¹ / ₄ "	4"	4"
SDR1006A-__	17 ⁷ / ₈ "	14 ⁷ / ₈ "	5"	20 ³ / ₈ "	17 ³ / ₈ "	3 ¹ / ₂ "	3 ¹ / ₂ "

- (1) Cylindrical sink, 12" inside diameter
(2) Cylindrical sink, 18" inside diameter

- (3) Triangular sink used in X92 and X94 assemblies
(4) Stepped sink, 4³/₄" deep at front

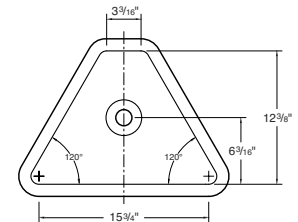

PLAN VIEW
SDR1048N-

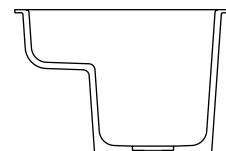

SECTION VIEW
SDR1049N-

Kemresin Colors

Kemresin® Drop-in Sinks are available in four standard colors: Black, Grey, Putty, and Slate.

To specify color, add color code to the end of the part number.

Example: SDR1005N-BK





Color Codes


BK = Black **PT** = Putty
GR = Grey **SL** = Slate

Laboratory Electrical Fixtures

Electrical Receptacles

All flush mount fixtures furnished complete with device, stainless steel face plate and electrical box. Pedestals are furnished with device, wrinkled black epoxy powder coated pedestal box, stainless steel face plate, tank nipple and locknut. All devices and boxes are UL listed. (Ship loose)

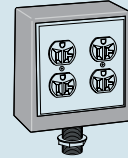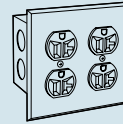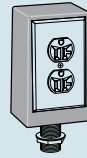
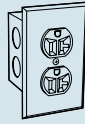




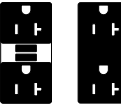
Type 2 and Type 3 pedestal boxes furnished with ¾" IPS tank nipple, Type 7 and Type 8 pedestal boxes furnished with 1" IPS tank nipple.			Flush mount single gang box	Single Wide Pedestal Type 2 – Single Faced Type 3 – Double Faced		Flush mount double gang box	Double Wide Pedestal Type 7 – Single Faced Type 8 – Double Faced	
Electrical Device		Color	Type 1	Type 2	Type 3	Type 6	Type 7	Type 8
Specification Grade 120 Volt AC, 20 amp, 3-wire polarized duplex receptacle NEMA 5-20		Black Ivory White Grey Red	S581-1K S581-1V S581-1W S581-1G S581-1R	S581-2K S581-2V S581-2W S581-2G S581-2R	S581-3K S581-3V S581-3W S581-3G S581-3R	S581-6K S581-6V S581-6W S581-6G S581-6R	S581-7K S581-7V S581-7W S581-7G S581-7R	S581-8K S581-8V S581-8W S581-8G S581-8R
Specification Grade 120 Volt AC, 20 amp, 3-wire polarized duplex receptacle NEMA 5-20		Black Ivory White Grey Red	S681-1K S681-1V S681-1W S681-1G S681-1R	S681-2K S681-2V S681-2W S681-2G S681-2R	S681-3K S681-3V S681-3W S681-3G S681-3R	S681-6K S681-6V S681-6W S681-6G S681-6R	S681-7K S681-7V S681-7W S681-7G S681-7R	S681-8K S681-8V S681-8W S681-8G S681-8R
Hospital Grade 120 Volt AC, 20 amp, 3-wire polarized duplex receptacle NEMA 5-20		Black Ivory White Grey Red	H581-1K H581-1V H581-1W H581-1G H581-1R	H581-2K H581-2V H581-2W H581-2G H581-2R	H581-3K H581-3V H581-3W H581-3G H581-3R	H581-6K H581-6V H581-6W H581-6G H581-6R	H581-7K H581-7V H581-7W H581-7G H581-7R	H581-8K H581-8V H581-8W H581-8G H581-8R
Specification Grade 240 Volt AC, 20 amp, 3-wire polarized duplex receptacle NEMA 6-20		Black Ivory White Grey Red	S582-1K S582-1V S582-1W S582-1G S582-1R	S582-2K S582-2V S582-2W S582-2G S582-2R	S582-3K S582-3V S582-3W S582-3G S582-3R	S582-6K S582-6V S582-6W S582-6G S582-6R	S582-7K S582-7V S582-7W S582-7G S582-7R	S582-8K S582-8V S582-8W S582-8G S582-8R
Ground Fault Protected 120 Volt AC, 20 amp, specification grade, duplex receptacle NEMA 5-20		Black Ivory White Grey Red	S656-1K S656-1V S656-1W S656-1G S656-1R	S656-2K S656-2V S656-2W S656-2G S656-2R	S656-3K S656-3V S656-3W S656-3G S656-3R			
Ground Fault Protected 120 Volt AC, 20 amp, duplex receptacle and 3-wire polarized duplex receptacle NEMA 5-20		Black Ivory White Grey Red				S657-6K S657-6V S657-6W S657-6G S657-6R	S657-7K S657-7V S657-7W S657-7G S657-7R	S657-8K S657-8V S657-8W S657-8G S657-8R
Ground Fault Protected 120 Volt AC, 20 amp, duplex receptacle and 3-wire polarized duplex receptacle NEMA 5-20		Black Ivory White Grey Red				S658-6K S658-6V S658-6W S658-6G S658-6R	S658-7K S658-7V S658-7W S658-7G S658-7R	S658-8K S658-8V S658-8W S658-8G S658-8R

Laboratory Electrical Fixtures

Electrical Switches and Data Ports


All flush mount fixtures furnished complete with device, stainless steel face plate and electrical box. All devices and boxes are UL listed. (Ship loose)

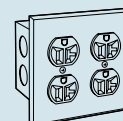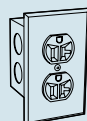

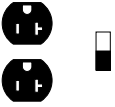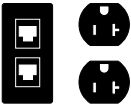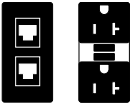
		Flush mount single gang box			Flush mount double gang box	
Electrical Device	Color		Cat-6	Cat-6a	Cat-6	Cat-6a
Specification Grade Receptacle, 20 amp, duplex, and single pole, 120/240 volt AC, 20 amp toggle switch. 	Black Ivory White Grey Red	S589-1K S589-1V S589-1W S589-1G S589-1R				
Specification Grade Receptacle, 20 amp; single pole, 20 amp. 120/240 VAC, switch. 	Black Ivory White Grey Red	S593-1K S593-1V S593-1W S593-1G S593-1R				
Fan Switch Motor-rated starter switch with pilot light, mounted in a single-gang receptacle box with faceplate. Includes 120-volt pilot light and double-pole toggle switch with thermal overload protection for up to 1 HP, single-phase, 60 Hz, 120/240-volt AC motors. (Thermal unit not included.) 	Black Ivory White Grey Red	S695-1K S695-1V S695-1W S695-1G S695-1R				
Specification Grade Single pole, 120/240 Volt AC, 20 amp toggle switch. 	Black Ivory White Grey Red	S655-1K S655-1V S655-1W S655-1G S655-1R				
RJ-45 Data Port Duplex Cat-6 Data Port 	Black Ivory White Grey Red		0720-1K-6 0720-1V-6 0720-1W-6 0720-1G-6 0720-1R-6	0720-1K-6a 0720-1V-6a 0720-1W-6a 0720-1G-6a 0720-1R-6a		
RJ-45 Data Port Duplex Cat-6 Data Port with 120 Volt AC, 20 amp, specification grade, duplex receptacle. NEMA 5-20 	Black Ivory White Grey Red				S721-6K-6 S721-6V-6 S721-6W-6 S721-6G-6 S721-6R-6	S721-6K-6a S721-6V-6a S721-6W-6a S721-6G-6a S721-6R-6a
RJ-45 Data Port Duplex Cat-6 Data Port with 120 Volt AC, 20 amp, specification grade, GFI duplex receptacle. NEMA 5-20 	Black Ivory White Grey Red				S751-6K-6 S751-6V-6 S751-6W-6 S751-6G-6 S751-6R-6	S751-6K-6a S751-6V-6a S751-6W-6a S751-6G-6a S751-6R-6a

Steel Cabinet Technical Information

Inset Style Wall Bench End Views

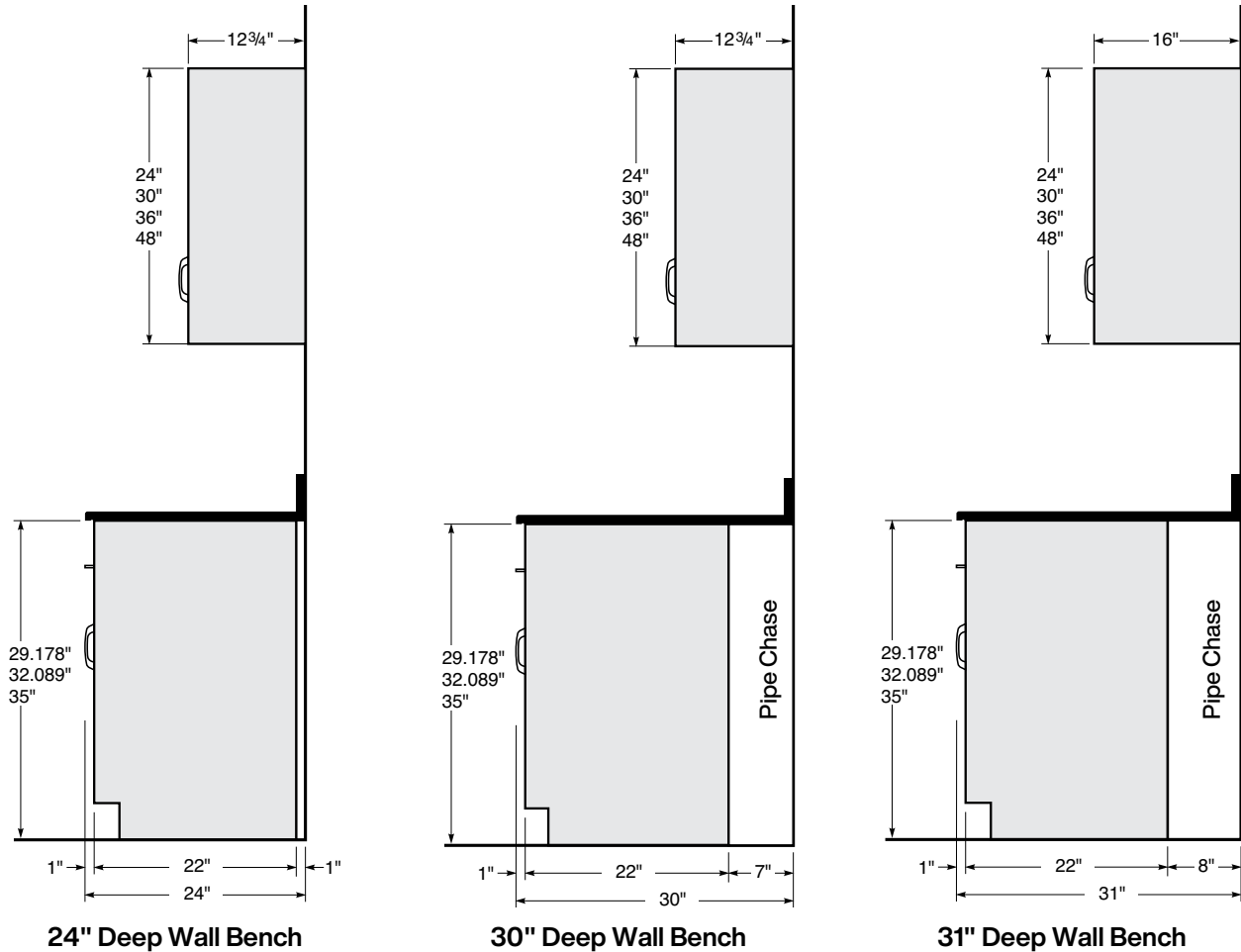


Inset Style Center Bench End Views

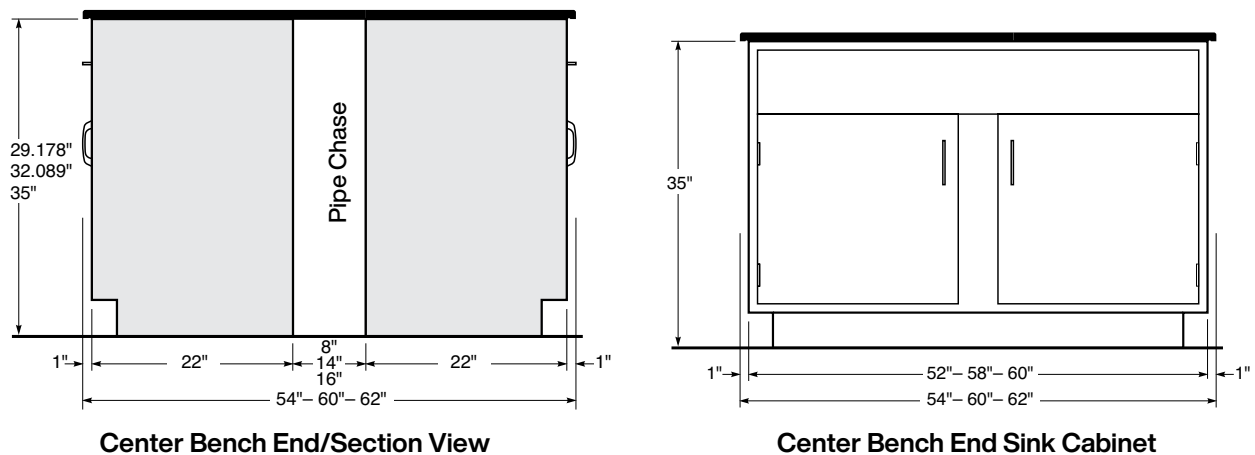

Steel Cabinet Technical Information

Inset Style Full Height Storage Cabinet End Views

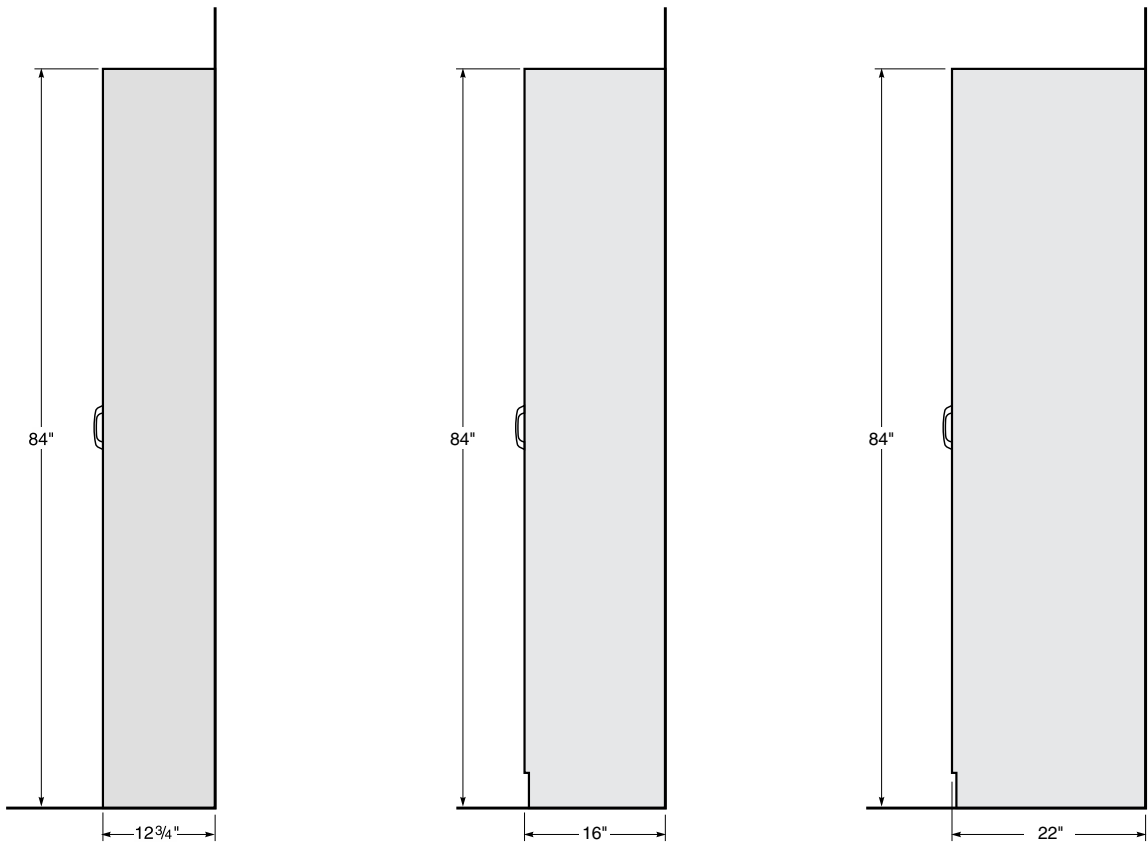


12 3/4" Deep Full Height Cabinet

16" Deep Full Height Cabinet

22" Deep Full Height Cabinet


Inset Style Center Bench Plan View

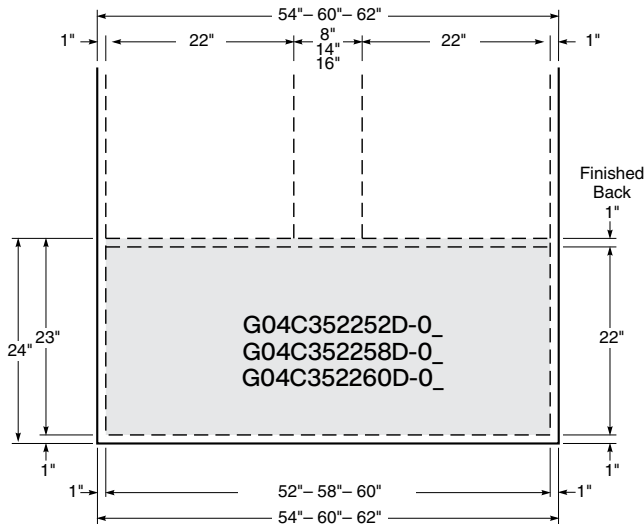


Plan View
Center Bench End Sink Cabinet

Steel Cabinet Technical Information


Full Overlay Style Wall Bench End Views

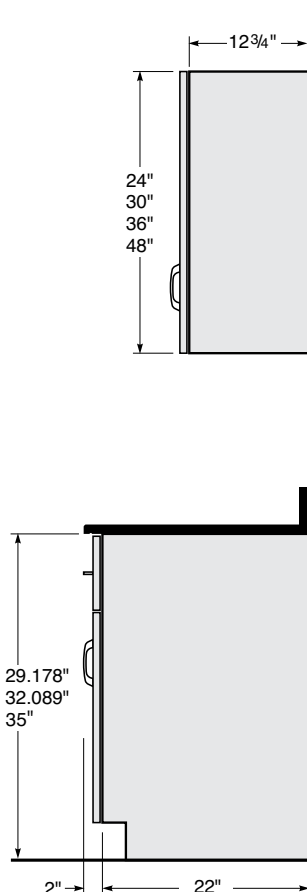

24" Deep Wall Bench

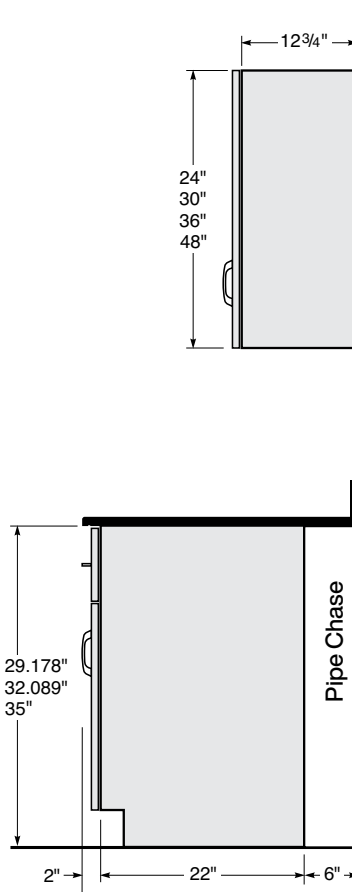


30" Deep Wall Bench

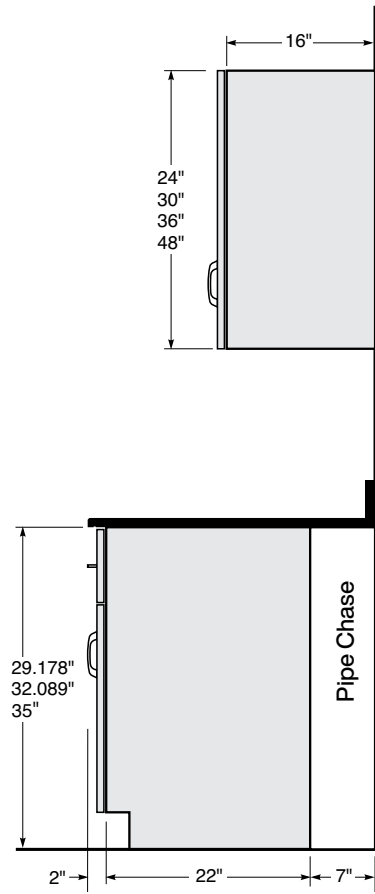


31" Deep Wall Bench

Full Overlay Style Center Bench End Views

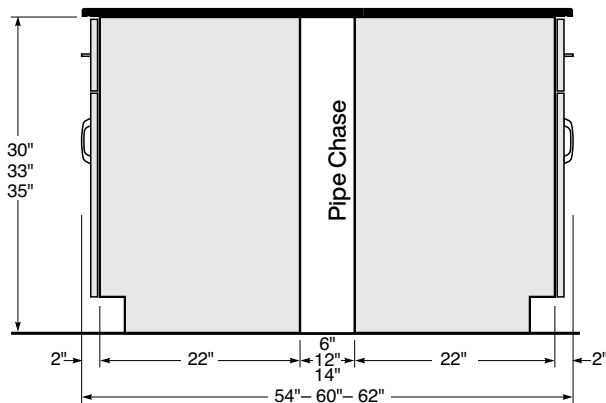


Center Bench End/Section View

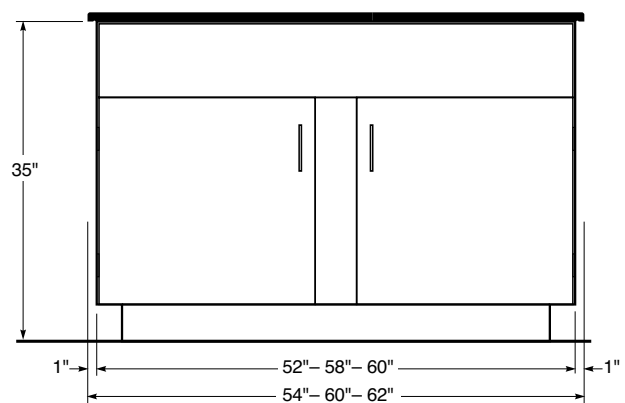

Center Bench End Sink Cabinet

Steel Cabinet Technical Information

Full Overlay Style Full Height Storage Cabinet End Views

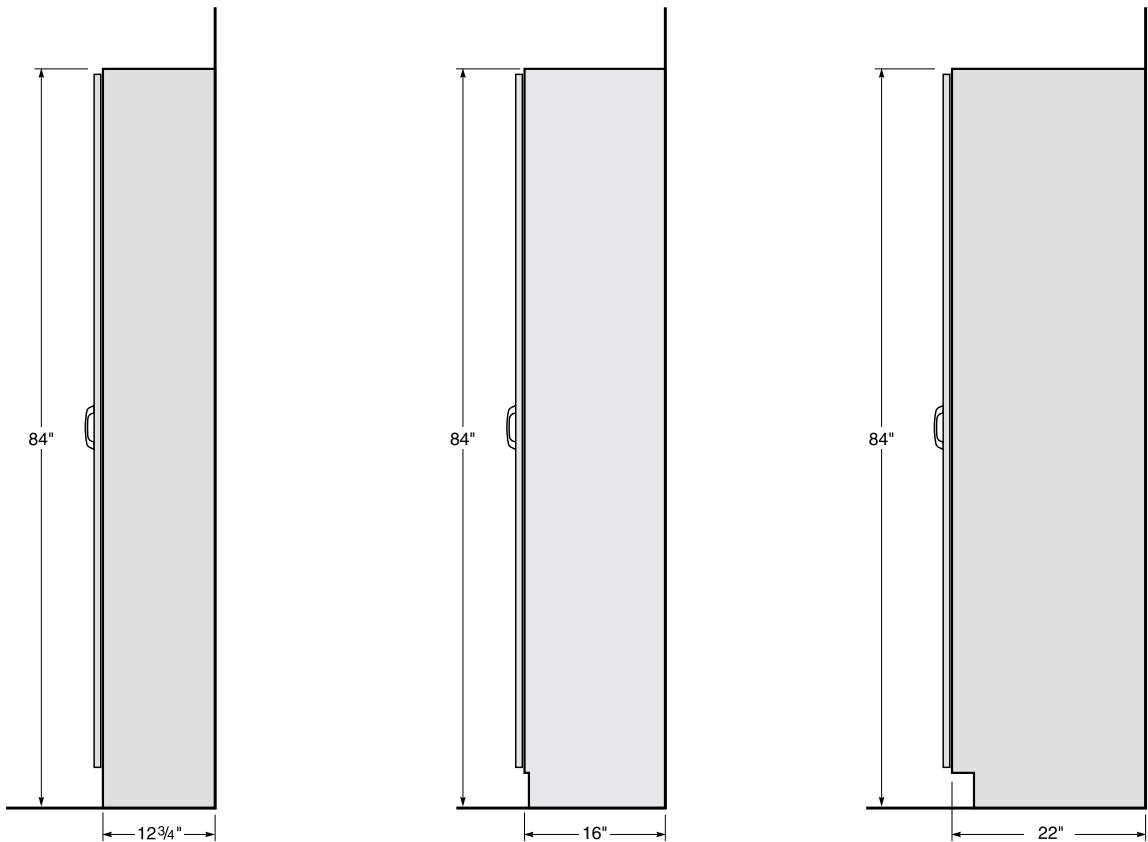


12 3/4" Deep Full Height Cabinet

16" Deep Full Height Cabinet

22" Deep Full Height Cabinet


Full Overlay Style Center Bench Plan View

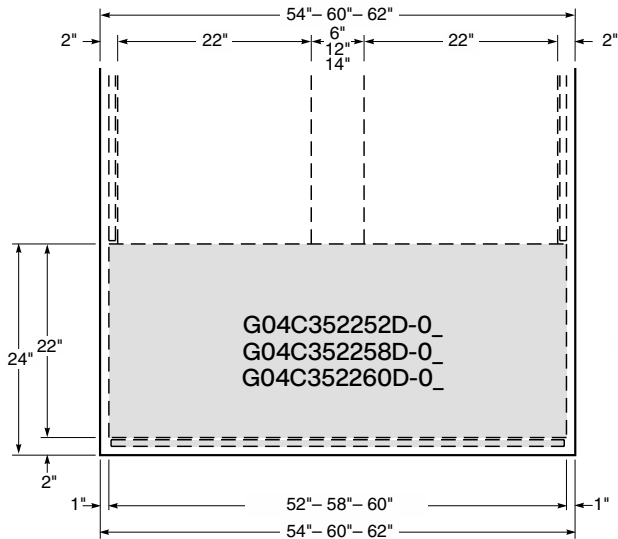

Plan View
Center Bench End Sink Cabinet

Steel Cabinet Technical Information

Flammable & Acid Safety Cabinet Vent Locations

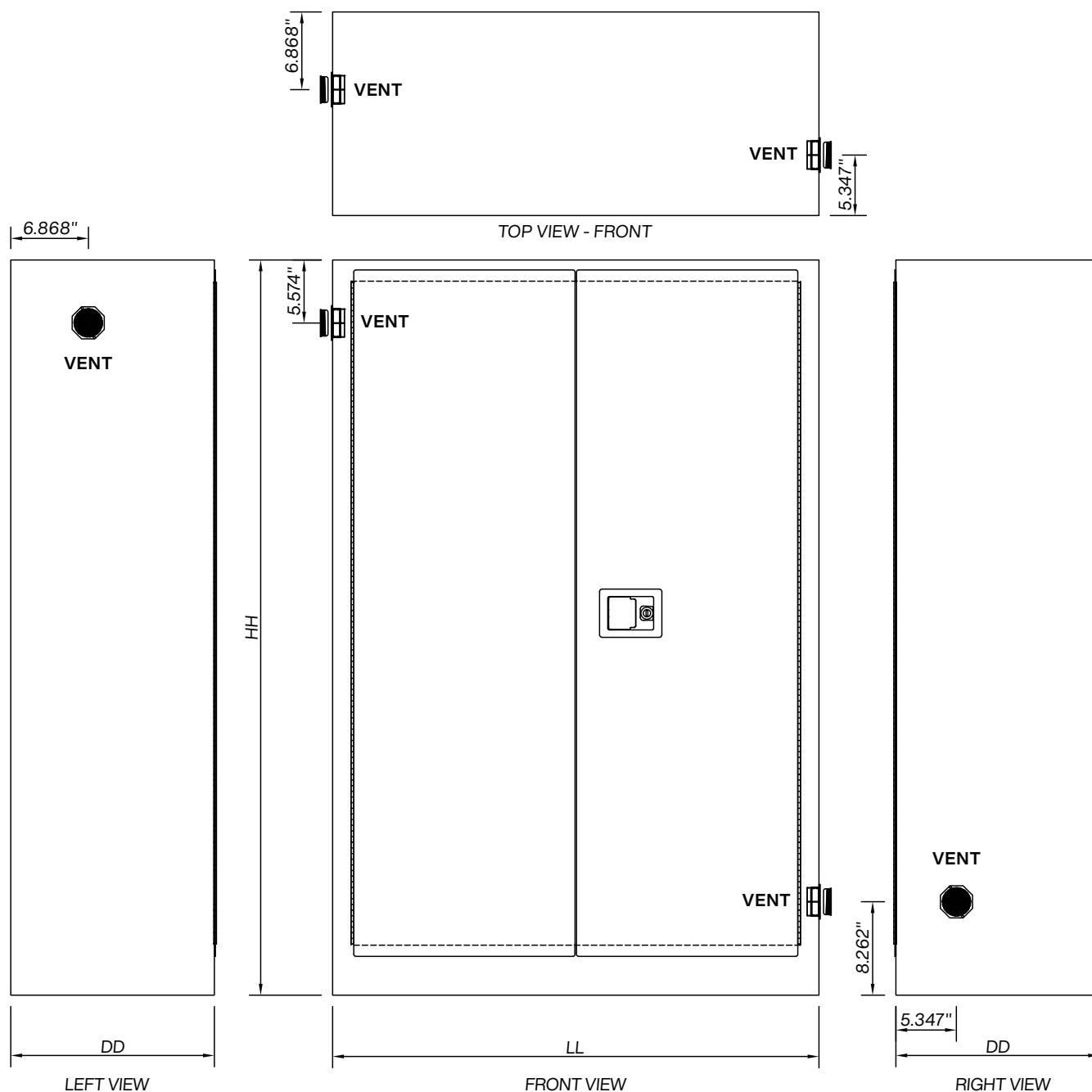

Lock Schedule Sample and Examples

Lock Schedule Examples

Keyed Different & Masterkeyed

An example of a simple lock schedule with disk tumbler locks. All locks on the project are keyed different with a single master key that can open all locks. Either disk tumbler or pin tumbler locks could be used for this arrangement.

LOCK SCHEDULE

Sheet ____ of ____

Project Name: SAMPLE PROJECT				Project/Quote Number		Date:		for Order Dept. use	
Lock Type	☒ Disk Tumbler ☐ Pin Tumbler	Number of Keys per Lock	2	Number of Master Keys per change	4	Number of Grand-Master Keys per change			

Room Number	Item Number	Cabinet (all cabinets with locks must be listed)	Number of Locks	Keyed Different	Keyed Alike	Key Group	Masterkey Group	Grand-Masterkey Group	
100	1	D30C352218N-0101J1	4	☒	☐		A		
100	1	D32C352230N-0101J1	6	☒	☐		A		
100	2	E40C352221L-0101J1	2	☒	☐		A		
105	6	E40C352221R-0101J1	4	☒	☐		A		
105	6	G00C352236D-0101J1	3	☒	☐		A		
105	6	G00C352242D-0101J1	1	☒	☐		A		

Keyed Alike in Sets

An example of another simple lock schedule with disk tumbler locks. On this project all the locks in each room are keyed alike, but each room is different from the other. A master key could be used in this type of keying arrangement, but in this example it is not. Either disk tumbler or pin tumbler locks could be used for this arrangement.

LOCK SCHEDULE

Sheet ____ of ____

Project Name: SAMPLE PROJECT				Project/Quote Number:		Date:		for Order Dept. use	
Lock Type	☒ Disk Tumbler ☐ Pin Tumbler	Number of Keys per Lock	2	Number of Master Keys per change		Number of Grand-Master Keys per change			

Room Number	Item Number	Cabinet (all cabinets with locks must be listed)	Number of Locks	Keyed Different	Keyed Alike	Key Group	Masterkey Group	Grand-Masterkey Group	
115	5	D30C352218N-0101K1	4	☐	☒	A			
115	5	D32C352230N-0101K1	4	☐	☒	A			
115	5	E40C352221L-0101K1	3	☐	☒	A			
126	16	E40C352221R-0101K1	2	☐	☒	B			
126	16	G00C352236D-0101K1	5	☐	☒	B			
126	17	G00C352242D-0101K1	1	☐	☒	B			
126	17		4	☐	☒	B			

Keyed Different – Masterkeyed & Grand-Masterkeyed

An example of a more complex lock schedule with pin tumbler locks. All locks are keyed different, with a separate master key for each room, and a grand-master key for the entire project. This keying arrangement could not be used with disk tumbler locks.

LOCK SCHEDULE

Sheet ____ of ____

Project Name: SAMPLE PROJECT				Project/Quote Number:		Date:		for Order Dept. use	
Lock Type	☐ Disk Tumbler ☒ Pin Tumbler	Number of Keys per Lock	2	Number of Master Keys per change	5	Number of Grand-Master Keys per change	5		

Room Number	Item Number	Cabinet (all cabinets with locks must be listed)	Number of Locks	Keyed Different	Keyed Alike	Key Group	Masterkey Group	Grand-Masterkey Group	
102	1	D30C352218N-0101J2	6	☒	☐		A	A	
102	1	D32C352230N-0101J2	2	☒	☐		A	A	
102	1	E40C352221L-0101J2	4	☒	☐		A	A	
105	3	E40C352221R-0101J2	2	☒	☐		B	A	
105	3	G00C352236D-0101J2	4	☒	☐		B	A	
1113	9	G00C352242D-0101J2	8	☒	☐		C	A	
1113	9		1	☒	☐		C	A	
1113	10		4	☒	☐		C	A	



Sheet _____ of _____

[illegible]

Typical Installations



Trademarks

Kewaunee, Alpha, Basic Bench, EarthLine, Enterprise, Evolution, Interceptor, Kemresin, Kewaunee Matrix, Research Collection, Signature, Sturdilite, and Supreme Air are registered trademarks of Kewaunee Scientific Corporation.

Patents:

Products in this catalog are protected by one or more of the following patents:

4,856,420	5,241,788	5,302,030	5,407,389
5,582,472	5,728,626	D388,986	DES. 346,437

Additional patents are pending.

Part numbers, illustrations, descriptions, specifications and other information contained in this catalog are based on the most current design and production information available

at the time of publication. Kewaunee Scientific Corporation reserves the right to make changes at any time without notice.

Typical Installations



We are pleased to have been selected to supply laboratory furniture for the facilities of some of the best companies, organizations and institutions in the world. We believe the repeat opportunities, over time, to participate in their construction projects is good evidence of our product quality and experience in successfully completing their laboratory projects.

Whether your laboratory projects are large or small, new construction or remodeling, Kewaunee has the qualified people and quality products to transform your visions into successful reality.

From pre-planning to final completion, experienced Kewaunee representatives are available to help you through the intricacies of laboratory design and layout, product selection, ordering, scheduling and installation.

RESEARCH COLLECTION

Steel Laboratory Furniture

Kewaunee Scientific Corporation is dedicated to manufacturing high quality products for the laboratory marketplace. We offer **Total Laboratory Solutions, Innovative Designs, and Technological Expertise** with a worldwide distribution network to provide laboratories that are truly World Class.

Laboratory Products from Kewaunee®:

ALPHA® – Adaptable Casework System
ALPHA Columns – Expandable Workstations
ALPHA Overhead Carrier Systems
BasikBench® – Simple 4-Leg Benches or Tables
Element Bench – Efficient Benches with Plumbed Posts
Enterprise® – Free-Standing Workstations
Evolution® – Column Based Workstations
Kemresin® – Epoxy Resin Counters & Sinks
Kewaunee Matrix® – Educational Stations
Research Collection® – Steel Casework
Signature® Series – Wood Casework
Supreme Air® – Fume Hoods

For our complete product offering, specifications, catalogs, brochures, case studies, and dealer list, go to kewaunee.com, or contact your local sales representative directly.



The mark of responsible forestry



Kewaunee Scientific Corporation

P.O. Box 1842 • Statesville, NC 28687-1842
704-873-7202 • kscmarketing@kewaunee.com • kewaunee.com

Kewaunee Labway India Pvt. Ltd.

No. CA-9A, 2nd Floor, Jigani Link Road
Jigani, Bangalore 562 106, India
Phone: 9180-27826725 • Fax: 9180-27826724
Web Site: kewaunee.in

Kewaunee Labway Asia Pte. Ltd.

194 Pandan Loop, #06-22
Pantech Business Hub, Singapore 128 383
Phone: 65-6773-0288 • Fax: 65-6773-2322



KEWAUNEE®

Kewaunee, Alpha, BasikBench, Enterprise, Evolution, Kemresin, Kewaunee Matrix, Research Collection, Signature, and Supreme Air are registered trademarks of Kewaunee Scientific Corporation